

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM² -> SINGLE-MOLECULE STORM

8


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

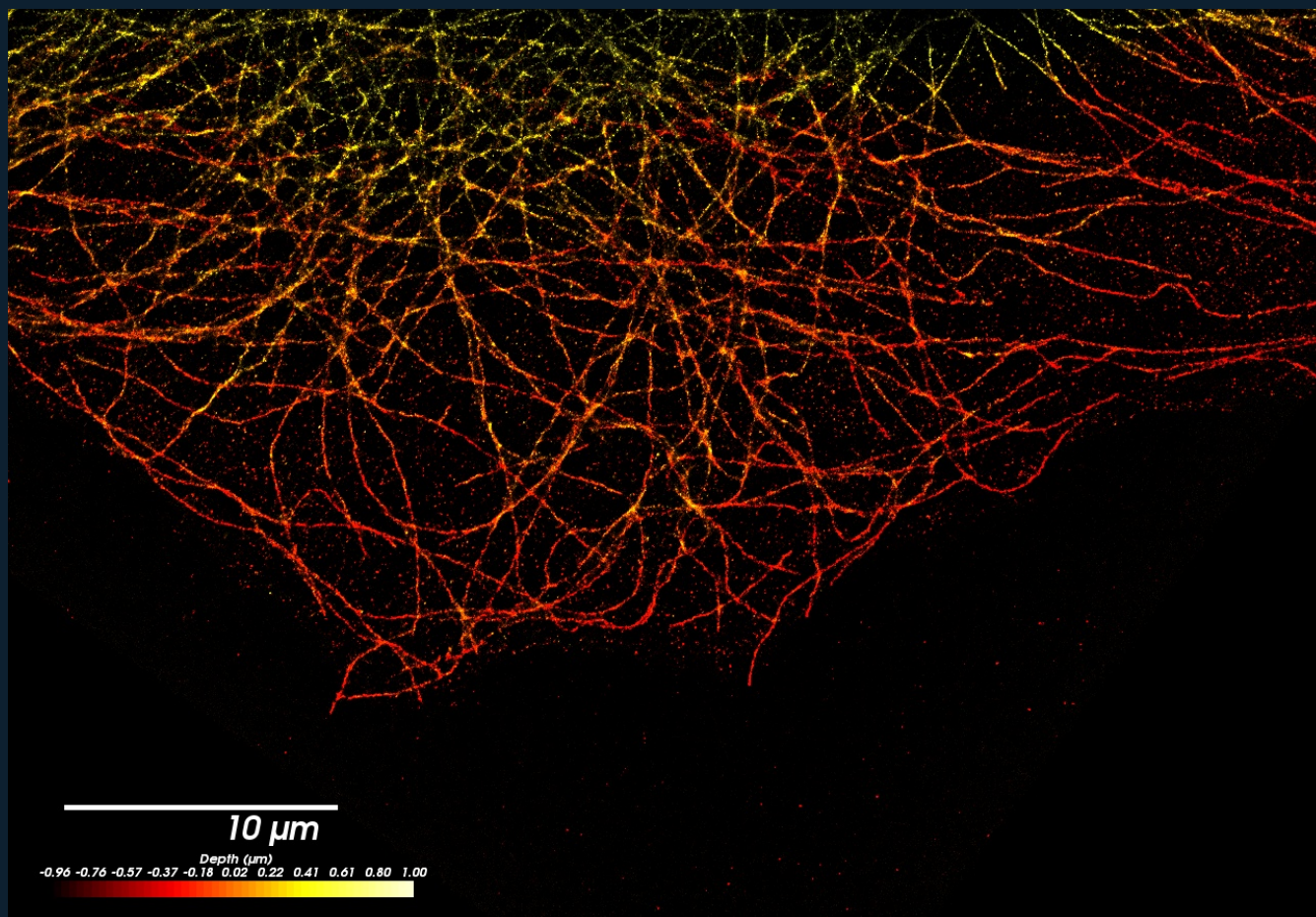

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SINGLE-MOLECULE STORM

SA-000089 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Single-molecule STORM presentation workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Single-molecule STORM
-> Control

BENNETT ASSESSMENT

The preserved SA-000089 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / FULL STEPDOWN / SOURCE-BOUNDED REVIEW

The preserved SA-000089 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

SETTINGS ASSESSMENT

The record preserves 8 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	555	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 555; 594	2; 45 Texas Red; 49 DAPI; 540 - 580; 335 - 383; 593 - 668; 420 - 470; 553
Widefield SIM — Apotome	405/DAPI; 555; 594	2; 45 Texas Red; 49 DAPI; 540 - 580; 335-383; 593 - 668; 420-470; 553
Confocal	405/DAPI; 488; 555	2; None; 561nm @0.6%; 405nm @0.3%; 553; 353; 568; 465
Airyscan	405/DAPI; 488; 555	2; None; 561nm @0.5%; 405nm @1.0%; 553; 353; 568; 465
SIM	405/DAPI; 488; 555	2; 561nm @1.80%; 561; 570; 561.000000
SIM ²	405/DAPI; 488; 555	2; 561nm @1.80%; 405nm @1.80%; 561; 405; 570; 515; 561.000000
Single-molecule STORM	555	Both Channels
Control	405/DAPI; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335 - 383; 555 - 595; 420 - 470; 553

MODALITY ASSESSMENT / PART 1 OF 9

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 9

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.102µm X 0.102µm; Image size (Pixels): 1388 X 1040; Image Size (Pixels): 1388 X 1040; Bit Depth: 12 Bit. Detector/camera: Imaging Device: Axiocam MR R3. Timing: Exposure Time: 100ms; 80ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @20%; X-Cite Xylis Lamp Intensity @15%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 30.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 555; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 9

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1165 X 992; Image Size (Pixels): 1165 X 992; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807. Timing: Exposure Time: 80ms; 25ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @29%; X-Cite Xylis Lamp Intensity @8%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 555; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 9

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 90 Slices (6.23µm); Channels: 2; Scaling (Per Pixel): 71.00nm X 71.00nm X 70.00µm; Image size (Pixels): 753 X 818; Image Size (Pixels): 753 X 818; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.70µs; Frame Time: 8.05s. Illumination: Laser Wavelength: 561nm @0.6%; 405nm @0.3%. Optical/scan: Pinhole: 1.00AU / 49µm; 1.00 AU / 37µm; Scan Zoom: 1.3; LSM Scan Speed: 8; Scan Direction: Bidirectional; Averaging: 8; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 44.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 9

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 94 Slices (4.65µm); Channels: 2; Scaling (Per Pixel): 35.00nm X 35.00nm X 50.00nm; Image size (Pixels): 2186 X 2186; Image Size (Pixels): 2186 X 2186; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.96µs; Frame Time: 10.94s. Illumination: Laser Wavelength: 561nm @0.5%; 405nm @1.0%. Optical/scan: Pinhole: 5.00 AU / 244µm; 5.00 AU / 189µm; Scan Zoom: 1.3; LSM Scan Speed: 6; Scan Direction: Bidirectional; Averaging: 2; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 8.7 (3D, Auto); Super Resolution 8.3 (3D, Auto). Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 9

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 31.30nm X 31.30nm; Image size (Pixels): 2128 X 2128; Image Size (Pixels): 2128 X 2128; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 50ms; Frame Time: 26.51s. Illumination: Laser Wavelength: 561nm @1.80%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 16; Effective NA: 1.46. Structured source metadata values: 43.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 9

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 31.30nm X 31.30nm; Image size (Pixels): 2073 X 1827; Image Size (Pixels): 2073 X 1827; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 50ms; Frame Time: 13.26s. Illumination: Laser Wavelength: 561nm @1.80%; 405nm @1.80%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 8; Effective NA: 1.46. Processing: Process Sampling: 2. Structured source metadata values: 56.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 9

Single-molecule STORM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Bruker Vutara VXL. Objective: Olympus 60X/1.30 silicon oil objective. Platform: Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:.

Detector/camera: Camera: Both Cameras. Timing: Exposure Time: 20ms Selected; Super Resolution Camera Exposure Time: 20ms; Widefield Camera Exposure Time: 100ms. Illumination: Photoactivation/Imaging Laser: 0% Selected.

Structured source metadata values: 71.

VISIBLE OBSERVATION

A preserved presentation image is available for Single-molecule STORM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Single-molecule STORM, documents platform context as Bruker Vutara VXL, and records detected dye/channel descriptors as 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Without localization tables, uncertainty, blinking, drift, and independent cluster evidence, the presentation cannot establish molecule count, cluster truth, exact precision, or molecular architecture.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Single-molecule STORM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 9 OF 9

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1600 X 1104; Image Size (Pixels): 1600 X 1104; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000089 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Confocal / Scaling (Per Pixel): 71.00nm X 71.00nm X 70.00µm -> 0.07100 µm X 0.07100 µm X 0.07000 µm. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=753 X 818; Image Size (Scaled)=53.16µm X 57.75µm; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.57%.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.00nm X 35.00nm X 50.00nm -> 0.03500 µm X 0.03500 µm X 0.05000 µm. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2186 X 2186; Image Size (Scaled)=77.16µm X 77.16µm; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.84%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 31.30nm X 31.30nm -> 0.03130 µm X 0.03130 µm. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2128 X 2128; Image Size (Scaled)=66.61µm X 66.61µm; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 31.30nm X 31.30nm -> 0.03130 µm X 0.03130 µm. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2073 X 1827; Image Size (Scaled)=64.89µm X 57.19µm; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked Stepdown order. Scientific wording is preserved; missing modalities are not invented. A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000089. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra 7	PRESERVED
7	SIM ²	ZEISS Elyra 7	PRESERVED
8	Single-molecule STORM	Bruker Vutara VXL	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

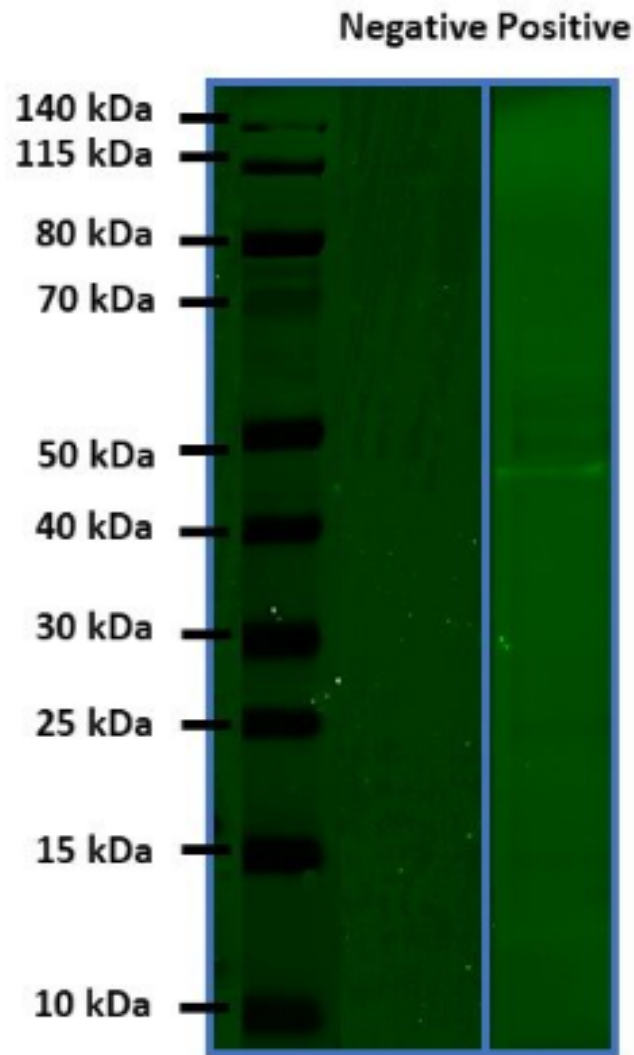

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	555
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	532nm
Emission Filter	572±28nm
Detector	Avalanche Photodiode
Intensity	5
Pixel size	100µm
Scan Speed	High
Quality	High
Focus	Membrane (0µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Donkey Anti-Mouse IgG (H + L)

Cat ID# - SA 000089

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Alpha Tubulin Monoclonal Antibody (YL1/2) for 2 hours, followed by 5X PBS washes. Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat anti-Rat IgG (H+L) was incubated for 1 hour, followed by washing 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 532nm

Emission Filter = 572±28nm

Detector = Avalanche Photodiode

Intensity = 5

Pixel size = 100µm

Scan Speed = High

Quality = High

Focus = Membrane (0µm in Z-axis)

Post Processing

None

Image: Identifyn® LLC

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000089
Stage	Western blot
Instrument	Azure Sapphire FL
Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

SOURCE PROVENANCE

Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	56EB19AF95AA6818CE401C5A56F11DAB6CC202B1ED4D1C5557D2A7135F7BE058
Method PDF SHA-256	FA28BE5BAD83D95F272F9C1E86AFCE43E2E10F9826B0D8E3A0516E0957FC693D

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

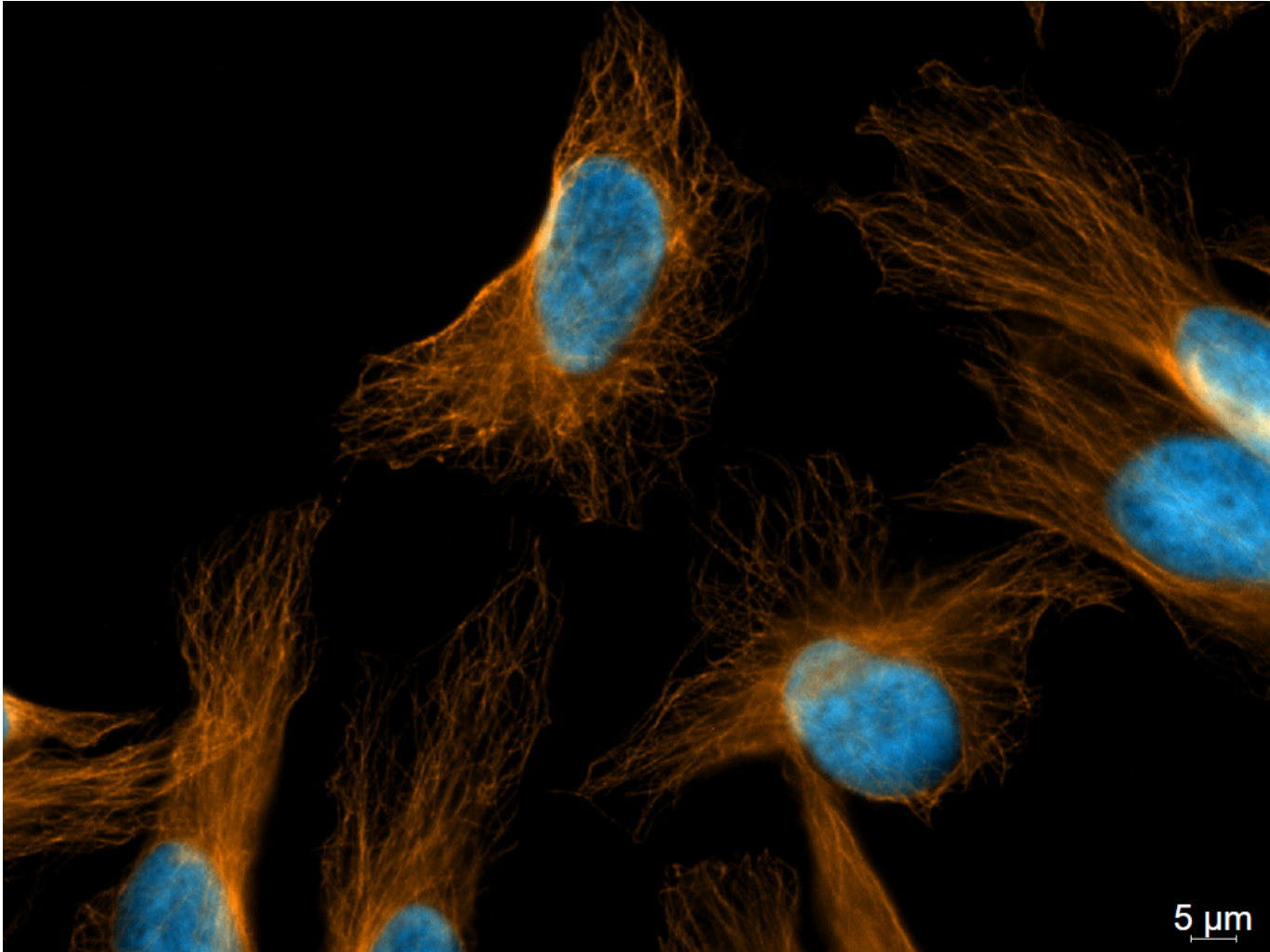

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 555; 594
METADATA VALUES	30

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.102µm X 0.102µm
Image size (Pixels)	1388 X 1040
Image Size (Scaled)	142.10µm X 106.48µm
Bit Depth	12 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	45 Texas Red 49 DAPI
Beam Splitter	585 395
Filter Excitation Wavelength	540 - 580 335 - 383
Filter Emission Wavelength	593 - 668 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @20% X-Cite Xylis Lamp Intensity @15%
Exposure Time	100ms 80ms
Excitation Wavelength	553 353
Emission Wavelength	568 465
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	0.88µm 0.72µm
Binning Mode	1,1

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)
Cat ID# - SA 000089

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, alpha Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2
Scaling (Per Pixel) = 0.102µm X 0.102µm
Image size (Pixels) = 1388 X 1040
Image Size (Scaled) = 142.10µm X 106.48µm
Bit Depth = 12 Bit
Image Center Position = X: 0.00µm, Y: 0.00µm
ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflector = 45 Texas Red
Beam Splitter = 585
Filter Excitation Wavelength = 540 - 580
Filter Emission Wavelength = 593 - 668
Contrast Method = Fluorescence
Light source = X-Cite Xylis Lamp Intensity @20%
Exposure Time = 100ms
Excitation Wavelength = 553
Emission Wavelength = 568
Effective NA = 1.4
Imaging Device = Axiocam MR R3
Camera Adapter = 1X
Depth of Focus = 0.88µm
Binning Mode = 1,1

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @15%
 Exposure Time = 80ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Imaging Device = AxioCam MR R3
 Camera Adapter = 1X
 Depth of Focus = 0.72µm
 Binning Mode = 1,1

Post Processing

None

Image Correction


None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000089
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	30
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	F6A46DE562488E1A8B39A0A5F2A337DF38E9FD79E616095FB60376810A0582EF
Method PDF SHA-256	895663EFE66E04552912A608FF6D1A13AF6522DAE8F3F9C3D70A271FE72D99B2

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

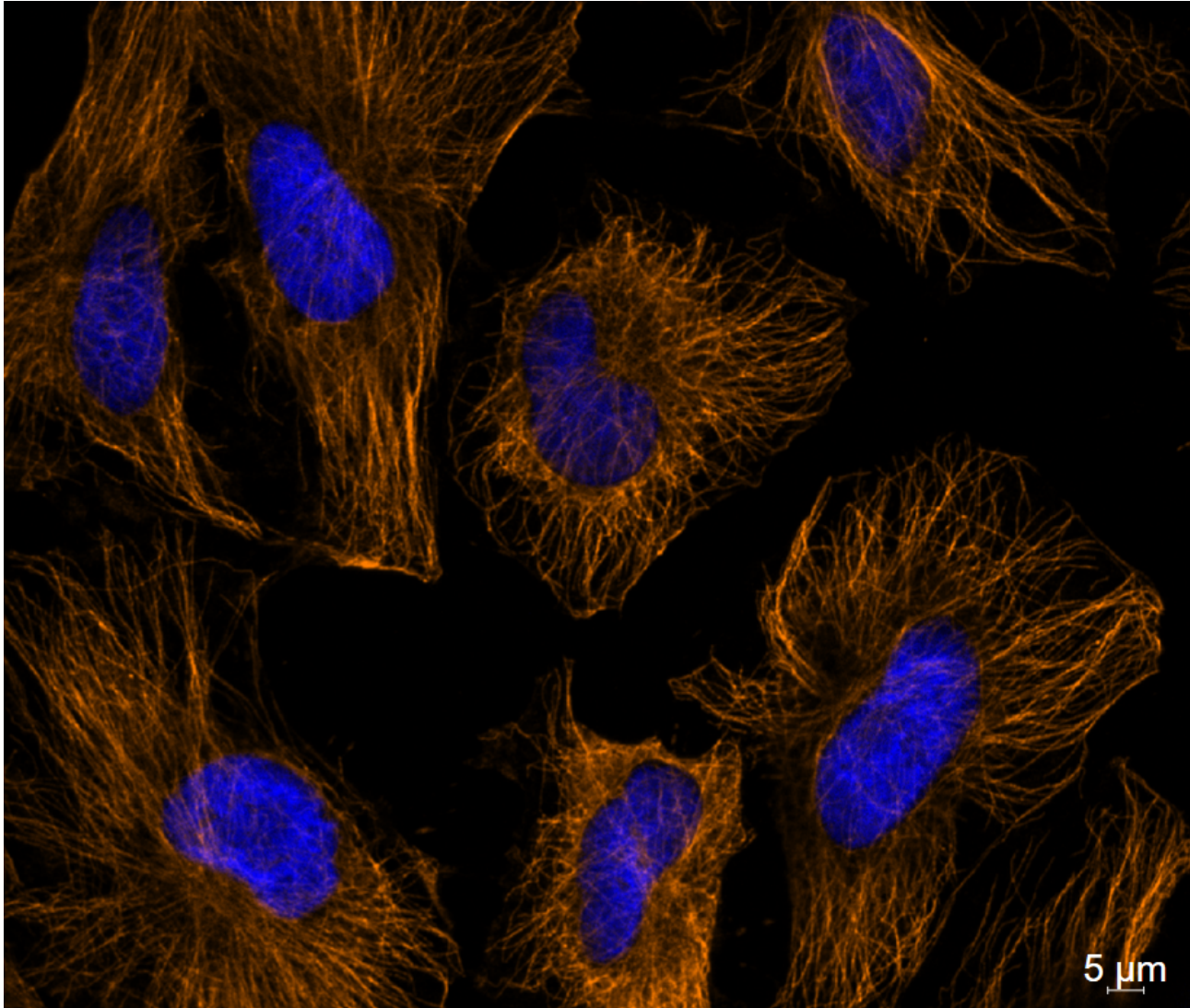

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 555; 594
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1165 X 992
Image Size (Scaled)	166.43µm X 141.73µm
Bit Depth	14 Bit
Image Center Position	X: 76.67mm, Y:49.23mm
Stage Position	X: X: 76.67mm, Y:49.23mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	45 Texas Red 49 DAPI
Beam Splitter	585 395
Filter Excitation Wavelength	540 - 580 335-383
Filter Emission Wavelength	593 - 668 420-470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @29% X-Cite Xylis Lamp Intensity @8%
Exposure Time	80ms 25ms
Excitation Wavelength	553 353
Emission Wavelength	568 465
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.10µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat IgG (H + L)

Cat ID# - SA 000089

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, alpha Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1165 X 992

Image Size (Scaled) = 166.43µm X 141.73µm

Bit Depth = 14 Bit

Image Center Position = X: 76.67mm, Y:49.23mm

Stage Position = X: X: 76.67mm, Y:49.23mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

594 -

Reflector = 45 Texas Red

Beam Splitter = 585

Filter Excitation Wavelength = 540 - 580

Filter Emission Wavelength = 593 - 668

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @29%

Exposure Time = 80ms

Excitation Wavelength = 553

Emission Wavelength = 568

Effective NA = 1.25

Imaging Device = Axiocam 807

Camera Adapter = 1X

Depth of Focus = 1.10µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335-383
 Filter Emission Wavelength = 420-470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @8%
 Exposure Time = 25ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = AxioCam 807
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections


None

Image by B.T. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000089
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	07DD7BAF5A332AC8717DCEE9DE5E155ADC30C8AF1BB65A217B994106EF88EC5D
Method PDF SHA-256	B17D6F4AFF18D989E2C63D5F3FE3B6C0CD7C1AAE418393E4DEB50246386EBD4C

ORIGINAL IMAGE EVIDENCE / CONFOCAL

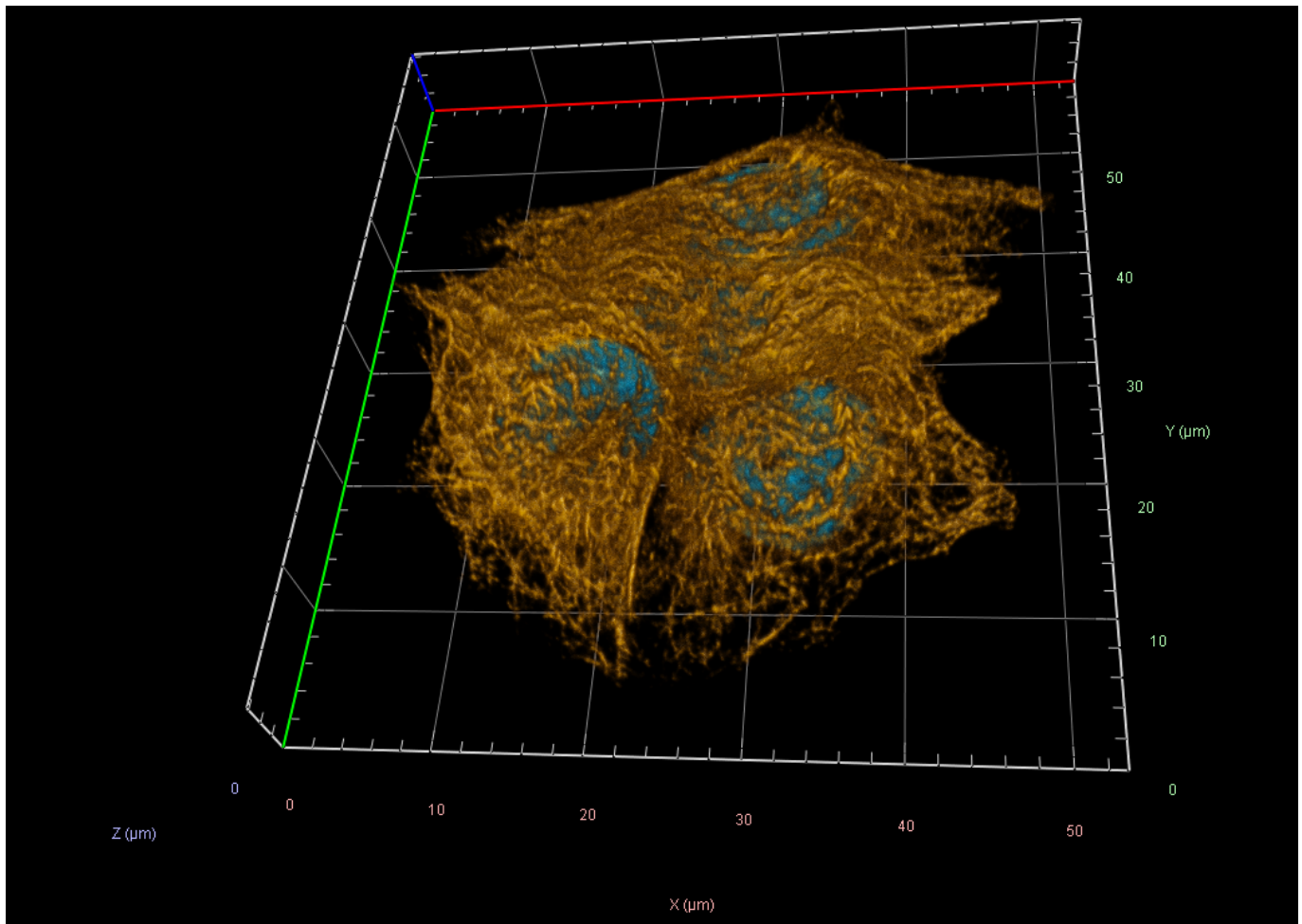

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HCT116
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	44

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	90 Slices (6.23µm)
Channels	2
Scaling (Per Pixel)	0.07100 µm X 0.07100 µm X 0.07000 µm
Image Size (Pixels)	753 X 818
Image Size (Scaled)	53.16µm X 57.75µm
Bit Depth	16 Bit
Image Center Position	X: 70.20mm, Y: 51.12mm
Stage Position	X: 70.20mm, Y: 51.12mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 49µm 1.00 AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	561nm @0.6% 405nm @0.3%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.3
Rotation	0°
Sampling	1.00
Pixel Time	0.70µs
Frame Time	8.05s
LSM Scan Speed	8
Scan Direction	Bidirectional
Line Step	1
Excitation Wavelength	553 353
Emission Wavelength	568 465
Effective NA	1.4
Detection Wavelength	540 - 615 410-470
Imaging Device	Airyscan
Detector Type	GaAsp-PMT
Detector Gain	650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
Averaging	8
3D Volume Projection	Precise
Appearance	Threshold 3% all channels, Ramp97% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness80%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

Cat ID# - SA 000089

HCT116 Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, alpha Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 90 Slices (6.23µm)

Channels = 2

Scaling (Per Pixel) = 71.00nm X 71.00nm X 70.00µm

Image Size (Pixels) = 753 X 818

Image Size (Scaled) = 53.16µm X 57.75µm

Bit Depth = 16 Bit

Image Center Position = X: 70.20mm, Y: 51.12mm

Stage Position = X: 70.20mm, Y: 51.12mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

555 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 49µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @0.6%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.70µs
 Frame Time = 8.05s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 8
 Excitation Wavelength = 553
 Emission Wavelength = 568
 Effective NA = 1.4
 Detection Wavelength = 540 - 615
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector - None
 Contrast Method - Fluorescence
 Pinhole = 1.00 AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.3%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.70µs
 Frame Time = 8.05s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 8
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410-470
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 3% all channels, Ramp97% all channels, Maximum 100% all channels.
 Background = Black
 Light = Brightness80%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.
 Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000089

SOURCE PROVENANCE

Stage	Confocal
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	44
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HCT116
Image SHA-256	2E845B70D0CB828582A0BDEEDB591E88D45C9C0081AAA3343D28B90C413E39DE
Method PDF SHA-256	D1440FCA5240E8D611C13946C289448F46593621338BA85E03F343B1A5A321A7

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

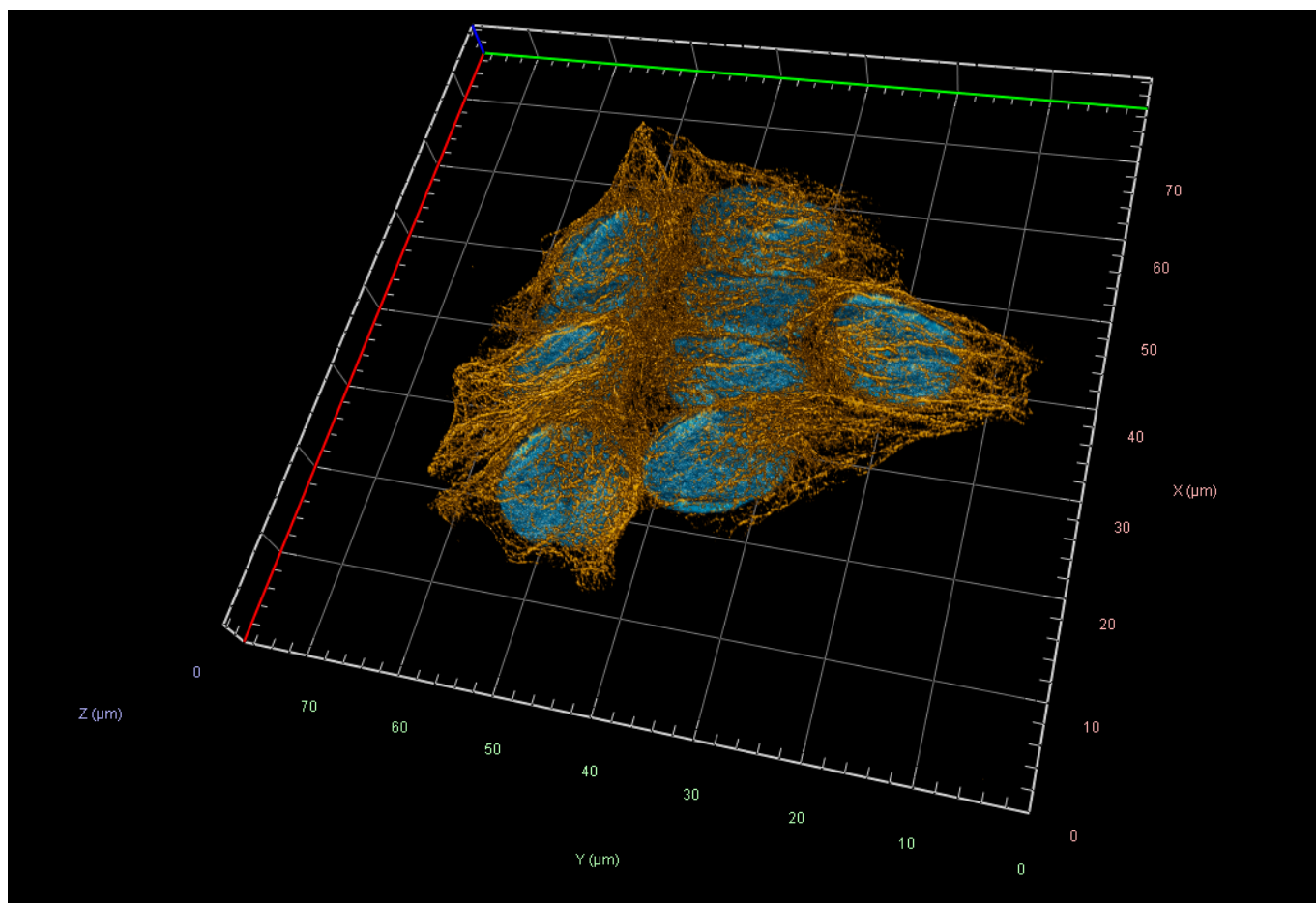

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HCT116
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	94 Slices (4.65µm)
Channels	2
Scaling (Per Pixel)	0.03500 µm X 0.03500 µm X 0.05000 µm
Image Size (Pixels)	2186 X 2186
Image Size (Scaled)	77.16µm X 77.16µm
Bit Depth	16 Bit
Image Center Position	X: 65.32mm, Y: 53.19mm
Stage Position	X: 65.32mm, Y: 513.19mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00 AU / 244µm 5.00 AU / 189µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	561nm @0.5% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.3
Rotation	0°
Sampling	2.00
Pixel Time	0.96µs
Frame Time	10.94s
LSM Scan Speed	6
Scan Direction	Bidirectional
Line Step	1
Averaging	2
Excitation Wavelength	553 353
Emission Wavelength	568 465
Effective NA	1.4
Detection Wavelength	540-604 400-500
Imaging Device	Airyscan
Detector Type	GaAsp-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	800 V
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 8.7 (3D, Auto) Super Resolution 8.3 (3D, Auto)
3D Volume Projection	Precise
Appearance	Threshold 8% all channels, Ramp 92% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness 128%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat IgG (H + L)

Cat ID# - SA 000089

HCT116 Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 94 Slices (4.65µm)

Channels = 2

Scaling (Per Pixel) = 35.00nm X 35.00nm X 50.00nm

Image Size (Pixels) = 2186 X 2186

Image Size (Scaled) = 77.16µm X 77.16µm

Bit Depth = 16 Bit

Image Center Position = X: 65.32mm, Y: 53.19mm

Stage Position = X: 65.32mm, Y: 513.19mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

555 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 244µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @0.5%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 10.94s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 553
 Emission Wavelength = 568
 Effective NA = 1.4
 Detection Wavelength = 540-604
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.7 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 189µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 10.94s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400-500
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.3 (3D, Auto)

Post Processing - AiryScan 3D Processing

Display Mode "SR"
 Super Resolution Auto Filter Selected
 Process 3D Current Image
 Created

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 8% all channels, Ramp 92% all channels, Maximum 100% all channels.
 Background = Black
 Light = Brightness 128%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections -

None


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000089
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HCT116
Image SHA-256	F90BC98506EF3C81C13C96EDBB94F6A7860F39FBB3DC47BCB9126CED44E9B304
Method PDF SHA-256	8E345AF686F6B59A01F80EDACE79354EA87A8EC051B647903036B088E4EBC372

ORIGINAL IMAGE EVIDENCE / SIM

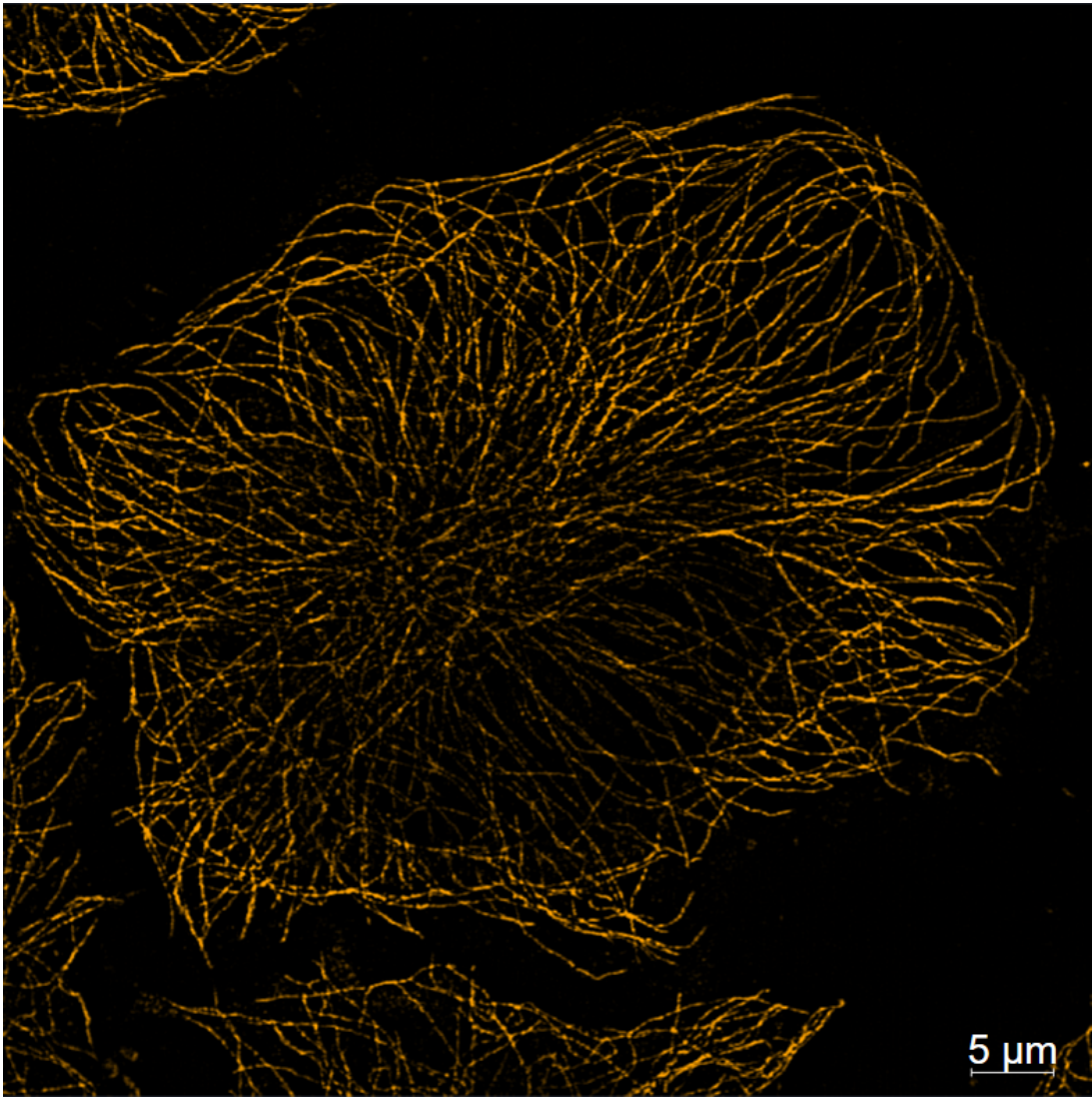

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	43

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 555	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 560
Beam Splitter 405	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Channels	2
Scaling (Per Pixel)	0.03130 µm X 0.03130 µm
Image Size (Pixels)	2128 X 2128
Image Size (Scaled)	66.61µm X 66.61µm
Bit Depth	16 Bit
Image Center Position	X: -928.10mm, Y: -1.26mm
Stage Position	X: -928.10mm, Y: -1.26mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	561nm @1.80%
Scan Mode	Frame
Scan Zoom	1.0
Rotation	0°
Frame Time	26.51s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	16
Excitation Wavelength	561
Emission Wavelength	570
Effective NA	1.46
Detection Wavelength	570-1000
Imaging Device	Edge 4.2 SIM
Exposure Time	50ms
Depth of Focus	0.81µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	561.000000
Grating Period (μm)	27.500000
Processing	2D
Input SNR	Yes
Sharpness	8.611391
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

Cat ID# - SA 000089

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 555 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Beam Splitter 405 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 31.30nm X 31.30nm

Image Size (Pixels) = 2128 X 2128

Image Size (Scaled) = 66.61µm X 66.61µm

Bit Depth = 16 Bit

Image Center Position = X: -928.10mm, Y: -1.26mm

Stage Position = X: -928.10mm, Y: -1.26mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

555 -

Contrast Method = Fluorescence
 Laser Wavelength = 561nm @1.80%
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 26.51s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 16
 Excitation Wavelength = 561
 Emission Wavelength = 570
 Effective NA = 1.46
 Detection Wavelength = 570-1000
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.81µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM Processing

Channel = 561.000000
Grating Period (µm): 27.500000

SIM

Processing: 2D
 Input SNR: Yes
 Sharpness: 8.611391
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Detrend: No
 Experimental PSF: No

Image Corrections

None


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000089
Stage	SIM
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	43
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	9779FD1FFBFD1EE00E669F48CDE8931061EDEF9652F868E5D7778598291FB768
Method PDF SHA-256	4C2364F7754B9159C7E1C20C20EC9631011E99C310415810A84A3A444B0B02AB

ORIGINAL IMAGE EVIDENCE / SIM²

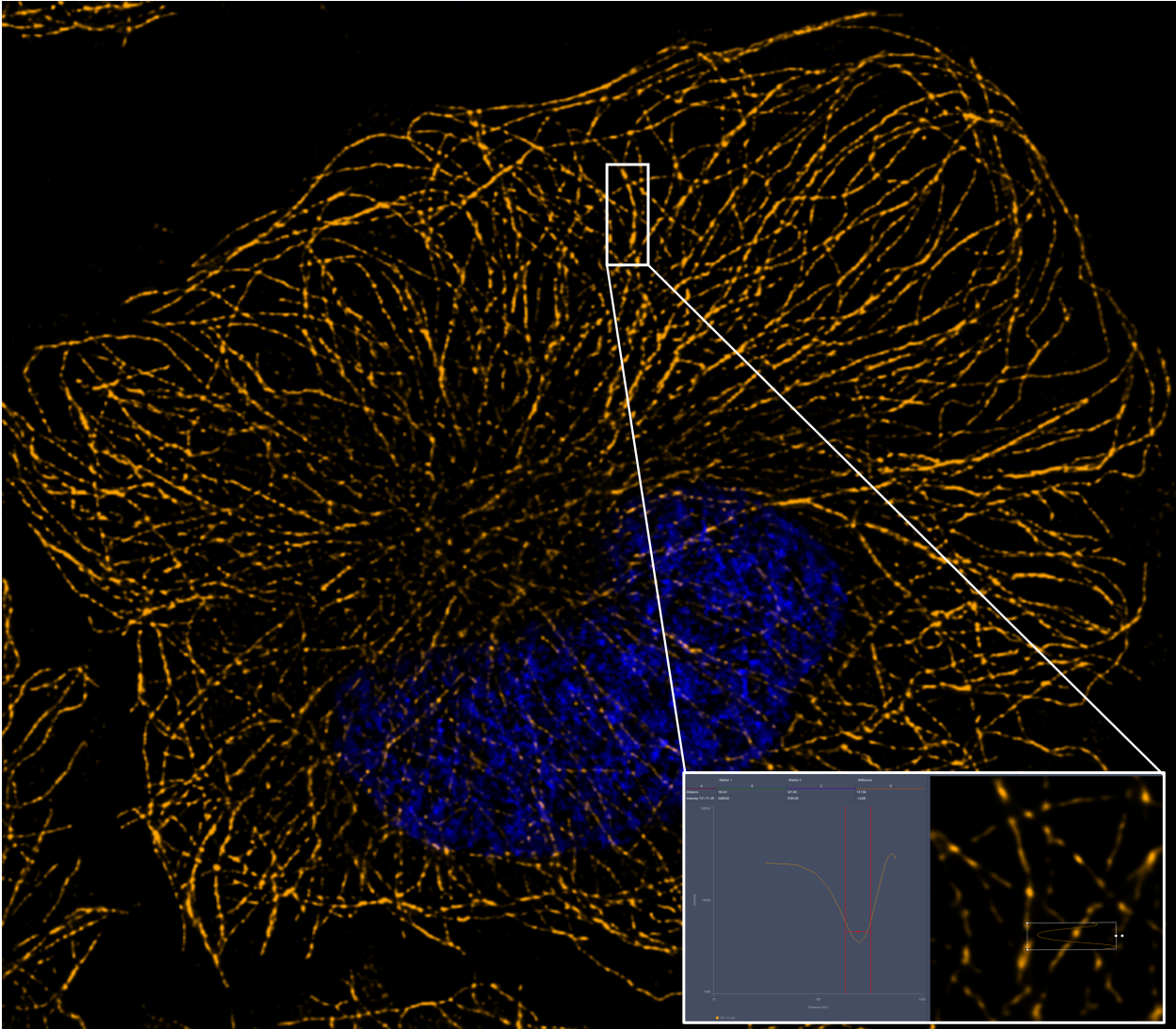

EVIDENCE TIER	SIM²
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	56

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 555	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 560
Beam Splitter 405	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Channels	2
Scaling (Per Pixel)	0.03130 µm X 0.03130 µm
Image Size (Pixels)	2073 X 1827
Image Size (Scaled)	64.89µm X 57.19µm
Bit Depth	16 Bit
Image Center Position	X: -928.10mm, Y: -1.26mm
Stage Position	X: -928.10mm, Y: -1.26mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	561nm @1.80% 405nm @1.80%
Scan Mode	Frame
Scan Zoom	1.0
Rotation	0°
Frame Time	13.26s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	8
Excitation Wavelength	561 405
Emission Wavelength	570 515
Effective NA	1.46
Detection Wavelength	570-1000 495-550
Imaging Device	Edge 4.2 SIM
Exposure Time	50ms
Depth of Focus	0.81µm 0.73µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	561.000000 405.000000
Grating Period (µm)	23.000000 36.500000
Processing	2D, Leap
Input SNR	Low
Iterations	16
Regularization Weight	0.065000
Process Sampling	2
Output Sampling	2
Filter	Median
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No
Profile Definition (left Inset)	Stroke thickness 1, showing profile in graphics was selected, and the Caliper X was employed in measurement between objects as shown in the picture. Grid distance = 1
Profile View	Image was selected, Profile Table and Int. Measurements were NOT selected. Log X Scale was chosen, but Log Y Scale was NOT selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

Cat ID# - SA 000089

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (YL1/2) was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 555 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Beam Splitter 405 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 31.30nm X 31.30nm

Image Size (Pixels) = 2073 X 1827

Image Size (Scaled) = 64.89µm X 57.19µm

Bit Depth = 16 Bit

Image Center Position = X: -928.10mm, Y: -1.26mm

Stage Position = X: -928.10mm, Y: -1.26mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

555 -

Contrast Method = Fluorescence
 Laser Wavelength = 561nm @1.80%
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 13.26s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 8
 Excitation Wavelength = 561
 Emission Wavelength = 570
 Effective NA = 1.46
 Detection Wavelength = 570-1000
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.81µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 561.000000
Grating Period (µm): 23.000000

SIM2

Processing: 2D, Leap
 Input SNR: Low
 Iterations: 16
 Regularization Weight: 0.065000
 Process Sampling: 2
 Output Sampling: 2
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

405 -

Contrast Method = Fluorescence
 Laser Wavelength = 405nm @1.80%
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 13.26s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 8
 Excitation Wavelength = 405
 Emission Wavelength = 515
 Effective NA = 1.46
 Detection Wavelength = 495-550
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.73µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 405.000000
Grating Period (µm): 36.500000

SIM2

Processing: 2D, Leap
 Input SNR: Low
 Iterations: 16
 Regularization Weight: 0.065000
 Process Sampling: 2
 Output Sampling: 2
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

Profile

Profile Definition (left Inset) = Stroke thickness 1, showing profile in graphics was selected, and the Caliper X was employed in measurement between objects as shown in the picture. Grid distance = 1

Profile View = Image was selected, Profile Table and Int. Measurements were NOT selected. Log X Scale was chosen, but Log Y Scale was NOT selected.

The graphical representation approximates the distance between signals using the calipers (red) against the signal curve (orange). Right is a profile view of the line between the signals covering the area without a signal. While this is a subjective measurement, it approximates ~137.58nm between the structure and within the gap.

Summary

The super-resolution SIM2 Image is meant to demonstrate the power of imaging below the diffraction limit and, importantly, that a known structure, in this case, tubulin, with the Identifyn® SRM Alexa Fluor™ 555 is bright and target specific. With increased resolution, off-target binding becomes increasingly visible; in this case, the detail demonstrates the product's specificity and stability. The image measurements were not meant to imply biological relevance but to demonstrate the detail achieved at such resolutions.

Image Corrections


None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000089
Stage	SIM²
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	56
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	445F77A4325C83E9902A6E08EF09F63A690961D41D5C7C46E450FBCE1898DEBA
Method PDF SHA-256	CF790A5D864FA828747AFBA3DB4D4C46B539D321E686E4AB92BC4FE564D9FB8A

ORIGINAL IMAGE EVIDENCE / SINGLE-MOLECULE STORM

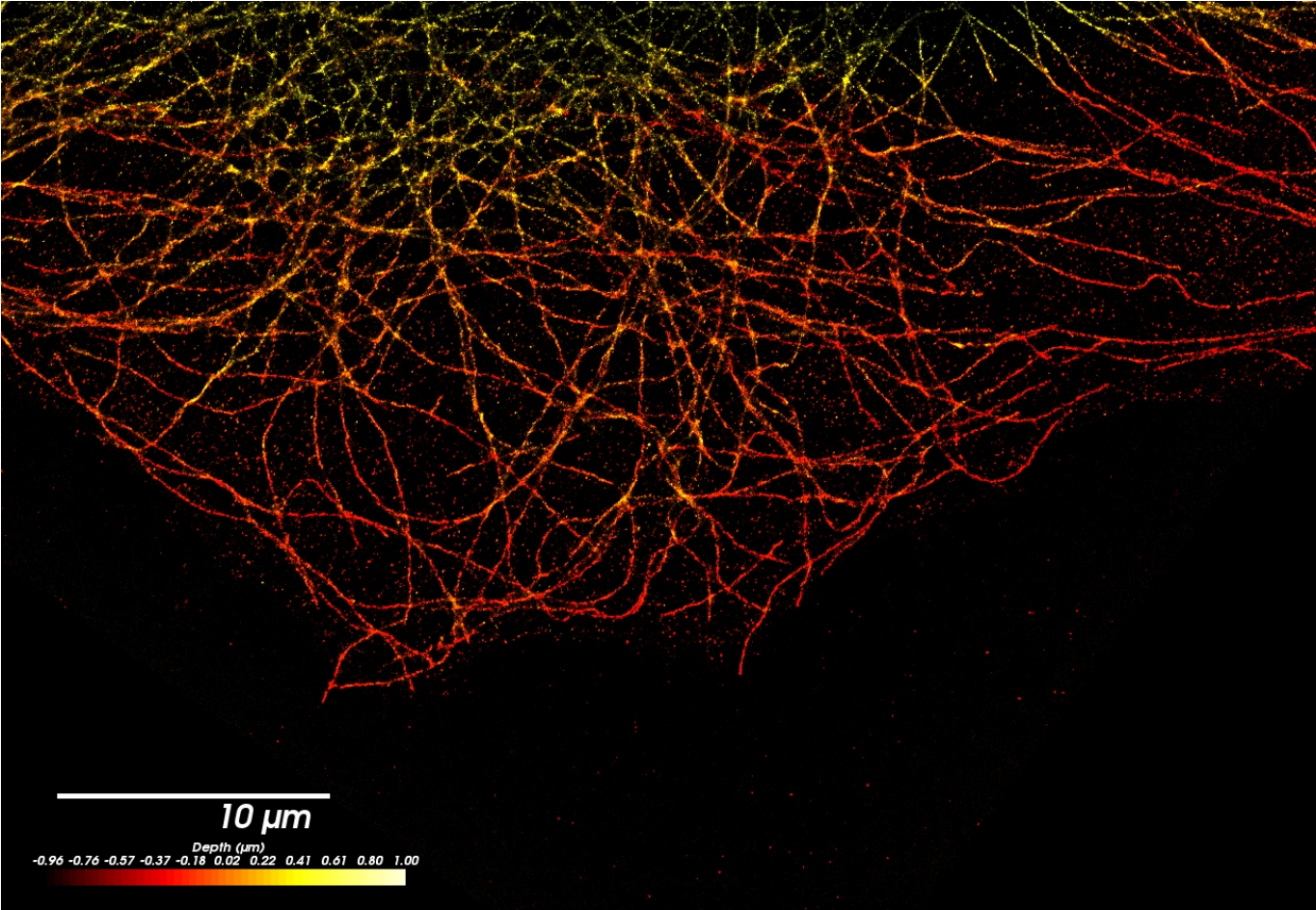

EVIDENCE TIER	Single-molecule STORM
INSTRUMENT	Bruker Vutara VXL
CELL MODEL	HeLa
DETECTED DYES	555
METADATA VALUES	71

STRUCTURED ACQUISITION METADATA / SINGLE-MOLECULE STORM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Objective	Olympus 60X/1.30 silicon oil objective
Number of Probes	1
Number of Simultaneous Probes	1
Number of Reference Probes	1
Camera	Both Cameras
Color Channel	Both Channels
Probe Capture Mode	Interleaved
Z -Stack Capture Mode	Interleaved
Multiple Location Capture	Not selected
Post Record Reference	Not Selected
Fluidics	Not Selected
Autofocus Drift Correction	On
Biplane Module	635nm Dichroic
Objective	Silicon Selected
Stage Insert	Rectangular Selected
Course Focus Position	8-Well Chamber Selected
Dichroic	SuperRes
PSF	(Red, Orange)
Heater	Current: 26
Recording Vertical Field of View	100%
Widefield Camera Exposure Time	100ms 20ms
Super Resolution Camera Exposure Time	20ms
Widefield Camera Binning	1
Per-Probe Exposure Time	Selected
Predefined Probe	First probe: 555 nm (635 nm Dichroic); First reference probe: 555 nm (635 nm Dichroic)
Photoactivation/Imaging Laser	0% Selected
Exposure Time	20ms Selected
Record Workflow	Selected Options
View Mode	Widefield 200µm
Capture Mode	Live
Piezo Current	178.5 µm

FIELD	SOURCE-DOCUMENTED VALUE
Widefield 555nm	0.1X Neutral Density Filter: Selected, @0.2% Laser Input
Reference Probe 1 Laser control (555 nm (635 nm Dichroic))	Not Selected
Widefield camera Exposure Time	100ms
Widefield Laser	"Widefield 555nm" was deselected, and "calibrate" under Autofocus was selected.
Current Intensity	4.5
Calibration Value	-5070 nm/unit
Autofocus Mode	Pre-Capture
SuperRes	50µm is Selected
Piezo Record	Begin: 178.5; End: 178.5
SuperRes 555nm	0.1X Neutral Density Filter: Selected, @0.02% Laser Input
Record Workflow - Record Tab	Capture Parameters
Frames	Probe 1: 40000
Z Positions	1 using Steps: 0.1µm (100nm)
Cycles	1
Time Points	1
Pause Between Cycles	Not Selected
Automatic Photoactivation Laser Control	Not Selected
Laser Power	90%
Cutout Width	16px
Max Frame Accumulation	1 frame
Max Accumulation Offset	1px
Fiducial at Max Accumulation	Not selected
Max Particles Per Frame	No Max
Background SLPE Removal	Not Selected
Cutout/Background Removal	Not Selected
Preview Cutout Positions	Selected
Probe1	BG threshold: 20
Sizing Mode	Radial precision
Particle size	10nm
Color by	Probe
Color Table	Single
Opacity Mode	Constant; Opacity 1: 0.02
Front-to-back Attenuation	Not Selected
Method	Particle (direct)
Time Intervals	10 Intervals

FIELD	SOURCE-DOCUMENTED VALUE
Pixel Size	50nm
Smoothing	8.00
PSF Z offset	-450 to +450
Radial precision	30
Axial precision	60

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

Cat ID# - SA 000089

Hela cells were grown on μ -Slide 8 Well high Glass Bottom Ibidi chamber, followed by -20 methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Tubulin Monoclonal Antibody (YL1/2), was incubated for 1 hour at room temperature at a dilution of 1:100. Cells were washed 3X in PBS, and Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:100. Cells were washed 5X in PBS followed by the addition of Identifyn's™ proprietary STORM Buffer to the sample chamber enabling image acquisition.

Microscope

Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:

Calibration -

Calibration follows Identifyn's Calibration Standard Operating Procedure, which can be found here.

<https://identifyn.com/storm-calibration>

Image Acquisition -

SRX Software - Under the 'Single Molecule Localization' tab, 'Record' is selected and contains the following workflow tabs and parameters.

Configuration Tab - Capture Configuration

Number of Probes: 1

Number of Simultaneous Probes: 1

Number of Reference Probes: 1

Camera: Both Cameras

Color Channel: Both Channels

Probe Capture Mode: Interleaved

Z -Stack Capture Mode: Interleaved

Multiple Location Capture: Not selected

Post Record Reference: Not Selected

Fluidics: Not Selected

Autofocus Drift Correction: On

Microscope configuration

Biplane Module: 635nm Dichroic
 Objective: Silicon Selected
 Stage Insert: Rectangular Selected
 Course Focus Position: 8-Well Chamber Selected
 Dichroic: SuperRes
 PSF: (Red, Orange)
 Heater: Current: 26

Camera Configuration

Recording Vertical Field of View: 100%
 Widefield Camera Exposure Time: 100ms
 Super Resolution Camera Exposure Time: 20ms
 Widefield Camera Binning: 1
 Per-Probe Exposure Time: Selected

Probe Configuration

Predefined Probe: First probe: 555 nm (635 nm Dichroic); First reference probe: 555 nm (635 nm Dichroic)
 Photoactivation/Imaging Laser: 0% Selected
 Exposure Time: 20ms Selected

Experiment Configuration

Record Workflow: Selected Options
 View Mode: Widefield 200µm
 Capture Mode: Live

Stage and Piezo Selection:

Examination of Sample:

Cell(s) or Cellular Structure Localization - The red-bordered box identifies the area of the SuperRes view and defines a suitable cell(s) or cellular structure. The cell(s) or structures are optimally focused.

Piezo Current: 178.5 µm

Probe 1 Laser control (555 nm (635 nm Dichroic))

Widefield 555nm: 0.1X Neutral Density Filter: Selected, @0.2% Laser Input
 Reference Probe 1 Laser control (555 nm (635 nm Dichroic): Not Selected

Advanced Tools Option

Widefield camera Exposure Time: 100ms
 Default values are maintained unless otherwise noted.

Widefield Laser: "Widefield 555nm" was deselected, and "calibrate" under Autofocus was selected.

Autofocus - Calibration Values

Current Position" -0.979

Current Intensity: 4.5

Calibration Value: -5070 nm/unit

Autofocus Mode: Pre-Capture

View Mode

SuperRes: 50µm is Selected

Capture Mode

Live

Stage and Piezo Selection

Cell(s) or Cellular Structure Localization - Positioning of the target Cell or Subsection of the Target Cell is optimized and optimally focused.

Piezo Record: Begin: 178.5; End: 178.5

Probe 1 Laser control (555 nm (635 nm Dichroic)

SuperRes 555nm: 0.1X Neutral Density Filter: Selected, @0.02% Laser Input

Reference Probe 1 Laser control (555 nm (635 nm Dichroic)

Not Selected

Advanced Tools

Widefield Camera Exposure Time: 20ms

Default values are maintained unless otherwise noted.

Record Workflow - Record Tab: Capture Parameters

Frames: Probe 1: 40000

Z Positions: 1 using Steps: 0.1µm (100nm)

Cycles: 1

Time Points: 1

Pause Between Cycles: Not Selected

Automatic Photoactivation Laser Control: Not Selected

Probe 1 Laser control (555 nm (635 nm Dichroic)

Laser Power: 90%

Expert Configuration Options - The Localization Tab

Beyond the values shown, all other parameters were default values.

Localization

Cutout Width: 16px
 Max Frame Accumulation: 1 frame
 Max Accumulation Offset: 1px
 Fiducial at Max Accumulation: Not selected
 Max Particles Per Frame: No Max
 Background SLPE Removal: Not Selected
 Cutout/Background Removal: Not Selected
 Preview Cutout Positions: Selected

Record Workflow - Localization Tab

Data Collection

The folder is created automatically with the localization file for this image acquisition.

Recording Summary

All parameters were carried over from previous tabs and are kept except:

Probe1: BG threshold: 20

Region of Interest

Not selected

Display

Was deployed to inspect if the BG threshold applied chooses localizations in an expected manner

Record Workflow - View Tab

Region of Interest

Default values

Visualization Parameters

Point Splatting
 Sizing Mode: Radial precision
 Particle size: 10nm
 Color by: Probe
 Color Table: Single
 Opacity Mode: Constant; Opacity 1: 0.02
 Front-to-back Attenuation: Not Selected

Post Acquisition analysis - View Tab

Drift Correction - Computed with the following parameters:

Method: Particle (direct)
 Time Intervals: 10 Intervals
 Pixel Size: 50nm
 Smoothing: 8.00

Advanced Particle Filters - Applied to the drift-corrected data with the following parameters

PSF Z offset: -450 to +450

Radial precision: 30

Axial precision: 60

Statistical analysis - N/A


Image by S. Thakur

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000089
Stage	Single-molecule STORM
Instrument	Bruker Vutara VXL
Microscope configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Structured source metadata values	71
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	C74F572020E7053F552E05750E163FE997625E1ED6F613CFABC35C08440CD6D7
Method PDF SHA-256	2D923F3772450FE1E1A69684503FBCDAF1384F3C60B1552620E1B60E8FCCE9F8

ORIGINAL IMAGE EVIDENCE / CONTROL

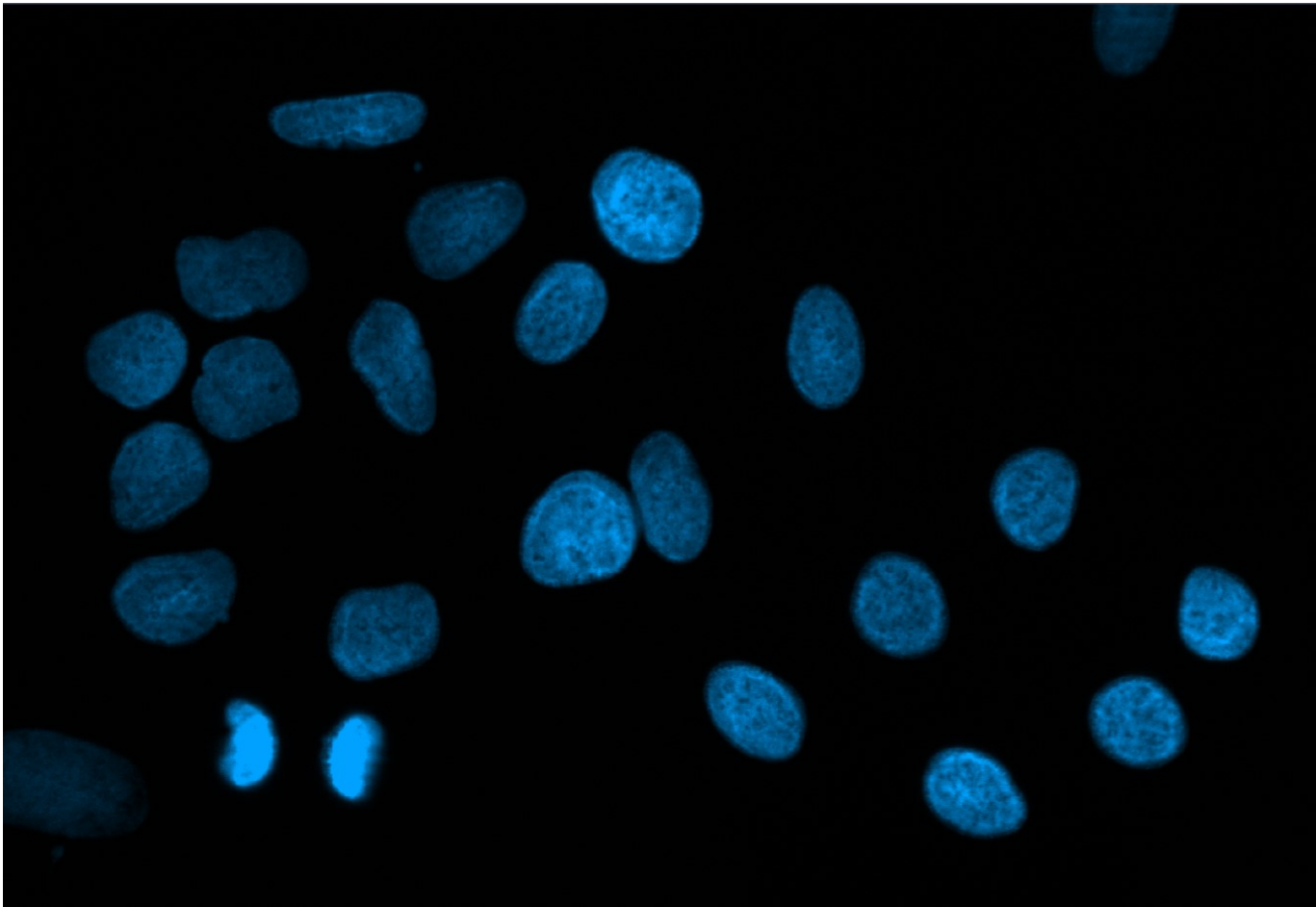

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 555
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1600 X 1104
Image Size (Scaled)	228.57µm X 157.71µm
Bit Depth	14 Bit
Image Center Position	X: 51.06mm, Y: 52.69mm
Stage Position	X: 51.06mm, Y: 52.69mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335 - 383
Filter Emission Wavelength	555 - 595 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	553 353
Emission Wavelength	568 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.10µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Goat Anti-Rat, IgG (H + L)

Cat ID# - SA 000089

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1600 X 1104

Image Size (Scaled) = 228.57µm X 157.71µm

Bit Depth = 14 Bit

Image Center Position = X: 51.06mm, Y: 52.69mm

Stage Position = X: 51.06mm, Y: 52.69mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

555 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 553

Emission Wavelength = 568

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.10µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000089
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	B02432BA821F7559BD800B3825683C1080526842530F61EBB9A70CB99AE0D03B
Method PDF SHA-256	DE3D72912E6BC4DF0E19D47A0AF0DD3FC4AD4E3DD69B8AC7944F7064B1249031