

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit, F(ab')₂ (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM² -> SINGLE-MOLECULE STORM

8


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

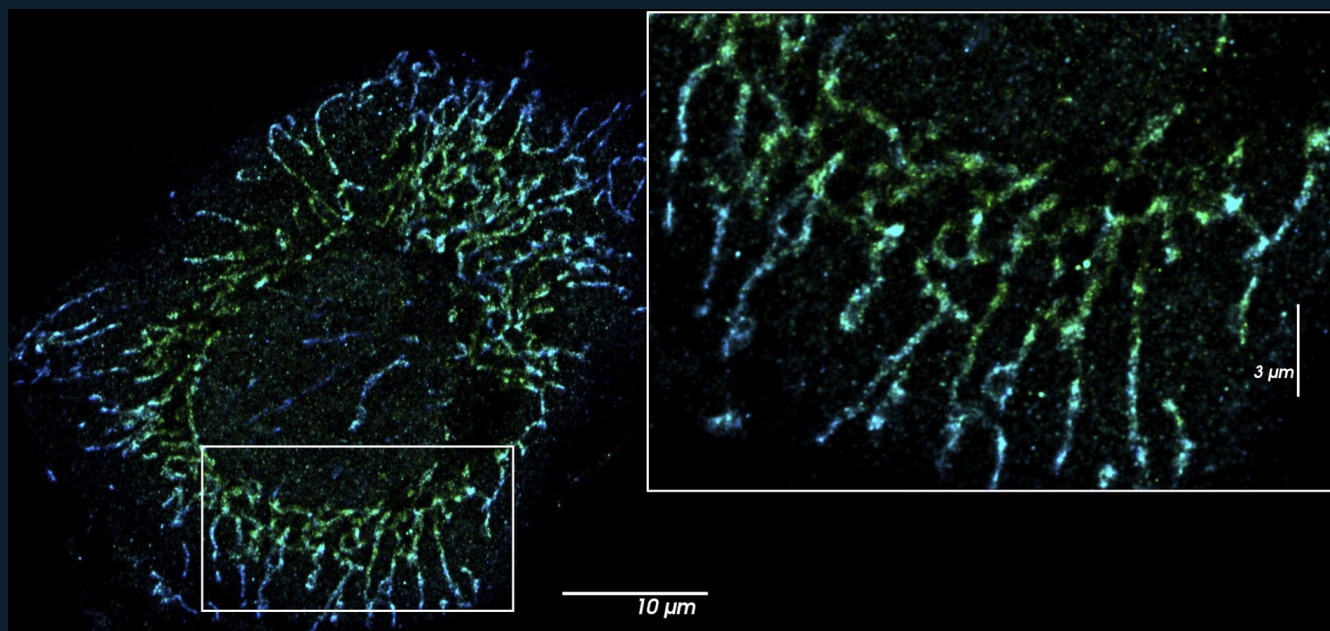

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SINGLE-MOLECULE STORM

SA-000018 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Single-molecule STORM presentation workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit, F(ab')₂ (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Single-molecule STORM
-> Control

BENNETT ASSESSMENT

The preserved SA-000018 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / FULL STEPDOWN / SOURCE-BOUNDED REVIEW

The preserved SA-000018 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

SETTINGS ASSESSMENT

The record preserves 8 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	488	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 488; 555	3; 38 HE eGFP; AF532-49014; 450 - 490; 515 - 545; 500 - 550; 555 - 595; 493
Widefield SIM — Apotome	405/DAPI; 488; 555	3; 38 HE eGFP; AF532-49014; 49 DAPI; 450 - 490; 515 - 545; 335-383; 500 - 550
Confocal	405/DAPI; 488	2; None; 488nm @0.2%; 405nm @0.3%; 493; 353; 517; 465
Airyscan	405/DAPI; 488	2; None; 488nm @0.09%; 405nm @1.0%; 493; 353; 517; 465
SIM	405/DAPI; 488; 594	2; 561nm @5.00%; 488nm @1.50%; 561; 488; 515; 655; 561.000000 (594)
SIM ²	405/DAPI; 488; 594	2; 488nm @1.50%; 561nm @5.0%; 488; 561; 655; 515; 488.000000
Single-molecule STORM	488	Both Channels
Control	405/DAPI; 488	2; 38 HE eGFP; 49 DAPI; 450 - 490; 335-383; 500 -550; 420-470; 493

MODALITY ASSESSMENT / PART 1 OF 9

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 488.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 9

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 3; Scaling (Per Pixel): 0.102µm X 0.102µm; Image size (Pixels): 968 X 793; Image Size (Pixels): 968 X 793; Bit Depth: 12 Bit. Detector/camera: Imaging Device: AxioCam MR R3. Timing: Exposure Time: 100ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @20%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 28.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 9

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Z-Stack: 14 Slices (2.6µm); Channels: 3; Scaling (Per Pixel): 0.143µm X 0.143µm X 0.200µm; Image size (Pixels): 723 X 628; Image Size (Pixels): 723 X 628; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807. Timing: Exposure Time: 100ms; 75ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @20%; X-Cite Xylis Lamp Intensity @5%. Optical/scan: Effective NA: 1.25. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 488; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 9

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 56 Slices (3.3µm); Channels: 2; Scaling (Per Pixel): 0.071µm X 0.071µm X 0.060µm; Image size (Pixels): 684 X 734; Image Size (Pixels): 684 X 734; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; GaAsP-PMT 1; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.65µs; Frame Time: 17.44s. Illumination: Laser Wavelength: 488nm @0.2%; 405nm @0.3%. Optical/scan: Pinhole: 1.00AU / 44µm; 1.00 AU / 37µm; Scan Zoom: 1.2; LSM Scan Speed: 8; Scan Direction: Bidirectional; Averaging: 16; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 45.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 9

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 39 Slices (1.9µm); Channels: 2; Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm; Image size (Pixels): 449 X 591; Image Size (Pixels): 449 X 591; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.96µs; Frame Time: 21.88s. Illumination: Laser Wavelength: 488nm @0.09%; 405nm @1.0%. Optical/scan: Pinhole: 5.00AU / 220µm; 5.00 AU / 181µm; Scan Zoom: 1.3; LSM Scan Speed: 6; Scan Direction: Bidirectional; Averaging: 4; Effective NA: 1.4. Processing: Projection: View Angle 57°, Scale Z 114%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 10.0 (3D, Auto); Super Resolution 8.6 (3D, Auto). Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 9

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 22 Slices (1.47µm); Channels: 2; Scaling (Per Pixel): 31.30nm X 31.30nm X 70.00nm; Image size (Pixels): 1429 X 2055; Image Size (Pixels): 1429 X 2055; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 100ms; 50ms; Frame Time: 1.00s. Illumination: Laser Wavelength: 561nm @5.00%; 488nm @1.50%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 4; Effective NA: 1.46. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 56.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 9

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 20 Slices (1.14µm); Channels: 2; Scaling (Per Pixel): 15.65nm X 15.65nm X 60.00nm; Image size (Pixels): 1506 X 1301; Image Size (Pixels): 1506 X 1301; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 50ms; 100ms; Frame Time: 1.0s. Illumination: Laser Wavelength: 488nm @1.50%; 561nm @5.0%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 4; Effective NA: 1.46. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 59.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 9

Single-molecule STORM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Bruker Vutara VXL. Objective: Olympus 60X/1.30 silicon oil objective. Platform: Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:.

Detector/camera: Camera: Both Cameras. Timing: Exposure Time: 20ms Selected; Super Resolution Camera Exposure Time: 20ms; Widefield Camera Exposure Time: 100ms. Illumination: Photoactivation/Imaging Laser: 0% Selected.

Structured source metadata values: 71.

VISIBLE OBSERVATION

A preserved presentation image is available for Single-molecule STORM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Single-molecule STORM, documents platform context as Bruker Vutara VXL, and records detected dye/channel descriptors as 488.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Without localization tables, uncertainty, blinking, drift, and independent cluster evidence, the presentation cannot establish molecule count, cluster truth, exact precision, or molecular architecture.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Single-molecule STORM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 9 OF 9

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1; ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1366 X 1080; Image Size (Pixels): 1366 X 1080; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 30.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Imager.M1; ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 488.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000018 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm -> 0.03530 μ m X 0.03530 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=449 X 591; Image Size (Scaled)=15.85 μ m X 20.86 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 31.30nm X 31.30nm X 70.00nm -> 0.03130 μ m X 0.03130 μ m X 0.07000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1429 X 2055; Image Size (Scaled)=44.73 μ m X 64.33 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 15.65nm X 15.65nm X 60.00nm -> 0.01565 μ m X 0.01565 μ m X 0.06000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1506 X 1301; Image Size (Scaled)=23.57 μ m X 20.36 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.00%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked
Stepdown order. Scientific wording is preserved; missing modalities are not invented.
A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000018. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra 7	PRESERVED
7	SIM ²	ZEISS Elyra 7	PRESERVED
8	Single-molecule STORM	Bruker Vutara VXL	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

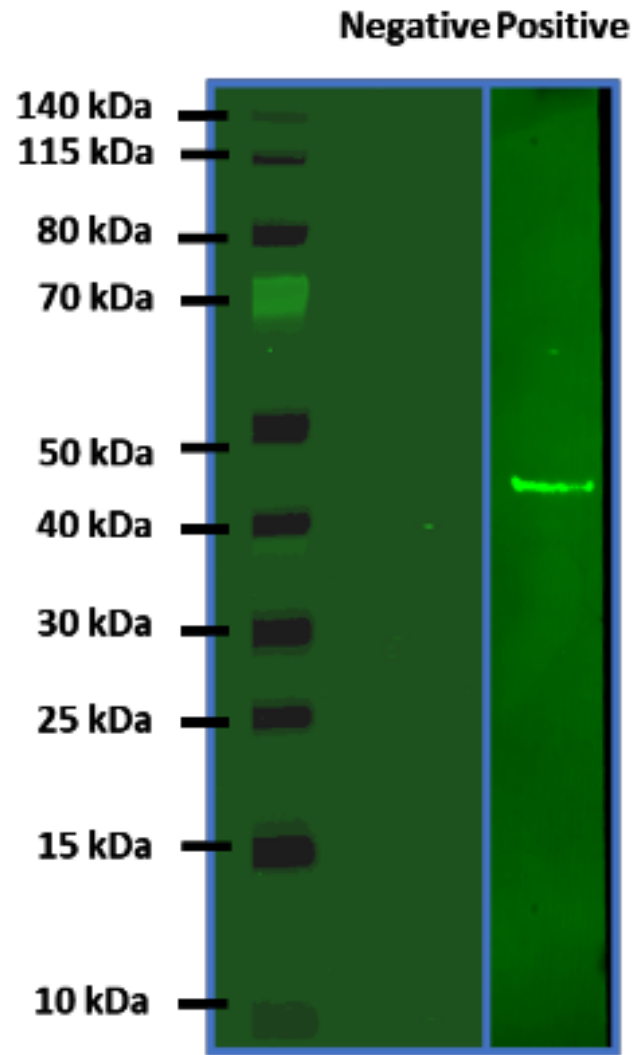

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	488
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	488nm
Emission Filter	513±17nm
Detector	PMT
Intensity	L2
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Actin Recombinant Monoclonal Antibody (JF47-01) (MA5-32530) for 2 hours, followed by 5X PBS washes. Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat anti-Rabbit F(ab)₂ was incubated for 1 hour and washed 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 488nm

Emission Filter = 513±17nm

Detector = PMT

Intensity = L2

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: Identifyn® LLC

Copyright 2026 GMD12 LLC.

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000018

SOURCE PROVENANCE

Stage	Western blot
Instrument	Azure Sapphire FL
Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	3E1A2D3ECC3BC959200F5D8F6DD6C3EE9945A72229E39B7C112FD388ABCC1E29
Method PDF SHA-256	0DD5805D2BFAF780745E4601455CE2B0CAF5505D2FCB295E3B86A3AE2C866A35

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

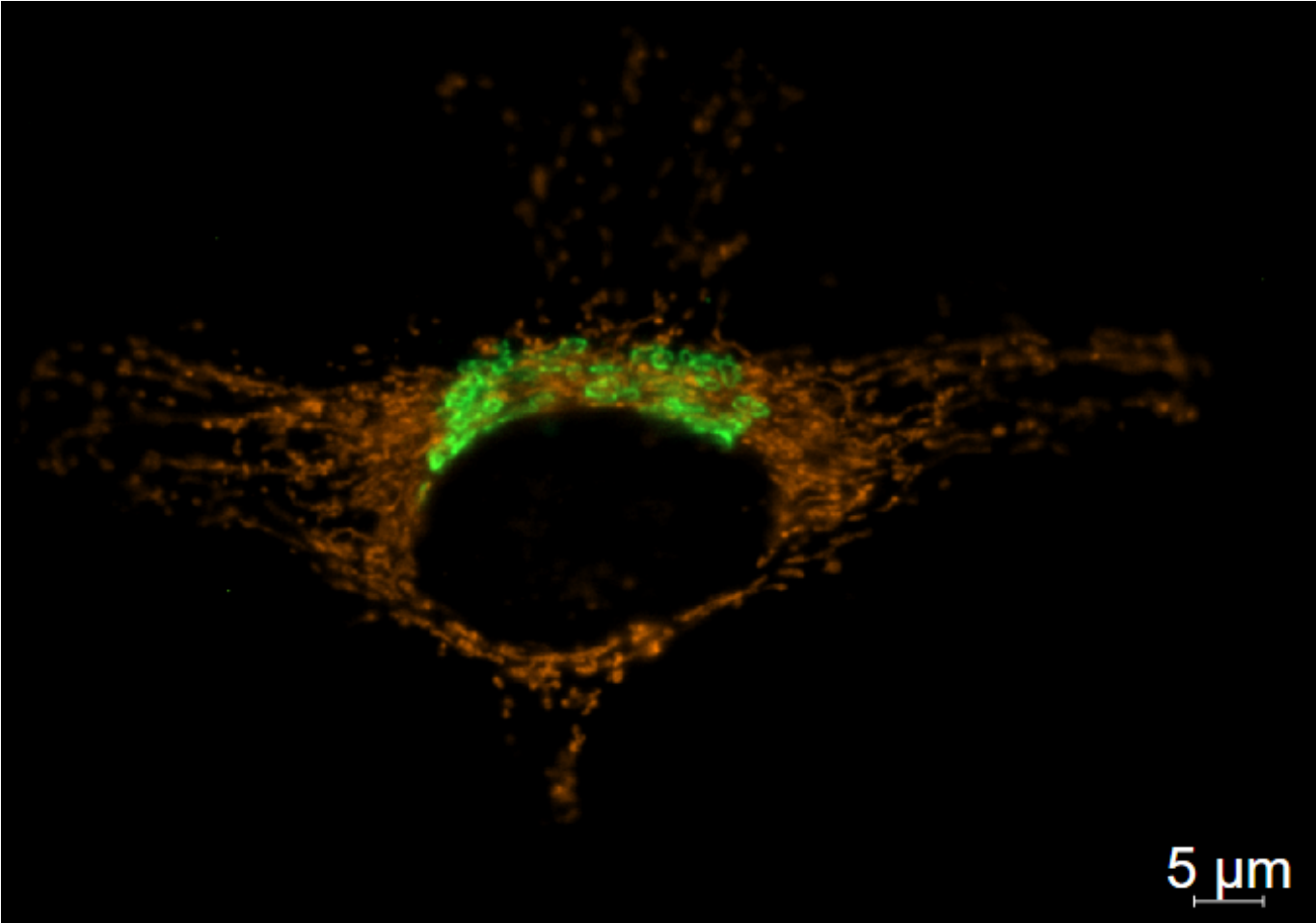

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	28

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	3
Scaling (Per Pixel)	0.102µm X 0.102µm
Image size (Pixels)	968 X 793
Image Size (Scaled)	142.10µm X 106.48µm
Bit Depth	12 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	38 HE eGFP AF532-49014
Beam Splitter	495 550
Filter Excitation Wavelength	450 - 490 515 - 545
Filter Emission Wavelength	500 - 550 555 - 595
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @20%
Exposure Time	100ms
Excitation Wavelength	493 553
Emission Wavelength	517 568
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	0.80µm 0.88µm
Binning Mode	1,1

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Donkey Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000042

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017)), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 3

Scaling (Per Pixel) = 0.102µm X 0.102µm

Image size (Pixels) = 968 X 793

Image Size (Scaled) = 142.10µm X 106.48µm

Bit Depth = 12 Bit

Image Center Position = X: 0.00µm, Y: 0.00µm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

488 -

Reflector = 38 HE eGFP

Beam Splitter = 495

Filter Excitation Wavelength = 450 - 490

Filter Emission Wavelength = 500 - 550

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @20%

Exposure Time = 100ms

Excitation Wavelength = 493

Emission Wavelength = 517

Effective NA = 1.4

Imaging Device = Axiocam MR R3

Camera Adapter = 1X

Depth of Focus = 0.80µm

Binning Mode = 1,1

555 -

Reflector = AF532-49014
 Beam Splitter = 550
 Filter Excitation Wavelength = 515 - 545
 Filter Emission Wavelength = 555 - 595
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @20%
 Exposure Time = 100ms
 Excitation Wavelength = 553
 Emission Wavelength = 568
 Effective NA = 1.4
 Imaging Device = AxioCam MR R3
 Camera Adapter = 1X
 Depth of Focus = 0.88µm
 Binning Mode = 1,1

Post Processing

None

Image Correction


None

Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000018
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	28
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	09E31FEE37ED6E71B5E27FAD051D03CC6384B3AD40C30BBB8629CAC7246B312E
Method PDF SHA-256	BCDC1E01BD7854D0262B9016F45145583ED1F15FFF40620CE6CCBF800E26281A

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

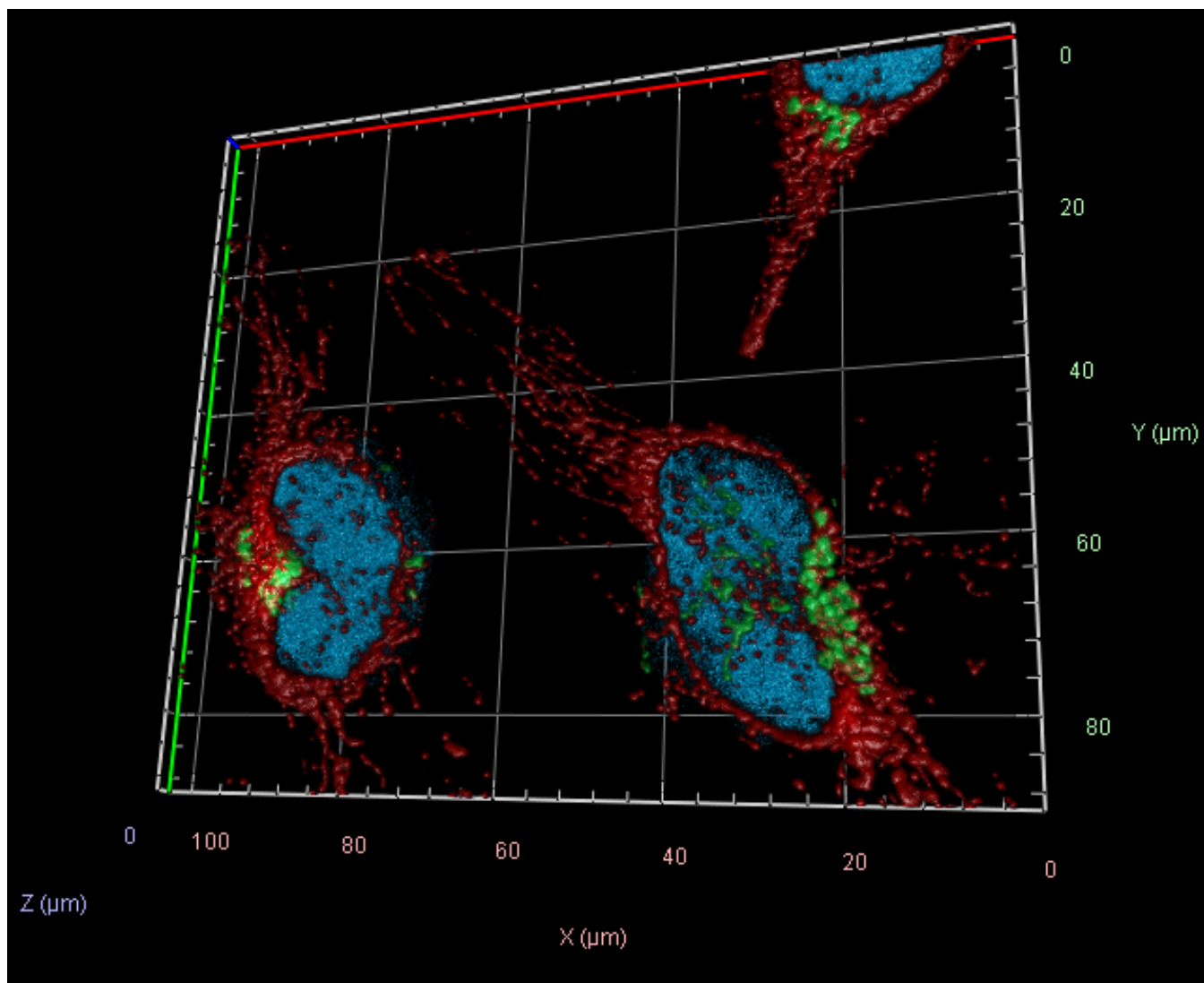

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 555
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Z-Stack	14 Slices (2.6µm)
Channels	3
Scaling (Per Pixel)	0.143µm X 0.143µm X 0.200µm
Image Size (Pixels)	723 X 628
Image Size (Scaled)	103.29µm X 89.71µm
Bit Depth	14 Bit
Image Center Position	X: 50.11mm, Y: 50.32mm
Stage Position	X: 50.11mm, Y: 50.32mm
ROI Center Offset	X: 1.14, Y: 0.00µm
Reflector	38 HE eGFP AF532-49014 49 DAPI
Beam Splitter	495 550 395
Filter Excitation Wavelength	450 - 490 515 - 545 335-383
Filter Emission Wavelength	500 - 550 555 - 595 420-470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @20% X-Cite Xylis Lamp Intensity @5%
Exposure Time	100ms 75ms 10ms
Excitation Wavelength	493 553 353
Emission Wavelength	517 568 465
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.00µm 1.10µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.
3D Volume Projection	Precise
Appearance	Threshold 1% (595), Threshold 3% (488), Threshold 14% (DAPI), Ramp 50% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness 100%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Donkey Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000042

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 555 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 14 Slices (2.6µm)

Channels = 3

Scaling (Per Pixel) = 0.143µm X 0.143µm X 0.200µm

Image Size (Pixels) = 723 X 628

Image Size (Scaled) = 103.29µm X 89.71µm

Bit Depth = 14 Bit

Image Center Position = X: 50.11mm, Y: 50.32mm

Stage Position = X: 50.11mm, Y: 50.32mm

ROI Center Offset = X: 1.14, Y: 0.00µm

488 -

Reflector = 38 HE eGFP
 Beam Splitter = 495
 Filter Excitation Wavelength = 450 - 490
 Filter Emission Wavelength = 500 - 550
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @20%
 Exposure Time = 100ms
 Excitation Wavelength = 493
 Emission Wavelength = 517
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 1.00µm
 Binning Mode = 2,2

555 -

Reflector = AF532-49014
 Beam Splitter = 550
 Filter Excitation Wavelength = 515 - 545
 Filter Emission Wavelength = 555 - 595
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @20%
 Exposure Time = 75ms
 Excitation Wavelength = 553
 Emission Wavelength = 568
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 1.10µm
 Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335-383
 Filter Emission Wavelength = 420-470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @5%
 Exposure Time = 10ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Correction, Fourier Filter, and Deconvolution were NOT used.

3D Image Processing

3D Volume Projection = Precise

Appearance = Threshold 1% (595), Threshold 3% (488), Threshold 14% (DAPI), Ramp 50% all channels, Maximum 100% all channels.

Background = Black

Light = Brightness 100%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping

Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	EAAF94A97310EEDA75AC6328565DA65C607DEC96B6F9B3C658B03748BF716163
Method PDF SHA-256	81542F89EA2BCE44C6DF197D5CBC9EDEC71AEDDC296FB4879F4E9B4A7D77742B

ORIGINAL IMAGE EVIDENCE / CONFOCAL

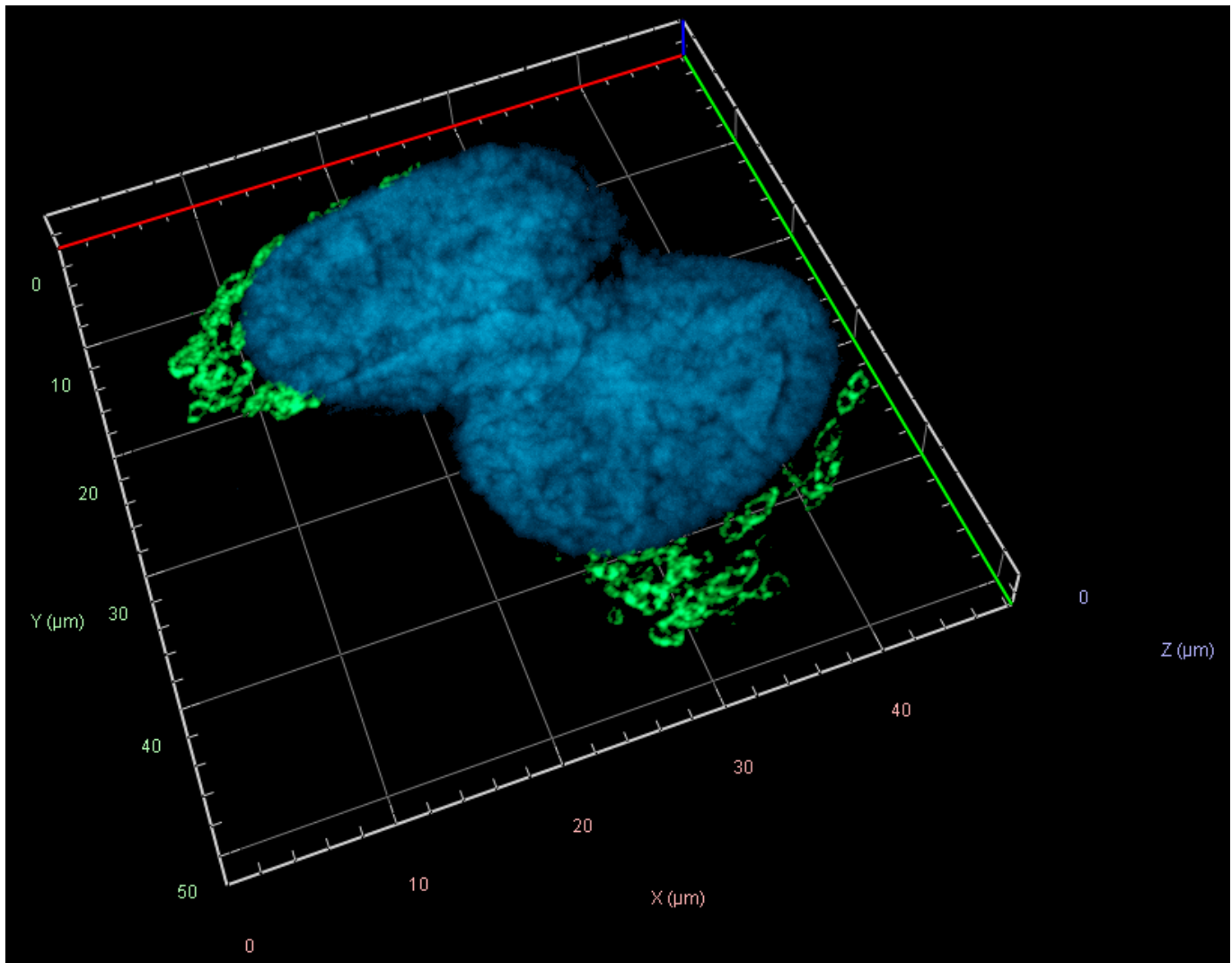

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488
METADATA VALUES	45

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	56 Slices (3.3µm)
Channels	2
Scaling (Per Pixel)	0.071µm X 0.071µm X 0.060µm
Image Size (Pixels)	684 X 734
Image Size (Scaled)	48.29µm X 51.82µm
Bit Depth	16 Bit
Image Center Position	X: 64.70mm, Y: 46.85mm
Stage Position	X: 64.70mm, Y: 46.85mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 44µm 1.00 AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @0.2% 405nm @0.3%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.2
Rotation	0°
Sampling	1.00
Pixel Time	0.65µs
Frame Time	17.44s
LSM Scan Speed	8
Scan Direction	Bidirectional
Line Step	1
Excitation Wavelength	493 353
Emission Wavelength	517 465
Effective NA	1.4
Detection Wavelength	498 - 556 410-470
Imaging Device	Airyscan GaAsp-PMT 1
Detector Type	GaAsp-PMT
Detector Gain	650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
Averaging	16
3D Maximum Projection	Precise
Appearance	Threshold 8% (488), Threshold 4% (DAPI), Ramp0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness110%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack =56 Slices (3.3µm)

Channels = 2

Scaling (Per Pixel) = 0.071µm X 0.071µm X 0.060µm

Image Size (Pixels) = 684 X 734

Image Size (Scaled) = 48.29µm X 51.82µm

Bit Depth = 16 Bit

Image Center Position = X: 64.70mm, Y: 46.85mm

Stage Position =X: 64.70mm, Y: 46.85mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

488 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 44µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @0.2%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.2
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.65µs
 Frame Time = 17.44s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 16
 Excitation Wavelength = 493
 Emission Wavelength = 517
 Effective NA = 1.4
 Detection Wavelength = 498 - 556
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector - None
 Contrast Method - Fluorescence
 Pinhole = 1.00 AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.3%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.2
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.65µs
 Frame Time = 17.44s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 16
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410-470
 Imaging Device = GaAsP-PMT 1
 Detector Type = GaAsP-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 8% (488), Threshold 4% (DAPI), Ramp0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness110%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000018
Stage	Confocal

SOURCE PROVENANCE

Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	45
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	99BD3D22E3100FE649C1A78AA544F8932EC3CFDA0F07018D59B61E5B8717E4AE
Method PDF SHA-256	B85944A3B80DB59EC63CF6CD80F0C30A87BC15080D530F66274702BEB1608C05

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

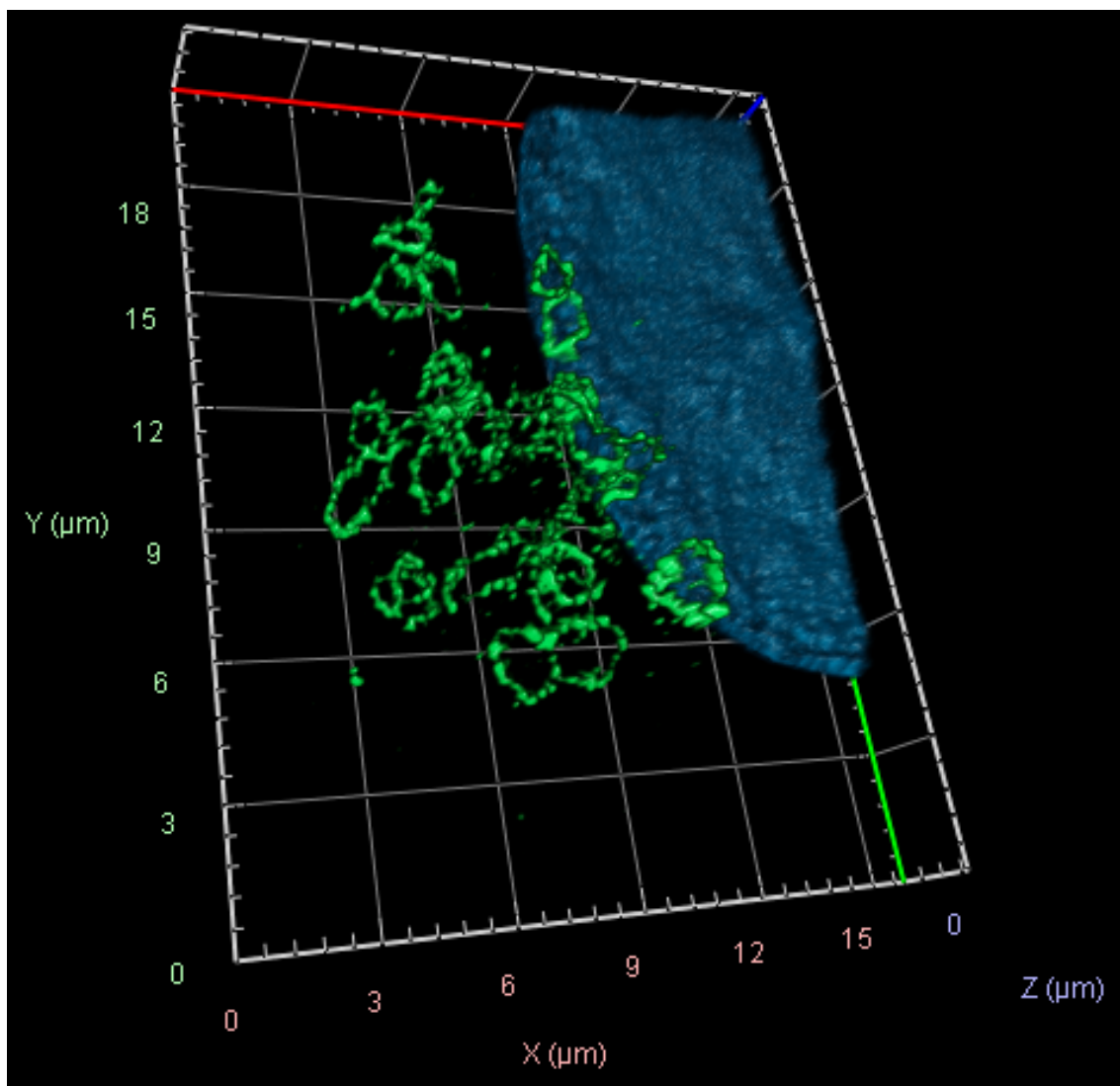

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	39 Slices (1.9µm)
Channels	2
Scaling (Per Pixel)	0.03530 µm X 0.03530 µm X 0.05000 µm
Image Size (Pixels)	449 X 591
Image Size (Scaled)	15.85µm X 20.86µm
Bit Depth	16 Bit
Image Center Position	X: 72.18mm, Y: 54.41mm
Stage Position	X: 72.18mm, Y: 54.41mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00AU / 220µm 5.00 AU / 181µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @0.09% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.3
Rotation	0°
Sampling	2.00
Pixel Time	0.96µs
Frame Time	21.88s
LSM Scan Speed	6
Scan Direction	Bidirectional
Line Step	1
Averaging	4
Excitation Wavelength	493 353
Emission Wavelength	517 465
Effective NA	1.4
Detection Wavelength	490 - 560 400-460
Imaging Device	Airyscan
Detector Type	GaAsp-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	800 V
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 10.0 (3D, Auto) Super Resolution 8.6 (3D, Auto)
3D Volume Projection	Precise
Appearance	Threshold 1% (488), Threshold 5% (DAPI), Ramp 98% (488), Ramp 95% (DAPI), Maximum 100% all channels
Background	Black
Light	Brightness 105%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Tone Mapping NOT Enabled
Projection	View Angle 57°, Scale Z 114%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods™

Identifyfyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyfyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 39 Slices (1.9µm)

Channels = 2

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 449 X 591

Image Size (Scaled) = 15.85µm X 20.86µm

Bit Depth = 16 Bit

Image Center Position = X: 72.18mm, Y: 54.41mm

Stage Position = X: 72.18mm, Y: 54.41mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

488 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00AU / 220µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @0.09%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 21.88s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 493
 Emission Wavelength = 517
 Effective NA = 1.4
 Detection Wavelength = 490 - 560
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 10.0 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 181µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 21.88s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400-460
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.6 (3D, Auto)

Post Processing - AiryScan 3D Processing

Display Mode "SR"
 Super Resolution Auto Filter Selected
 Process 3D Current Image
 Created

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 1% (488), Threshold 5% (DAPI), Ramp 98% (488), Ramp 95% (DAPI), Maximum 100% all channels
 Background = Black
 Light = Brightness 105%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Tone Mapping NOT Enabled
 Projection = View Angle 57°, Scale Z 114%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections -

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	BE39AA0A0A46D4D245D0F81C67564E07353CCA2D26266C141C894AE15DA1B9D3
Method PDF SHA-256	E1D8D8CDE63610F95B337F545F5E962309D471DEC551EE8C2132B6FCAF86DA60

ORIGINAL IMAGE EVIDENCE / SIM

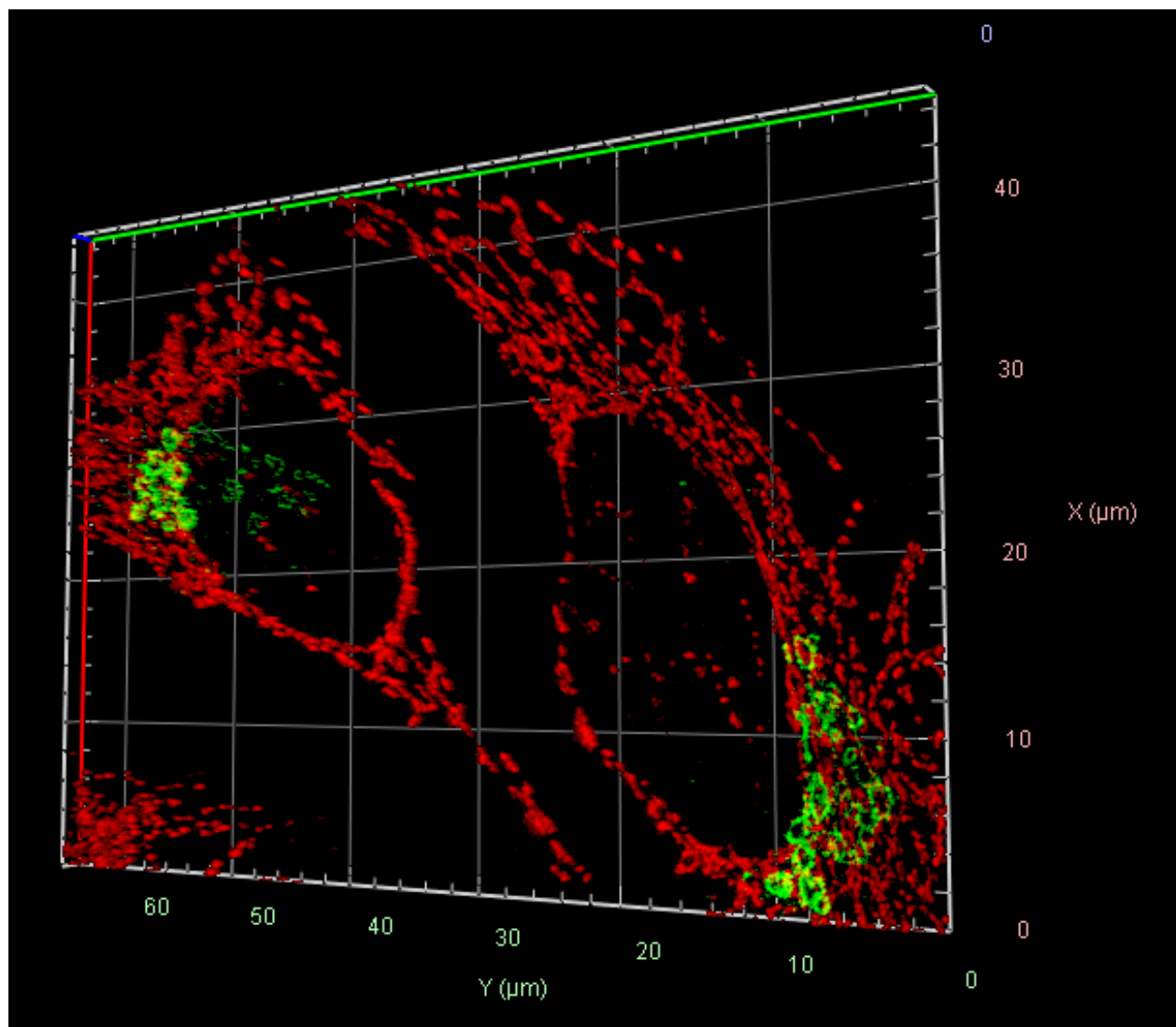

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	56

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 594	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 560 SBS LP 640
Beam Splitter 488	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Z-Stack	22 Slices (1.47µm)
Channels	2
Scaling (Per Pixel)	0.03130 µm X 0.03130 µm X 0.07000 µm
Image Size (Pixels)	1429 X 2055
Image Size (Scaled)	44.73µm X 64.33µm
Bit Depth	16 Bit
Image Center Position	X: 2.36mm, Y: -3.70mm
Stage Position	X: 2.36mm, Y: -3.70mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	561nm @5.00% 488nm @1.50%
Scan Mode	Other
Scan Zoom	1.0
Rotation	0°
Frame Time	1.00s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	4
Excitation Wavelength	561 488
Emission Wavelength	515 655
Effective NA	1.46
Imaging Device	Edge 4.2 SIM
Exposure Time	100ms 50ms
Depth of Focus	0.73µm 0.93µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	561.000000 (594) 488.000000
Grating Period (µm)	27.500000
Processing	3D
Auto Sharpness	Yes
Sharpness	10.612758 11.268022
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No
3D Maximum Projection	Precise
Appearance	Threshold 6% (594), Threshold 8% (488), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 100%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Donkey Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000053

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 594= Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Beam Splitter 488 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 640

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Z Stack - (3D)

Z-Stack = 22 Slices (1.47µm)

Channels = 2

Scaling (Per Pixel) = 31.30nm X 31.30nm X 70.00nm

Image Size (Pixels) = 1429 X 2055

Image Size (Scaled) = 44.73µm X 64.33µm

Bit Depth = 16 Bit

Image Center Position = X: 2.36mm, Y: -3.70mm

Stage Position = X: 2.36mm, Y: -3.70mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Contrast Method = Fluorescence
 Laser Wavelength = 561nm @5.00%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.00s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 561
 Emission Wavelength = 515
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 100ms
 Depth of Focus = 0.73µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM Processing

Channel = 561.000000 (594)

Grating Period (µm): 27.500000

SIM

Processing: 3D
 Auto Sharpness: Yes
 Sharpness: 10.612758
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

488 -

Contrast Method = Fluorescence
 Laser Wavelength = 488nm @1.50%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.00s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 488
 Emission Wavelength = 655
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.93µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM Processing

Channel = 488.000000

Grating Period (µm): 27.500000

SIM

Processing: 3D
 Auto Sharpness: Yes
 Sharpness: 11.268022
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 6% (594), Threshold 8% (488), Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 100%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	SIM
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	56
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	EC69BF58B9AB094E300CE17DFF6175906BB458E56F7DCC806867521293D17697
Method PDF SHA-256	93E61184E8268A77BBE2C0EF8C42F69CFAF1A70BBC2F5D58C428E1668347FE4F

ORIGINAL IMAGE EVIDENCE / SIM²

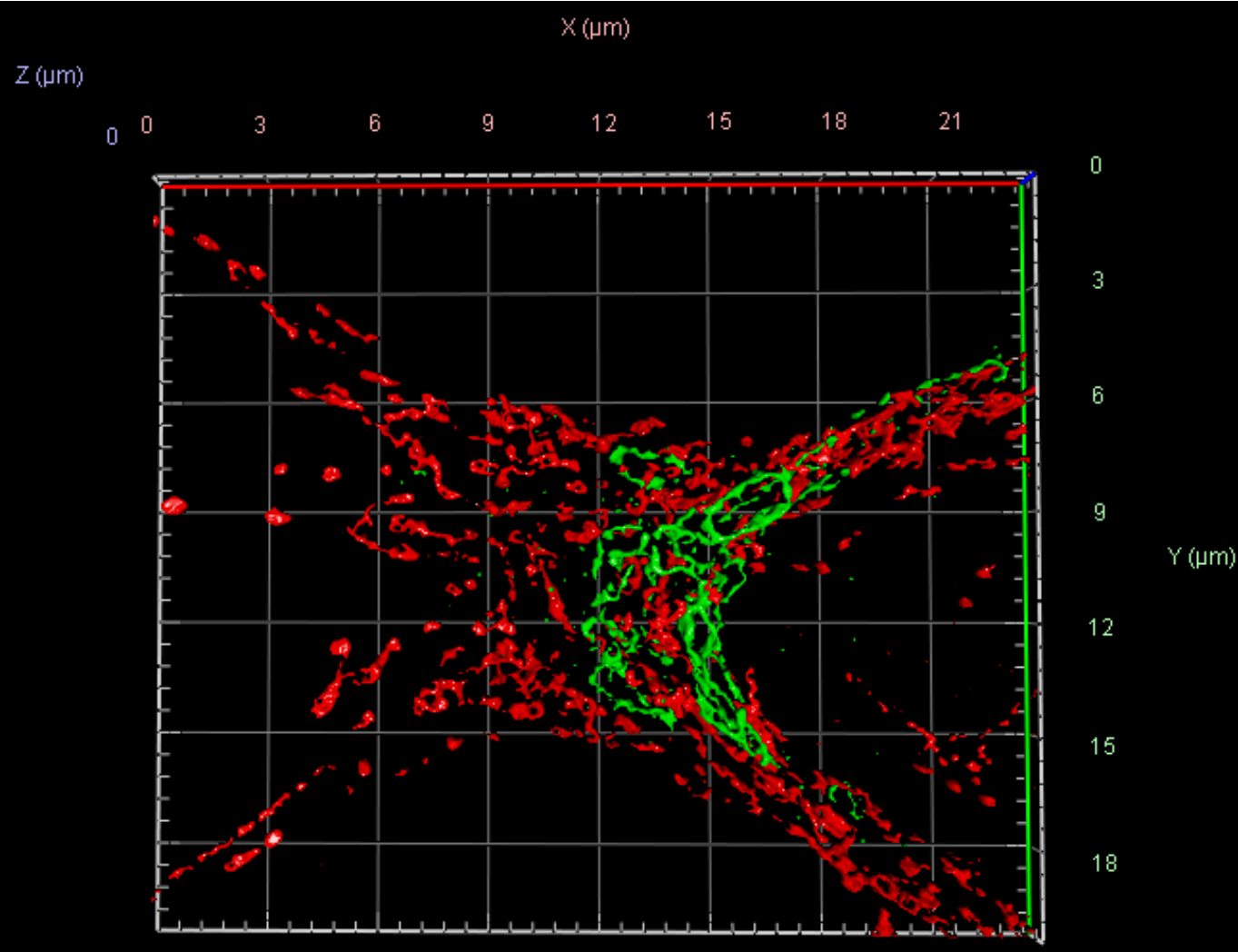

EVIDENCE TIER	SIM ²
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	59

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 594	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 560
Beam Splitter 488	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Z-Stack	20 Slices (1.14µm)
Channels	2
Scaling (Per Pixel)	0.01565 µm X 0.01565 µm X 0.06000 µm
Image Size (Pixels)	1506 X 1301
Image Size (Scaled)	23.57µm X 20.36µm
Bit Depth	16 Bit
Image Center Position	X: 4.68mm, Y: 2.07mm
Stage Position	X: 4.68mm, Y: 2.07mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	488nm @1.50% 561nm @5.0%
Scan Mode	Other
Scan Zoom	1.0
Rotation	0°
Frame Time	1.0s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	4
Excitation Wavelength	488 561
Emission Wavelength	655 515
Effective NA	1.46
Imaging Device	Edge 4.2 SIM
Exposure Time	50ms 100ms
Depth of Focus	0.93µm 0.73µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	488.000000 561.000000 (594)
Grating Period (μm)	27.500000 32.000000
Processing	3D
Input SNR	Medium
Iterations	16
Regularization Weight	0.065000
Processing Sampling	4
Output Sampling	4
Filter	Median
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No
3D Surface Projection	Precise
Appearance	Threshold 20% (594), Threshold 20% (488), Ambient Light 0% all channels, Specular Light 100% all channels, Shininess 30% all channels
Background	Black
Light	Azimuth 80°, Elongation -130°
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000018

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Donkey Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000053

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 594 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Beam Splitter 488 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 560

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Z Stack - (3D)

Z-Stack = 20 Slices (1.14µm)

Channels = 2

Scaling (Per Pixel) = 15.65nm X 15.65nm X 60.00nm

Image Size (Pixels) = 1506 X 1301

Image Size (Scaled) = 23.57µm X 20.36µm

Bit Depth = 16 Bit

Image Center Position = X: 4.68mm, Y: 2.07mm

Stage Position = X: 4.68mm, Y: 2.07mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

488 -

Contrast Method = Fluorescence
 Laser Wavelength = 488nm @1.50%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 488
 Emission Wavelength = 655
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.93µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 488.000000
Grating Period (µm): 27.500000

SIM2

Processing: 3D
 Input SNR: Medium
 Iterations: 16
 Regularization Weight: 0.065000
 Processing Sampling: 4
 Output Sampling: 4
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

594 -

Contrast Method = Fluorescence
 Laser Wavelength = 561nm @5.0%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 561
 Emission Wavelength = 515
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 100ms
 Depth of Focus = 0.73µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 561.000000 (594)

Grating Period (µm): 32.000000

SIM2
 Processing: 3D
 Input SNR: Medium
 Iterations: 16
 Regularization Weight: 0.065000
 Processing Sampling: 4
 Output Sampling: 4
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

3D Image Processing

3D Surface Projection = Precise
 Appearance = Threshold 20% (594), Threshold 20% (488), Ambient Light 0% all channels, Specular Light 100% all channels, Shininess 30% all channels
 Background = Black
 Light = Azimuth 80°, Elongation -130°
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	SIM²
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	59
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	2B5EED1A95198A688B80191D59CC9557B36ABFB966E397BA3DC8EB59C96D8A4D
Method PDF SHA-256	E790A51D08D7B078E40F75C784B5CB64CCE389E424A78599323038E821AE4097

ORIGINAL IMAGE EVIDENCE / SINGLE-MOLECULE STORM

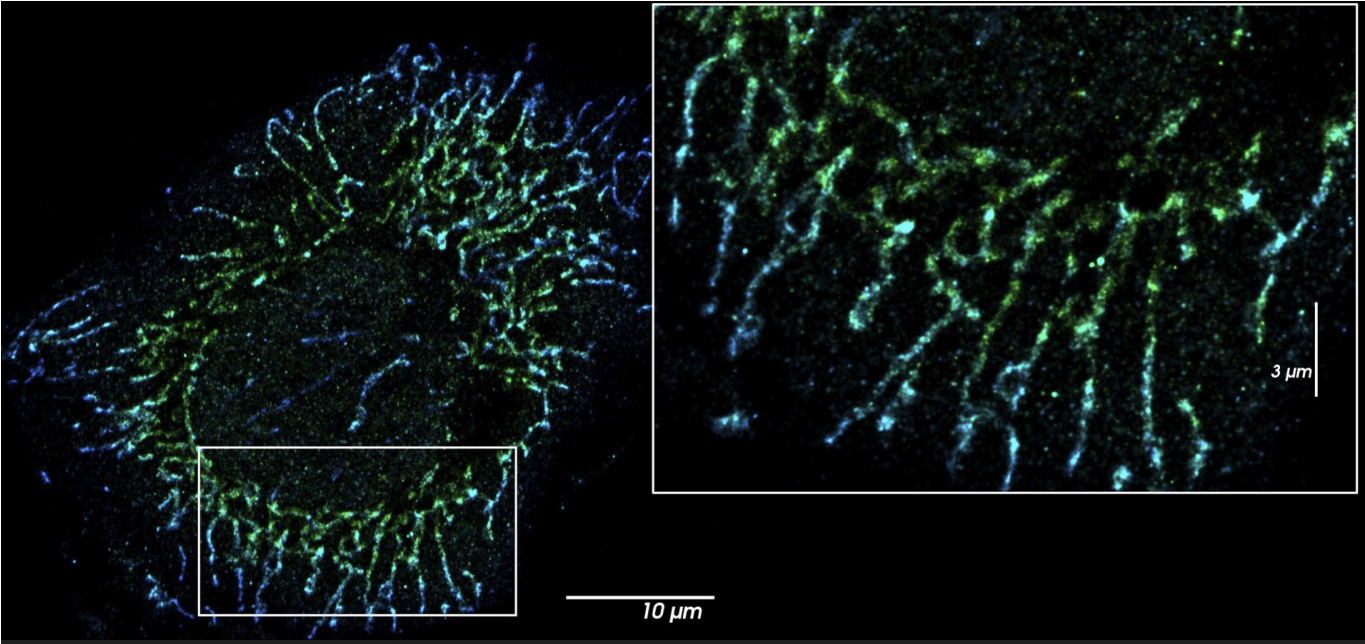

EVIDENCE TIER	Single-molecule STORM
INSTRUMENT	Bruker Vutara VXL
CELL MODEL	HeLa
DETECTED DYES	488
METADATA VALUES	71

STRUCTURED ACQUISITION METADATA / SINGLE-MOLECULE STORM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Objective	Olympus 60X/1.30 silicon oil objective
Number of Probes	1
Number of Simultaneous Probes	1
Number of Reference Probes	1
Camera	Both Cameras
Color Channel	Both Channels
Probe Capture Mode	Interleaved
Z -Stack Capture Mode	Interleaved
Multiple Location Capture	Not selected
Post Record Reference	Not Selected
Fluidics	Not Selected
Autofocus Drift Correction	On
Biplane Module	635nm Dichroic
Objective	Silicon Selected
Stage Insert	Rectangular Selected
Course Focus Position	8-Well Chamber Selected
Dichroic	SuperRes
PSF	(Red, Orange, Green)
Heater	Current: 26
Recording Vertical Field of View	100%
Widefield Camera Exposure Time	100ms 20ms
Super Resolution Camera Exposure Time	20ms
Widefield Camera Binning	1
Per-Probe Exposure Time	Selected
Predefined Probe	First probe: 488 nm (635 nm Dichroic); First reference probe: 488 nm (635 nm Dichroic)
Photoactivation/Imaging Laser	0% Selected
Exposure Time	20ms Selected
Record Workflow	Selected Options
View Mode	Widefield 200µm
Capture Mode	Live
Piezo Current	180 µm

FIELD	SOURCE-DOCUMENTED VALUE
Widefield 488nm	0.1X Neutral Density Filter: Selected, @0.02% Laser Input
Reference Probe 1 Laser control (488 nm (635 nm Dichroic))	Not Selected
Widefield camera Exposure Time	100ms
Widefield Laser	"Widefield 488nm" was deselected, and "calibrate" under Autofocus was selected.
Current Intensity	4.5
Calibration Value	-5084 nm/unit
Autofocus Mode	Pre-Capture
SuperRes	50µm is Selected
Piezo Record	Begin: 179; End: 179
SuperRes 488nm	0.1X Neutral Density Filter: Selected, @0.01% Laser Input
Record Workflow - Record Tab	Capture Parameters
Frames	Probe 1: 20000
Z Positions	1 using Steps: 0.2µm (200nm)
Cycles	1
Time Points	1
Pause Between Cycles	Not Selected
Automatic Photoactivation Laser Control	Not Selected
Laser Power	9%
Cutout Width	16px
Max Frame Accumulation	1 frame
Max Accumulation Offset	1px
Fiducial at Max Accumulation	Not selected
Max Particles Per Frame	No Max
Background SLPE Removal	Not Selected
Cutout/Background Removal	Not Selected
Preview Cutout Positions	Selected
Probe1	BG threshold: 7
Sizing Mode	Radial Precision
Particle size	4nm
Color by	Depth
Color Table	Single
Opacity Mode	Constant; Opacity 1: 05
Front-to-back Attenuation	Not Selected
Method	Particle (direct)
Time Intervals	5 Intervals

FIELD	SOURCE-DOCUMENTED VALUE
Pixel Size	20nm
Smoothing	8.00
PSF Z offset	-500 to +500
Radial precision	35
Axial precision	65

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit, F(ab')₂ (H + L)

Cat ID# - SA 000018

Hela cells were grown on μ-Slide 8 Well high Glass Bottom Ibidi chamber, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Mitochondria Rabbit Recombinant Antibody (MTC02,2860R), was incubated for 1 hour at room temperature at a dilution of 1:100. Cells were washed 3X in Identifyn's™ proprietary Wash Buffer, and Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:100. Cells were washed 5X in Identifyn's™ proprietary Wash Buffer followed by the addition of Identifyn's™ proprietary STORM Buffer to the sample chamber enabling image acquisition.

Microscope

Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:

Calibration -

Calibration follows Identifyn's Calibration Standard Operating Procedure, which can be found here.

<https://identifyn.com/storm-calibration>

Image Acquisition -

SRX Software - Under the 'Single Molecule Localization' tab, 'Record' is selected and contains the following workflow tabs and parameters.

Configuration Tab - Capture Configuration

Number of Probes: 1

Number of Simultaneous Probes: 1

Number of Reference Probes: 1

Camera: Both Cameras

Color Channel: Both Channels

Probe Capture Mode: Interleaved

Z -Stack Capture Mode: Interleaved

Multiple Location Capture: Not selected

Post Record Reference: Not Selected

Fluidics: Not Selected

Autofocus Drift Correction: On

Microscope configuration

Biplane Module: 635nm Dichroic
 Objective: Silicon Selected
 Stage Insert: Rectangular Selected
 Course Focus Position: 8-Well Chamber Selected
 Dichroic: SuperRes
 PSF: (Red, Orange, Green)
 Heater: Current: 26

Camera Configuration

Recording Vertical Field of View: 100%
 Widefield Camera Exposure Time: 100ms
 Super Resolution Camera Exposure Time: 20ms
 Widefield Camera Binning: 1
 Per-Probe Exposure Time: Selected

Probe Configuration

Predefined Probe: First probe: 488 nm (635 nm Dichroic); First reference probe: 488 nm (635 nm Dichroic)
 Photoactivation/Imaging Laser: 0% Selected
 Exposure Time: 20ms Selected

Experiment Configuration

Record Workflow: Selected Options
 View Mode: Widefield 200µm
 Capture Mode: Live

Stage and Piezo Selection

Examination of Sample

Cell(s) or Cellular Structure Localization - The red-bordered box identifies the area of the SuperRes view and defines a suitable cell(s) or cellular structure. The cell(s) or structures are optimally focused.

Piezo Current: 180 µm

Probe 1 Laser control (488 nm (635 nm Dichroic))

Widefield 488nm: 0.1X Neutral Density Filter: Selected, @0.02% Laser Input
 Reference Probe 1 Laser control (488 nm (635 nm Dichroic): Not Selected

Advanced Tools Option

Widefield camera Exposure Time: 100ms
 Default values are maintained unless otherwise noted.

Widefield Laser: "Widefield 488nm" was deselected, and "calibrate" under Autofocus was selected.

Autofocus - Calibration Values

Current Position" -0.979

Current Intensity: 4.5

Calibration Value: -5084 nm/unit

Autofocus Mode: Pre-Capture

View Mode

SuperRes: 50µm is Selected

Capture Mode

Live

Stage and Piezo Selection

Cell(s) or Cellular Structure Localization - Positioning of the target Cell or Subsection of the Target Cell is optimized and optimally focused.

Piezo Record: Begin: 179; End: 179

Probe 1 Laser control (488 nm (635 nm Dichroic)

SuperRes 488nm: 0.1X Neutral Density Filter: Selected, @0.01% Laser Input

Reference Probe 1 Laser control (488 nm (635 nm Dichroic)

Not Selected

Advanced Tools

Widefield Camera Exposure Time: 20ms

Default values are maintained unless otherwise noted.

Record Workflow - Record Tab: Capture Parameters

Frames: Probe 1: 20000

Z Positions: 1 using Steps: 0.2µm (200nm)

Cycles: 1

Time Points: 1

Pause Between Cycles: Not Selected

Automatic Photoactivation Laser Control: Not Selected

Probe 1 Laser control (488 nm (635 nm Dichroic)

Laser Power: 9%

Expert Configuration Options - The Localization Tab

Beyond the values shown, all other parameters were default values.

Localization

Cutout Width: 16px
 Max Frame Accumulation: 1 frame
 Max Accumulation Offset: 1px
 Fiducial at Max Accumulation: Not selected
 Max Particles Per Frame: No Max
 Background SLPE Removal: Not Selected
 Cutout/Background Removal: Not Selected
 Preview Cutout Positions: Selected

Record Workflow - Localization Tab

Data Collection

The folder is created automatically with the localization file for this image acquisition.

Recording Summary

All parameters were carried over from previous tabs and are kept except:

Probe1: BG threshold: 7

Region of Interest -

Not selected

Display

Was deployed to inspect if the BG threshold applied chooses localizations in an expected manner

Record Workflow - View Tab

Region of Interest

Default values

Visualization Parameters

Point Splatting
 Sizing Mode: Radial Precision
 Particle size: 4nm
 Color by: Depth
 Color Table: Single
 Opacity Mode: Constant; Opacity 1: 05
 Front-to-back Attenuation: Not Selected

Post Acquisition analysis - View Tab

Drift Correction - Computed with the following parameters:

Method: Particle (direct)
 Time Intervals: 5 Intervals
 Pixel Size: 20nm
 Smoothing: 8.00

Advanced Particle Filters - Applied to the drift-corrected data with the following parameters

PSF Z offset: -500 to +500

Radial precision: 35

Axial precision: 65

The other parameters were not applied.

Statistical analysis - N/A


Image by S. Thakur

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	Single-molecule STORM
Instrument	Bruker Vutara VXL
Microscope configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Structured source metadata values	71
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	F49938C2C7D177B9CEC63D59EA2CDE5C2C6C6FA3CFDDF9E52053523B99EFDDF5
Method PDF SHA-256	AD36193EAE580A9820EDFB18C64F0DF35BF5338D1127A1D584D9242D10D87B2D

ORIGINAL IMAGE EVIDENCE / CONTROL

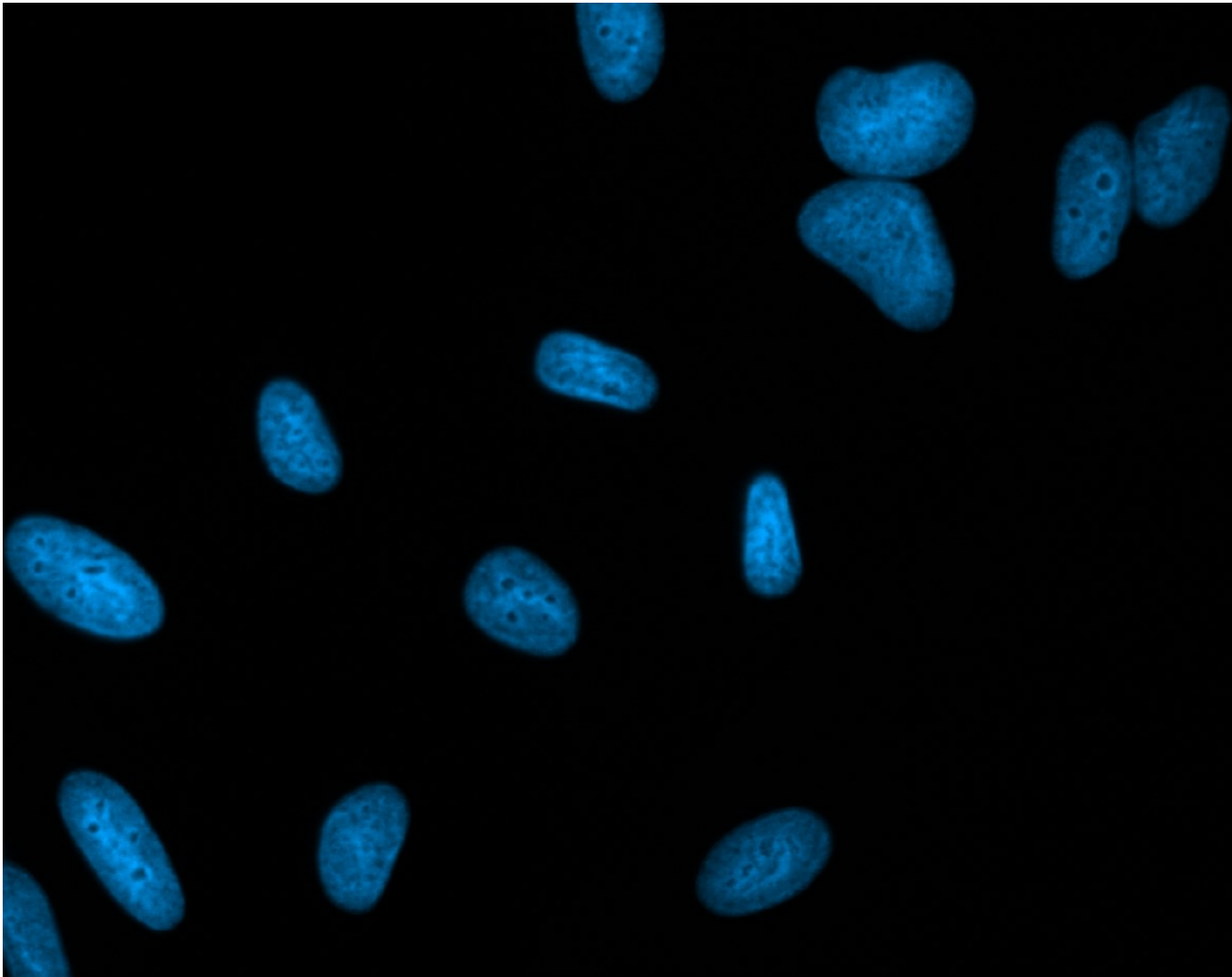

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Imager.M1; ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488
METADATA VALUES	30

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1366 X 1080
Image Size (Scaled)	195.14µm X 154.29µm
Bit Depth	14 Bit
Image Center Position	X: 48.83mm, Y: 50.26mm
ROI Center Offset	X: 48.83µm, Y: 50.26µm
Reflector	38 HE eGFP 49 DAPI
Beam Splitter	495 395
Filter Excitation Wavelength	450 - 490 335-383
Filter Emission Wavelength	500 -550 420-470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	493 353
Emission Wavelength	517 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.0µm 0.90µm
Binning Mode	2,2

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Goat Anti-Rabbit, F(ab')₂ (H + L)

Cat ID# - SA 000018

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1366 X 1080

Image Size (Scaled) = 195.14µm X 154.29µm

Bit Depth = 14 Bit

Image Center Position = X: 48.83mm, Y: 50.26mm

ROI Center Offset = X: 48.83µm, Y: 50.26µm

488 -

Reflector = 38 HE eGFP

Beam Splitter = 495

Filter Excitation Wavelength = 450 - 490

Filter Emission Wavelength = 500 -550

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 493

Emission Wavelength = 517

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.0µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335-383
 Filter Emission Wavelength = 420-470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing -

None

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Axio Imager.M1 Widefield. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The data suggests the secondary antibody is highly specific and free of unconjugated (free) dye.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000018
Stage	Control
Instrument	ZEISS Axio Imager.M1; ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	30
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa

SOURCE PROVENANCE

Image SHA-256	0F04B0CD25D42FEAAC42A9AF8A46296EA4CF35E77C2633728A87146411B45F74
Method PDF SHA-256	EA1254C5422833900E017EDC8FF406BADBA6DFA28A3D60A633A533C264962488