

IDENTIFYN RESOLUTION STEPDOWN

ATR.P

MICROSCOPY EVIDENCE ASSESSMENT

Rabbit anti-Human ATR.p (Threonine 1985) Recombinant Monoclonal - Identifyn® SRM Alexa Fluor™ 647

WESTERN BLOT

1


ORDERED EVIDENCE TIERS

0

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

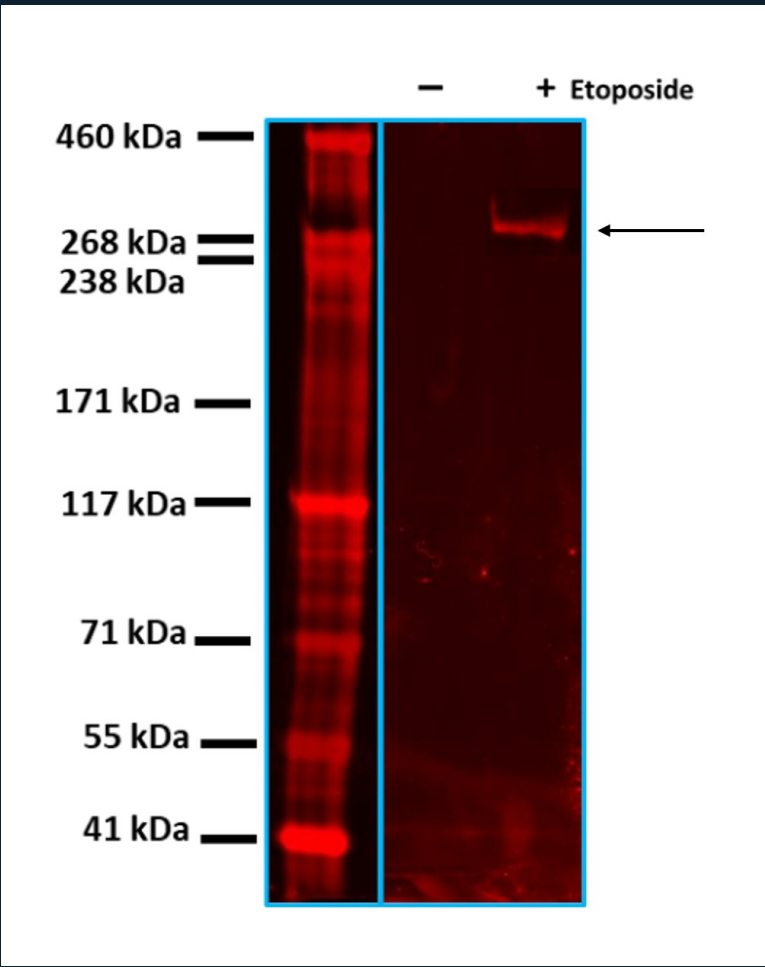

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / WESTERN BLOT

DC-000081 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Treatment-response experimental workflow

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for ATR.p, and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot

BENNETT ASSESSMENT

The preserved DC-000081 record contains one available modality (Western blot). No cross-modality synthesis is possible, and missing stages are not inferred.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / LIMITED SINGLE-MODALITY EVIDENCE

The preserved DC-000081 record contains one available modality (Western blot). No cross-modality synthesis is possible, and missing stages are not inferred.

SETTINGS ASSESSMENT

The record preserves 1 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

No separate control record is preserved; none is inferred from method or template language.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	647	Not deterministically stated in the acquisition block

MODALITY ASSESSMENT / PART 1 OF 1

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

No cross-modality synthesis is possible because one evidence modality is preserved.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - SCIENTIST-AUTHORIZED CORRECTION

Cross-tier / treatment concentration: 200 μ M -> Etoposide = 2 μ M. The governing scientist identified every 200 μ M occurrence as a technician transcription error. The canonical historical concentration is 2 μ M. Supporting evidence: Scientist-authorized correction supplied for the Identifyn historical archive; original technician text remains preserved unchanged. Authority: Dr. Brian T. Bennett. Preservation: Original source text and hashes remain unchanged; only the derivative canonical field is corrected.

HIGH CONFIDENCE - SOURCE-DERIVED METADATA CANONICALIZATION

Product identity / phosphosite: Blank structured canonical field -> Threonine 1985. The phosphorylation site is explicitly stated in the preserved source product identity. Supporting evidence: Rabbit anti-Human ATR.p (Threonine 1985) Recombinant Monoclonal - Identifyn® SRM Alexa Fluor™ 647

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked
Stepdown order. Scientific wording is preserved; missing modalities are not invented.
No separate control record is present in the preserved product record.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for DC-000081. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. No separate control record is present in the preserved product record, and none is inferred.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

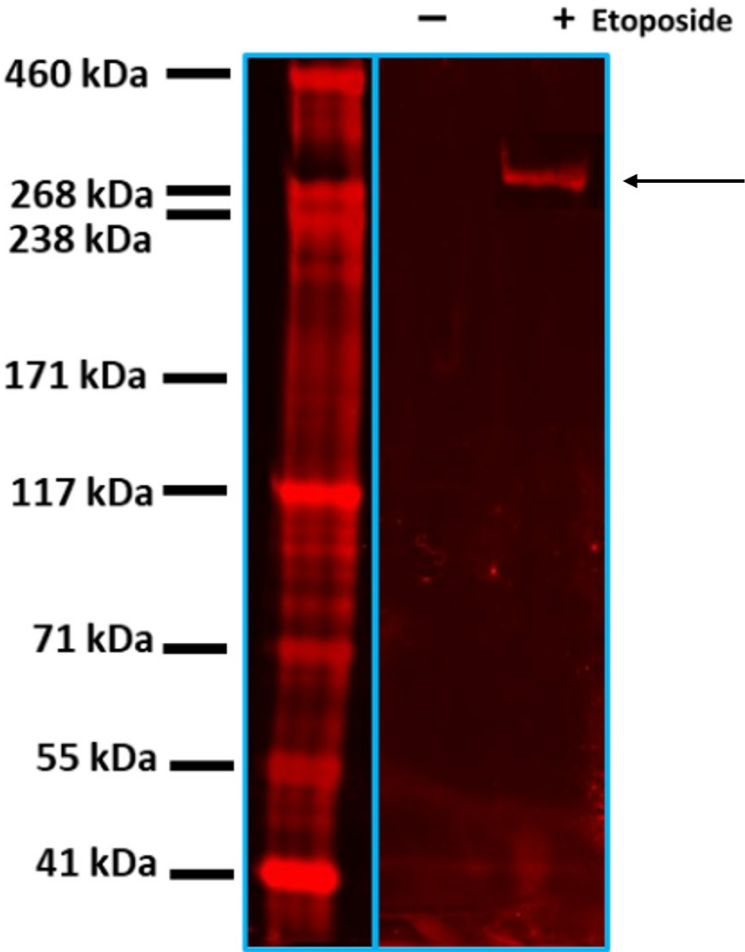

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	647
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	638nm
Emission Filter	676±37nm
Detector	PMT
Intensity	3
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	S. Khanal

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® Primary Antibody

Rabbit anti-Human Phospho- ATR Recombinant Monoclonal-Identifyn® SRM Alexa Fluor™ 647 IgG - (BLR032P)
Cat ID# - DC 000081

Expected MW: 301 kD

HeLa cells were treated with 200 µM etoposide for 20 hours to induce DNA damage. Cells were homogenized using an ultrasonic probe and a nuclear lysis buffer, then centrifuged to pellet the cell debris. The lysate supernatant (100 µg) was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80 °C for 10 min, and subjected to electrophoresis using NuPAGE™ Tris-Acetate Mini Protein Gels, 3 to 8%, 1.0 mm in a Mini Gel Tank. Thermo Fisher's HiMark™ Pre-stained Protein Standard was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using the iBlot™ 3 Starter Kit PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour, then washed 5 times with PBS. The blocked membrane was incubated with Rabbit anti-Human Phospho-ATR Recombinant Monoclonal-Identifyn® SRM Alexa Fluor™ 647 IgG (BLR032P) for 20 hours, followed by five washes in 1X PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 638nm

Emission Filter = 676±37nm

Detector = PMT

Intensity = 3

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: S. Khanal

© Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	DC-000081
Stage	Western blot
Instrument	Azure Sapphire FL
Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	97086330D835F9381F4E44B858E53AF8D3666424ED5193DD36E2CB19BD8B59B3
Method PDF SHA-256	6FDF8077CDD18A9A6F12FE3EFDBA0482AACD1BE14CAEE5B835170F2593397BA2