

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse, F(ab')₂ (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM²

7


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

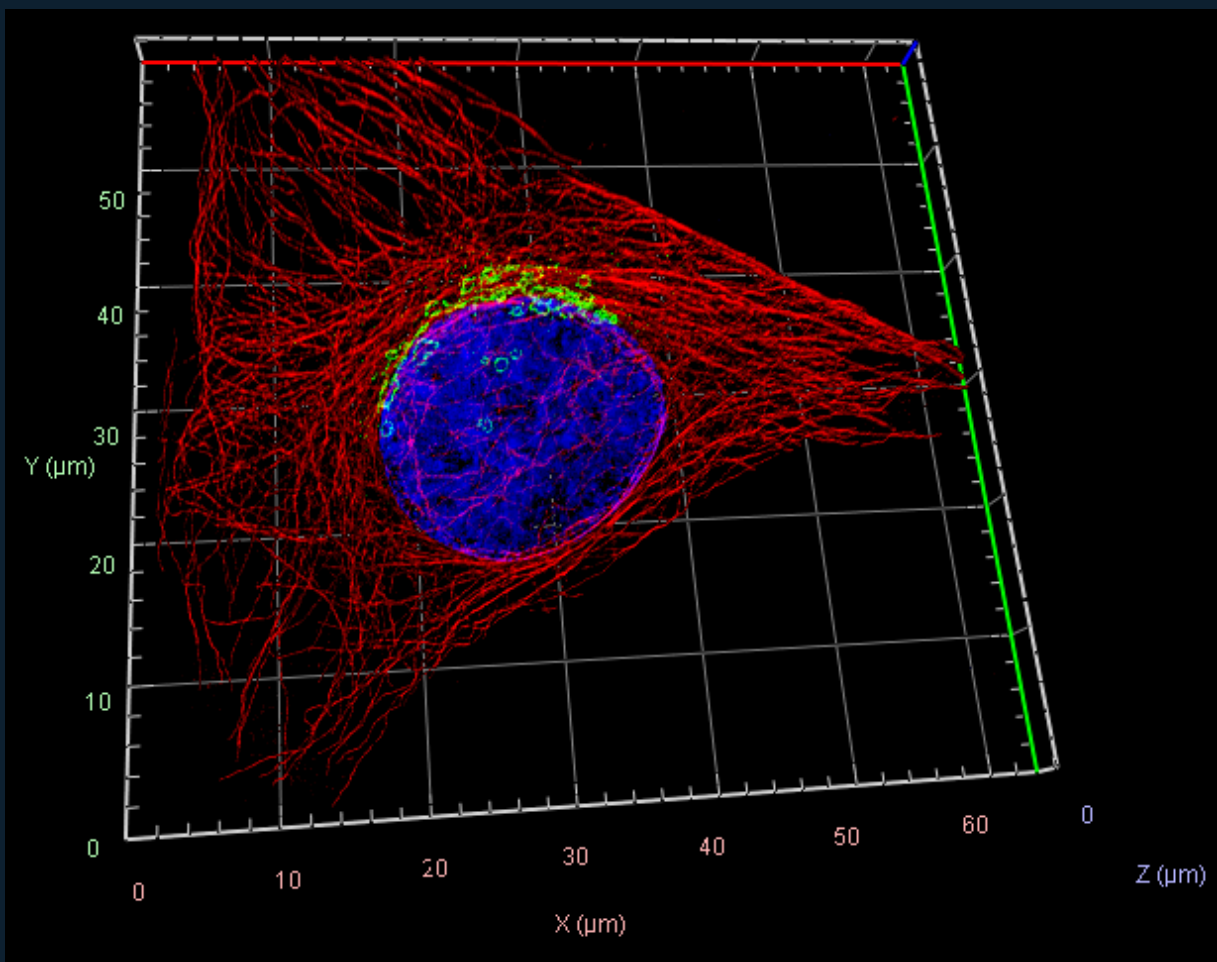

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SIM²

SA-000055 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse, F(ab')₂ (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Control

BENNETT ASSESSMENT

The preserved SA-000055 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / PARTIAL STEPDOWN / SOURCE-BOUNDED REVIEW

The preserved SA-000055 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

SETTINGS ASSESSMENT

The record preserves 7 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	594	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 488; 594	2; 45 Texas Red; 38 eGFP; 540 - 580; 450 - 490; 593 - 668; 500 - 550; 590
Widefield SIM — Apotome	405/DAPI; 594	2; 45 Texas Red; 49 DAPI; 540 - 580; 335 - 383; 593 - 668; 420 - 470; 590
Confocal	405/DAPI; 488; 594	2; None; 561nm @0.09%; 405nm @0.7%; 590; 353; 618; 465
Airyscan	405/DAPI; 488; 594	2; None; 561nm @1.70%; 405nm @1.0%; 590; 353; 618; 465
SIM	405/DAPI; 488; 594	3; 1; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T2-Axiocam 820 #2-SR; T2-Axiocam 820 -SR; 561
SIM²	405/DAPI; 488; 594	3; 1; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T2-Axiocam 820 #2-SR; T2-Axiocam 820 -SR; 561
Control	405/DAPI; 594	2; 45 Texas Red; 49 DAPI; 540 - 580; 335 - 383; 593 - 668; 420 - 470; 590

MODALITY ASSESSMENT / PART 1 OF 8

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 8

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.102µm X 0.102µm; Image size (Pixels): 1388 X 1040; Image Size (Pixels): 1388 X 1040; Bit Depth: 12 Bit. Detector/camera: Imaging Device: Axiocam MR R3. Timing: Exposure Time: 100ms; 200ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @40%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 29.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 8

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Z-Stack: 63 Slices (6.20µm); Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm X .0100µm; Image size (Pixels): 346 X 376; Image Size (Pixels): 346 X 376; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807. Timing: Exposure Time: 200ms; 15ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @18%; X-Cite Xylis Lamp Intensity @8%. Optical/scan: Effective NA: 1.25. Processing: Projection: View Angle 34°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: White. Structured source metadata values: 38.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 8

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.071m X 0.071µm X 0.060µm; Image size (Pixels): 1037 X 1431; Image Size (Pixels): 1037 X 1431; Bit Depth: 16 Bit. Detector/camera: Imaging Device: GaAsP-Pmt2; GaAsP-Pmt1; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.37µs; Frame Time: 3.55s. Illumination: Laser Wavelength: 561nm @0.09%; 405nm @0.7%. Optical/scan: Pinhole: 1.00AU / 54µm; 1.00AU / 37µm; Scan Zoom: .080; 1.0; LSM Scan Speed: 9; Scan Direction: Bidirectional; Averaging: 4; Effective NA: 1.4. Processing: Projection: View Angle 69°, Scale Z 150%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 8

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 109 Slices (5.40µm); Channels: 2; Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm; Image size (Pixels): 1382 X 1410; Image Size (Pixels): 1382 X 1410; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.96µs; Frame Time: 21.88s. Illumination: Laser Wavelength: 561nm @1.70%; 405nm @1.0%. Optical/scan: Pinhole: 5.00 AU / 251µm; 5.00 AU / 188µm; Scan Zoom: 1.3; LSM Scan Speed: 6; Scan Direction: Bidirectional; Averaging: 4; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 9.4 (3D, Auto); Super Resolution 8.0 (3D, Auto). Structured source metadata values: 47.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 8

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 70 Slices (3.45µm); Channels: 3; 1; Scaling (Per Pixel): 21.11nm X 21.11nm X 50.00nm; Image size (Pixels): 3744 X 3779; Image Size (Pixels): 3744 X 3779; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820 #2; AxioCam 820. Timing: Exposure Time: 120ms; 100ms; Frame Time: 1.33s; 676.00ms. Illumination: Laser Wavelength: 561nm @3.50%; 488nm @1.0%; Illumination Wavelength: 561; 488. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM; Process Sampling: 2; Result Sampling: 4; Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 63.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 8

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 45 Slices (2.20µm); Channels: 3; 1; Scaling (Per Pixel): 16.89nm X 16.89nm X 50.00nm; Image size (Pixels): 3761 X 3542; Image Size (Pixels): 3761 X 3542; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Axiocam 820 #2; Axiocam 820. Timing: Exposure Time: 120ms; 100ms; Frame Time: 1.33s; 676.00ms. Illumination: Laser Wavelength: 561nm @3.50%; 488nm @1.0%; Illumination Wavelength: 561; 488. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM2; Process Sampling: 4; Result Sampling: 4; Projection: View Angle 45°, Scale Z 141%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 63.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 8

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1078 X 823; Image Size (Pixels): 1078 X 823; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 594.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000055 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm -> 0.03530 μ m X 0.03530 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1382 X 1410; Image Size (Scaled)=48.78 μ m X 49.77 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 21.11nm X 21.11nm X 50.00nm -> 0.02111 μ m X 0.02111 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=3744 X 3779; Image Size (Scaled)=79.05 μ m X 79.78 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.02%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 16.89nm X 16.89nm X 50.00nm -> 0.01689 μ m X 0.01689 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=3761 X 3542; Image Size (Scaled)=63.52 μ m X 59.82 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked Stepdown order. Scientific wording is preserved; missing modalities are not invented. A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000055. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra	PRESERVED
7	SIM ²	ZEISS Elyra	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

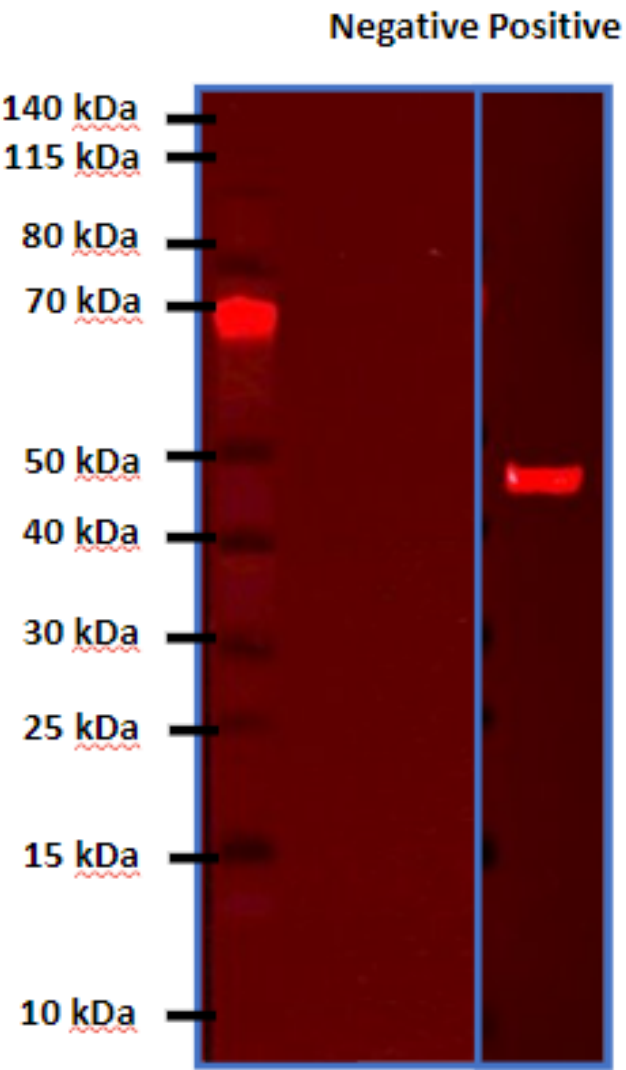

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	594
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	532nm
Emission Filter	572±28nm
Detector	Avalanche Photodiode
Intensity	6
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods™

Identifyfyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Actin Monoclonal Antibody (mAbGEa) (MA1-744) for 2 hours, followed by 5X PBS washes. Identifyfyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit anti-Mouse F(ab)₂ was incubated for 1 hour and washed 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 532nm

Emission Filter = 572±28nm

Detector = Avalanche Photodiode

Intensity = 6

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: Identifyfyn® LLC

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000055
Stage	Western blot
Instrument	Azure Sapphire FL

SOURCE PROVENANCE

Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	CFEA26375E9B5220582092201594FF9E068AA54C262FD13E872387BC7272CF9A
Method PDF SHA-256	2830556E7C4A067618AA287886F8092DB4F2E4CA4D50E3FD9E566F32630E512B

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

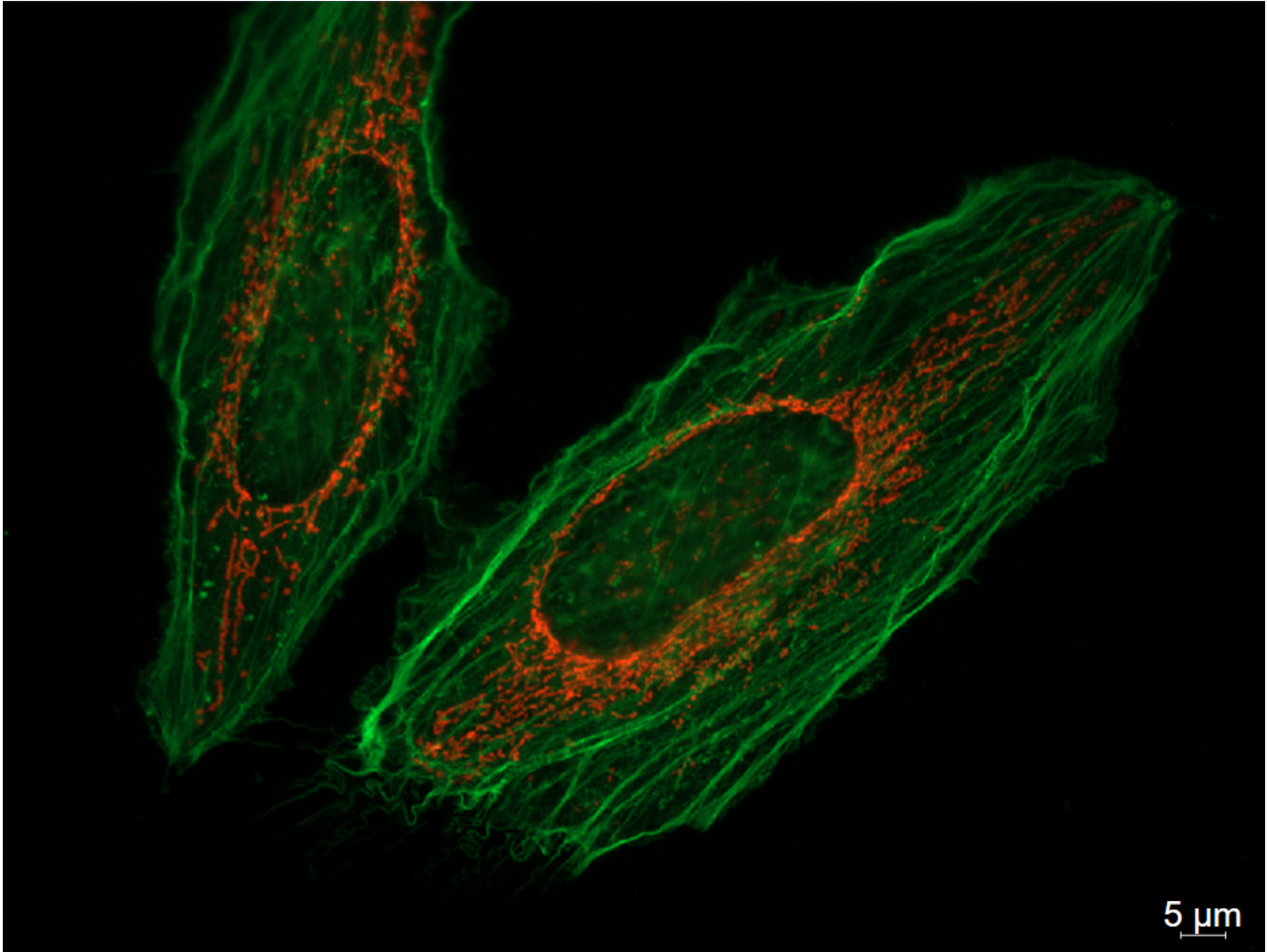

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	29

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.102µm X 0.102µm
Image size (Pixels)	1388 X 1040
Image Size (Scaled)	142.10µm X 106.48µm
Bit Depth	12 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	45 Texas Red 38 eGFP
Beam Splitter	585 495
Filter Excitation Wavelength	540 - 580 450 - 490
Filter Emission Wavelength	593 - 668 500 - 550
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @40%
Exposure Time	100ms 200ms
Excitation Wavelength	590 493
Emission Wavelength	618 517
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	0.96µm 0.80µm
Binning Mode	1,1

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Donkey anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000020

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Mitochondria Antibody (MA5-12017), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Actin Recombinant Rabbit Monoclonal Antibody (JF47-01), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.102µm X 0.102µm

Image size (Pixels) = 1388 X 1040

Image Size (Scaled) = 142.10µm X 106.48µm

Bit Depth = 12 Bit

Image Center Position = X: 0.00µm, Y: 0.00µm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflector = 45 Texas Red

Beam Splitter = 585

Filter Excitation Wavelength = 540 - 580

Filter Emission Wavelength = 593 - 668

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @40%

Exposure Time = 100ms

Excitation Wavelength = 590

Emission Wavelength = 618

Effective NA = 1.4

Imaging Device = Axiocam MR R3

Camera Adapter = 1X

Depth of Focus = 0.96µm

Binning Mode = 1,1

488 -

Reflector = 38 eGFP
 Beam Splitter = 495
 Filter Excitation Wavelength = 450 - 490
 Filter Emission Wavelength = 500 - 550
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @40%
 Exposure Time = 200ms
 Excitation Wavelength = 493
 Emission Wavelength = 517
 Effective NA = 1.4
 Imaging Device = AxioCam MR R3
 Camera Adapter = 1X
 Depth of Focus = 0.80µm
 Binning Mode = 1,1

Post Processing

None

Image Correction


Identifyn® SRM Alexa Fluor™ 594 was pseudo-colored in Zen Software from red, the default color, to orange to provide visual contrast between the 594/488 channels.

Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000055
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	29
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	FDF1BE0FDE5914116D1B94BDD51134C1D8164DD59673CFD81D96D61F036DDBF
Method PDF SHA-256	CFC1BF537D2EA0C0E37AD26F3A970818EA9E762B1284811593F99FD5812539B3

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

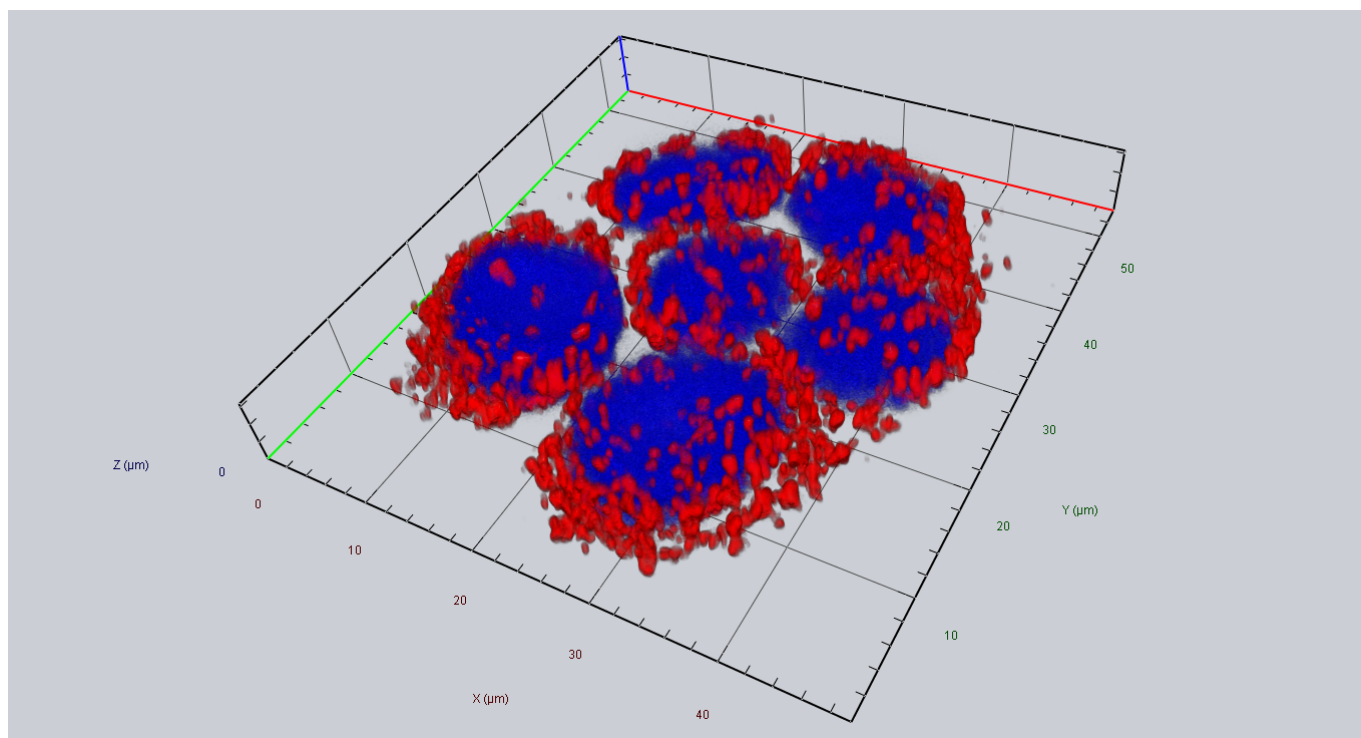

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	MCF-7
DETECTED DYES	405/DAPI; 594
METADATA VALUES	38

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Z-Stack	63 Slices (6.20µm)
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm X .0100µm
Image Size (Pixels)	346 X 376
Image Size (Scaled)	49.43µm X 53.71µm
Bit Depth	14 Bit
Image Center Position	X: 103.83mm, Y: 51.41mm
Stage Position	X: 103.83mm, Y: 51.41mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	45 Texas Red 49 DAPI
Beam Splitter	585 395
Filter Excitation Wavelength	540 - 580 335 - 383
Filter Emission Wavelength	593 - 668 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @18% X-Cite Xylis Lamp Intensity @8%
Exposure Time	200ms 15ms
Excitation Wavelength	590 353
Emission Wavelength	618 465
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.20µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.
3D Transparency Projection	Precise
Appearance	Threshold 10%, Ramp 20%, Maximum 100%
Background	White
Light	Brightness 103%, Enabled Tone Mapping
Projection	View Angle 34°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

MCF7 were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Mitochondria Antibody (MA5-12017), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 63 Slices (6.20µm)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm X .0100µm

Image Size (Pixels) = 346 X 376

Image Size (Scaled) = 49.43µm X 53.71µm

Bit Depth = 14 Bit

Image Center Position = X: 103.83mm, Y: 51.41mm

Stage Position = X: 103.83mm, Y: 51.41mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

594 -

Reflector = 45 Texas Red

Beam Splitter = 585

Filter Excitation Wavelength = 540 - 580

Filter Emission Wavelength = 593 - 668

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @18%

Exposure Time = 200ms

Excitation Wavelength = 590

Emission Wavelength = 618

Effective NA = 1.25

Imaging Device = Axiocam 807

Camera Adapter = 1X

Depth of Focus = 1.20µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @8%
 Exposure Time = 15ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Correction, Fourier Filter, and Deconvolution were NOT used.

3D Image Processing

3D Transparency Projection = Precise
 Appearance = Threshold 10%, Ramp 20%, Maximum 100%
 Background = White
 Light = Brightness 103%, Enabled Tone Mapping
 Projection = View Angle 34°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000055
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	38
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT

SOURCE PROVENANCE

Cell model	MCF-7
Image SHA-256	1E59D3AF30FCE64F2E1A38C2BE27047CEDA0368D83DC025BD62B40373E1AEF91
Method PDF SHA-256	D686BDDBE3BD2E31479334CA98C4C3EF7189F119372D5E092C16182F0AD046BF

ORIGINAL IMAGE EVIDENCE / CONFOCAL

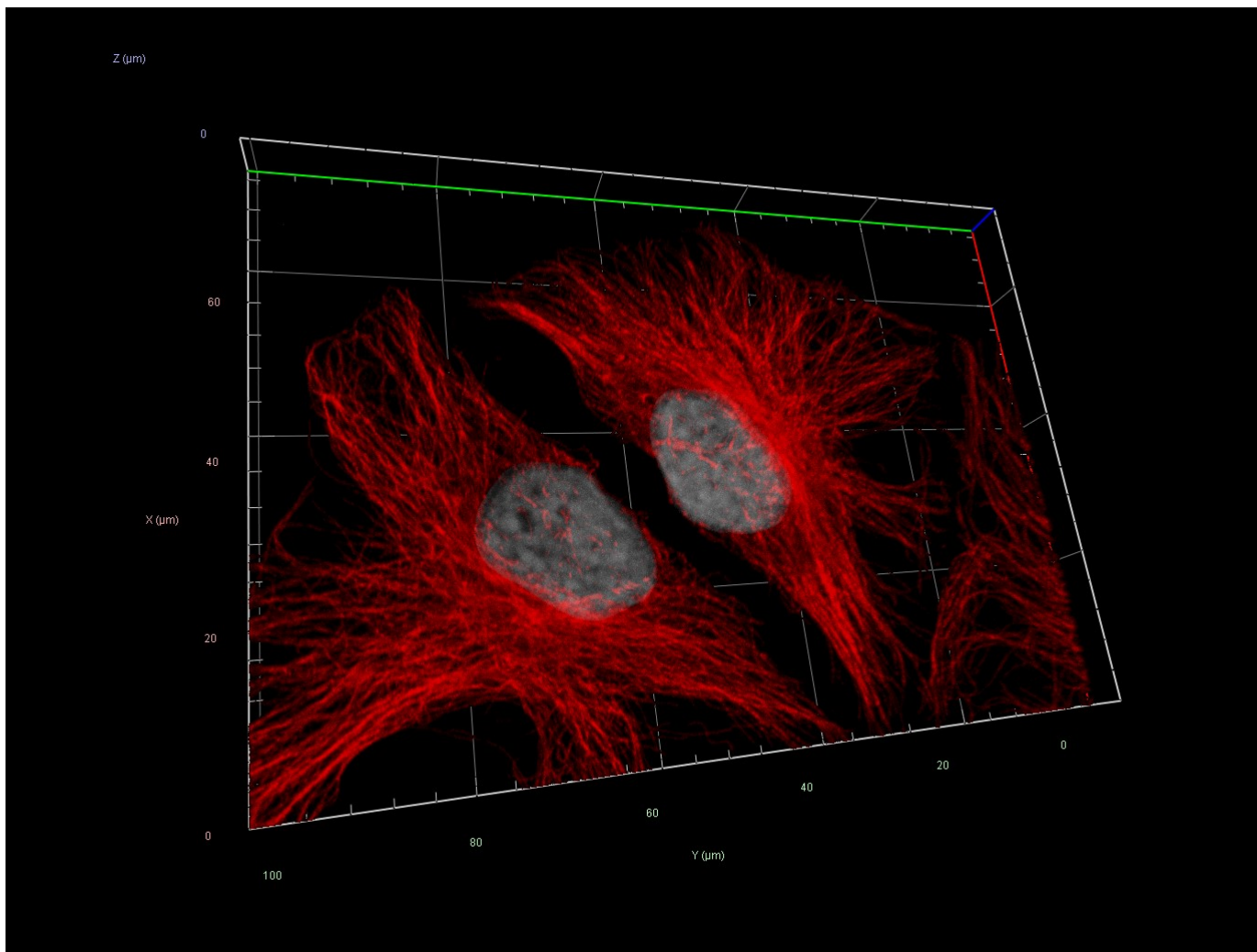

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Channels	2
Scaling (Per Pixel)	0.071m X 0.071µm X 0.060µm
Image size (Pixels)	1037 X 1431
Image Size (Scaled)	73.18µm X 100.99µm
Bit Depth	16 Bit
Image Center Position	X: 88.63mm, Y: 54.52mm
Stage Position	X: 88.63mm, Y: 54.52mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 54µm 1.00AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	561nm @0.09% 405nm @0.7%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	.080 1.0
Rotation	0°
Sampling	1.00
Pixel Time	0.37µs
Frame Time	3.55s
LSM Scan Speed	9
Scan Direction	Bidirectional
Line Step	1
Averaging	4
Excitation Wavelength	590 353
Emission Wavelength	618 465
Effective NA	1.4
Detection Wavelength	590 - 700 410 - 470
Imaging Device	GaAsP-Pmt2 GaAsP-Pmt1
Detector Type	GaAsP-PMT
Detector Gain	700 V 650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
3D Maximum Projection	Precise
Appearance	Threshold 4% all channels, Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 65%, Enabled Tone Mapping
Projection	View Angle 69°, Scale Z 150%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

HeLa Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.071m X 0.071µm X 0.060µm

Image size (Pixels) = 1037 X 1431

Image Size (Scaled) = 73.18µm X 100.99µm

Bit Depth = 16 Bit

Image Center Position = X: 88.63mm, Y: 54.52mm

Stage Position = X: 88.63mm, Y: 54.52mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 54µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @0.09%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = .080
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.37µs
 Frame Time = 3.55s
 LSM Scan Speed = 9
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 590
 Emission Wavelength = 618
 Effective NA = 1.4
 Detection Wavelength = 590 - 700
 Imaging Device = GaAsP-Pmt2
 Detector Type = GaAsP-PMT
 Detector Gain = 700 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.7%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.37µs
 Frame Time = 3.55s
 LSM Scan Speed = 9
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410 - 470
 Imaging Device = GaAsP-Pmt1
 Detector Type = GaAsP-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 4% all channels, Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 65%, Enabled Tone Mapping
 Projection = View Angle 69°, Scale Z 150%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000055
Stage	Confocal

SOURCE PROVENANCE

Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	7735D5B59A9F78496CEC76CCFFCFECD749DC062C44C57B086ED02A3BB5E3897B
Method PDF SHA-256	28FE1FF917FFB223FCC414936B2F7842C5AB00572C6A48B34A60C9147D5CE49A

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

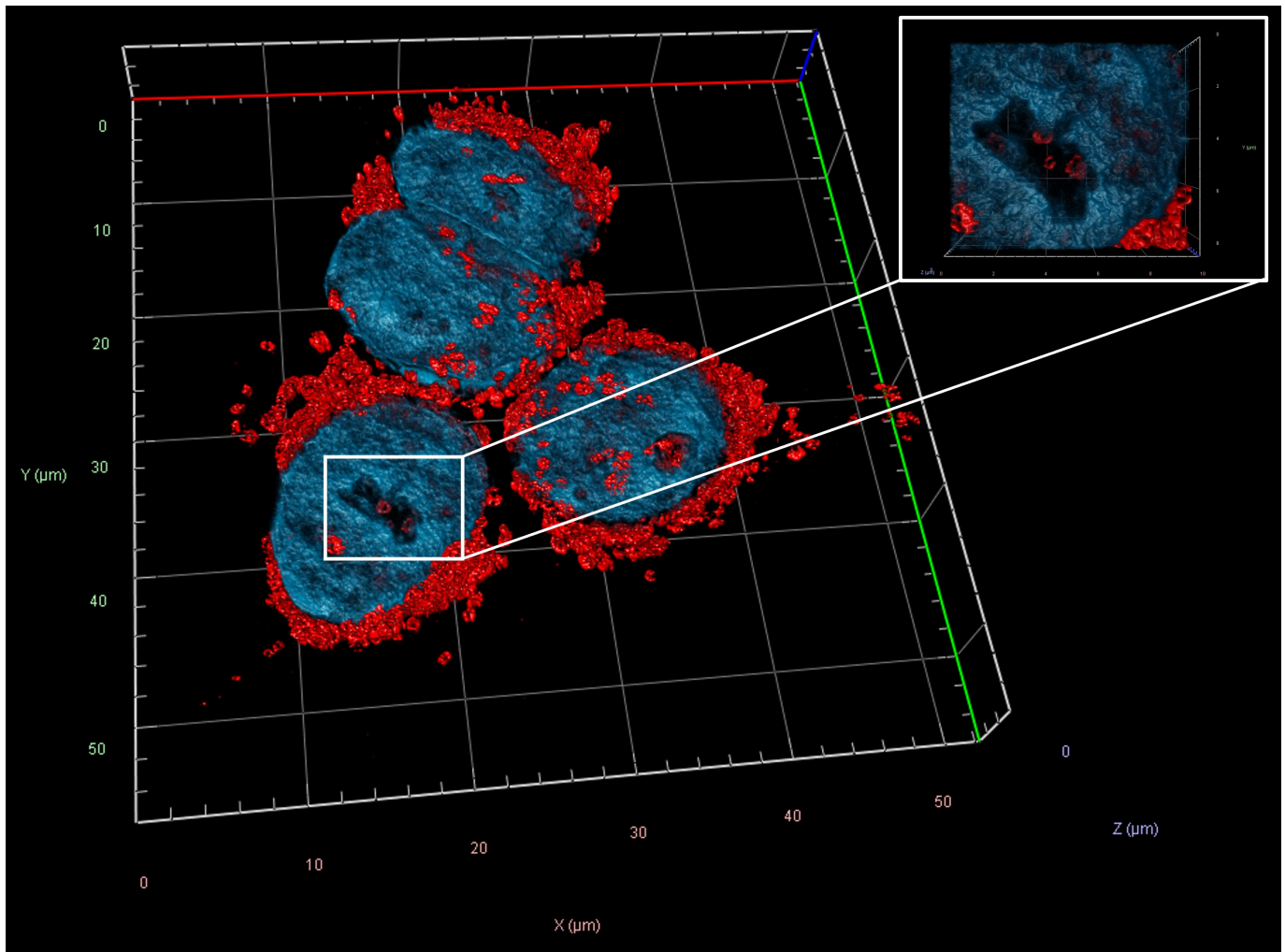

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	MCF-7
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	47

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	109 Slices (5.40µm)
Channels	2
Scaling (Per Pixel)	0.03530 µm X 0.03530 µm X 0.05000 µm
Image Size (Pixels)	1382 X 1410
Image Size (Scaled)	48.78µm X 49.77µm
Bit Depth	16 Bit
Image Center Position	X: 52.74mm, Y: 54.33mm
Stage Position	X: 52.74mm, Y: 54.33mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00 AU / 251µm 5.00 AU / 188µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	561nm @1.70% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.3
Rotation	0°
Sampling	2.00
Pixel Time	0.96µs
Frame Time	21.88s
LSM Scan Speed	6
Scan Direction	Bidirectional
Line Step	1
Averaging	4
Excitation Wavelength	590 353
Emission Wavelength	618 465
Effective NA	1.4
Detection Wavelength	560-640 400-495
Imaging Device	Airyscan
Detector Type	GaAsp-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	850 V 800 V
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 9.4 (3D, Auto) Super Resolution 8.0 (3D, Auto)
3D Volume Projection	Precise
Appearance	Threshold 2%, Ramp 98%, Maximum 100%
Background	Black
Light	Brightness 182%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods

Identifyfyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

MCF7 were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Mitochondria Antibody (MA5-12017), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyfyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 109 Slices (5.40µm)

Channels = 2

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 1382 X 1410

Image Size (Scaled) = 48.78µm X 49.77µm

Bit Depth = 16 Bit

Image Center Position = X: 52.74mm, Y: 54.33mm

Stage Position = X: 52.74mm, Y: 54.33mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 251µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @1.70%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 21.88s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 590
 Emission Wavelength = 618
 Effective NA = 1.4
 Detection Wavelength = 560-640
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 9.4 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 188µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 21.88s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400-495
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.0 (3D, Auto)

Post Processing - AiryScan 3D Processing

Display Mode "SR"
 Super Resolution Auto Filter Selected
 Process 3D Current Image
 Created

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 2%, Ramp 98%, Maximum 100%
 Background = Black
 Light = Brightness 182%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections -

None


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000055
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	47
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	MCF-7
Image SHA-256	E7DAF441C01CF23A06207ED768F8C37E8B5A51F55E20BFFB4E17706CD09E8CEA
Method PDF SHA-256	BD11BD76D33300E3BB95C7FC3FFF75ED44965C9E609F694A02599D5DC9E1C196

ORIGINAL IMAGE EVIDENCE / SIM

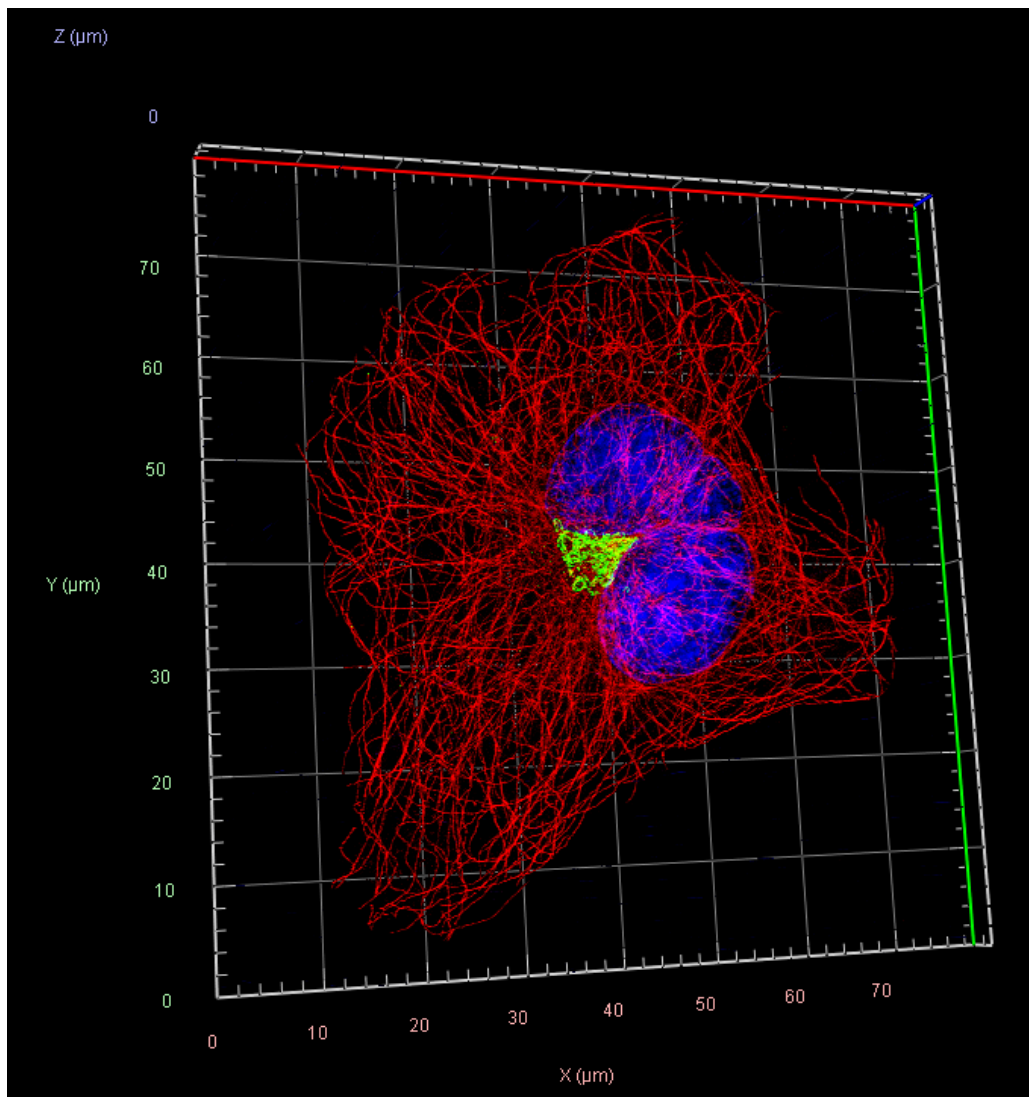

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	63

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	70 Slices (3.45µm)
Channels	3 1
Scaling (Per Pixel)	0.02111 µm X 0.02111 µm X 0.05000 µm
Image Size (Pixels)	3744 X 3779
Image Size (Scaled)	79.05µm X 79.78µm
Bit Depth	16 Bit
Image Center Position	X: 56.00mm, Y: 36.24mm
Stage Position	X: 56.00mm, Y: 36.24mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	561 488 405
Channel Name	T2-Axiocam 820 #2-SR T2-Axiocam 820 -SR
Excitation Wavelength	561 488 405
Emission Wavelength	583 729 525
Effective NA	1.4
Detection Wavelength	503 - 546, 576 - 590 657 - 800
Imaging Device	Axiocam 820 #2 Axiocam 820
Camera Adapter	1X
Exposure Time	120ms 100ms 50ms
Depth of Focus	0.90µm 1.13µm 0.81µm
Binning Mode	2,2
Laser Wavelength	561nm @3.50% 488nm @1.0% 405nm @5.0%
Frame Time	1.33s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	2
Algorithm	SIM

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	11.6303
Input SNR	N/A
Regularization Weight	N/A
Iterations	N/A
Filter	N/A
Fitting	N/A
Normalization	Automatic
Experimental PSF	N/A
3D Maximum Projection	Precise
Appearance	Threshold 1% all channels, Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 128%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Donkey Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000020

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody that followed this, Golgi GOLPH2 Polyclonal Antibody), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant catalog number P36961.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 70 Slices (3.45µm)

Channels = 3

Scaling (Per Pixel) = 21.11nm X 21.11nm X 50.00nm

Image Size (Pixels) = 3744 X 3779

Image Size (Scaled) = 79.05µm X 79.78µm

Bit Depth = 16 Bit

Image Center Position = X: 56.00mm, Y: 36.24mm

Stage Position = X: 56.00mm, Y: 36.24mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 561
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 561
 Emission Wavelength = 583
 Effective NA = 1.4
 Detection Wavelength = 503 - 546, 576 - 590
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.90µm
 Binning Mode = 2,2
 Laser Wavelength = 561nm @3.50%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 32.0µm grating

488 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820 -SR
 Excitation Wavelength = 488
 Emission Wavelength = 729
 Effective NA = 1.4
 Detection Wavelength = 657 - 800
 Imaging Device = Axiocam 820
 Camera Adapter = 1X
 Exposure Time = 100ms
 Depth of Focus = 1.13µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @1.0%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 27.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 503 - 546, 576 - 590
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @5.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM Processing

Processing Dimension - 3D
 Result Sampling = 4
 Process Sampling = 2
 Channels = 1
 Algorithm = SIM
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 11.6303
 Input SNR = N/A
 Regularization Weight = N/A
 Iterations = N/A
 Filter = N/A
 Fitting = N/A
 Normalization = Automatic
 Experimental PSF = N/A

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 1% all channels, Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 128%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

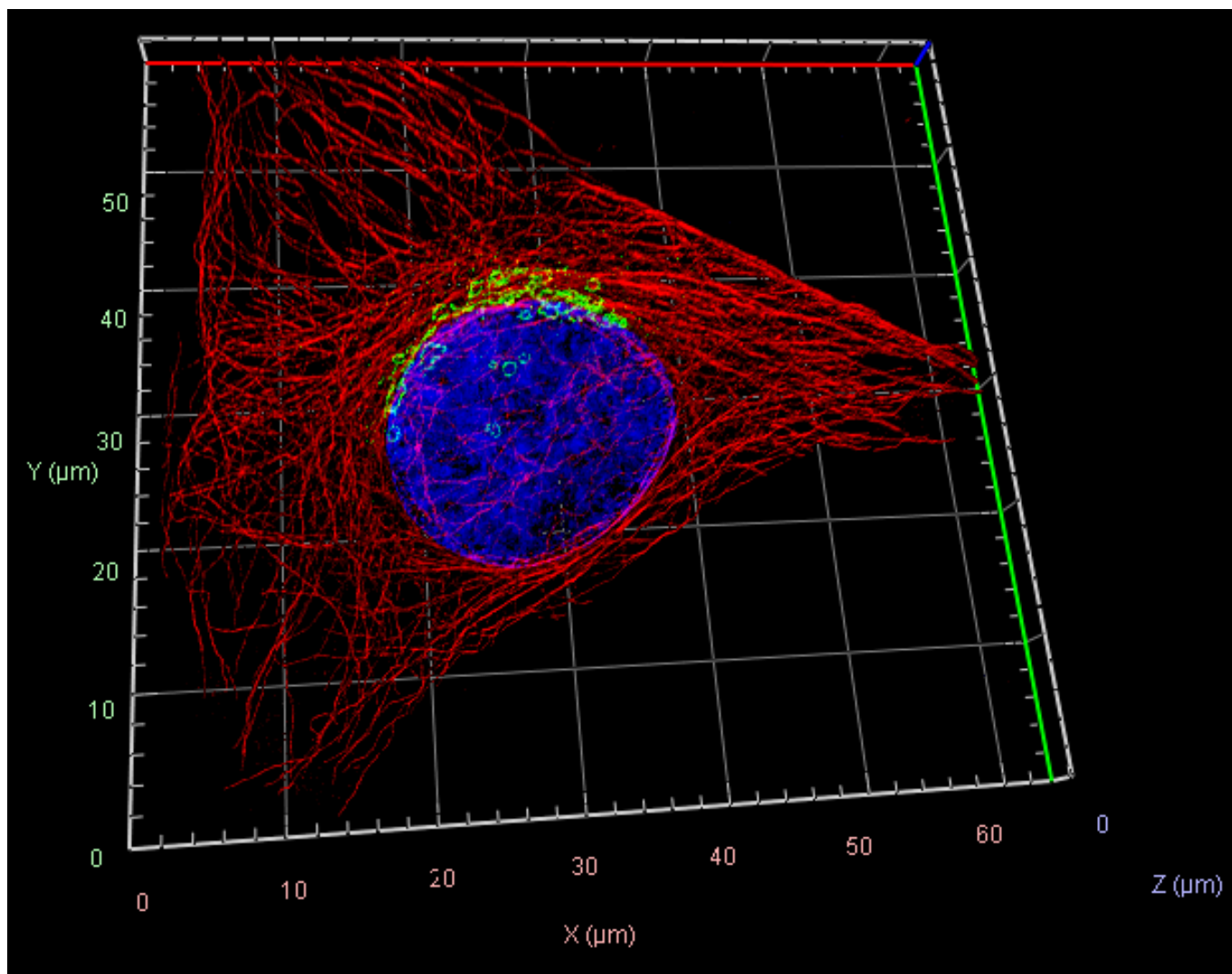
None

Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000055
Stage	SIM
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	63
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	E1A4B1FB09AF3663067B079116C11A6CA45DEBD994098DAE179AFA6ECB895DD9
Method PDF SHA-256	4520771DDD05E45509B82B40D0F3495D3CC9D4AA560D98949E0825EEC35D94DB

ORIGINAL IMAGE EVIDENCE / SIM²

EVIDENCE TIER	SIM ²
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 594
METADATA VALUES	63

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	45 Slices (2.20µm)
Channels	3 1
Scaling (Per Pixel)	0.01689 µm X 0.01689 µm X 0.05000 µm
Image Size (Pixels)	3761 X 3542
Image Size (Scaled)	63.52µm X 59.82µm
Bit Depth	16 Bit
Image Center Position	X: 56.20mm, Y: 33.93mm
Stage Position	X: 56.20mm, Y: 33.93mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	561 488 405
Channel Name	T2-Axiocam 820 #2-SR T2-Axiocam 820 -SR
Excitation Wavelength	561 488 405
Emission Wavelength	583 729 525
Effective NA	1.4
Detection Wavelength	503 - 546, 576 - 590 657 - 800
Imaging Device	Axiocam 820 #2 Axiocam 820
Camera Adapter	1X
Exposure Time	120ms 100ms 50ms
Depth of Focus	0.90µm 1.13µm 0.81µm
Binning Mode	2,2
Laser Wavelength	561nm @3.50% 488nm @1.0% 405nm @5.0%
Frame Time	1.33s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	4
Algorithm	SIM2

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	10.3460
Input SNR	N/A
Regularization Weight	N/A
Iterations	N/A
Filter	N/A
Fitting	N/A
Normalization	Automatic
Experimental PSF	N/A
3D Maximum Projection	Precise
Appearance	Threshold 2% all channels, Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 110%, Enabled Tone Mapping
Projection	View Angle 45°, Scale Z 141%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

Identifyn® SRM Alexa Fluor™ 488 Affinity Purified Donkey Anti-Rabbit F(ab')₂ (H + L)

Cat ID# - SA 000020

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody that followed this, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 488 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant catalog number P36961.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 45 Slices (2.20µm)

Channels = 3

Scaling (Per Pixel) = 16.89nm X 16.89nm X 50.00nm

Image Size (Pixels) = 3761 X 3542

Image Size (Scaled) = 63.52µm X 59.82µm

Bit Depth = 16 Bit

Image Center Position = X: 56.20mm, Y: 33.93mm

Stage Position = X: 56.20mm, Y: 33.93mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

594 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 561
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 561
 Emission Wavelength = 583
 Effective NA = 1.4
 Detection Wavelength = 503 - 546, 576 - 590
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.90µm
 Binning Mode = 2,2
 Laser Wavelength = 561nm @3.50%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 32.0µm grating

488 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820 -SR
 Excitation Wavelength = 488
 Emission Wavelength = 729
 Effective NA = 1.4
 Detection Wavelength = 657 - 800
 Imaging Device = Axiocam 820
 Camera Adapter = 1X
 Exposure Time = 100ms
 Depth of Focus = 1.13µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @1.0%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 27.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 503 - 546, 576 - 590
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @5.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM2 Processing

Processing Dimension - 3D
 Result Sampling = 4
 Process Sampling = 4
 Channels = 1
 Algorithm = SIM2
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 10.3460
 Input SNR = N/A
 Regularization Weight = N/A
 Iterations = N/A
 Filter = N/A
 Fitting = N/A
 Normalization = Automatic
 Experimental PSF = N/A

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 2% all channels, Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 110%, Enabled Tone Mapping
 Projection = View Angle 45°, Scale Z 141%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000055
Stage	SIM ²
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	63
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	EBCFDBFF5C4E4DB0AFC32B9039A2A8911BF053E0912F86CF8F94C577A38D2053
Method PDF SHA-256	DCF1846D81136ED553BF321332C4F4110BBAAD9C0616FE154FC247E2BBE168B1

ORIGINAL IMAGE EVIDENCE / CONTROL

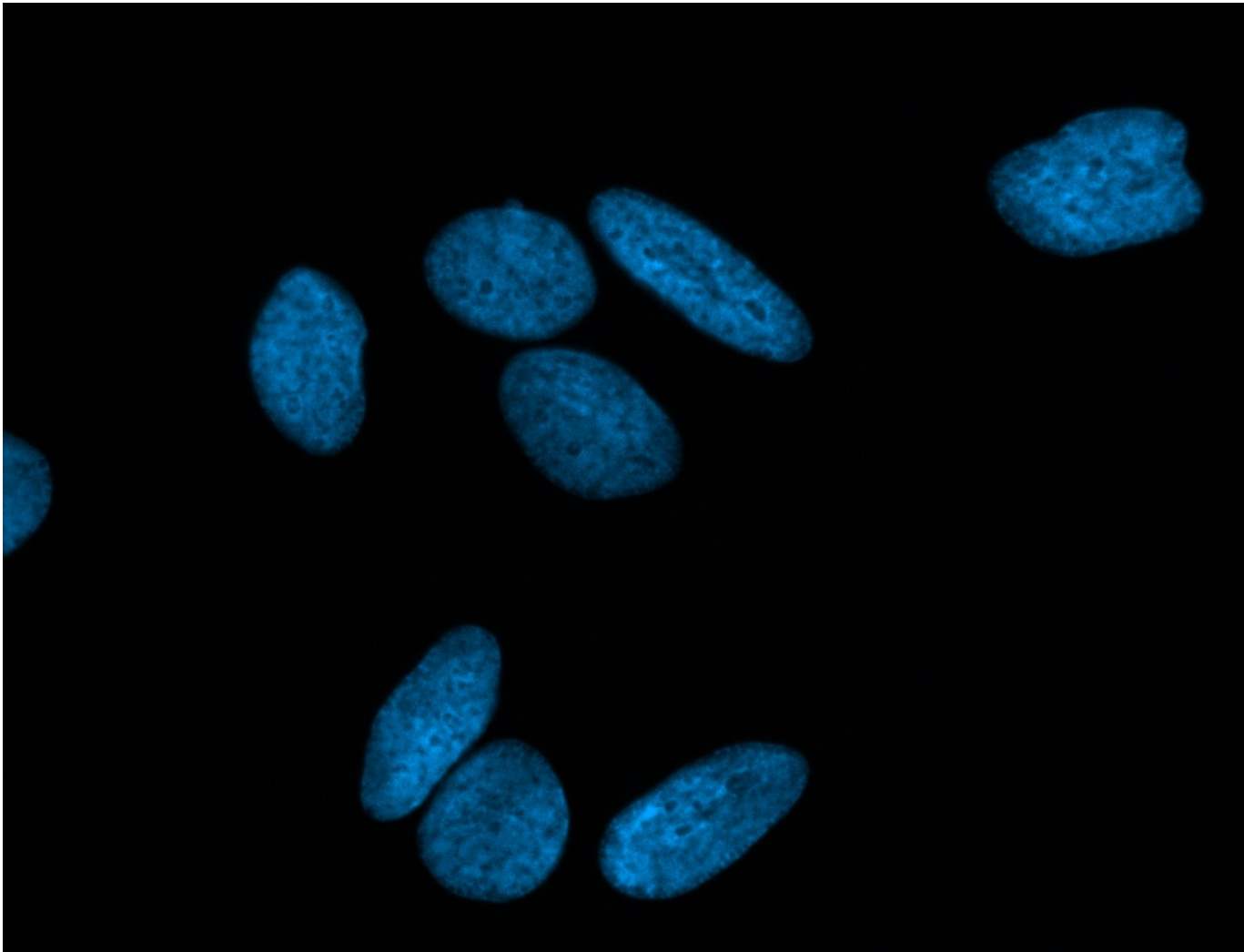

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 594
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1078 X 823
Image Size (Scaled)	154.00µm X 117.57µm
Bit Depth	14 Bit
Image Center Position	X: 57.36mm, Y: 46.56mm
Stage Position	X: 57.36mm, Y: 46.56mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	45 Texas Red 49 DAPI
Beam Splitter	585 395
Filter Excitation Wavelength	540 - 580 335 - 383
Filter Emission Wavelength	593 - 668 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	590 353
Emission Wavelength	618 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.20µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods

Identifyn® SRM Alexa Fluor™ 594 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000055

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 594 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1078 X 823

Image Size (Scaled) = 154.00µm X 117.57µm

Bit Depth = 14 Bit

Image Center Position = X: 57.36mm, Y: 46.56mm

Stage Position = X: 57.36mm, Y: 46.56mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

594 -

Reflector = 45 Texas Red

Beam Splitter = 585

Filter Excitation Wavelength = 540 - 580

Filter Emission Wavelength = 593 - 668

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 590

Emission Wavelength = 618

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.20µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000055
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	57B40ED2BBB4BF6F0DFD4FC27CFAB11392ECD9879ED01B93096F8A3C01AD08D7
Method PDF SHA-256	FC9C8BC2539238FD4C7F169321F37DB0F6D2D959BE366F26E3C0211EFF25E923