

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse, F(ab')₂ (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM²

7


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

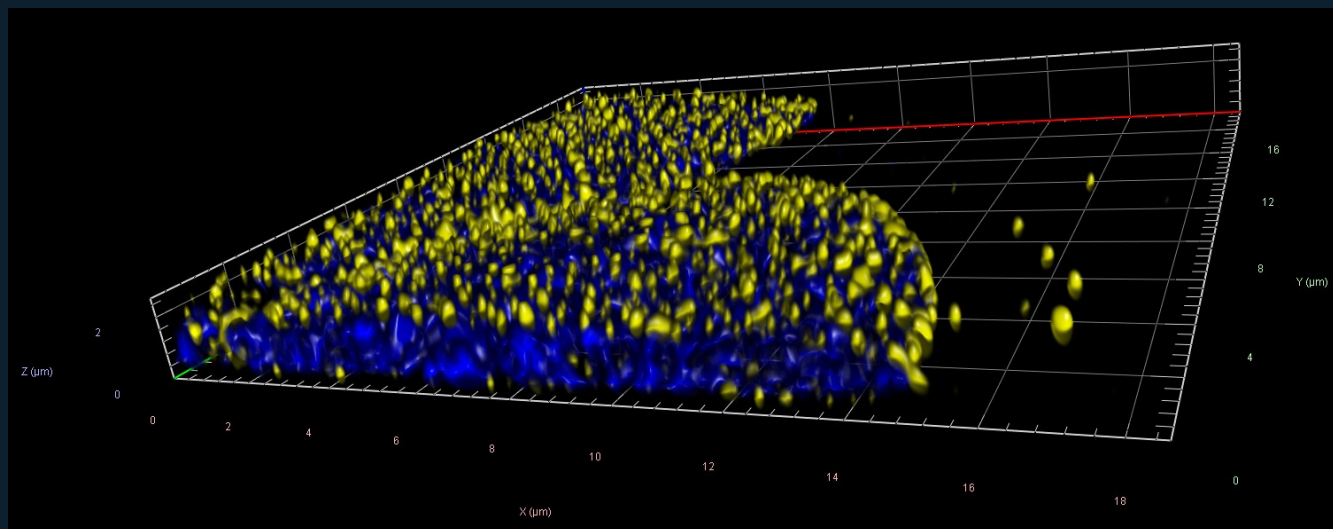

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SIM²

SA-000033 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse, F(ab')₂ (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Control

BENNETT ASSESSMENT

The preserved SA-000033 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / PARTIAL STEPDOWN / SOURCE-BOUNDED REVIEW

The preserved SA-000033 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

SETTINGS ASSESSMENT

The record preserves 7 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	532	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 532; 555; 647	2; AF532-49014; 50 Cy5; 515 - 545; 625 - 655; 555 - 595; 665 - 715; 528
Widefield SIM — Apotome	405/DAPI; 532; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335-383; 555 - 595; 420-470; 528
Confocal	405/DAPI; 488; 532	2; None; 561nm @3.50%; 405nm @0.2%; 528; 353; 553; 465
Airyscan	405/DAPI; 488; 532	2; None; 488nm @0.7%; 405nm @1.0%; 528; 353; 553; 465
SIM	405/DAPI; 488; 532	2; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T2-Axiocam 820 #2-SR; T3-Axiocam 820 #2-SR; 488; 405
SIM ²	405/DAPI; 488; 532	2; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T2-Axiocam 820 #2-SR; T3-Axiocam 820 #2-SR; 488; 405
Control	405/DAPI; 532; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335 - 383; 555 - 595; 420 - 470; 528

MODALITY ASSESSMENT / PART 1 OF 8

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 8

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.102µm X 0.102µm; Image size (Pixels): 1388 X 1040; Image Size (Pixels): 1388 X 1040; Bit Depth: 12 Bit. Detector/camera: Imaging Device: Axiocam MR R3. Timing: Exposure Time: 125ms; 150ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @30%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 29.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 532; 555; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 8

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 833 X 763; Image Size (Pixels): 833 X 763; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807. Timing: Exposure Time: 250ms; 50ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @30%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 8

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.071m X 0.071µm; Image size (Pixels): 1437 X 1437; Image Size (Pixels): 1437 X 1437; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.54µs; Frame Time: 20.90s. Illumination: Laser Wavelength: 561nm @3.50%; 405nm @0.2%. Optical/scan: Pinhole: 1.00AU / 48µm; 1.00AU / 37µm; Scan Zoom: 1.0; LSM Scan Speed: 8; Scan Direction: Bidirectional; Averaging: 16; Effective NA: 1.4. Structured source metadata values: 39.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 8

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 32 Slices (1.55µm); Channels: 2; Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm; Image size (Pixels): 2028 X 2028; Image Size (Pixels): 2028 X 2028; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.51µs; Frame Time: 5.07s. Illumination: Laser Wavelength: 488nm @0.7%; 405nm @1.0%. Optical/scan: Pinhole: 5.00AU / 229µm; 5.00AU / 181µm; Scan Zoom: 1.4; LSM Scan Speed: 7; Scan Direction: Bidirectional; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 9.0 (3D, Auto); Super Resolution 8.2 (3D, Auto). Structured source metadata values: 47.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 8

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 62 Slices (3.05µm); Channels: 2; Scaling (Per Pixel): 16.89nm X 16.89nm X 50.00nm; Image size (Pixels): 2488 X 1700; Image Size (Pixels): 2488 X 1700; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820 #2. Timing: Exposure Time: 120ms; 50ms; Frame Time: 1.59s; 676.00ms. Illumination: Laser Wavelength: 488nm @6.00%; 405nm @6.0%; Illumination Wavelength: 488; 405. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM; Process Sampling: N/A; Result Sampling: 4. Structured source metadata values: 49.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 8

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 52 Slices (2.55µm); Channels: 2; Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm; Image size (Pixels): 1408 X 1338; Image Size (Pixels): 1408 X 1338; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820 #2. Timing: Exposure Time: 120ms; 50ms; Frame Time: 1.59s; 676.00ms. Illumination: Laser Wavelength: 488nm @8.00%; 405nm @6.0%; Illumination Wavelength: 488; 405. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM2; Process Sampling: 4; Result Sampling: 4; Projection: View Angle 35°, Scale Z 66%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 55.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 8

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1492 X 776; Image Size (Pixels): 1492 X 776; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000033 record contains Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². Missing modalities are not represented or inferred.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm -> 0.03530 μ m X 0.03530 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2028 X 2028; Image Size (Scaled)=71.59 μ m X 71.11 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.67%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 16.89nm X 16.89nm X 50.00nm -> 0.01689 μ m X 0.01689 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2488 X 1700; Image Size (Scaled)=42.02 μ m X 28.71 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm -> 0.01320 μ m X 0.01320 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1408 X 1338; Image Size (Scaled)=18.58 μ m X 17.66 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.03%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked
Stepdown order. Scientific wording is preserved; missing modalities are not invented.
A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000033. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra	PRESERVED
7	SIM ²	ZEISS Elyra	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

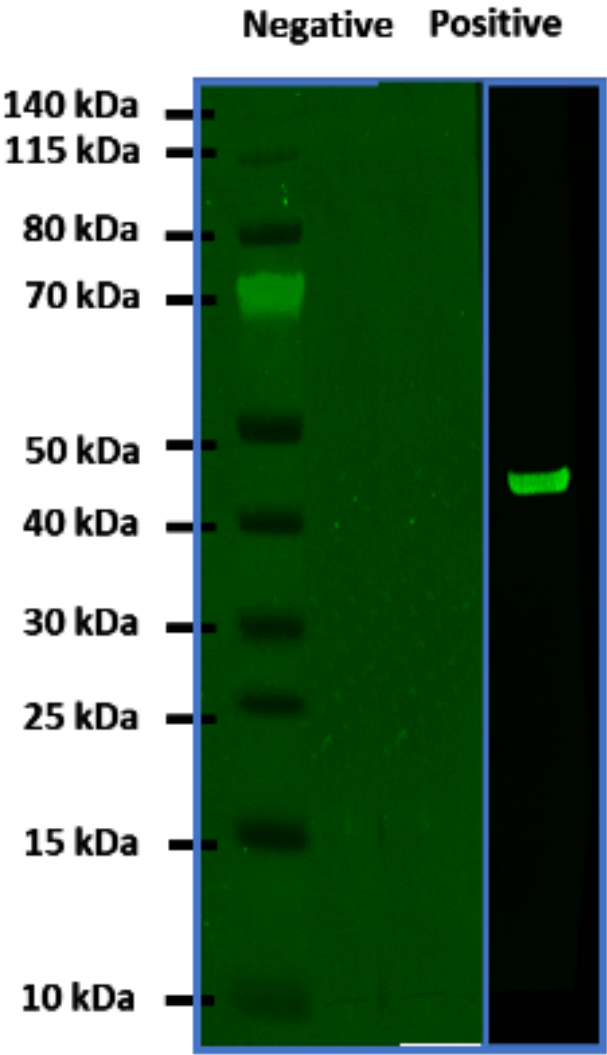

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	532
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	532nm
Emission Filter	572±28nm
Detector	Avalanche Photodiode
Intensity	2
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Actin Monoclonal Antibody (mAbGEa) (MA1-744) for 2 hours, followed by 5X PBS washes. Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit anti-Mouse F(ab)₂ was incubated for 1 hour and washed 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 532nm

Emission Filter = 572±28nm

Detector = Avalanche Photodiode

Intensity = 2

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: Identifyn® LLC

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000033
Stage	Western blot
Instrument	Azure Sapphire FL
Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

SOURCE PROVENANCE

Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	75CF7B98C61E16F224DC35BF537CF324805F1DA5FE0E427720E4E0FA787A75A6
Method PDF SHA-256	3C097412172A25C67350189D0F605501E6ECC78495359077D3A556E9970FA8AD

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

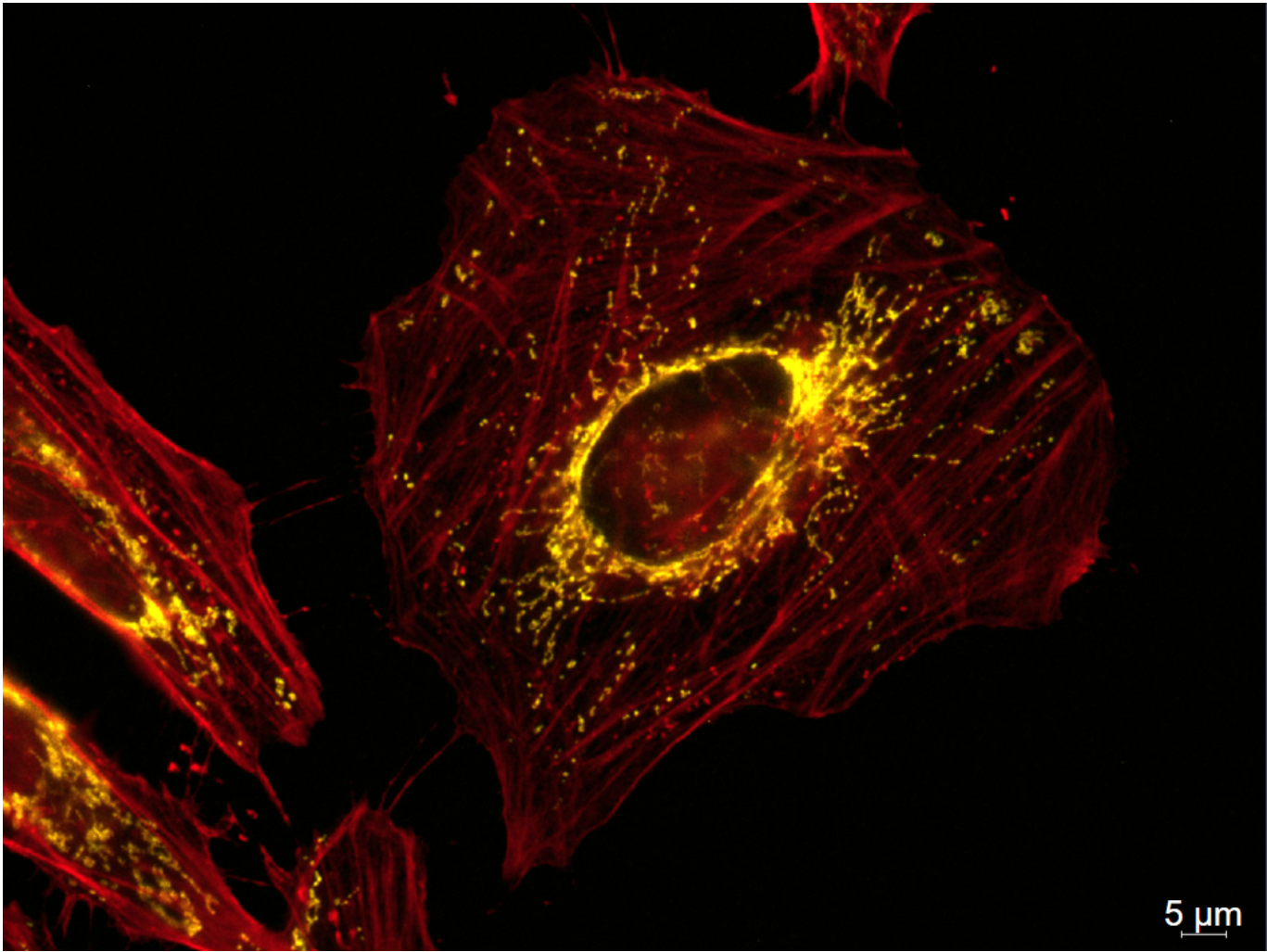

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555; 647
METADATA VALUES	29

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.102µm X 0.102µm
Image size (Pixels)	1388 X 1040
Image Size (Scaled)	142.10µm X 106.48µm
Bit Depth	12 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.10µm, Y: 0.00µm
Reflector	AF532-49014 50 Cy5
Beam Splitter	550 660
Filter Excitation Wavelength	515 - 545 625 - 655
Filter Emission Wavelength	555 - 595 665 - 715
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @30%
Exposure Time	125ms 150ms
Excitation Wavelength	528 653
Emission Wavelength	553 668
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	0.86µm 1.03µm
Binning Mode	1,1

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Donkey Anti-Rabbit, F(ab')₂ (H + L)

Cat ID# - SA 000065

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Mitochondria Antibody (MA5-12017), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Actin Recombinant Rabbit Monoclonal Antibody (JF47-01), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.102µm X 0.102µm

Image size (Pixels) = 1388 X 1040

Image Size (Scaled) = 142.10µm X 106.48µm

Bit Depth = 12 Bit

Image Center Position = X: 0.00µm, Y: 0.00µm

ROI Center Offset = X: 0.10µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @30%

Exposure Time = 125ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.4

Imaging Device = Axiocam MR R3

Camera Adapter = 1X

Depth of Focus = 0.86µm

Binning Mode = 1,1

647 -

Reflector = 50 Cy5
 Beam Splitter = 660
 Filter Excitation Wavelength = 625 - 655
 Filter Emission Wavelength = 665 - 715
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @30%
 Exposure Time = 150ms
 Excitation Wavelength = 653
 Emission Wavelength = 668
 Effective NA = 1.4
 Imaging Device = AxioCam MR R3
 Camera Adapter = 1X
 Depth of Focus = 1.03µm
 Binning Mode = 1,1

Post Processing

None

Image Correction


None

Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000033
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	29
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	BDBF54F302BDB54A84FE45E998EE127A029857C6A2A7F48FDED201D3DD51AE61
Method PDF SHA-256	94307F89B272172BC9104B98EE3C4FE7F551185DC13AB0D4C41C41BF07EFB788

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

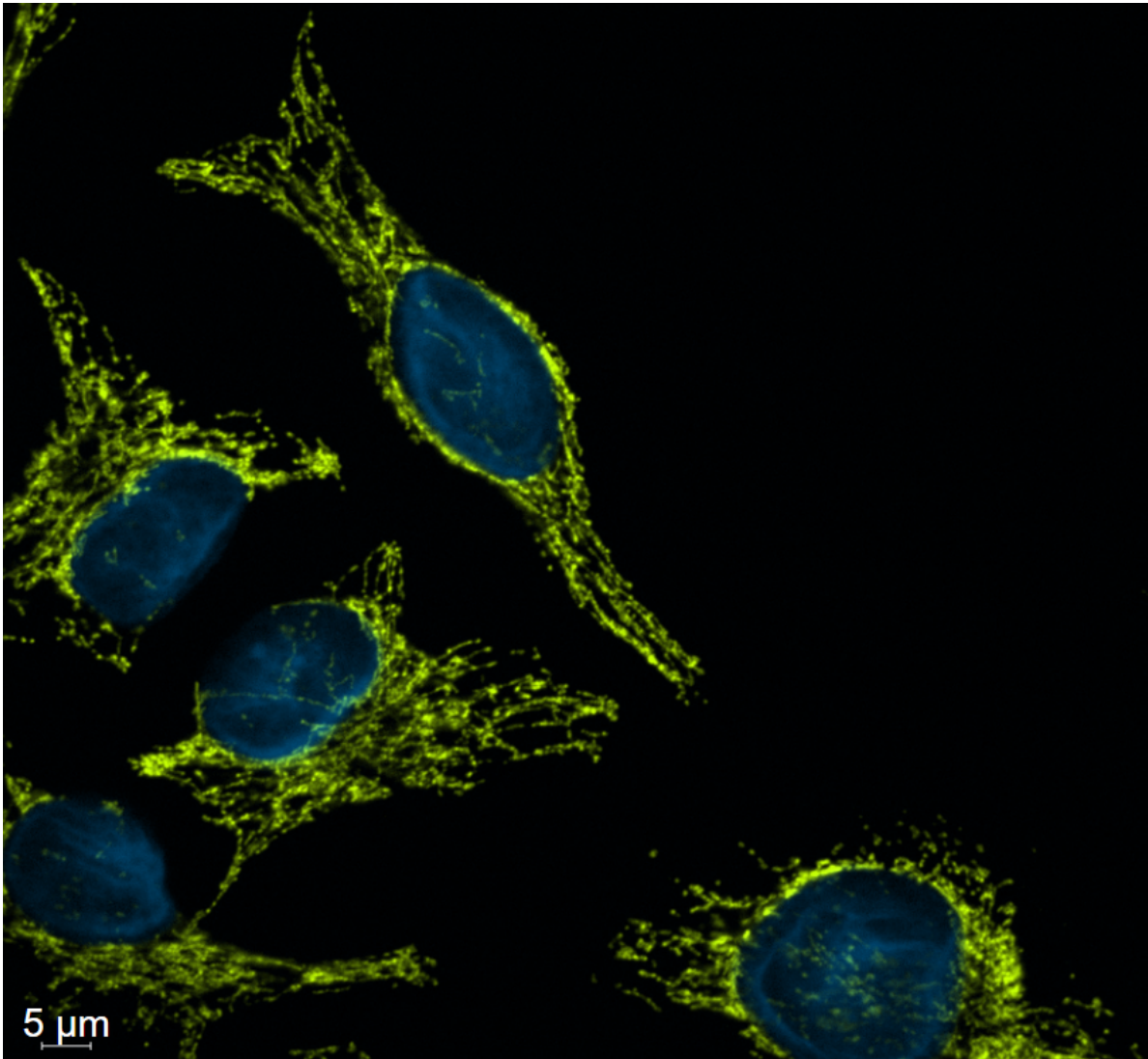

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	833 X 763
Image Size (Scaled)	119.00µm X 109.00µm
Bit Depth	14 Bit
Image Center Position	X: 75.57mm, Y:55.04mm
Stage Position	X: 75.57mm, Y:55.04mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335-383
Filter Emission Wavelength	555 - 595 420-470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @30% X-Cite Xylis Lamp Intensity @10%
Exposure Time	250ms 50ms
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.07µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Mitochondria Antibody (MA5-12017), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 833 X 763

Image Size (Scaled) = 119.00µm X 109.00µm

Bit Depth = 14 Bit

Image Center Position = X: 75.57mm, Y:55.04mm

Stage Position = X: 75.57mm, Y:55.04mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @30%

Exposure Time = 250ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.25

Imaging Device = Axiocam 807

Camera Adapter = 1X

Depth of Focus = 1.07µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335-383
 Filter Emission Wavelength = 420-470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 50ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = AxioCam 807
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None


Image by E.F Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000033
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	AC523FE14A984FF805A64E11AB7B9CF180B631E1545F84D3861C4245933EC295
Method PDF SHA-256	5F357E0D4AF018C9D90C04B532C47C588D97005739F7F7323DAA7BA4F3B4B238

ORIGINAL IMAGE EVIDENCE / CONFOCAL

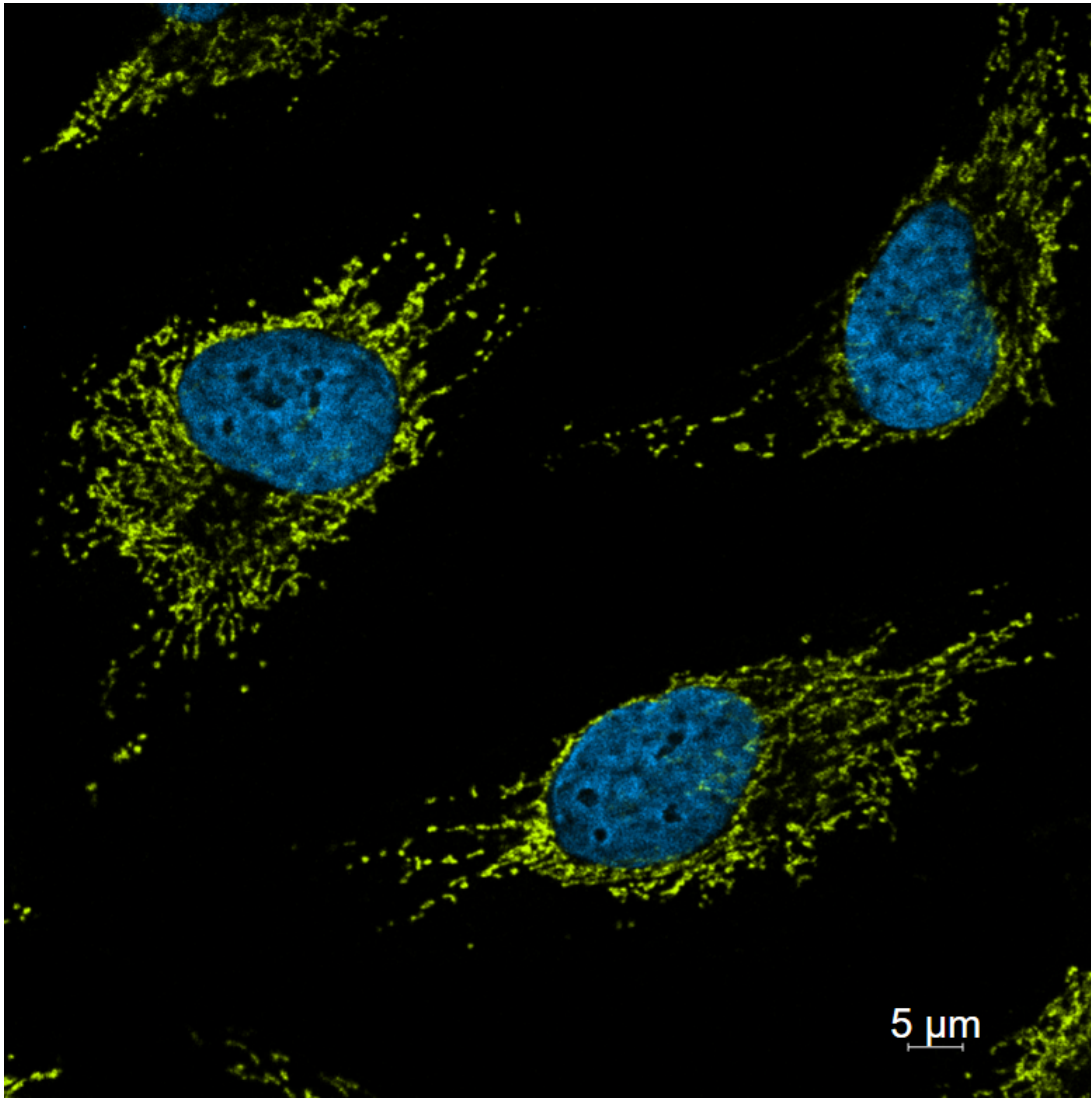

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	39

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Channels	2
Scaling (Per Pixel)	0.071m X 0.071µm
Image size (Pixels)	1437 X 1437
Image Size (Scaled)	101.41µm X 101.41µm
Bit Depth	16 Bit
Image Center Position	X: 51.29mm, Y: 50.79mm
Stage Position	X: 51.29mm, Y: 50.79mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 48µm 1.00AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	561nm @3.50% 405nm @0.2%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.0
Rotation	0°
Sampling	1.00
Pixel Time	0.54µs
Frame Time	20.90s
LSM Scan Speed	8
Scan Direction	Bidirectional
Line Step	1
Averaging	16
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Detection Wavelength	510 - 590 410 - 470
Imaging Device	Airyscan
Detector Type	GaAsp-PMT
Detector Gain	800 V 650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Mitochondria Antibody (MA5-12017), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.071m X 0.071µm

Image size (Pixels) = 1437 X 1437

Image Size (Scaled) = 101.41µm X 101.41µm

Bit Depth = 16 Bit

Image Center Position = X: 51.29mm, Y: 50.79mm

Stage Position = X: 51.29mm, Y: 50.79mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 48µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @3.50%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.54µs
 Frame Time = 20.90s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 16
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 510 - 590
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.2%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.54µs
 Frame Time = 20.90s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 16
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410 - 470
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

Post Processing

N/A

Image Corrections

None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000033
Stage	Confocal
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

SOURCE PROVENANCE

Structured source metadata values	39
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	B7D8AE365C82EF0FC5E4BE50AF55AD97EBD0E0B10235B8F0C28681BFDB240F5C
Method PDF SHA-256	28E90D3C790E56272B6231D7AB7421EB4FAA5FC20F88D671D90FE1492998C24A

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

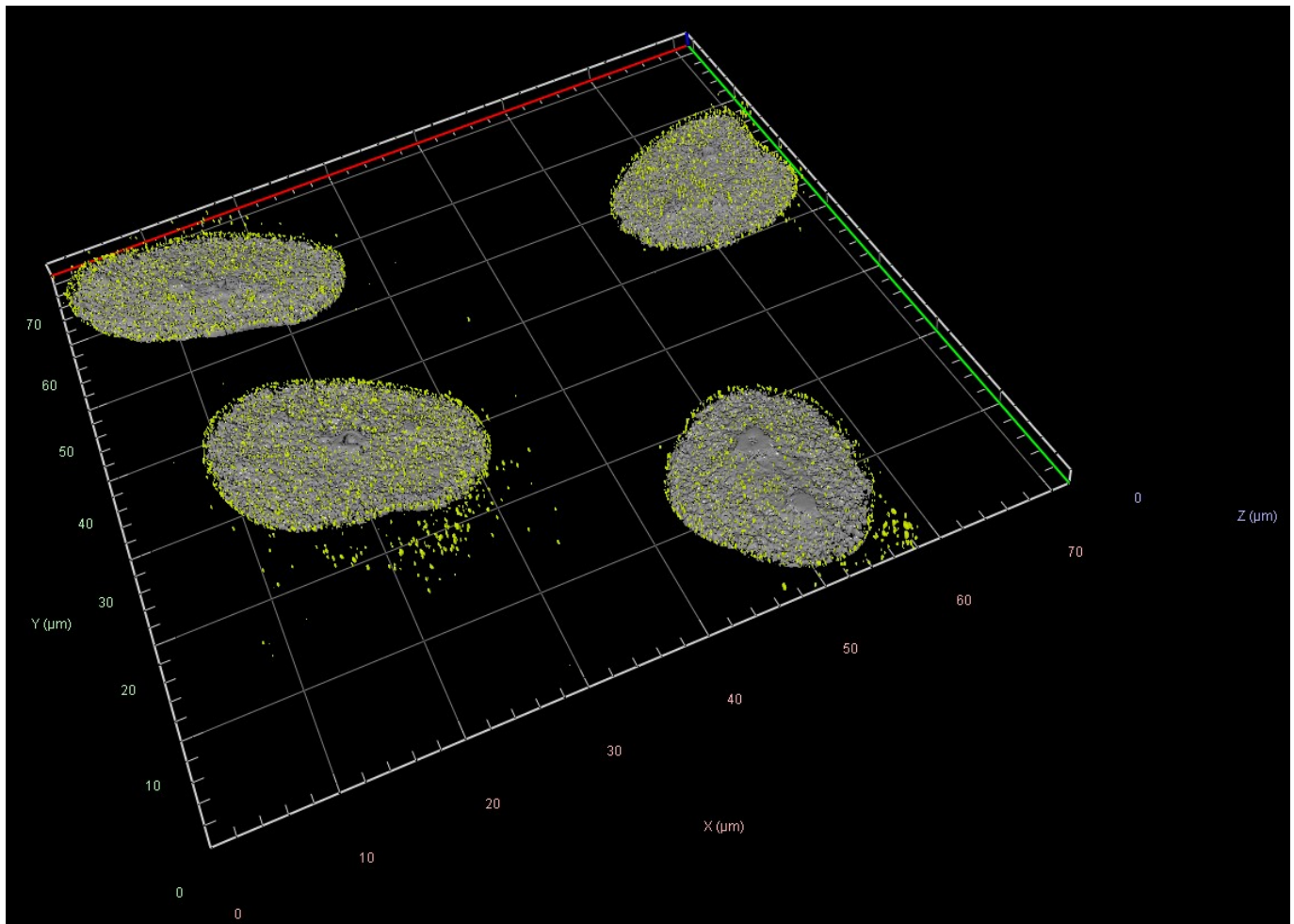

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	47

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	32 Slices (1.55µm)
Channels	2
Scaling (Per Pixel)	0.03530 µm X 0.03530 µm X 0.05000 µm
Image Size (Pixels)	2028 X 2028
Image Size (Scaled)	71.59µm X 71.11µm
Bit Depth	16 Bit
Image Center Position	X: 70.69mm, Y: 48.05mm
Stage Position	X: 70.69mm, Y: 48.05mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00AU / 229µm 5.00AU / 181µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @0.7% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.4
Rotation	32.73° 0°
Sampling	2.00
Pixel Time	0.51µs
Frame Time	5.07s
LSM Scan Speed	7
Scan Direction	Bidirectional
Line Step	1
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Detection Wavelength	515 - 580 400 - 465
Imaging Device	Airyscan
Detector Type	GaAsp-PMT
Detector Gain	850 V 800 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 9.0 (3D, Auto) Super Resolution 8.2 (3D, Auto)
3D Modeling Projection	Precise
Appearance	Threshold 20% all channels, Ambient Light 0% all channels, Specular Light 100% all channels, Shininess 30% all channels
Background	Black
Light	Azimuth 80°, Elongation -130°, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Nuclear Pore Complex Monoclonal Antibody (39C7), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 32 Slices (1.55µm)

Channels = 2

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 2028 X 2028

Image Size (Scaled) = 71.59µm X 71.11µm

Bit Depth = 16 Bit

Image Center Position = X: 70.69mm, Y: 48.05mm

Stage Position = X: 70.69mm, Y: 48.05mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00AU / 229µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @0.7%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.4
 Rotation = 32.73°
 Sampling = 2.00
 Pixel Time = 0.51µs
 Frame Time = 5.07s
 LSM Scan Speed = 7
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 2
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 515 - 580
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 9.0 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00AU / 181µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.4
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.51µs
 Frame Time = 5.07s
 LSM Scan Speed = 7
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 2
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400 - 465
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.2 (3D, Auto)

3D Image Processing

3D Modeling Projection = Precise
 Appearance = Threshold 20% all channels, Ambient Light 0% all channels, Specular Light 100% all channels,
 Shininess 30% all channels
 Background = Black
 Light = Azimuth 80°, Elongation -130°, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera
 Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

DAPI was pseudo-colored dark grey instead of its default blue coloring to demonstrate contrast and the specificity and brightness of the Identifyn® SRM Alexa Fluor™ 532 dye product.


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.
 Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000033
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	47
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	1D27F63F5DEA90F38F4A9D03540C508F1B65712F19E2C6774EC4783410AB4B18
Method PDF SHA-256	EF5210EAC67578D68E816B4D7954F5B78E649B856A4052321946C1C3D96FFF8B

ORIGINAL IMAGE EVIDENCE / SIM

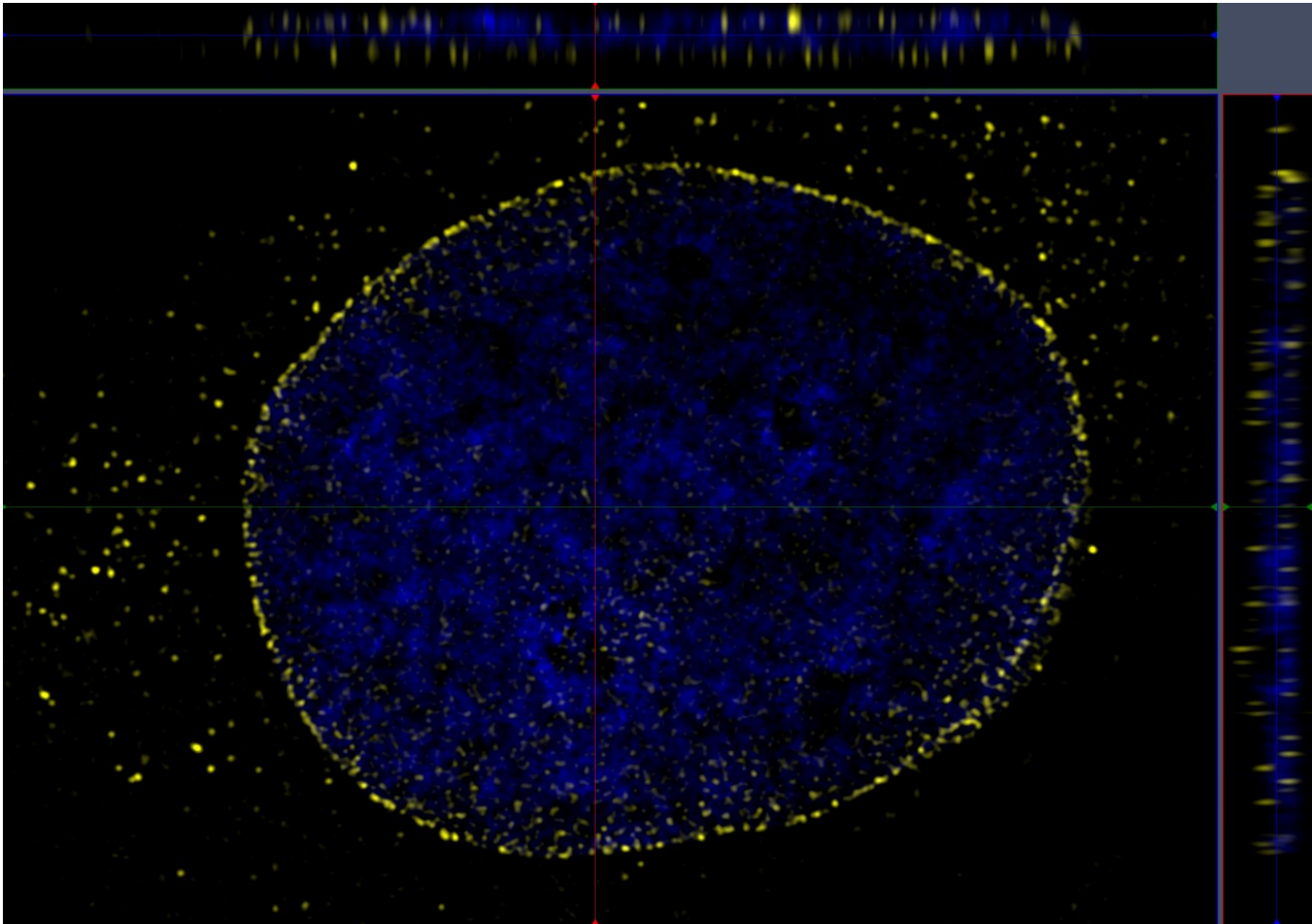

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	49

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	62 Slices (3.05µm)
Channels	2
Scaling (Per Pixel)	0.01689 µm X 0.01689 µm X 0.05000 µm
Image Size (Pixels)	2488 X 1700
Image Size (Scaled)	42.02µm X 28.71µm
Bit Depth	16 Bit
Image Center Position	X: 45.02mm, Y: 39.94mm
Stage Position	X: 45.02mm, Y: 39.94mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	488 405
Channel Name	T2-Axiocam 820 #2-SR T3-Axiocam 820 #2-SR
Excitation Wavelength	488 405
Emission Wavelength	525 445
Effective NA	1.4
Detection Wavelength	419-470, 503-546, 576-617, 657-800 419-470, 503 - 546, 576 - 617, 657-800
Imaging Device	Axiocam 820 #2
Camera Adapter	1X
Exposure Time	120ms 50ms
Depth of Focus	0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	488nm @6.00% 405nm @6.0%
Frame Time	1.59s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	N/A
Algorithm	SIM

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	9.6639 (532), 8.8597 (DAPI)
Fitting	Fast Fit
Histogram	Scale to Min/Max
Baseline	Cut
Cut Lines	X 1218, Y 842, Z 25
Line Width	1 all axis
Cit Line Opacity	50

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Nuclear Pore Complex Monoclonal Antibody (39C7), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 62 Slices (3.05µm)

Channels = 2

Scaling (Per Pixel) = 16.89nm X 16.89nm X 50.00nm

Image Size (Pixels) = 2488 X 1700

Image Size (Scaled) = 42.02µm X 28.71µm

Bit Depth = 16 Bit

Image Center Position = X: 45.02mm, Y: 39.94mm

Stage Position = X: 45.02mm, Y: 39.94mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503-546, 576-617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @6.00%
 Frame Time = 1.59s
 Scan Direction = Unidirectional
 Grating 36.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T3-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503 - 546, 576 - 617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @6.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM Processing

Processing Dimension - 3D

Result Sampling = 4

Process Sampling = N/A

Channels = 2

Algorithm = SIM

SIM Settings = Standard

Auto Sharpness = On

Sharpness = 9.6639 (532), 8.8597 (DAPI)

Fitting = Fast Fit

Histogram = Scale to Min/Max

Baseline = Cut

Orthogonal View

Cut Lines = X 1218, Y 842, Z 25

Line Width = 1 all axis

Cut Line Opacity = 50

Maximum Projection and Measure 3D Distance NOT enabled

Image Corrections

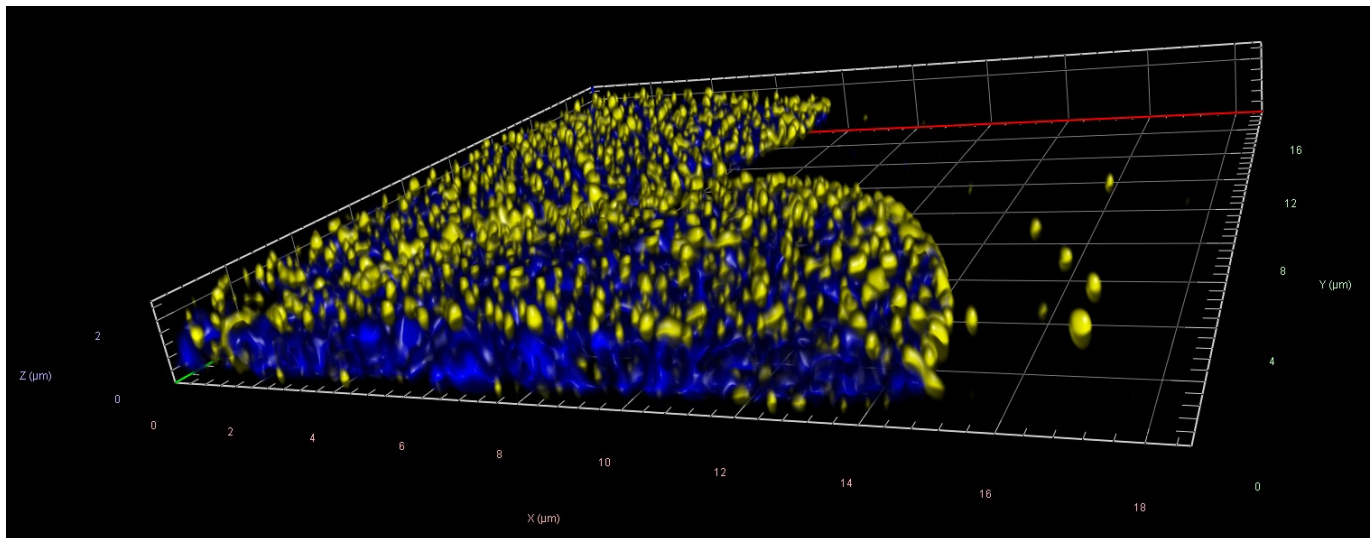
None

Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000033
Stage	SIM
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	49
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	5CA6DE9F1A4E19EDDFF4813711236D00EF2C6679A670CC2B7C86FBAEC92029AD
Method PDF SHA-256	F420D20F765C717427A120F77AF4F0C2757248E4A5DEE509F3858E066BE0E0D1

ORIGINAL IMAGE EVIDENCE / SIM²

EVIDENCE TIER	SIM ²
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	55

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	52 Slices (2.55µm)
Channels	2
Scaling (Per Pixel)	0.01320 µm X 0.01320 µm X 0.05000 µm
Image Size (Pixels)	1408 X 1338
Image Size (Scaled)	18.58µm X 17.66µm
Bit Depth	16 Bit
Image Center Position	X: 48.45mm, Y: 36.69mm
Stage Position	X: 48.45mm, Y: 36.69mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	488 405
Channel Name	T2-Axiocam 820 #2-SR T3-Axiocam 820 #2-SR
Excitation Wavelength	488 405
Emission Wavelength	525 445
Effective NA	1.4
Detection Wavelength	419-470, 503-546, 576-617, 657-800 419-470, 503 - 546, 576 - 617, 657-800
Imaging Device	Axiocam 820 #2
Camera Adapter	1X
Exposure Time	120ms 50ms
Depth of Focus	0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	488nm @8.00% 405nm @6.0%
Frame Time	1.59s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	4
Algorithm	SIM2

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	8.8315 (532), 7.7892 (DAPI)
Order Combination	Weiner filter
Regularization	None
Regularization Weight	0.0000
Iterations	3
Filter	No Filter
Fitting	Fast Fit
Histogram	Scale to Min/Max
3D Volume Projection	Precise
Appearance	Threshold 4% (532), Threshold 4% (DAPI) Ramp 96% (532), Ramp 96% (DAPI), Maximum 100% both channels
Background	Black
Light	Brightness 186%, Azimuth 73°, Elongation 29°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 66%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Nuclear Pore Complex Monoclonal Antibody (39C7), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 52 Slices (2.55µm)

Channels = 2

Scaling (Per Pixel) = 13.20nm X 13.20nm X 50.00nm

Image Size (Pixels) = 1408 X 1338

Image Size (Scaled) = 18.58µm X 17.66µm

Bit Depth = 16 Bit

Image Center Position = X: 48.45mm, Y: 36.69mm

Stage Position = X: 48.45mm, Y: 36.69mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503-546, 576-617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @8.00%
 Frame Time = 1.59s
 Scan Direction = Unidirectional
 Grating 36.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T3-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503 - 546, 576 - 617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @6.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM2 Processing

Processing Dimension - 3D
 Result Sampling = 4
 Process Sampling = 4
 Channels = 2
 Algorithm = SIM2
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 8.8315 (532), 7.7892 (DAPI)
 Order Combination = Weiner filter
 Regularization = None
 Regularization Weight = 0.0000
 Iterations = 3
 Filter = No Filter
 Fitting = Fast Fit
 Histogram = Scale to Min/Max

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 4% (532), Threshold 4% (DAPI) Ramp 96% (532), Ramp 96% (DAPI), Maximum 100% both channels
 Background = Black
 Light = Brightness 186%, Azimuth 73°, Elongation 29°, Enabled Directional Lighting, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 66%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE	
Catalog	SA-000033
Stage	SIM²
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	55
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	29E0F27250F196AF44BCB3B7A1D62E81E0D9EB906F79070885DB737717C2E674

SOURCE PROVENANCE

Method PDF SHA-256

2CD35B31143FC737213FB7D8B9EBD6E2CE6306B5566AFFCA1C2EDBD536354949

ORIGINAL IMAGE EVIDENCE / CONTROL

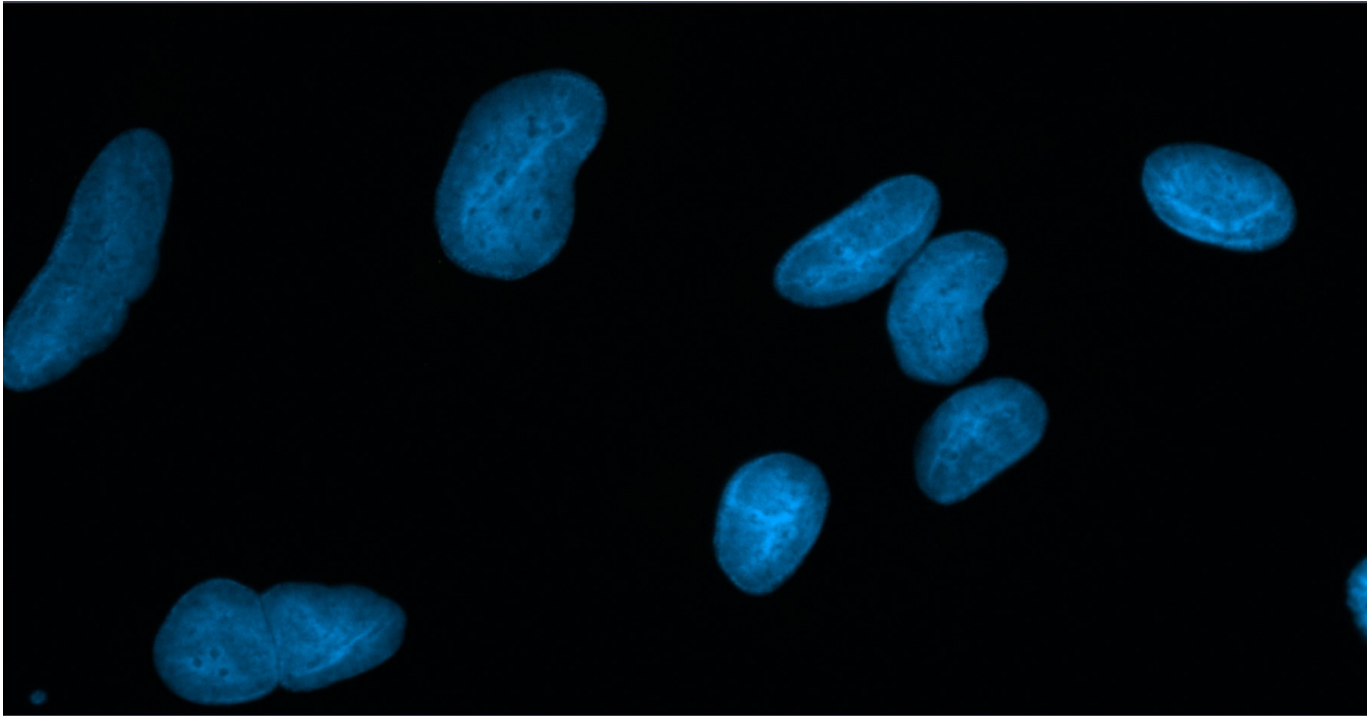

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1492 X 776
Image Size (Scaled)	213.14µm X 110.86µm
Bit Depth	14 Bit
Image Center Position	X: 50.92mm, Y: 49.47mm
Stage Position	X: 50.92mm, Y: 49.47mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335 - 383
Filter Emission Wavelength	555 - 595 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.07µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Rabbit Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000033

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1492 X 776

Image Size (Scaled) = 213.14µm X 110.86µm

Bit Depth = 14 Bit

Image Center Position = X: 50.92mm, Y: 49.47mm

Stage Position = X: 50.92mm, Y: 49.47mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.07µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000033
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	58FC70EA272DB175253ABDE2D7329D86F094F8FA0C8FB7287626A4F035209785
Method PDF SHA-256	5BEEB230BD50B4C02964FEDBB13DC6448525FE37AFA11F1000EEE52938E6CFF1