

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse, F(ab')₂ (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM²

7


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

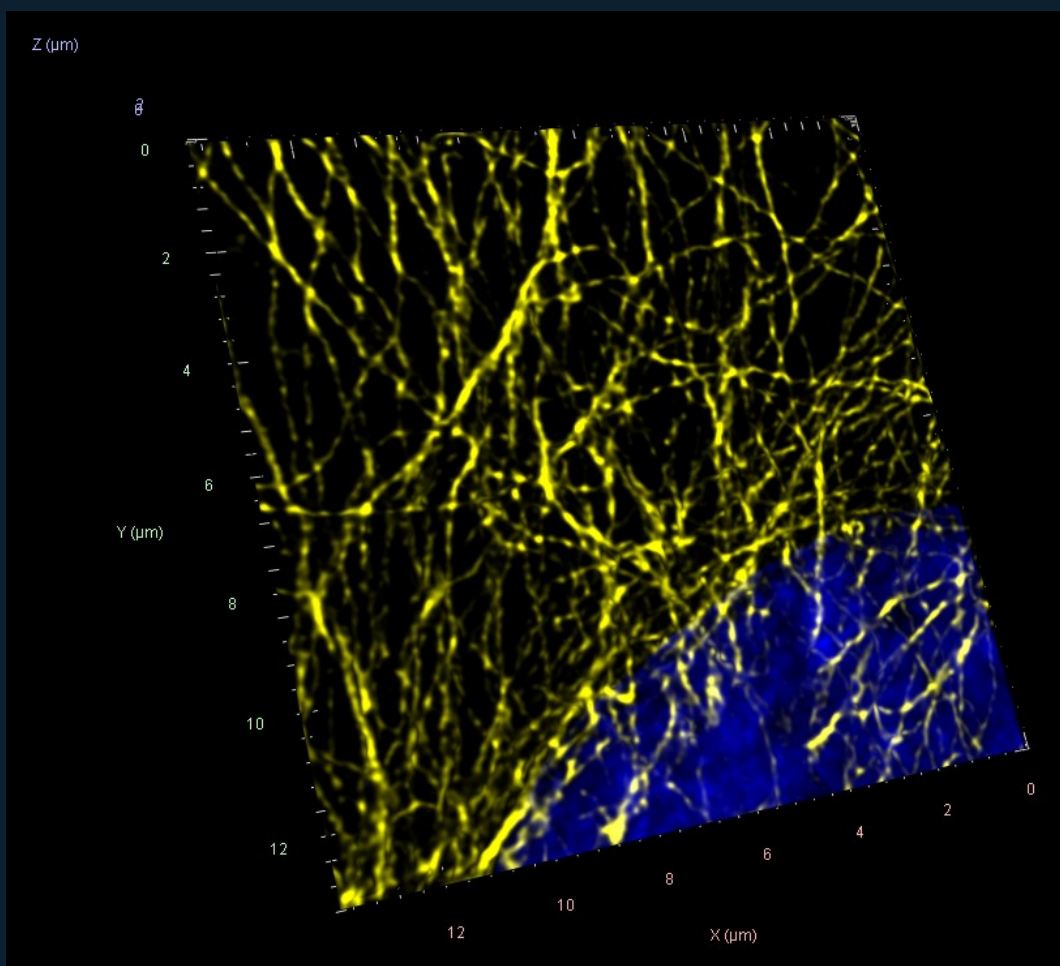

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SIM²

SA-000029 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse, F(ab')₂ (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Control

BENNETT ASSESSMENT

The preserved SA-000029 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / NIR-BOUNDED PARTIAL STEPDOWN

The preserved SA-000029 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

SETTINGS ASSESSMENT

The record preserves 7 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	532	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 532; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335 - 383; 555 - 595; 420 - 470; 528
Widefield SIM — Apotome	405/DAPI; 532; 555; 647	3; AF532-49014; 50 Cy5; 49 DAPI; 515 - 545; 625 - 655; 335-383; 555 - 695
Confocal	405/DAPI; 488; 532; 750	2; None; 488nm @0.9%; 405nm @1.0%; 528; 353; 553; 465
Airyscan	405/DAPI; 488; 532; 647	3; None; 488nm @1.0%; 640nm @0.5%; 405nm @1.0%; 528; 653; 353
SIM	405/DAPI; 488; 532	2; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T1-Axiocam 820-SR; T2-Axiocam 820 #2-SR; 488; 405
SIM ²	405/DAPI; 488; 532	2; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T1-Axiocam 820-SR; T2-Axiocam 820 #2-SR; 488; 405
Control	405/DAPI; 532; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335 - 383; 555 - 595; 420 - 470; 528

MODALITY ASSESSMENT / PART 1 OF 8

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 8

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.102µm X 0.102µm; Image size (Pixels): 1216 X 1028; Image Size (Pixels): 1216 X 1028; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Axiocam 705. Timing: Exposure Time: 50ms; 20ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @6.0%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 30.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 8

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Z-Stack: 24 Slices (4.60µm); Channels: 3; Scaling (Per Pixel): 0.143µm X 0.143µm X 0.200µm; Image size (Pixels): 636 X 563; Image Size (Pixels): 636 X 563; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807. Timing: Exposure Time: 50ms; 200ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @60%. Optical/scan: Effective NA: 1.25. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 47.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 8

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.071m X 0.071µm X 0.070µm; Image size (Pixels): 1306 X 1306; Image Size (Pixels): 1306 X 1306; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; GaAsP-Pmt1; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.59µs; Frame Time: 4.75s. Illumination: Laser Wavelength: 488nm @0.9%; 405nm @1.0%. Optical/scan: Pinhole: 1.00AU / 45µm; 1.00AU / 37µm; Scan Zoom: 1.1; LSM Scan Speed: 8; Scan Direction: Bidirectional; Averaging: 4; Effective NA: 1.4. Structured source metadata values: 40.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 750.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 8

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 61 Slices (3.0µm); Channels: 3; Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm; Image size (Pixels): 2588 X 2588; 460 X 480; Image Size (Pixels): 2588 X 2588; 460 X 480; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.81µs; Frame Time: 12.92s. Illumination: Laser Wavelength: 488nm @1.0%; 640nm @0.5%. Optical/scan: Pinhole: 5.00 AU / 229µm; 5.00 AU / 280µm; Scan Zoom: 1.1; LSM Scan Speed: 6; Scan Direction: Bidirectional; Averaging: 2; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 8.7 (3D, Auto); Super Resolution 9.3 (3D, Auto). Structured source metadata values: 54.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 8

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 60 Slices (2.95µm); Channels: 2; Scaling (Per Pixel): 21.11nm X 21.11nm X 50.00nm; Image size (Pixels): 3429 X 3429; Image Size (Pixels): 3429 X 3429; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820; AxioCam 820 #2. Timing: Exposure Time: 100ms; 50ms; Frame Time: 1.33s; 676.00ms. Illumination: Laser Wavelength: 488nm @4.00%; 405nm @7.0%; Illumination Wavelength: 488; 405. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM; Result Sampling: 4; Projection: View Angle 45°, Scale Z 93%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 51.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 8

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 47 Slices (2.30µm); Channels: 2; Scaling (Per Pixel): 21.11nm X 21.11nm X 50.00nm; Image size (Pixels): 658 X 626; Image Size (Pixels): 658 X 626; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820; AxioCam 820 #2. Timing: Exposure Time: 100ms; 50ms; Frame Time: 1.33s; 676.00ms. Illumination: Laser Wavelength: 488nm @4.00%; 405nm @7.0%; Illumination Wavelength: 488; 405. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM2; Result Sampling: 4; Projection: View Angle 35°, Scale Z 10%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 56.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 8

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1298 X 1067; Image Size (Pixels): 1298 X 1067; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000029 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm -> 0.03530 μ m X 0.03530 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2588 X 2588; Image Size (Scaled)=91.35 μ m X 91.35 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 21.11nm X 21.11nm X 50.00nm -> 0.02111 μ m X 0.02111 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=658 X 626; Image Size (Scaled)=13.89 μ m X 13.22 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.04%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked
Stepdown order. Scientific wording is preserved; missing modalities are not invented.
A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000029. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra	PRESERVED
7	SIM ²	ZEISS Elyra	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

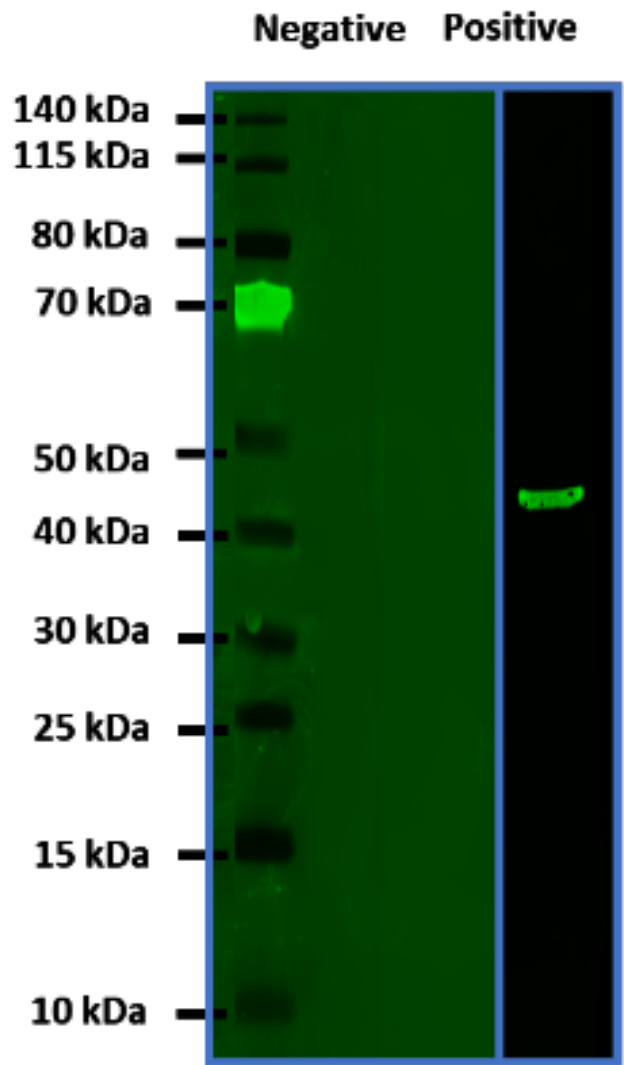

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	532
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	532nm
Emission Filter	572±28nm
Detector	Avalanche Photodiode
Intensity	L2
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods™

Identifyfyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Actin Monoclonal Antibody (mAbGEa) (MA1-744) for 2 hours, followed by 5X PBS washes. Identifyfyn® SRM Alexa Fluor™ 532 Affinity Purified Goat anti-Mouse F(ab)₂ was incubated for 1 hour and washed 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 532nm

Emission Filter = 572±28nm

Detector = Avalanche Photodiode

Intensity = L2

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: Identifyfyn® LLC

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000029
Stage	Western blot
Instrument	Azure Sapphire FL

SOURCE PROVENANCE

Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	B25E304F694CD7E771767E1C6227F3F37FE0128C28A12FDDB069C13192CF1DBA
Method PDF SHA-256	0C587DC43EB0BFD008FAF79F07945E24D9F060735E0D40756FBC3598B4D84DDA

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

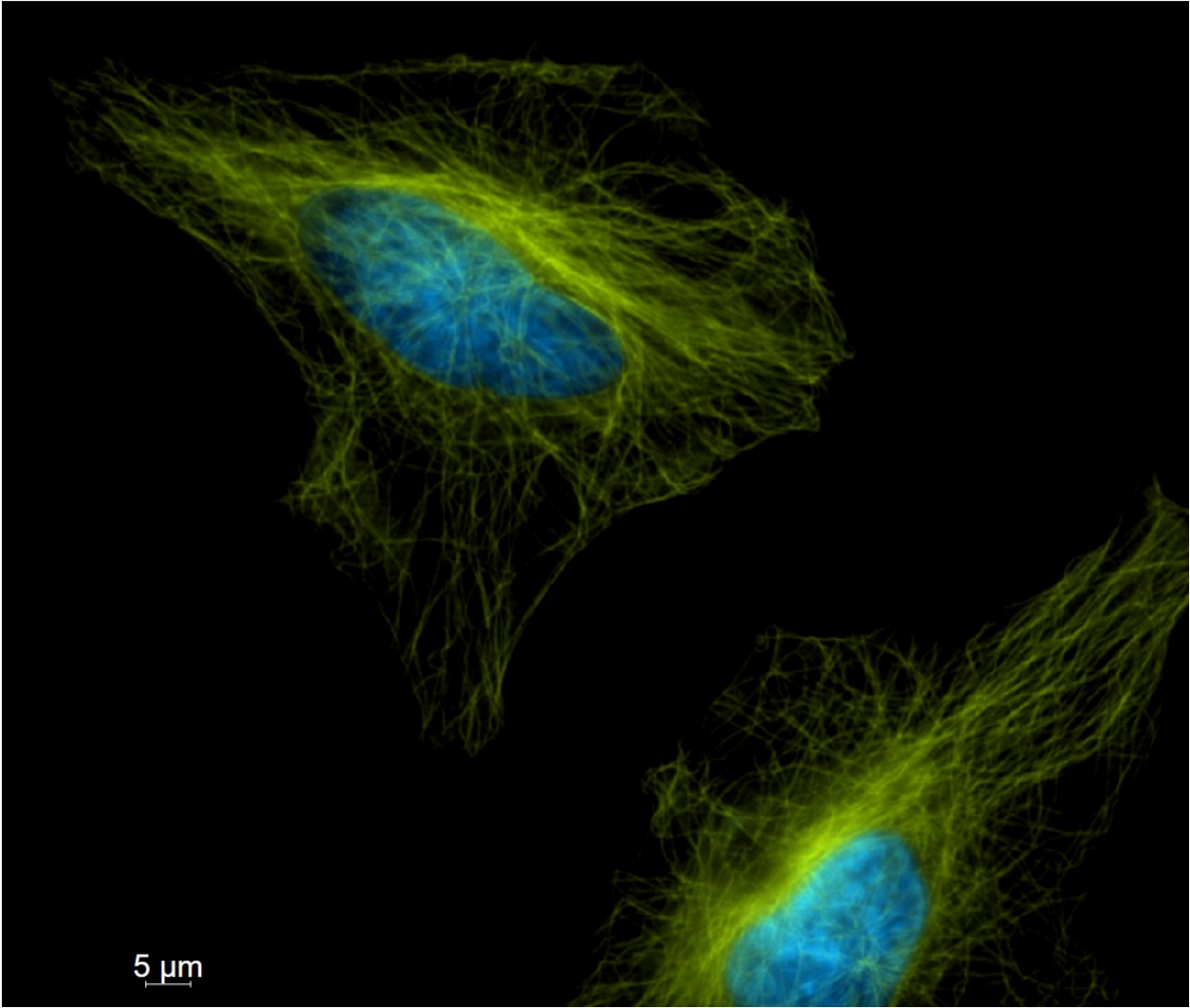

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	30

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.102µm X 0.102µm
Image size (Pixels)	1216 X 1028
Image Size (Scaled)	113.18µm X 112.59µm
Bit Depth	16 Bit
Image Center Position	X: 0.00mm, Y: 0.00mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335 - 383
Filter Emission Wavelength	555 - 595 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @6.0% X-Cite Xylis Lamp Intensity @10%
Exposure Time	50ms 20ms
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Imaging Device	Axiocam 705
Camera Adapter	1X
Depth of Focus	0.86µm 0.72µm
Binning Mode	2,2

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.102µm X 0.102µm

Image size (Pixels) = 1216 X 1028

Image Size (Scaled) = 113.18µm X 112.59µm

Bit Depth = 16 Bit

Image Center Position = X: 0.00mm, Y: 0.00mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @6.0%

Exposure Time = 50ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.4

Imaging Device = Axiocam 705

Camera Adapter = 1X

Depth of Focus = 0.86µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 20ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Imaging Device = AxioCam 705
 Camera Adapter = 1X
 Depth of Focus = 0.72µm
 Binning Mode = 2,2

Post Processing

None

Image Correction

None


Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	30
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	F2FD064F3473487F01945E578668164FDA1C5E527BBF53818A95FA91590EEB77
Method PDF SHA-256	3C85F0666C640459D00B2C02A6DD50FD2B47D89D9E7D741BDCBC1066DBD52460

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

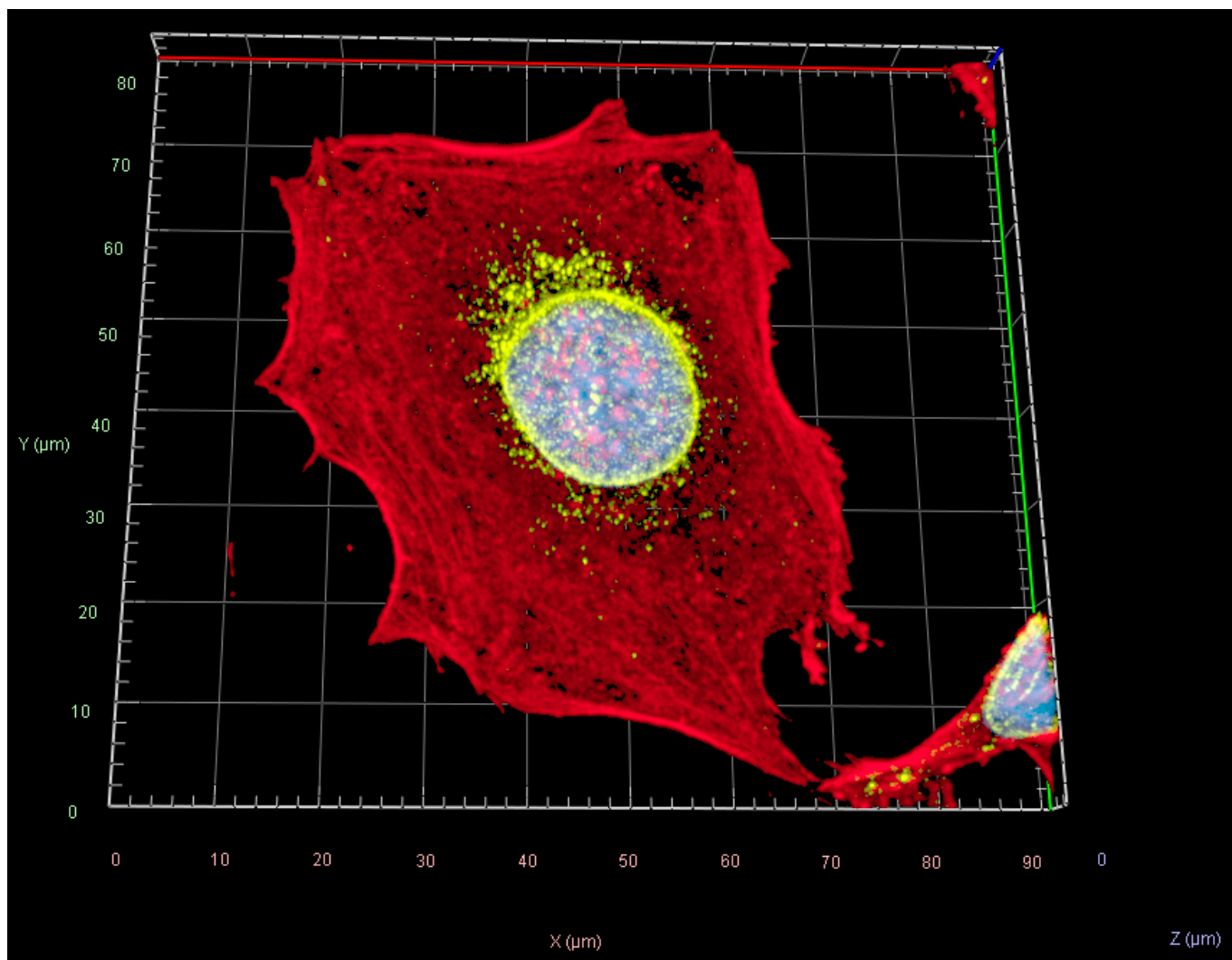

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555; 647
METADATA VALUES	47

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Z-Stack	24 Slices (4.60µm)
Channels	3
Scaling (Per Pixel)	0.143µm X 0.143µm X 0.200µm
Image Size (Pixels)	636 X 563
Image Size (Scaled)	90.86µm X 80.43µm
Bit Depth	14 Bit
Image Center Position	X: 50.20mm, Y: 49.90mm
Stage Position	X: 50.20mm, Y: 49.90mm
ROI Center Offset	X: 1.14, Y: 0.00µm
Reflector	AF532-49014 50 Cy5 49 DAPI
Beam Splitter	550 660 395
Filter Excitation Wavelength	515 - 545 625 - 655 335-383
Filter Emission Wavelength	555 - 695 665 - 715 420-470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @60% X-Cite Xylis Lamp Intensity @8%
Exposure Time	50ms 200ms 20ms
Excitation Wavelength	528 653 353
Emission Wavelength	553 668 465
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.07µm 1.30µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.
3D Maximum Projection	Precise
Appearance	Threshold 13% (532), Threshold 13% (647), Threshold 13% (DAPI), Ramp 0% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness 100%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Donkey Anti-Rabbit IgG (H + L)

Cat ID# - SA 000060

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Nuclear Pore Complex Monoclonal Antibody (39C7), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Actin Recombinant Rabbit Monoclonal Antibody (JF47-01), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 24 Slices (4.60µm)

Channels = 3

Scaling (Per Pixel) = 0.143µm X 0.143µm X 0.200µm

Image Size (Pixels) = 636 X 563

Image Size (Scaled) = 90.86µm X 80.43µm

Bit Depth = 14 Bit

Image Center Position = X: 50.20mm, Y: 49.90mm

Stage Position = X: 50.20mm, Y: 49.90mm

ROI Center Offset = X: 1.14, Y: 0.00µm

532 -

Reflector = AF532-49014
 Beam Splitter = 550
 Filter Excitation Wavelength = 515 - 545
 Filter Emission Wavelength = 555 - 695
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @25%
 Exposure Time = 50ms
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 1.07µm
 Binning Mode = 2,2

647 -

Reflector = 50 Cy5
 Beam Splitter = 660
 Filter Excitation Wavelength = 625 - 655
 Filter Emission Wavelength = 665 - 715
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @60%
 Exposure Time = 200ms
 Excitation Wavelength = 653
 Emission Wavelength = 668
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 1.30µm
 Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335-383
 Filter Emission Wavelength = 420-470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @8%
 Exposure Time = 20ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Correction, Fourier Filter, and Deconvolution were NOT used.

3D Image Processing

3D Maximum Projection = Precise

Appearance = Threshold 13% (532), Threshold 13% (647), Threshold 13% (DAPI), Ramp 0% all channels, Maximum 100% all channels.

Background = Black

Light = Brightness 100%, Enabled Tone Mapping

Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	47
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	1A4626ABEE0E75356C1BEF8FC28238AD793CFCEC16C0229469FC4162368DD006
Method PDF SHA-256	941DEEB628BFF6C41A068AACCE70175B5305516C2A62F798CBB70D901D345DA9

ORIGINAL IMAGE EVIDENCE / CONFOCAL

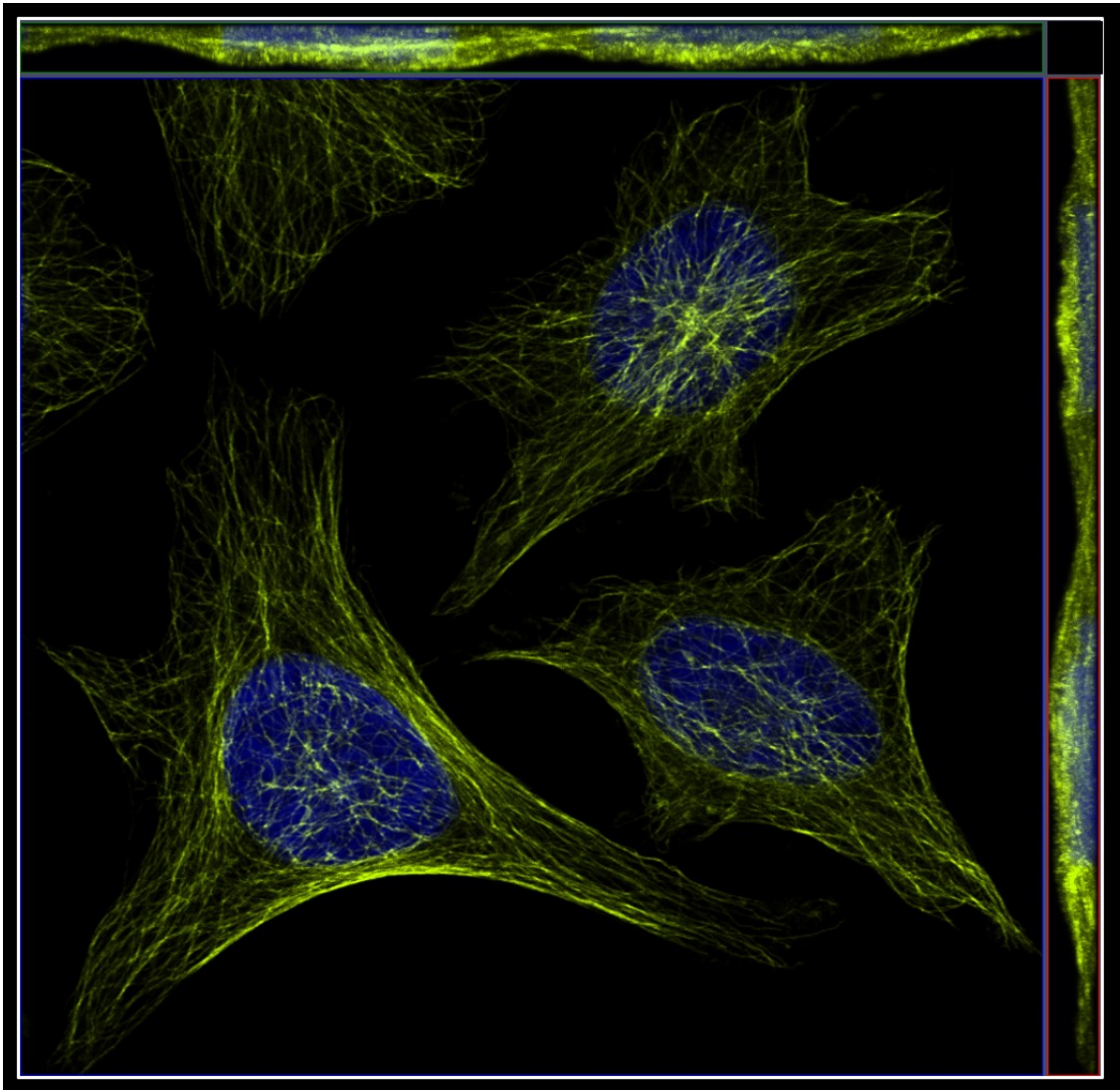

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 750
METADATA VALUES	40

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Channels	2
Scaling (Per Pixel)	0.071m X 0.071µm X 0.070µm
Image size (Pixels)	1306 X 1306
Image Size (Scaled)	92.19µm X 92.19µm
Bit Depth	16 Bit
Image Center Position	X: 86.48mm, Y: 51.12mm
Stage Position	X: 86.48mm, Y: 51.12mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 45µm 1.00AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @0.9% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.1
Rotation	0°
Sampling	1.00
Pixel Time	0.59µs
Frame Time	4.75s
LSM Scan Speed	8
Scan Direction	Bidirectional
Line Step	1
Averaging	4
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Detection Wavelength	500 - 575 410 - 467
Imaging Device	Airyscan GaAsP-Pmt1
Detector Type	GaAsP-PMT
Detector Gain	750 V 650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.071m X 0.071µm X 0.070µm

Image size (Pixels) = 1306 X 1306

Image Size (Scaled) = 92.19µm X 92.19µm

Bit Depth = 16 Bit

Image Center Position = X: 86.48mm, Y: 51.12mm

Stage Position = X: 86.48mm, Y: 51.12mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 45µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @0.9%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.1
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.59µs
 Frame Time = 4.75s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 500 - 575
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 750 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.1
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.59µs
 Frame Time = 4.75s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410 - 467
 Imaging Device = GaAsP-Pmt1
 Detector Type = GaAsP-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

Orthogonal View

Maximal Intensity Projection (MIP)

Image Corrections

None

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000029
Stage	Confocal
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

SOURCE PROVENANCE

Structured source metadata values	40
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	5C32C341105342FC3A658B1104CF3E30028C95A89D031188AB8D996319083804
Method PDF SHA-256	415B9934A373EB6206D2D88AFA1CFEC5C203C8C607D87CAA1139FA9949AD750C

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

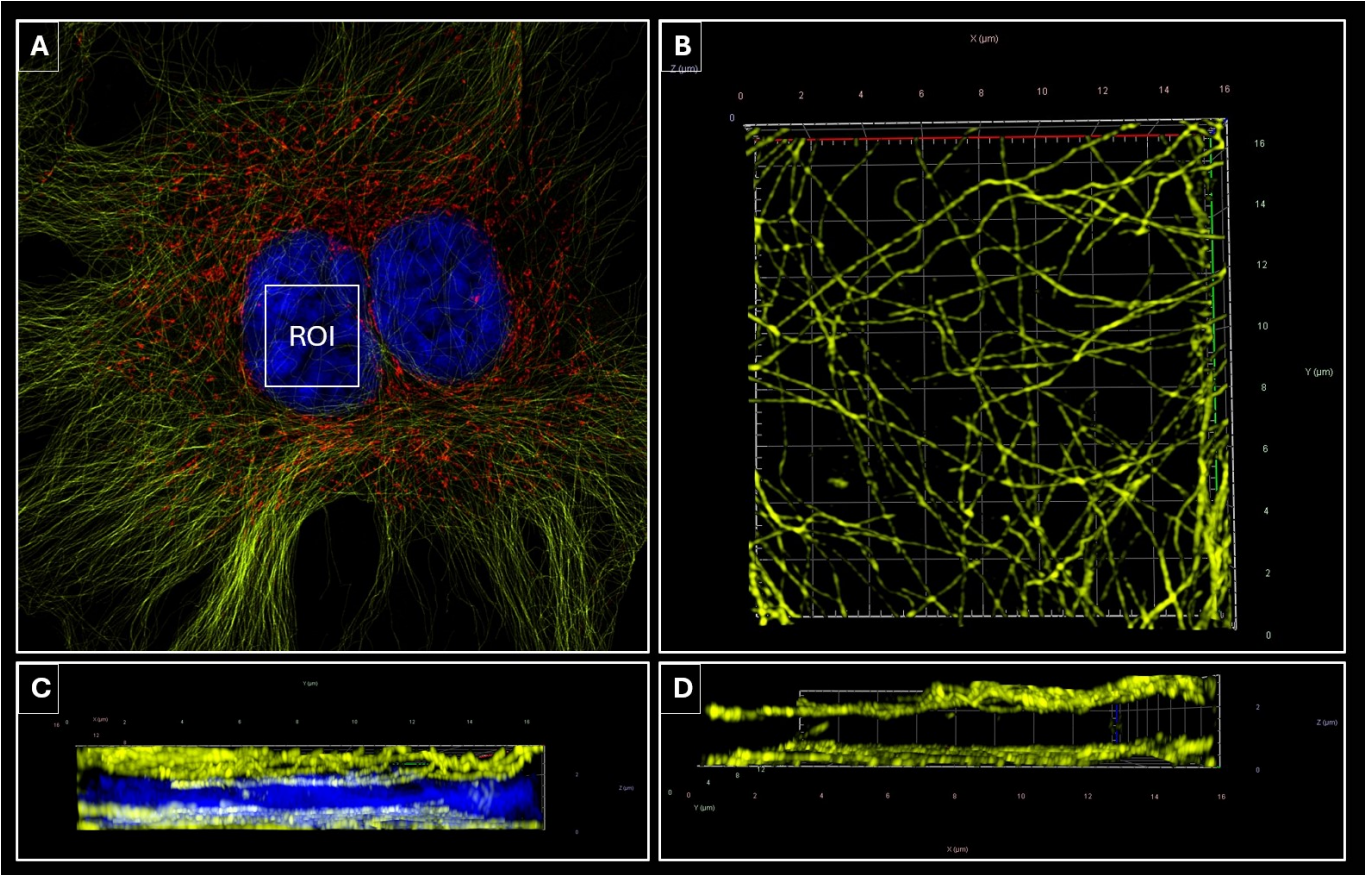

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 647
METADATA VALUES	54

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	61 Slices (3.0µm)
Channels	3
Scaling (Per Pixel)	0.03530 µm X 0.03530 µm X 0.05000 µm
Image Size (Pixels)	2588 X 2588 460 X 480
Image Size (Scaled)	91.35µm X 91.35µm 16.24µm X 16.94µm
Bit Depth	16 Bit
Image Center Position	X: 76.88mm, Y: 54.74mm
Stage Position	X: 76.88mm, Y: 54.74mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00 AU / 229µm 5.00 AU / 280µm 4.97 AU / 181µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @1.0% 640nm @0.5% 405nm @1.0%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.1
Rotation	0°
Sampling	2.00
Pixel Time	0.81µs
Frame Time	12.92s
LSM Scan Speed	6
Scan Direction	Bidirectional
Line Step	1
Averaging	2
Excitation Wavelength	528 653 353
Emission Wavelength	553 668 465
Effective NA	1.4
Detection Wavelength	515-580 635-700 400-465
Imaging Device	AiryScan
Detector Type	GaAsP-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	850 V
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 8.7 (3D, Auto) Super Resolution 9.3 (3D, Auto) Super Resolution 8.4 (3D, Auto)
3D Maximum Projection	Precise
Appearance	Threshold 1% (532), Threshold 1% (DAPI), Ramp 0% (532), Ramp 0% (DAPI) Maximum 100% all channels
Background	Black
Light	Brightness 100%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Donkey Anti-Rabbit, IgG (H + L)

Cat ID# - SA 000060

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Rabbit Recombinant Antibody (MTC02,2860R), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant catalog number P36961.

Microscope

Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D) - Panel A, Z-Plane 12 of 61 shown

Z-Stack = 61 Slices (3.0µm)

Channels = 3

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 2588 X 2588

Image Size (Scaled) = 91.35µm X 91.35µm

Bit Depth = 16 Bit

Image Center Position = X: 76.88mm, Y: 54.74mm

Stage Position = X: 76.88mm, Y: 54.74mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

Z Stack - (3D) - Panels B,C and D

Z-Stack = 61 Slices (3.0µm)

Channels = 3

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 460 X 480

Image Size (Scaled) = 16.24µm X 16.94µm

Bit Depth = 16 Bit

Image Center Position = X: 76.88mm, Y: 54.74mm

Stage Position = X: 76.88mm, Y: 54.74mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 229µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.1
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.81µs
 Frame Time = 12.92s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 515-580
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.7 (3D, Auto)

647 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 280µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 640nm @0.5%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.1
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.81µs
 Frame Time = 12.92s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 653
 Emission Wavelength = 668
 Effective NA = 1.4
 Detection Wavelength = 635-700
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 9.3 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 4.97 AU / 181µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.1
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.81µs
 Frame Time = 12.92s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400-465
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.4 (3D, Auto)

Post Processing - AiryScan 3D Processing

Display Mode "SR"
 Super Resolution Auto Filter Selected
 Process 3D Current Image
 Created

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 1% (532), Threshold 1% (DAPI), Ramp 0% (532), Ramp 0% (DAPI) Maximum 100% all channels
 Background = Black
 Light = Brightness 100%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Summary

Panel A displays Z-Plane 12 of 61 for all channels, 405 (DAPI), 532 (Alpha Tubulin), and 647 (Mito), noting the Region of Interest (ROI) chosen to be presented in Panels B, C, and D. Panel B is 3-Dimensional Maximum Projection of the 532 channel (Alpha tubulin). Panels C and D show the 3-dimensional Maximum Projection rotated to its Z-axis with Panel C, including DAPI and Panel D, excluding DAPI.

Image Corrections -

None


Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	54
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	28839CBCB7399E31C598DE083B634FD8EC68884F981A240CACA5EC857FEB5F40
Method PDF SHA-256	3DF6EB86BA7AE5E7E314CF42FA4882DA7D99F6A8213311204C908DF178C00F82

ORIGINAL IMAGE EVIDENCE / SIM

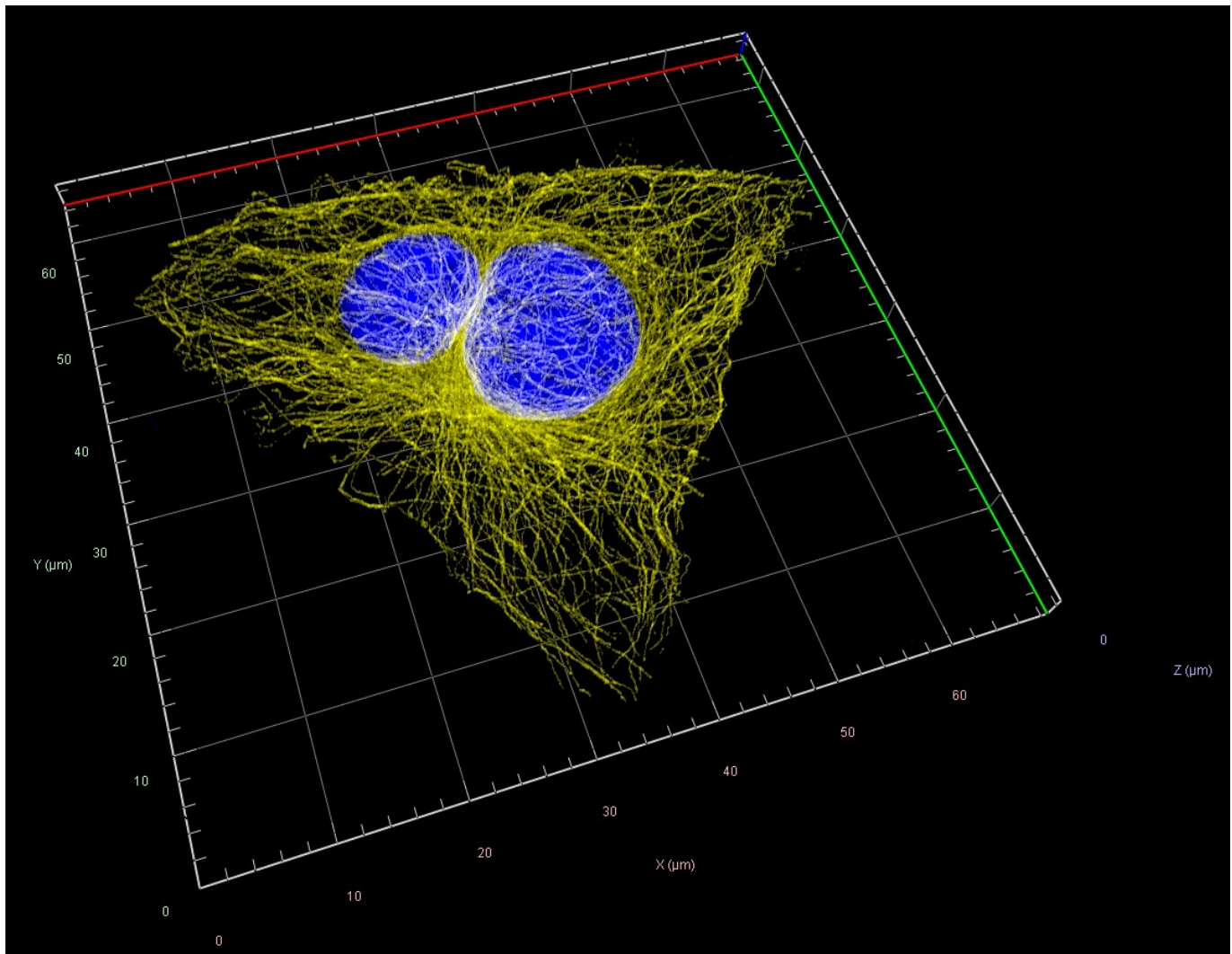

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	51

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	60 Slices (2.95µm)
Channels	2
Scaling (Per Pixel)	21.11nm X 21.11nm X 50.00nm
Image Size (Pixels)	3429 X 3429
Image Size (Scaled)	68.60µm X 64.73µm
Bit Depth	16 Bit
Image Center Position	X: 59.24mm, Y: 34.27mm
Stage Position	X: 59.24mm, Y: 34.27mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	488 405
Channel Name	T1-Axiocam 820-SR T2-Axiocam 820 #2-SR
Excitation Wavelength	488 405
Emission Wavelength	525 445
Effective NA	1.4
Detection Wavelength	419-470, 503-546, 576-617, 657-800 419-470, 503 - 546, 576 - 617, 657-800
Imaging Device	Axiocam 820 Axiocam 820 #2
Camera Adapter	1X
Exposure Time	100ms 50ms
Depth of Focus	0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	488nm @4.00% 405nm @7.0%
Frame Time	1.33s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Algorithm	SIM
SIM Settings	Standard

FIELD	SOURCE-DOCUMENTED VALUE
Auto Sharpness	On
Sharpness	11.3472 (532), 9.6738 (DAPI)
Fitting	Fast Fit
Histogram	Scale to Min/Max
Baseline	Cut
3D Maximum Projection	Precise
Appearance	Threshold 4% (532), Threshold 4% DAPI), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 115%, Enabled Tone Mapping
Projection	View Angle 45°, Scale Z 93%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

HeLa Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 60 Slices (2.95µm)

Channels = 2

Scaling (Per Pixel) = 21.11nm X 21.11nm X 50.00nm

Image Size (Pixels) = 3429 X 3429

Image Size (Scaled) = 68.60µm X 64.73µm

Bit Depth = 16 Bit

Image Center Position = X: 59.24mm, Y: 34.27mm

Stage Position = X: 59.24mm, Y: 34.27mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T1-Axiocam 820-SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503-546, 576-617, 657-800
 Imaging Device = Axiocam 820
 Camera Adapter = 1X
 Exposure Time = 100ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @4.00%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 32.0µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503 - 546, 576 - 617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @7.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM Processing

Processing Dimension - 3D

Result Sampling = 4

Channels = 2

Algorithm = SIM

SIM Settings = Standard

Auto Sharpness = On

Sharpness = 11.3472 (532), 9.6738 (DAPI)

Fitting = Fast Fit

Histogram = Scale to Min/Max

Baseline = Cut

3D Image Processing

3D Maximum Projection = Precise

Appearance = Threshold 4% (532), Threshold 4% DAPI, Ramp 0% all channels, Maximum 100% all channels

Background = Black

Light = Brightness 115%, Enabled Tone Mapping

Projection = View Angle 45°, Scale Z 93%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	SIM
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	51
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	F37B3AA0FAD07EA3682A507129F81F0D3E1036EDCEEE04644BB3F8185B60F6AE
Method PDF SHA-256	9CDC12D970E0E82DCF446834D835B39D6CC0DF55836F2EF2E0227D67EDFEBD79

ORIGINAL IMAGE EVIDENCE / SIM²

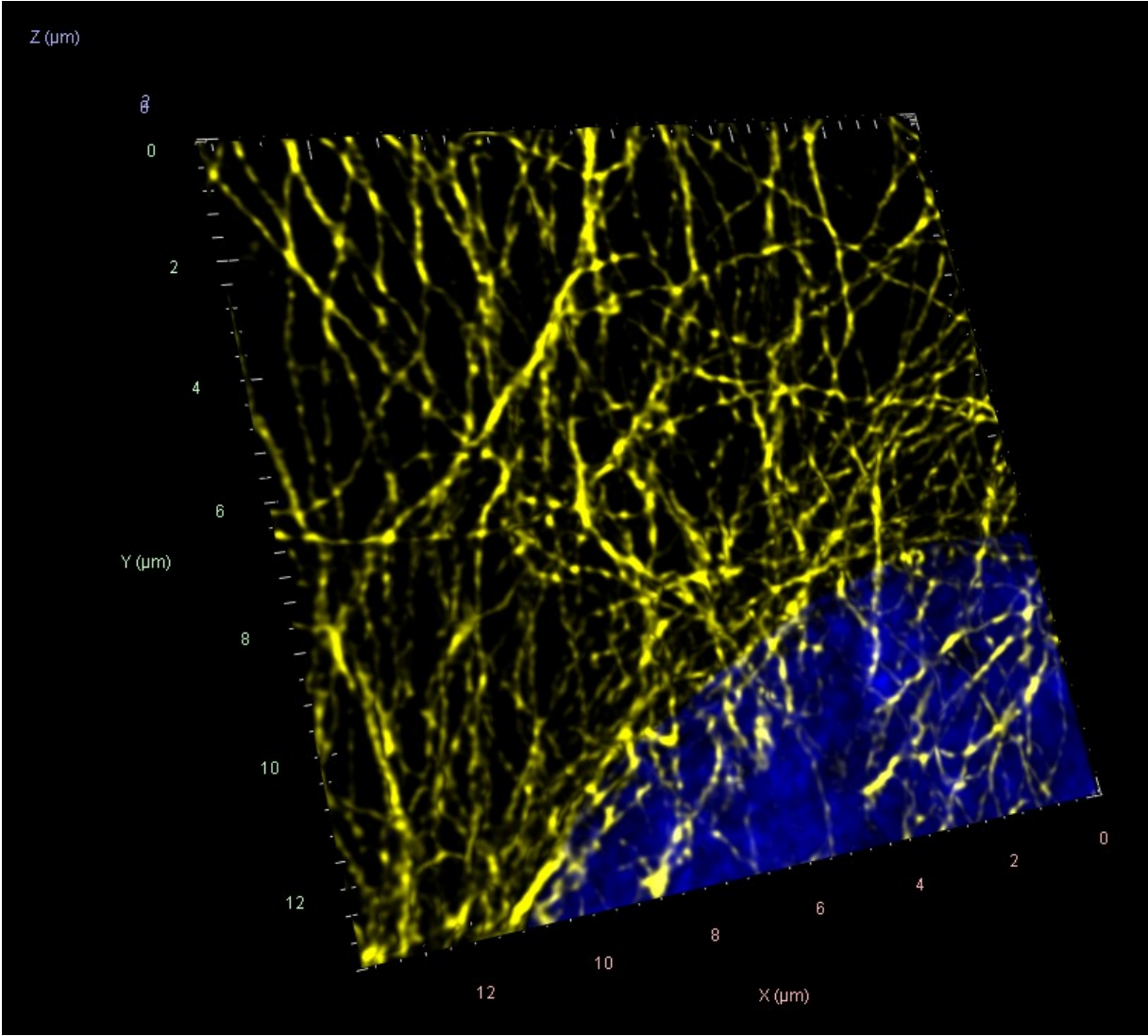

EVIDENCE TIER	SIM²
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	56

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	47 Slices (2.30µm)
Channels	2
Scaling (Per Pixel)	0.02111 µm X 0.02111 µm X 0.05000 µm
Image Size (Pixels)	658 X 626
Image Size (Scaled)	13.89µm X 13.22µm
Bit Depth	16 Bit
Image Center Position	X: 58.58mm, Y: 35.91mm
Stage Position	X: 58.58mm, Y: 35.91mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	488 405
Channel Name	T1-Axiocam 820-SR T2-Axiocam 820 #2-SR
Excitation Wavelength	488 405
Emission Wavelength	525 445
Effective NA	1.4
Detection Wavelength	419-470, 503-546, 576-617, 657-800 419-470, 503 - 546, 576 - 617, 657-800
Imaging Device	Axiocam 820 Axiocam 820 #2
Camera Adapter	1X
Exposure Time	100ms 50ms
Depth of Focus	0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	488nm @4.00% 405nm @7.0%
Frame Time	1.33s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Processing Sampling	4
Algorithm	SIM2

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	10.3999 (532), 8.9597 (DAPI)
Order Combination	Weiner Filter
Regularization	None
Regularization Weight	0.0000
Iterations	3
Filter	No Filter
Fitting	Fast Fit
Histogram	Scale to Min/Max
3D Maximum Projection	Precise
Appearance	Threshold 1% (532), Threshold 1% DAPI), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 114%, Tone Mapping is NOT Enabled
Projection	View Angle 35°, Scale Z 10%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

HeLa Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Alpha Tubulin Monoclonal Antibody (DM1A), was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 47 Slices (2.30µm)

Channels = 2

Scaling (Per Pixel) = 21.11nm X 21.11nm X 50.00nm

Image Size (Pixels) = 658 X 626

Image Size (Scaled) = 13.89µm X 13.22µm

Bit Depth = 16 Bit

Image Center Position = X: 58.58mm, Y: 35.91mm

Stage Position = X: 58.58mm, Y: 35.91mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T1-Axiocam 820-SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503-546, 576-617, 657-800
 Imaging Device = Axiocam 820
 Camera Adapter = 1X
 Exposure Time = 100ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @4.00%
 Frame Time = 1.33s
 Scan Direction = Unidirectional
 Grating 32.0µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T2-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419-470, 503 - 546, 576 - 617, 657-800
 Imaging Device = Axiocam 820 #2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @7.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM2 Processing

Processing Dimension - 3D
 Result Sampling = 4
 Processing Sampling = 4
 Channels = 2
 Algorithm = SIM2
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 10.3999 (532), 8.9597 (DAPI)
 Order Combination = Weiner Filter
 Regularization = None
 Regularization Weight = 0.0000
 Iterations = 3
 Filter = No Filter
 Fitting = Fast Fit
 Histogram = Scale to Min/Max

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 1% (532), Threshold 1% DAPI), Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 114%, Tone Mapping is NOT Enabled
 Projection = View Angle 35°, Scale Z 10%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	SIM ²
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	56
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	5761660B299544BDD049532123DB017A5FD86547BFADA21B9A0AAC5F124F824C
Method PDF SHA-256	B28039B18AA8EC801A95FB5CDEED18FAEA0DFF2A00260D7EAFE3B8A96F1F9A0C

ORIGINAL IMAGE EVIDENCE / CONTROL

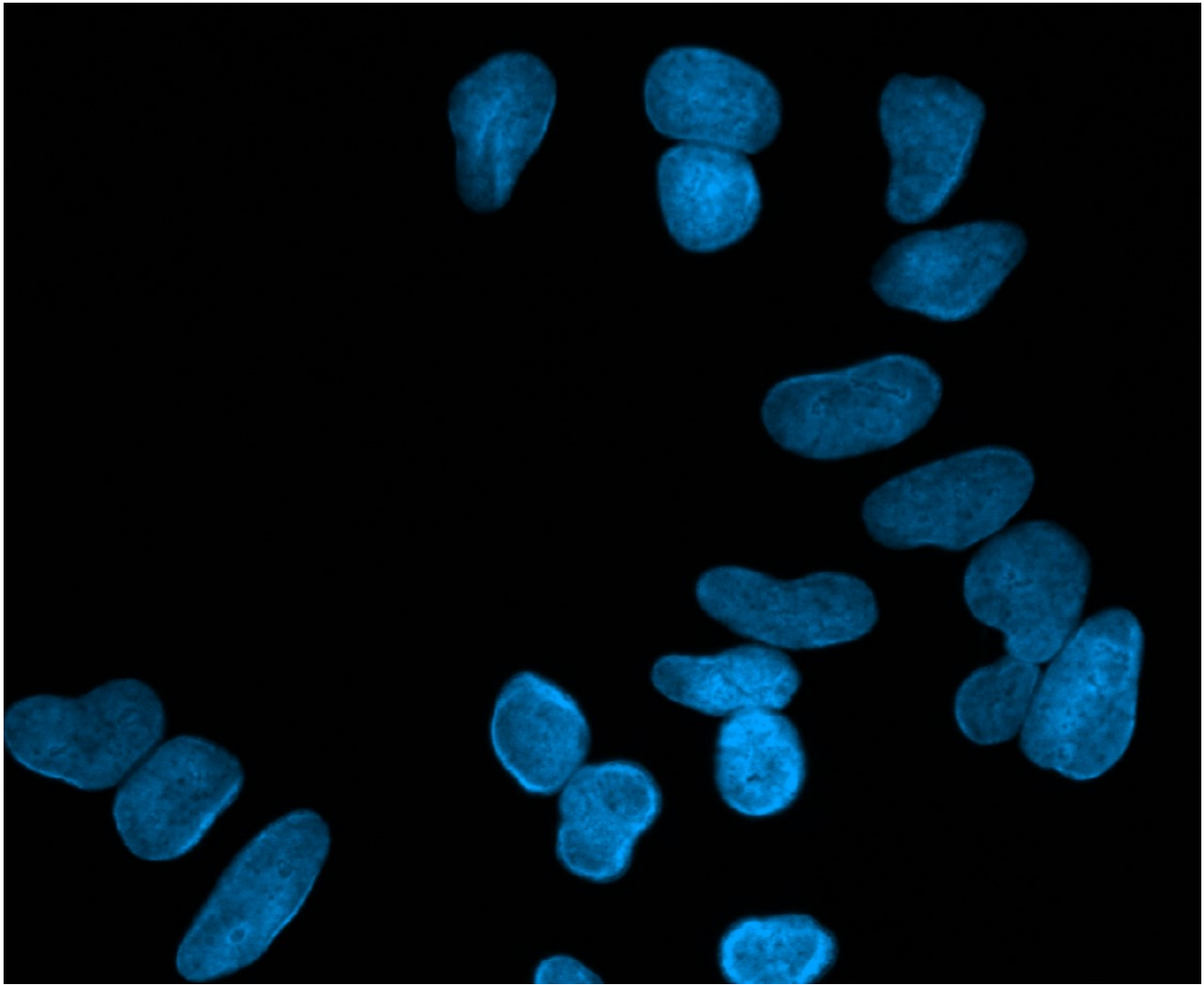

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1298 X 1067
Image Size (Scaled)	185.43µm X 152.43µm
Bit Depth	14 Bit
Image Center Position	X: 51.52mm, Y: 51.66mm
Stage Position	X: 51.52mm, Y: 51.66mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335 - 383
Filter Emission Wavelength	555 - 595 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.07µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Mouse F(ab')₂ (H + L)

Cat ID# - SA 000029

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1298 X 1067

Image Size (Scaled) = 185.43µm X 152.43µm

Bit Depth = 14 Bit

Image Center Position = X: 51.52mm, Y: 51.66mm

Stage Position = X: 51.52mm, Y: 51.66mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.07µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000029
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	BBDD0369F59FD8C819139A9614699CE934971E86E3F5C9CD57F9A716D165B880
Method PDF SHA-256	32310119029CA5A15AF151FD0F8DE6A66EC6FAB0E370AB7098138EB981A7E2F5