

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human, IgG (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM² -> SINGLE-MOLECULE STORM

8


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

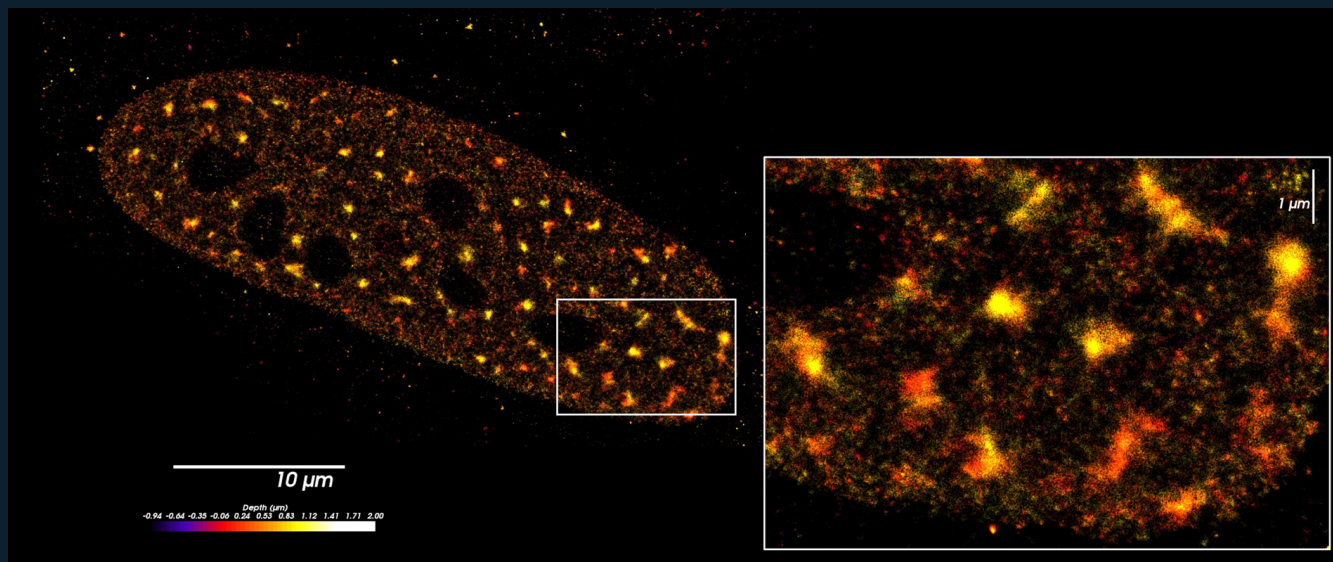

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SINGLE-MOLECULE STORM

SA-000058 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Single-molecule STORM presentation workflow
- Preserved control experiment
- Treatment-response experimental workflow
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human, IgG (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Single-molecule STORM -> Control

BENNETT ASSESSMENT

The preserved SA-000058 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / FULL STEPDOWN / SOURCE-BOUNDED REVIEW

The preserved SA-000058 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

SETTINGS ASSESSMENT

The record preserves 8 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	647	Not deterministically stated in the acquisition block
Widefield	405/DAPI; 532; 555; 647	2; 50 Cy5; AF532-49014; 625 - 655; 515 - 545; 665 - 715; 555 - 595; 653
Widefield SIM — Apotome	405/DAPI; 555; 647	2; 50 Cy5; AF532-49014; 625 - 655; 515 - 545; 665 - 715; 555 - 595; 653
Confocal	405/DAPI; 488; 532; 647	3; None; 640nm @0.2%; 488nm @1.0%; 405nm @0.2%; 653; 528; 353
Airyscan	405/DAPI; 488; 555; 647	3; None; 640nm @0.9%; 561nm @0.2%; 405nm @0.6%; 653; 553; 353
SIM	405/DAPI; 488; 532; 647	3; 1; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T1-Axiocam 820 -SR; T2-Axiocam 820#2 -SR; T3-Axiocam 820 #2-SR
SIM ²	405/DAPI; 488; 532; 647	3; 1; SR HE LBF 405/488/561/642; 200 - 405, 200 - 488, 200 - 561, 200 - 642; 419 - 470, 503 - 546, 576 - 617, 657 - 800; T1-Axiocam 820 -SR; T2-Axiocam 820#2 -SR; T3-Axiocam 820 #2-SR
Single-molecule STORM	647	Both Channels
Control	405/DAPI; 647	2; 50 Cy5; 49 DAPI; 625 - 655; 335 - 383; 665 - 715; 420 - 470; 653

MODALITY ASSESSMENT / PART 1 OF 9

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 9

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.110µm X 0.110µm; Image size (Pixels): 878 X 976; Image Size (Pixels): 878 X 976; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam MR R3. Timing: Exposure Time: 150ms; 80ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @20%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 29.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 405/DAPI; 532; 555; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 9

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Z-Stack: 113 Slices (11.2µm); Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm X 0.100µm; Image size (Pixels): 1283 X 864; Image Size (Pixels): 1283 X 864; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807. Timing: Exposure Time: 40ms; 35ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @6.0%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 555; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 9

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 11 Slices (700nm); Channels: 3; Scaling (Per Pixel): 0.071µm X 0.071µm X 0.050µm; Image size (Pixels): 1437 X 1437; Image Size (Pixels): 1437 X 1437; Bit Depth: 16 Bit. Detector/camera: Imaging Device: GaAsP-PMT 2; Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.54µs; Frame Time: 5.23s. Illumination: Laser Wavelength: 640nm @0.2%; 488nm @1.0%. Optical/scan: Pinhole: 1.00AU / 57µm; 1.00AU / 46µm; Scan Zoom: 1.0; LSM Scan Speed: 8; Scan Direction: Bidirectional; Averaging: 4; Effective NA: 1.4. Processing: Projection: View Angle 52°, Scale Z 244%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 52.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 9

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 61 Slices (3.0µm); Channels: 3; Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm; Image size (Pixels): 1657 X 1622; Image Size (Pixels): 1657 X 1622; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.96µs; Frame Time: 10.94s. Illumination: Laser Wavelength: 640nm @0.9%; 561nm @0.2%. Optical/scan: Pinhole: 5.00 AU / 282µm; 5.00 AU / 250µm; Scan Zoom: 1.3; LSM Scan Speed: 6; Scan Direction: Bidirectional; Averaging: 2; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 9.0 (3D, Auto); Super Resolution 8.5 (3D, Auto). Structured source metadata values: 56.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 555; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 9

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 29 Slices (1.40µm); Channels: 3; 1; Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm; Image size (Pixels): 3486 X 3954; Image Size (Pixels): 3486 X 3954; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820; AxioCam 820#2. Timing: Exposure Time: 100ms; 120ms; Frame Time: 1.33s; 1.59s. Illumination: Laser Wavelength: 640nm @2.0%; 488nm @8.0%; Illumination Wavelength: 640; 488. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM; Process Sampling: 2; Result Sampling: 4; Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 64.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 9

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra. Objective: Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:. Acquisition: Z-Stack: 40 Slices (1.95µm); Channels: 3; 1; Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm; Image size (Pixels): 2306 X 2107; Image Size (Pixels): 2306 X 2107; Bit Depth: 16 Bit. Detector/camera: Imaging Device: AxioCam 820; AxioCam 820#2. Timing: Exposure Time: 100ms; 120ms; Frame Time: 1.33s; 1.59s. Illumination: Laser Wavelength: 640nm @2.0%; 488nm @8.0%; Illumination Wavelength: 640; 488. Optical/scan: Scan Direction: Unidirectional; Effective NA: 1.4. Processing: Algorithm: SIM2; Process Sampling: 4; Result Sampling: 4; Projection: View Angle 35°, Scale Z 32%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 64.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 9

Single-molecule STORM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Bruker Vutara VXL. Objective: Olympus 60X/1.30 silicon oil objective. Platform: Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:.

Detector/camera: Camera: Both Cameras. Timing: Exposure Time: 20ms Selected; Super Resolution Camera Exposure Time: 20ms; Widefield Camera Exposure Time: 100ms. Illumination: Photoactivation/Imaging Laser: 0% Selected.

Structured source metadata values: 71.

VISIBLE OBSERVATION

A preserved presentation image is available for Single-molecule STORM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Single-molecule STORM, documents platform context as Bruker Vutara VXL, and records detected dye/channel descriptors as 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Without localization tables, uncertainty, blinking, drift, and independent cluster evidence, the presentation cannot establish molecule count, cluster truth, exact precision, or molecular architecture.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Single-molecule STORM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 9 OF 9

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, AxioCam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 998 X 1069; Image Size (Pixels): 998 X 1069; Bit Depth: 14 Bit. Detector/camera: Imaging Device: AxioCam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000058 record contains the full ordered evidence sequence from Western blot through single-molecule STORM. Each modality retains its own acquisition and processing boundary.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.30nm X 35.30nm X 50.00nm -> 0.03530 μ m X 0.03530 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1657 X 1622; Image Size (Scaled)=58.49 μ m X 57.25 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm -> 0.01320 μ m X 0.01320 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=3486 X 3954; Image Size (Scaled)=46.00 μ m X 52.17 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.04%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 13.20nm X 13.20nm X 50.00nm -> 0.01320 μ m X 0.01320 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=2306 X 2107; Image Size (Scaled)=30.43 μ m X 27.80 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.04%.

UNRESOLVED - PRESERVED AS CONFLICT

Cross-tier / treatment duration: {"Widefield": ["24 hours"], "Widefield SIM — Apotome": ["24 hours"], "Confocal": ["24 hours"], "Airyscan": ["24 hours"], "SIM": ["24 hours"], "SIM²": ["24 hours"], "Single-molecule STORM": ["6 hours"]} -> Preserved as modality-specific historical duration values. Genuinely different treatment durations are present across evidence stages, and no source or scientist correction resolves them. Supporting evidence: Only duration values from explicit etoposide-treatment sentences were classified; concentration, wavelength, resolution and pixel-size values were excluded.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked
Stepdown order. Scientific wording is preserved; missing modalities are not invented.
A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000058. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra	PRESERVED
7	SIM ²	ZEISS Elyra	PRESERVED
8	Single-molecule STORM	Bruker Vutara VXL	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

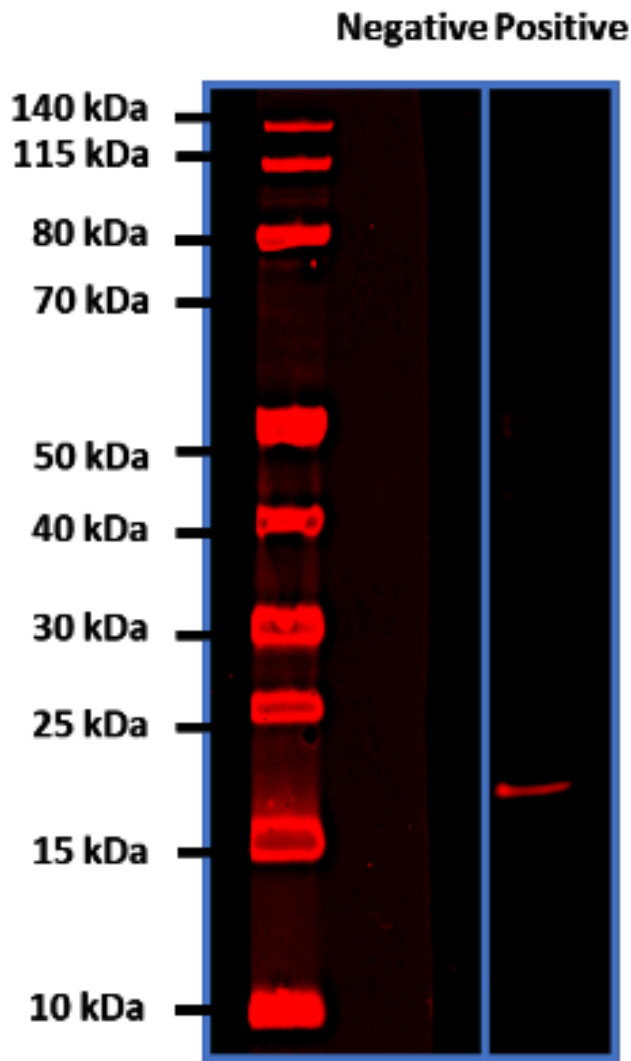

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	647
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	638nm
Emission Filter	676±37nm
Detector	Avalanche Photodiode
Intensity	L1
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)

Cat ID# - SA 000058

Etoposide-treated (2 μ M) HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Identifyn® Rabbit Anti-Human γ H2AX chimeric antibody (Human Fc) for 2 hours, followed by 5X PBS washes. Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat anti-Human IgG was incubated for 1 hour, followed by washing 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 638nm

Emission Filter = 676 $\pm$ 37nm

Detector = Avalanche Photodiode

Intensity = L1

Pixel size = 100 μ m

Scan Speed = High

Quality = High

Focus = Membrane (0 μ m in Z-axis)

Post Processing

None

Image: Identifyn® LLC

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000058

SOURCE PROVENANCE

Stage	Western blot
Instrument	Azure Sapphire FL
Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	F9A6828EB35D57B0A80EEC92C297AC19871C1DC7703452B987F114D6D62EF6E7
Method PDF SHA-256	5AE0D4166333CBF614FB34D6563A779050A3DD51C1BE44FFAAB59ACACB9C8681

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

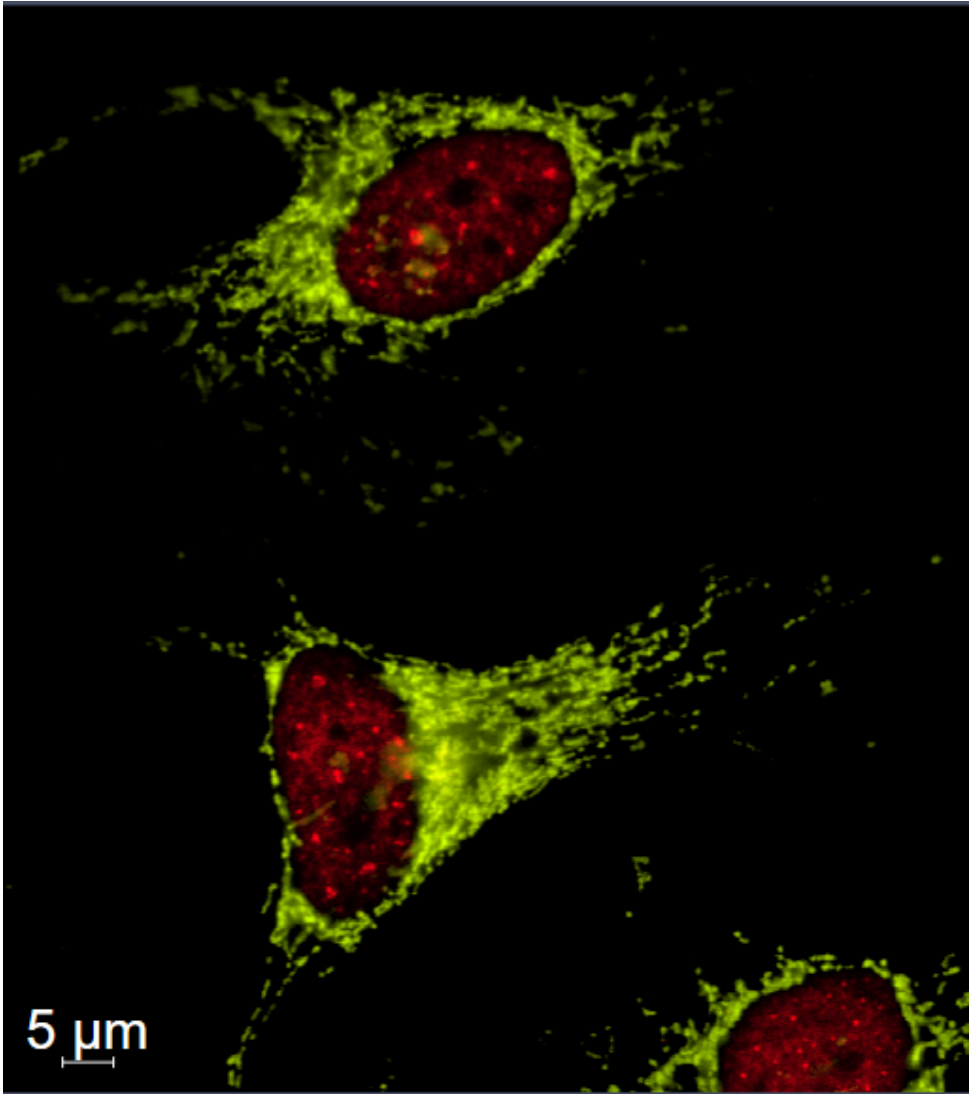

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555; 647
METADATA VALUES	29

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.110µm X 0.110µm
Image size (Pixels)	878 X 976
Image Size (Scaled)	96.16µm X 106.90µm
Bit Depth	14 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	50 Cy5 AF532-49014
Beam Splitter	660 550
Filter Excitation Wavelength	625 - 655 515 - 545
Filter Emission Wavelength	665 - 715 555 - 595
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @20%
Exposure Time	150ms 80ms
Excitation Wavelength	653 528
Emission Wavelength	668 553
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	1.03µm 0.86µm
Binning Mode	2,2

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000026

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam 705, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.110 μ m X 0.110 μ m

Image size (Pixels) = 878 X 976

Image Size (Scaled) = 96.16 μ m X 106.90 μ m

Bit Depth = 14 Bit

Image Center Position = X: 0.00 μ m, Y: 0.00 μ m

ROI Center Offset = X: 0.00 μ m, Y: 0.00 μ m

647 -

Reflector = 50 Cy5
 Beam Splitter = 660
 Filter Excitation Wavelength = 625 - 655
 Filter Emission Wavelength = 665 - 715
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @20%
 Exposure Time = 150ms
 Excitation Wavelength = 653
 Emission Wavelength = 668
 Effective NA = 1.4
 Imaging Device = Axiocam MR R3
 Camera Adapter = 1X
 Depth of Focus = 1.03µm
 Binning Mode = 2,2

532 -

Reflector = AF532-49014
 Beam Splitter = 550
 Filter Excitation Wavelength = 515 - 545
 Filter Emission Wavelength = 555 - 595
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @20%
 Exposure Time = 80ms
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Imaging Device = Axiocam MR R3
 Camera Adapter = 1X
 Depth of Focus = 0.86µm
 Binning Mode = 2,2

Post Processing

None

Image Correction

None

Image by K. Shivanna

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000058

SOURCE PROVENANCE

Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, AxioCam 705, and X-Cite Xylys Lamp was used with the following parameters for each dye:
Structured source metadata values	29
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	5395091AF40BDB4E016EA14F0B59E33C28C83EB410EEFB1E680A3D4C48D2F375
Method PDF SHA-256	B21371BA61D68B4C9D22EC125C8C9B15049D56CA16F8D5B8D78C365014A3BF1A

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

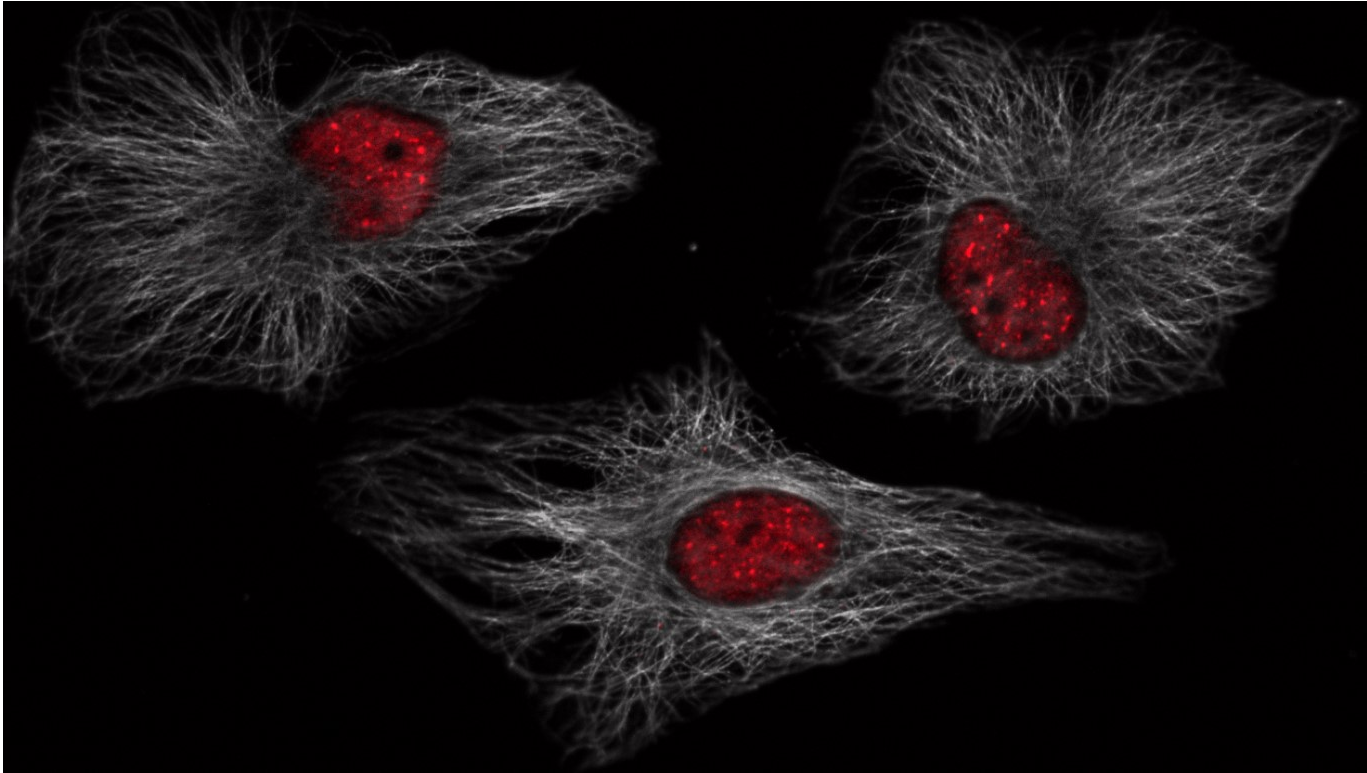

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 555; 647
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Z-Stack	113 Slices (11.2µm)
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm X 0.100µm
Image Size (Pixels)	1283 X 864
Image Size (Scaled)	183.29µm X 123.43µm
Bit Depth	14 Bit
Image Center Position	X: 103.37mm, Y: 46.98mm
Stage Position	X: 103.37mm, Y: 46.98mm
ROI Center Offset	X: 1.14, Y: 0.00µm
Reflector	50 Cy5 AF532-49014
Beam Splitter	660 550
Filter Excitation Wavelength	625 - 655 515 - 545
Filter Emission Wavelength	665 - 715 555 - 595
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @6.0%
Exposure Time	40ms 35ms
Excitation Wavelength	653 553
Emission Wavelength	668 568
Effective NA	1.25
Imaging Device	Axiocam 807
Camera Adapter	1X
Depth of Focus	1.30µm 1.10µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000037

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Alpha Tubulin Monoclonal Antibody (DM1A), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 555 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Z Stack - (3D) - Z-plane 42 of 113 is shown in a 2-Dimensional Image

Z-Stack = 113 Slices (11.2µm)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm X 0.100µm

Image Size (Pixels) = 1283 X 864

Image Size (Scaled) = 183.29µm X 123.43µm

Bit Depth = 14 Bit

Image Center Position = X: 103.37mm, Y: 46.98mm

Stage Position = X: 103.37mm, Y: 46.98mm

ROI Center Offset = X: 1.14, Y: 0.00µm

647 -

Reflector = 50 Cy5

Beam Splitter = 660

Filter Excitation Wavelength = 625 - 655

Filter Emission Wavelength = 665 - 715

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @6.0%

Exposure Time = 40ms

Excitation Wavelength = 653

Emission Wavelength = 668

Effective NA = 1.25

Imaging Device = AxioCam 807

Camera Adapter = 1X

Depth of Focus = 1.30µm

Binning Mode = 2,2

555 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @6.0%

Exposure Time = 35ms

Excitation Wavelength = 553

Emission Wavelength = 568

Effective NA = 1.25

Imaging Device = AxioCam 807

Camera Adapter = 1X

Depth of Focus = 1.10µm

Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Correction, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000058
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	12A585B9C29E6BA015E114BD8C0095551F99172331E00488F93FE77632CB63C1
Method PDF SHA-256	A4733A9DB3A3252BBE660CDCBD4224CC8BC5237FDA2B33FDF9368FEB727C46FC

ORIGINAL IMAGE EVIDENCE / CONFOCAL

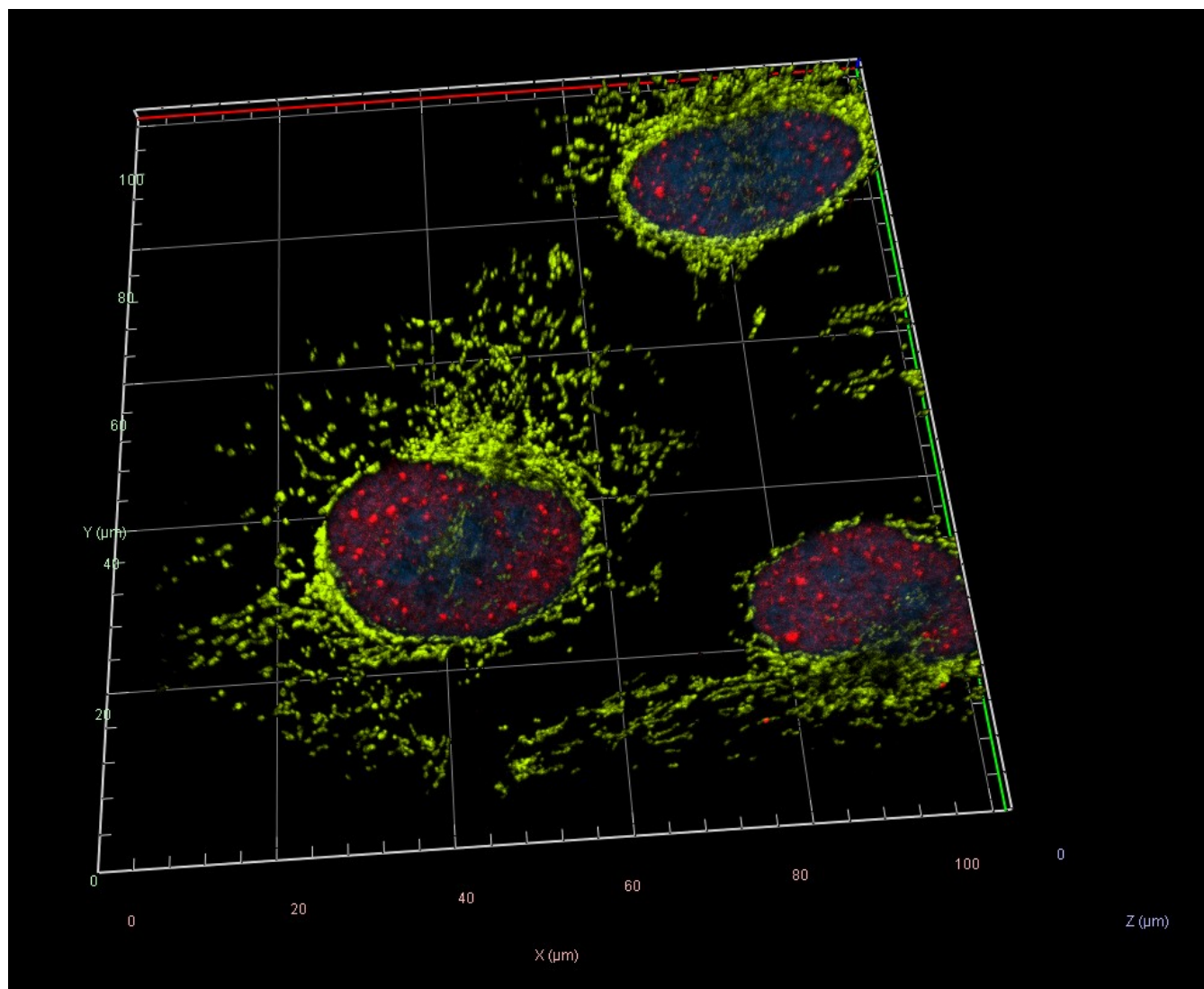

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 647
METADATA VALUES	52

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	11 Slices (700nm)
Channels	3
Scaling (Per Pixel)	0.071µm X 0.071µm X 0.050µm
Image Size (Pixels)	1437 X 1437
Image Size (Scaled)	101.41µm X 101.41µm
Bit Depth	16 Bit
Image Center Position	X: 80.87mm, Y: 55.45mm
Stage Position	X: 80.87mm, Y: 55.45mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 57µm 1.00AU / 46µm 1.00 AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	640nm @0.2% 488nm @1.0% 405nm @0.2%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.0
Rotation	0°
Sampling	1.00
Pixel Time	0.54µs
Frame Time	5.23s
LSM Scan Speed	8
Scan Direction	Bidirectional
Line Step	1
Averaging	4
Excitation Wavelength	653 528 353
Emission Wavelength	668 553 465
Effective NA	1.4
Detection Wavelength	656 - 700 517 - 580 410-470
Imaging Device	GaAsp-PMT 2 Airyscan GaAsp-PMT1
Detector Type	GaAsp-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	650 V 700 V
Detector Offset	0
Detector Digital Gain	1.0
3D Maximum Projection	Precise
Appearance	Threshold 2% (647), Threshold 2% (532), Threshold 2% (DAPI), Ramp 0% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness112%, Enabled Tone Mapping is NOT selected
Projection	View Angle 52°, Scale Z 244%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000026

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 11 Slices (700nm)

Channels = 3

Scaling (Per Pixel) = $0.071\mu\text{m} \times 0.071\mu\text{m} \times 0.050\mu\text{m}$

Image Size (Pixels) = 1437 X 1437

Image Size (Scaled) = $101.41\mu\text{m} \times 101.41\mu\text{m}$

Bit Depth = 16 Bit

Image Center Position = X: 80.87mm, Y: 55.45mm

Stage Position = X: 80.87mm, Y: 55.45mm

ROI Center Offset = X: $0.00\mu\text{m}$, Y: $0.00\mu\text{m}$

647 -

Reflector = None

Contrast Method = Fluorescence

Pinhole = 1.00AU / $57\mu\text{m}$

LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)

Laser Wavelength = 640nm @0.2%

Laser Blanking = Enabled

Scan Mode = Frame

Scan Zoom = 1.0

Rotation = 0°

Sampling = 1.00

Pixel Time = $0.54\mu\text{s}$

Frame Time = 5.23s

LSM Scan Speed = 8

Scan Direction = Bidirectional

Line Step = 1

Averaging = 4

Excitation Wavelength = 653

Emission Wavelength = 668

Effective NA = 1.4

Detection Wavelength = 656 - 700

Imaging Device = GaAsp-PMT 2

Detector Type = GaAsp-PMT

Detector Gain = 650 V

Detector Offset = 0

Detector Digital Gain = 1.0

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 46µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @1.0%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.54µs
 Frame Time = 5.23s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 517 - 580
 Imaging Device = Airyscan
 Detector Type = GaAsp-PMT
 Detector Gain = 700 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector - None
 Contrast Method - Fluorescence
 Pinhole = 1.00 AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.2%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.0
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.54µs
 Frame Time = 5.23s
 LSM Scan Speed = 8
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 4
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410-470
 Imaging Device = GaAsp-PMT1
 Detector Type = GaAsp-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 2% (647), Threshold 2% (532), Threshold 2% (DAPI), Ramp 0% all channels, Maximum 100% all channels.
 Background = Black
 Light = Brightness112%, Enabled Tone Mapping is NOT selected
 Projection = View Angle 52°, Scale Z 244%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000058

SOURCE PROVENANCE

Stage	Confocal
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	52
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	C7E163D9D8FCC8A36D38D71A956691BB30E4F404BE832B64779FB21C40C35FF8
Method PDF SHA-256	A622CE1DCFAB5E898CB03968A45F79CE6466F678EF43AE8E677156131F7EF498

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

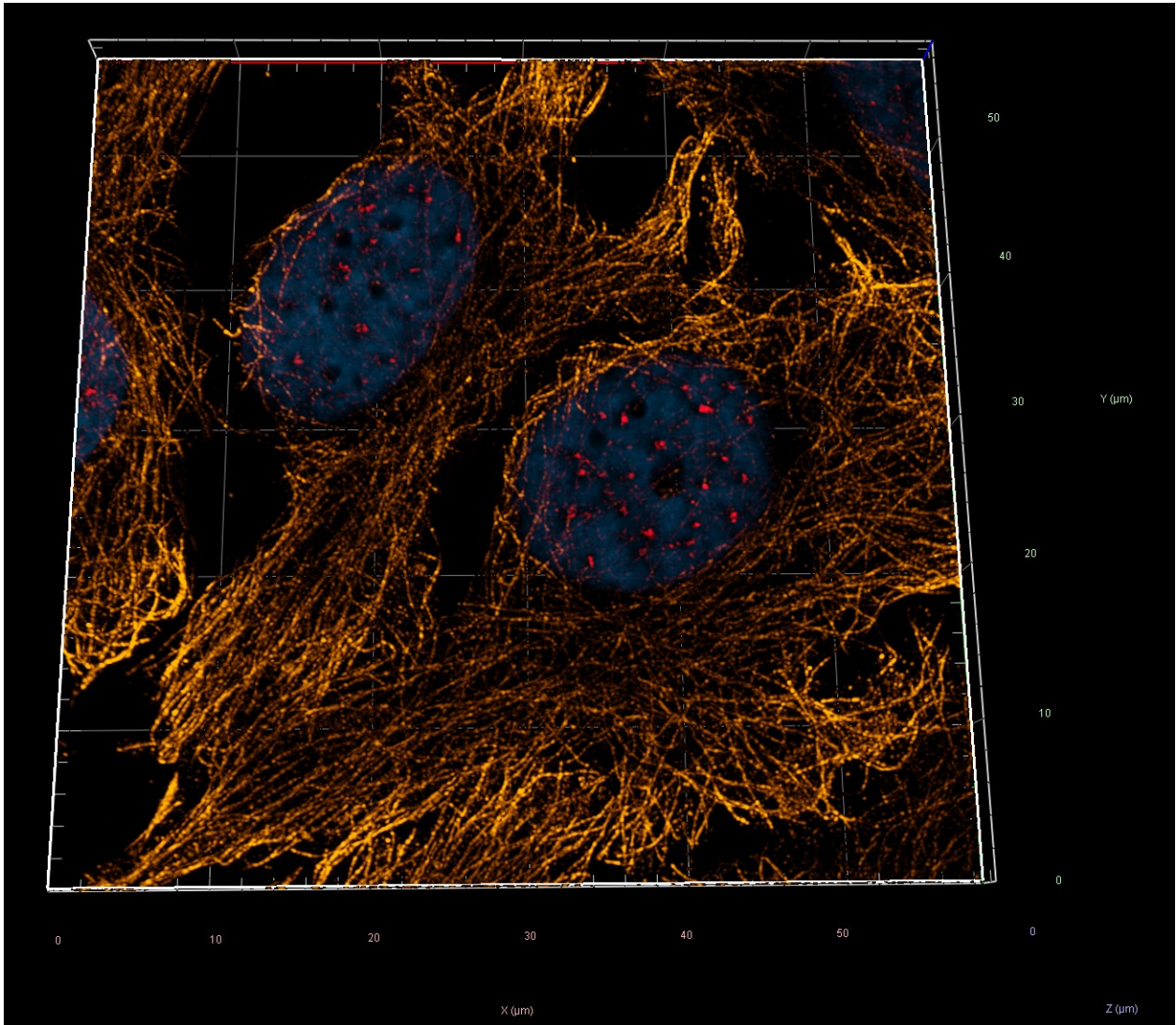

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 555; 647
METADATA VALUES	56

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	61 Slices (3.0µm)
Channels	3
Scaling (Per Pixel)	0.03530 µm X 0.03530 µm X 0.05000 µm
Image Size (Pixels)	1657 X 1622
Image Size (Scaled)	58.49µm X 57.25µm
Bit Depth	16 Bit
Image Center Position	X: 88.26mm, Y: 52.67mm
Stage Position	X: 88.26mm, Y: 52.67mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00 AU / 282µm 5.00 AU / 250µm 5.00 AU / 182µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	640nm @0.9% 561nm @0.2% 405nm @0.6%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.3
Rotation	-6.11°
Sampling	2.00
Pixel Time	0.96µs
Frame Time	10.94s
LSM Scan Speed	6
Scan Direction	Bidirectional
Line Step	1
Averaging	2
Excitation Wavelength	653 553 353
Emission Wavelength	668 568 465
Effective NA	1.4
Detection Wavelength	650-700 565 - 630 400-464
Imaging Device	AiryScan
Detector Type	GaAsP-PMT

FIELD	SOURCE-DOCUMENTED VALUE
Detector Gain	850 V
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 9.0 (3D, Auto) Super Resolution 8.5 (3D, Auto)
3D Maximum Projection	Precise
Appearance	Threshold 1% (647), Threshold 1% (555), Threshold 1% (DAPI), Ramp 0% all channels, Maximum 100% all channels.
Background	Black
Light	Brightness 00%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)
Clipping Image Process	The Clipping plane is introduced to demonstrate the specificity and localization of the Identifyn® SRM Alexa Fluor™ 647 and the Identifyn® Rabbit anti-Human gamma H2AX Recombinant Monoclonal, Human, Fc, to the DNA damage-induced gamma-H2AX foci within the nuclear compartment.
X - Y	Active was selected, and an invisible filter was chosen to remove the obstruction from the image-cut plane and highlight the nuclear compartment and gamma-H2AX foci. A white Outline of the clipped perimeter was shown, the Clip Back was Selected, and the Clip Front was NOT selected.
Position of Clip	-52%
Clip Y	-0°
Clip Z	0°

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 555 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000037

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Alpha Tubulin Monoclonal Antibody (DM1A), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 555 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 61 Slices (3.0µm)

Channels = 3

Scaling (Per Pixel) = 35.30nm X 35.30nm X 50.00nm

Image Size (Pixels) = 1657 X 1622

Image Size (Scaled) = 58.49µm X 57.25µm

Bit Depth = 16 Bit

Image Center Position = X: 88.26mm, Y: 52.67mm

Stage Position = X: 88.26mm, Y: 52.67mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

647 -

Reflector = None

Contrast Method = Fluorescence

Pinhole = 5.00 AU / 282µm

LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)

Laser Wavelength = 640nm @0.9%

Laser Blanking = Enabled

Scan Mode = Frame

Scan Zoom = 1.3

Rotation = -6.11°

Sampling = 2.00

Pixel Time = 0.96µs

Frame Time = 10.94s

LSM Scan Speed = 6

Scan Direction = Bidirectional

Line Step = 1

Averaging = 2

Excitation Wavelength = 653

Emission Wavelength = 668

Effective NA = 1.4

Detection Wavelength = 650-700

Imaging Device = Airyscan

Detector Type = GaAsP-PMT

Detector Gain = 850 V

Detector Offset = 0

Detector Digital Gain = 1.0

AiryScan Mode = Super Resolution 9.0 (3D, Auto)

555 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 250µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 561nm @0.2%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = -6.11°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 10.94s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 553
 Emission Wavelength = 568
 Effective NA = 1.4
 Detection Wavelength = 565 - 630
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.5 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00 AU / 182µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.6%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.3
 Rotation = -6.11°
 Sampling = 2.00
 Pixel Time = 0.96µs
 Frame Time = 10.94s
 LSM Scan Speed = 6
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 2
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400-464
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.5 (3D, Auto)

Post Processing - AiryScan 3D Processing

Display Mode "SR"
 Super Resolution Auto Filter Selected
 Process 3D Current Image
 Created

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 1% (647), Threshold 1% (555), Threshold 1% (DAPI), Ramp 0% all channels, Maximum 100% all channels.
 Background = Black
 Light = Brightness 00%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Clipping Image Process = The Clipping plane is introduced to demonstrate the specificity and localization of the Identifyn® SRM Alexa Fluor™ 647 and the Identifyn® Rabbit anti-Human gamma H2AX Recombinant Monoclonal, Human, Fc, to the DNA damage-induced gamma-H2AX foci within the nuclear compartment.

X - Y = Active was selected, and an invisible filter was chosen to remove the obstruction from the image-cut plane and highlight the nuclear compartment and gamma-H2AX foci. A white Outline of the clipped perimeter was shown, the Clip Back was Selected, and the Clip Front was NOT selected.

Position of Clip = -52%

Clip Y = -0°

Clip Z =0°

Image Corrections -


None

Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000058
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 AiryScan, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	56
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	23C4FE086A75BA08E8734999EEA109F957E85DB061F49512A38D1810E6C41585
Method PDF SHA-256	E02D1F6E17B22CA83799D32C11C615D5CE7BFD2F2F4898244D6D5CE7DE060A6D

ORIGINAL IMAGE EVIDENCE / SIM

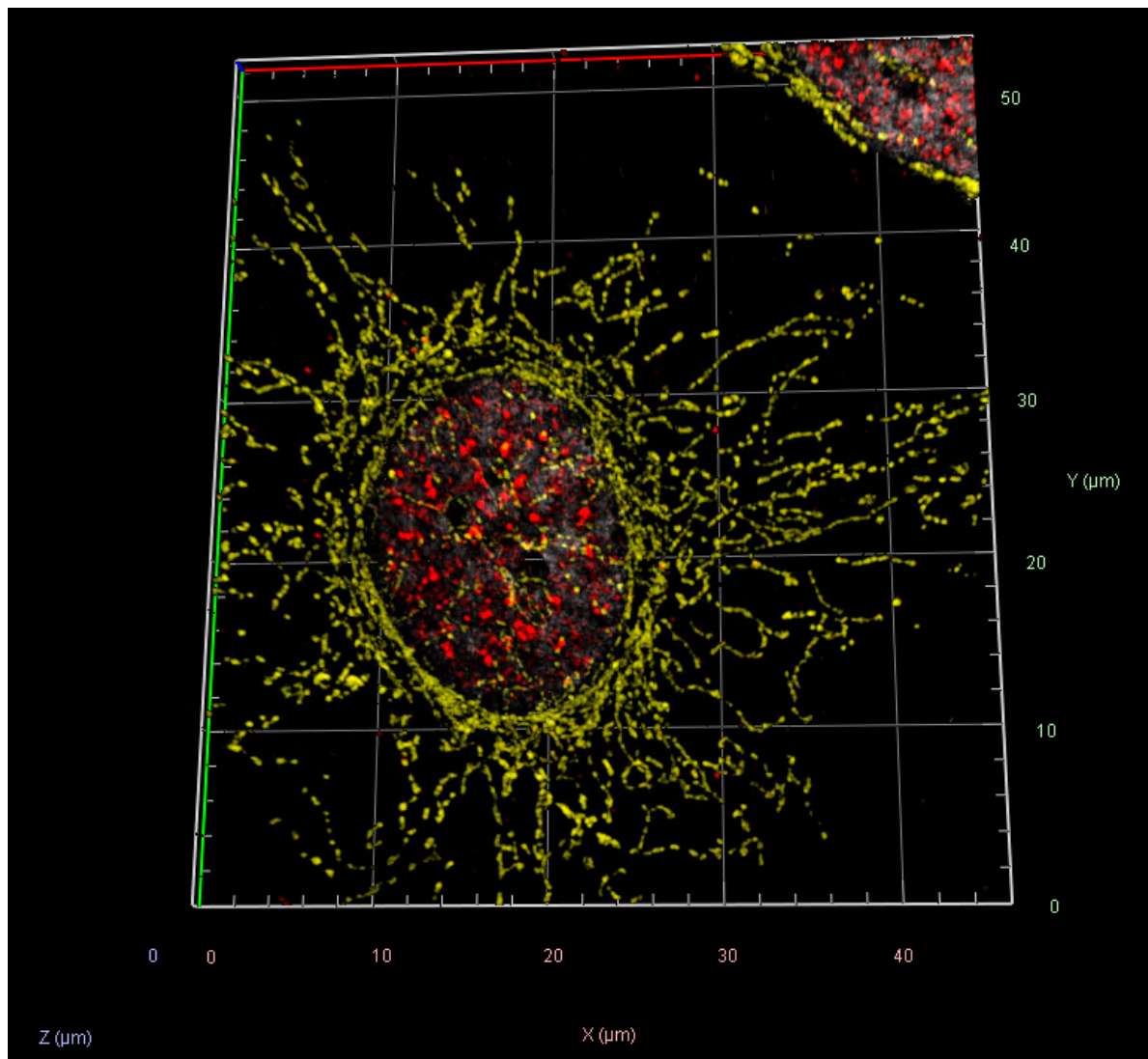

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 647
METADATA VALUES	64

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	29 Slices (1.40µm)
Channels	3 1
Scaling (Per Pixel)	0.01320 µm X 0.01320 µm X 0.05000 µm
Image Size (Pixels)	3486 X 3954
Image Size (Scaled)	46.00µm X 52.17µm
Bit Depth	16 Bit
Image Center Position	X: 64.97mm, Y: 35.85mm
Stage Position	X: 64.97mm, Y: 35.85mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	640 488 405
Channel Name	T1-Axiocam 820 -SR T2-Axiocam 820#2 -SR T3-Axiocam 820 #2-SR
Excitation Wavelength	640 488 405
Emission Wavelength	729 525 445
Effective NA	1.4
Detection Wavelength	419 - 470, 503 - 546, 576 - 617, 657 -800
Imaging Device	Axiocam 820 Axiocam 820#2
Camera Adapter	1X
Exposure Time	100ms 120ms 50ms
Depth of Focus	1.13µm 0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	640nm @2.0% 488nm @8.0% 405nm @8.0%
Frame Time	1.33s 1.59s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	2
Algorithm	SIM

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	9.3400
Input SNR	N/A
Regularization Weight	N/A
Iterations	N/A
Filter	N/A
Fitting	N/A
Normalization	Automatic
Experimental PSF	N/A
3D Maximum Projection	Precise
Appearance	Threshold 1% (647), Threshold 1% (532), Threshold 1% (DAPI), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 100%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000026

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 29 Slices (1.40µm)

Channels = 3

Scaling (Per Pixel) = 13.20nm X 13.20nm X 50.00nm

Image Size (Pixels) = 3486 X 3954

Image Size (Scaled) = 46.00µm X 52.17µm

Bit Depth = 16 Bit

Image Center Position = X: 64.97mm, Y: 35.85mm

Stage Position = X: 64.97mm, Y: 35.85mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

647 -

Reflectors = SR HE LBF 405/488/561/642

Beam Splitter = 405, 488, 561, 642

Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642

Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800

Contrast Method = Fluorescence

Illumination Wavelength = 640

Channel Name = T1-Axiocam 820 -SR

Excitation Wavelength = 640

Emission Wavelength = 729

Effective NA = 1.4

Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800

Imaging Device = Axiocam 820

Camera Adapter = 1X

Exposure Time = 100ms

Depth of Focus = 1.13µm

Binning Mode = 2,2

Laser Wavelength = 640nm @2.0%

Frame Time = 1.33s

Scan Direction = Unidirectional

Grating 36.5µm grating

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820#2 -SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 -800
 Imaging Device = Axiocam 820#2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @8.0%
 Frame Time = 1.59s
 Scan Direction = Unidirectional
 Grating 36.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T3-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 -800
 Imaging Device = Axiocam 820#2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @8.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM Processing

Processing Dimension - 3D
 Result Sampling = 4
 Process Sampling = 2
 Channels = 1
 Algorithm = SIM
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 9.3400
 Input SNR = N/A
 Regularization Weight = N/A
 Iterations = N/A
 Filter = N/A
 Fitting = N/A
 Normalization = Automatic
 Experimental PSF = N/A

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 1% (647), Threshold 1% (532), Threshold 1% (DAPI), Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 100%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

DAPI was pseudo-colored in Zen Software from blue, the default color, to dark grey to provide visual contrast between the 405/647 channels.

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000058
Stage	SIM
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	64
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa

SOURCE PROVENANCE

Image SHA-256	88ADC5D55CD2CA04BE1B234EFBDB770CC4A7EB944D002F358CD72BEEC68B813A
Method PDF SHA-256	8F218D16C32E1B8F14738B3DA452A911415F26027505555239589319414D0E34

ORIGINAL IMAGE EVIDENCE / SIM²

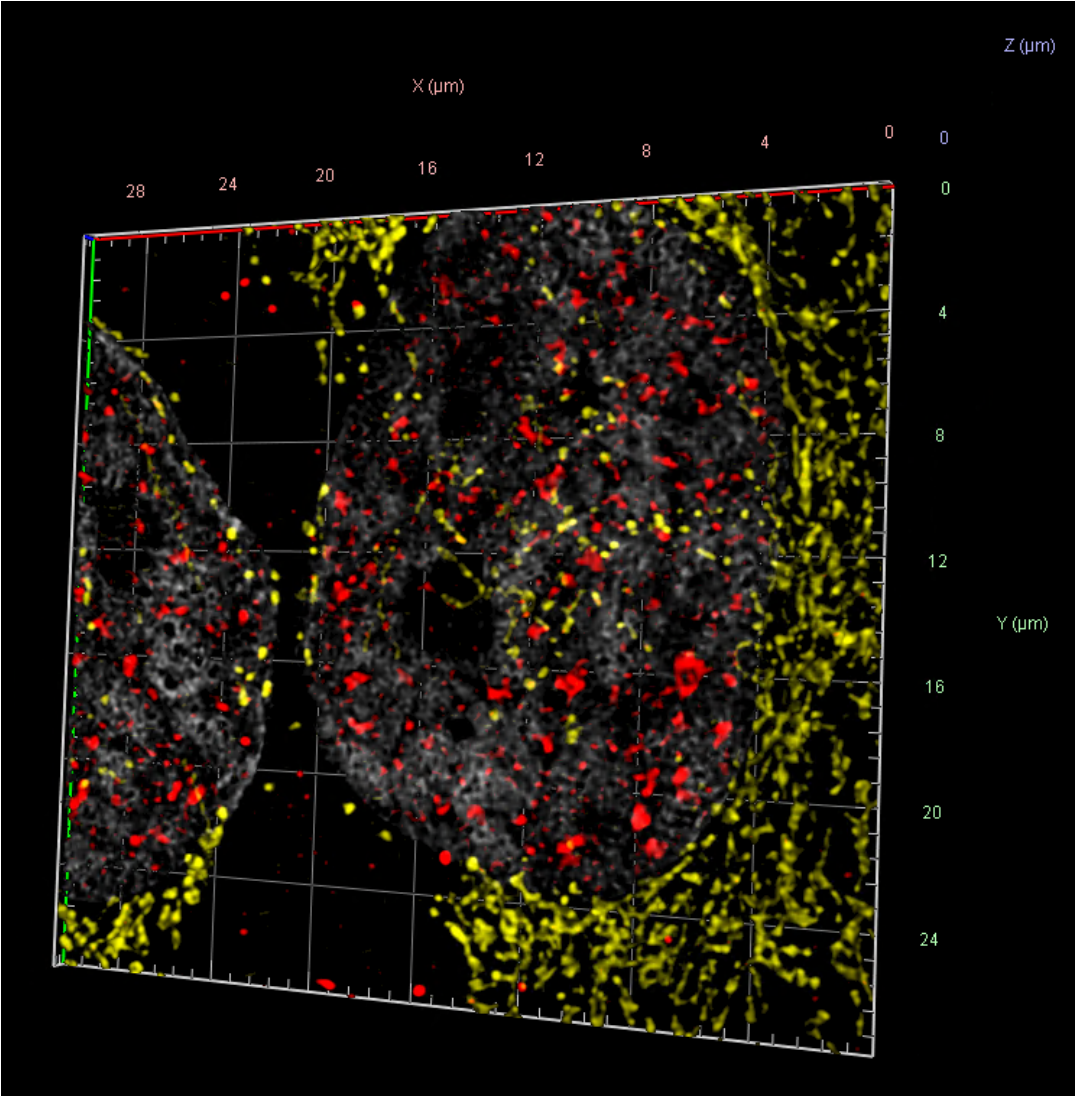

EVIDENCE TIER	SIM ²
INSTRUMENT	ZEISS Elyra
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 647
METADATA VALUES	64

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	40 Slices (1.95µm)
Channels	3 1
Scaling (Per Pixel)	0.01320 µm X 0.01320 µm X 0.05000 µm
Image Size (Pixels)	2306 X 2107
Image Size (Scaled)	30.43µm X 27.80µm
Bit Depth	16 Bit
Image Center Position	X: 67.73mm, Y: 37.84mm
Stage Position	X: 67.73mm, Y: 37.84mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflectors	SR HE LBF 405/488/561/642
Beam Splitter	405, 488, 561, 642
Filter Excitation Wavelength	200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength	419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method	Fluorescence
Illumination Wavelength	640 488 405
Channel Name	T1-Axiocam 820 -SR T2-Axiocam 820#2 -SR T3-Axiocam 820 #2-SR
Excitation Wavelength	640 488 405
Emission Wavelength	729 525 445
Effective NA	1.4
Detection Wavelength	419 - 470, 503 - 546, 576 - 617, 657 -800
Imaging Device	Axiocam 820 Axiocam 820#2
Camera Adapter	1X
Exposure Time	100ms 120ms 50ms
Depth of Focus	1.13µm 0.81µm 0.69µm
Binning Mode	2,2
Laser Wavelength	640nm @2.0% 488nm @8.0% 405nm @8.0%
Frame Time	1.33s 1.59s 676.00ms
Scan Direction	Unidirectional
Result Sampling	4
Process Sampling	4
Algorithm	SIM2

FIELD	SOURCE-DOCUMENTED VALUE
SIM Settings	Standard
Auto Sharpness	On
Sharpness	8.8593
Input SNR	N/A
Regularization Weight	N/A
Iterations	N/A
Filter	N/A
Fitting	N/A
Normalization	Automatic
Experimental PSF	N/A
3D Maximum Projection	Precise
Appearance	Threshold 5% (647), Threshold 3% (532), Threshold 4% (DAPI), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 98%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 32%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn Primary Antibody

Identifyn® Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG
Cat ID# - PA-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)
Cat ID# - SA 000058

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Donkey Anti-Mouse, IgG (H + L)
Cat ID# - SA 000026

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 24 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated for 1 hour at room temperature at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 40 Slices (1.95µm)

Channels = 3
Scaling (Per Pixel) = 13.20nm X 13.20nm X 50.00nm
Image Size (Pixels) = 2306 X 2107
Image Size (Scaled) = 30.43µm X 27.80µm
Bit Depth = 16 Bit
Image Center Position = X: 67.73mm, Y: 37.84mm
Stage Position = X: 67.73mm, Y: 37.84mm
ROI Center Offset = X: 0.00µm, Y: 0.00µm

647 -

Reflectors = SR HE LBF 405/488/561/642
Beam Splitter = 405, 488, 561, 642
Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
Contrast Method = Fluorescence
Illumination Wavelength = 640
Channel Name = T1-Axiocam 820 -SR
Excitation Wavelength = 640
Emission Wavelength = 729
Effective NA = 1.4
Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
Imaging Device = Axiocam 820
Camera Adapter = 1X
Exposure Time = 100ms
Depth of Focus = 1.13µm
Binning Mode = 2,2
Laser Wavelength = 640nm @2.0%
Frame Time = 1.33s
Scan Direction = Unidirectional
Grating 36.5µm grating

532 -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 488
 Channel Name = T2-Axiocam 820#2 -SR
 Excitation Wavelength = 488
 Emission Wavelength = 525
 Effective NA = 1.4
 Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 -800
 Imaging Device = Axiocam 820#2
 Camera Adapter = 1X
 Exposure Time = 120ms
 Depth of Focus = 0.81µm
 Binning Mode = 2,2
 Laser Wavelength = 488nm @8.0%
 Frame Time = 1.59s
 Scan Direction = Unidirectional
 Grating 36.5µm grating

DAPI -

Reflectors = SR HE LBF 405/488/561/642
 Beam Splitter = 405, 488, 561, 642
 Filter Excitation Wavelength = 200 - 405, 200 - 488, 200 - 561, 200 - 642
 Filter Emission Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 - 800
 Contrast Method = Fluorescence
 Illumination Wavelength = 405
 Channel Name = T3-Axiocam 820 #2-SR
 Excitation Wavelength = 405
 Emission Wavelength = 445
 Effective NA = 1.4
 Detection Wavelength = 419 - 470, 503 - 546, 576 - 617, 657 -800
 Imaging Device = Axiocam 820#2
 Camera Adapter = 1X
 Exposure Time = 50ms
 Depth of Focus = 0.69µm
 Binning Mode = 2,2
 Laser Wavelength = 405nm @8.0%
 Frame Time = 676.00ms
 Scan Direction = Unidirectional
 Grating 36.5µm grating

Post Processing - SIM2 Processing

Processing Dimension - 3D
 Result Sampling = 4
 Process Sampling = 4
 Channels = 1
 Algorithm = SIM2
 SIM Settings = Standard
 Auto Sharpness = On
 Sharpness = 8.8593
 Input SNR = N/A
 Regularization Weight = N/A
 Iterations = N/A
 Filter = N/A
 Fitting = N/A
 Normalization = Automatic
 Experimental PSF = N/A

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 5% (647), Threshold 3% (532), Threshold 4% (DAPI), Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 98%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 32%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

DAPI was pseudo-colored in Zen Software from blue, the default color, to dark grey to provide visual contrast between the 405/647 channels.

Image by P. Raut

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000058
Stage	SIM ²
Instrument	ZEISS Elyra
Microscope configuration	Zeiss Elyra 5 Lattice SIM, AxioObserver.Z1/ 7 using Plan-Apochromat 63X/1.40 oil DIC M27, 1X Tubelens, was used with the following parameters for each dye:
Structured source metadata values	64
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa

SOURCE PROVENANCE

Image SHA-256	4C37535B07AE8A5352D42A46BACD27DB4C49EEB6684C75B47C3806C288341891
Method PDF SHA-256	6DED0DEA8007ABB87363E6AEC3E4A3D52EFA1FDF6FE31C8C3D2541817BED9921

ORIGINAL IMAGE EVIDENCE / SINGLE-MOLECULE STORM

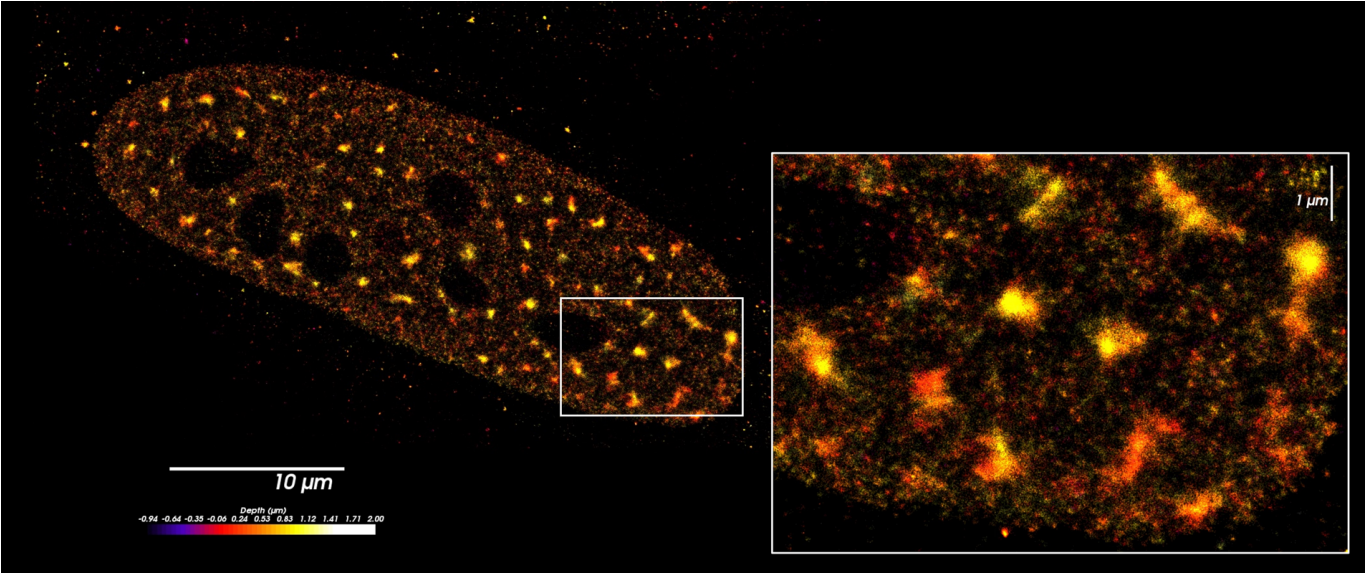

EVIDENCE TIER	Single-molecule STORM
INSTRUMENT	Bruker Vutara VXL
CELL MODEL	HeLa
DETECTED DYES	647
METADATA VALUES	71

STRUCTURED ACQUISITION METADATA / SINGLE-MOLECULE STORM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Objective	Olympus 60X/1.30 silicon oil objective
Number of Probes	1
Number of Simultaneous Probes	1
Number of Reference Probes	1
Camera	Both Cameras
Color Channel	Both Channels
Probe Capture Mode	Interleaved
Z -Stack Capture Mode	Interleaved
Multiple Location Capture	Not selected
Post Record Reference	Not Selected
Fluidics	Not Selected
Autofocus Drift Correction	On
Biplane Module	635nm Dichroic
Objective	Silicon Selected
Stage Insert	Rectangular Selected
Course Focus Position	8-Well Chamber Selected
Dichroic	SuperRes
PSF	(Red, Orange)
Heater	Current: 26
Recording Vertical Field of View	100%
Widefield Camera Exposure Time	100ms 20ms
Super Resolution Camera Exposure Time	20ms
Widefield Camera Binning	1
Per-Probe Exposure Time	Selected
Predefined Probe	First probe: 640 nm (635 nm Dichroic); First reference probe: 640 nm (635 nm Dichroic)
Photoactivation/Imaging Laser	0% Selected
Exposure Time	20ms Selected
Record Workflow	Selected Options
View Mode	Widefield 200µm
Capture Mode	Live
Piezo Current	182.5 µm

FIELD	SOURCE-DOCUMENTED VALUE
Widefield 640nm	0.1X Neutral Density Filter: Selected, @0.2% Laser Input
Reference Probe 1 Laser control (640 nm (640 nm Dichroic))	Not Selected
Widefield camera Exposure Time	100ms
Widefield Laser	"Widefield 640nm" was deselected, and "calibrate" under Autofocus was selected.
Current Intensity	4.5
Calibration Value	-5070 nm/unit
Autofocus Mode	Pre-Capture
SuperRes	50µm is Selected
Piezo Record	Begin: 180.0; End: 181.2
SuperRes 640nm	0.1X Neutral Density Filter: Selected, @0.02% Laser Input
Record Workflow - Record Tab	Capture Parameters
Frames	Probe 1: 500
Z Positions	6 using Steps: 0.2µm (200nm)
Cycles	10
Time Points	1
Pause Between Cycles	Not Selected
Automatic Photoactivation Laser Control	Not Selected
Laser Power	30%
Cutout Width	16px
Max Frame Accumulation	1 frame
Max Accumulation Offset	1px
Fiducial at Max Accumulation	Not selected
Max Particles Per Frame	No Max
Background SLPE Removal	Not Selected
Cutout/Background Removal	Not Selected
Preview Cutout Positions	Selected
Probe1	BG threshold: 20
Sizing Mode	Constant
Particle size	20nm
Color by	Depth
Color Table	Single
Opacity Mode	Constant; Opacity 1: 0.3
Front-to-back Attenuation	Not Selected
Method	Particle (direct)
Time Intervals	10 Intervals

FIELD	SOURCE-DOCUMENTED VALUE
Pixel Size	50nm
Smoothing	8.00
PSF Z offset	-450 to +450
Radial precision	25
Axial precision	60

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® Primary Antibody

Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG

Cat ID# - Pa-RR-000004

Identifyn® Secondary Antibodies

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human, IgG (H + L)

Cat ID# - SA 000058

Gamma H2AX (γ H2AX) is a variant of the histone H2A protein that plays a crucial role in the cellular response to DNA double-strand breaks (DSBs). Its primary function is as a marker of DNA damage, particularly in the context of DSBs. When DSBs occur in chromatin, γ H2AX is rapidly phosphorylated explicitly at the serine 139 residue (referred to as Ser139 phosphorylation or γ -phosphorylation) near the site of the break. This phosphorylation event forms γ H2AX foci, visualized via immunofluorescence microscopy techniques as an aggregate of fluorescently labeled protein at the site of breaks.

Method

Hela cells were grown on μ -Slide 8 Well high Glass Bottom Ibidi chamber. DNA damage was induced via treatment with 2 μ M Etoposide-spiked complete growth media for 6 hours, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Rabbit anti-Human γ H2Ax Recombinant Monoclonal, Human Fc, IgG, was incubated for 1 hour at room temperature at a dilution of 1:100. Cells were washed 3X in Identifyn's™ proprietary Wash Buffer and Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:50. Cells were washed 5X in Identifyn's™ proprietary Wash Buffer followed by the addition of Identifyn's™ proprietary STORM Buffer to the sample chamber enabling image acquisition.

Microscope

Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:

Calibration -

Calibration follows Identifyn's Calibration Standard Operating Procedure, which can be found here.

<https://identifyn.com/storm-calibration>

Image Acquisition -

SRX Software - Under the 'Single Molecule Localization' tab, 'Record' is selected and contains the following workflow tabs and parameters.

Configuration Tab - Capture Configuration

Number of Probes: 1
 Number of Simultaneous Probes: 1
 Number of Reference Probes: 1
 Camera: Both Cameras
 Color Channel: Both Channels
 Probe Capture Mode: Interleaved
 Z -Stack Capture Mode: Interleaved
 Multiple Location Capture: Not selected
 Post Record Reference: Not Selected
 Fluidics: Not Selected
 Autofocus Drift Correction: On

Microscope configuration

Biplane Module: 635nm Dichroic
 Objective: Silicon Selected
 Stage Insert: Rectangular Selected
 Course Focus Position: 8-Well Chamber Selected
 Dichroic: SuperRes
 PSF: (Red, Orange)
 Heater: Current: 26

Camera Configuration

Recording Vertical Field of View: 100%
 Widefield Camera Exposure Time: 100ms
 Super Resolution Camera Exposure Time: 20ms
 Widefield Camera Binning: 1
 Per-Probe Exposure Time: Selected

Probe Configuration

Predefined Probe: First probe: 640 nm (635 nm Dichroic); First reference probe: 640 nm (635 nm Dichroic)
 Photoactivation/Imaging Laser: 0% Selected
 Exposure Time: 20ms Selected

Experiment Configuration

Record Workflow: Selected Options
 View Mode: Widefield 200µm
 Capture Mode: Live

Stage and Piezo Selection

Examination of Sample

Cell(s) or Cellular Structure Localization - The red-bordered box identifies the area of the SuperRes view and defines a suitable cell(s) or cellular structure. The cell(s) or structures are optimally focused.

Piezo Current: 182.5 μ m

Probe 1 Laser control (640 nm (635 nm Dichroic)

Widefield 640nm: 0.1X Neutral Density Filter: Selected, @0.2% Laser Input
Reference Probe 1 Laser control (640 nm (640 nm Dichroic): Not Selected

Advanced Tools Option

Widefield camera Exposure Time: 100ms
Default values are maintained unless otherwise noted.

Widefield Laser: "Widefield 640nm" was deselected, and "calibrate" under Autofocus was selected.

Autofocus - Calibration Values

Current Position" -0.979
Current Intensity: 4.5
Calibration Value: -5070 nm/unit
Autofocus Mode: Pre-Capture

View Mode

SuperRes: 50 μ m is Selected

Capture Mode

Live

Stage and Piezo Selection

Cell(s) or Cellular Structure Localization - Positioning of the target Cell or Subsection of the Target Cell is optimized and optimally focused.

Piezo Record: Begin: 180.0; End: 181.2

Probe 1 Laser control (640 nm (635 nm Dichroic)

SuperRes 640nm: 0.1X Neutral Density Filter: Selected, @0.02% Laser Input

Reference Probe 1 Laser control (640 nm (635 nm Dichroic)

Not Selected

Advanced Tools

Widefield Camera Exposure Time: 20ms
Default values are maintained unless otherwise noted.

Record Workflow - Record Tab: Capture Parameters

Frames: Probe 1: 500

Z Positions: 6 using Steps: 0.2µm (200nm)

Cycles: 10

Time Points: 1

Pause Between Cycles: Not Selected

Automatic Photoactivation Laser Control: Not Selected

Probe 1 Laser control (640 nm (640 nm Dichroic)

Laser Power: 30%

Expert Configuration Options - The Localization Tab

Beyond the values shown, all other parameters were default values.

Localization

Cutout Width: 16px

Max Frame Accumulation: 1 frame

Max Accumulation Offset: 1px

Fiducial at Max Accumulation: Not selected

Max Particles Per Frame: No Max

Background SLPE Removal: Not Selected

Cutout/Background Removal: Not Selected

Preview Cutout Positions: Selected

Record Workflow - Localization Tab

Data Collection

The folder is created automatically with the localization file for this image acquisition.

Recording Summary

All parameters were carried over from previous tabs and are kept except:

Probe1: BG threshold: 20

Region of Interest -

Not selected

Display

Was deployed to inspect if the BG threshold applied chooses localizations in an expected manner

Record Workflow - View Tab

Region of Interest

Default values

Visualization Parameters

Point Splatting

Sizing Mode: Constant

Particle size: 20nm

Color by: Depth

Color Table: Single

Opacity Mode: Constant; Opacity 1: 0.3

Front-to-back Attenuation: Not Selected

Post Acquisition analysis - View Tab

Drift Correction - Computed with the following parameters:

Method: Particle (direct)

Time Intervals: 10 Intervals

Pixel Size: 50nm

Smoothing: 8.00

Advanced Particle Filters - Applied to the drift-corrected data with the following parameters

PSF Z offset: -450 to +450

Radial precision: 25

Axial precision: 60

The other parameters were not applied.

Statistical analysis - NA

Image by S. Thakur

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000058
Stage	Single-molecule STORM
Instrument	Bruker Vutara VXL
Microscope configuration	Bruker Vutara VXL Biplane Single Molecule Microscope, using an Olympus 60X/1.30 silicon oil objective, was used with the following parameters:
Structured source metadata values	71
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa

SOURCE PROVENANCE

Image SHA-256	12C375C5EF3D2FEAD8EFA4D7218B14D78A77BBE167883F7D1AD4505263C614B9
Method PDF SHA-256	9E9EFA07275BFA41014669659DABAFD31B0A607C807BAFED7473D10157489067

ORIGINAL IMAGE EVIDENCE / CONTROL

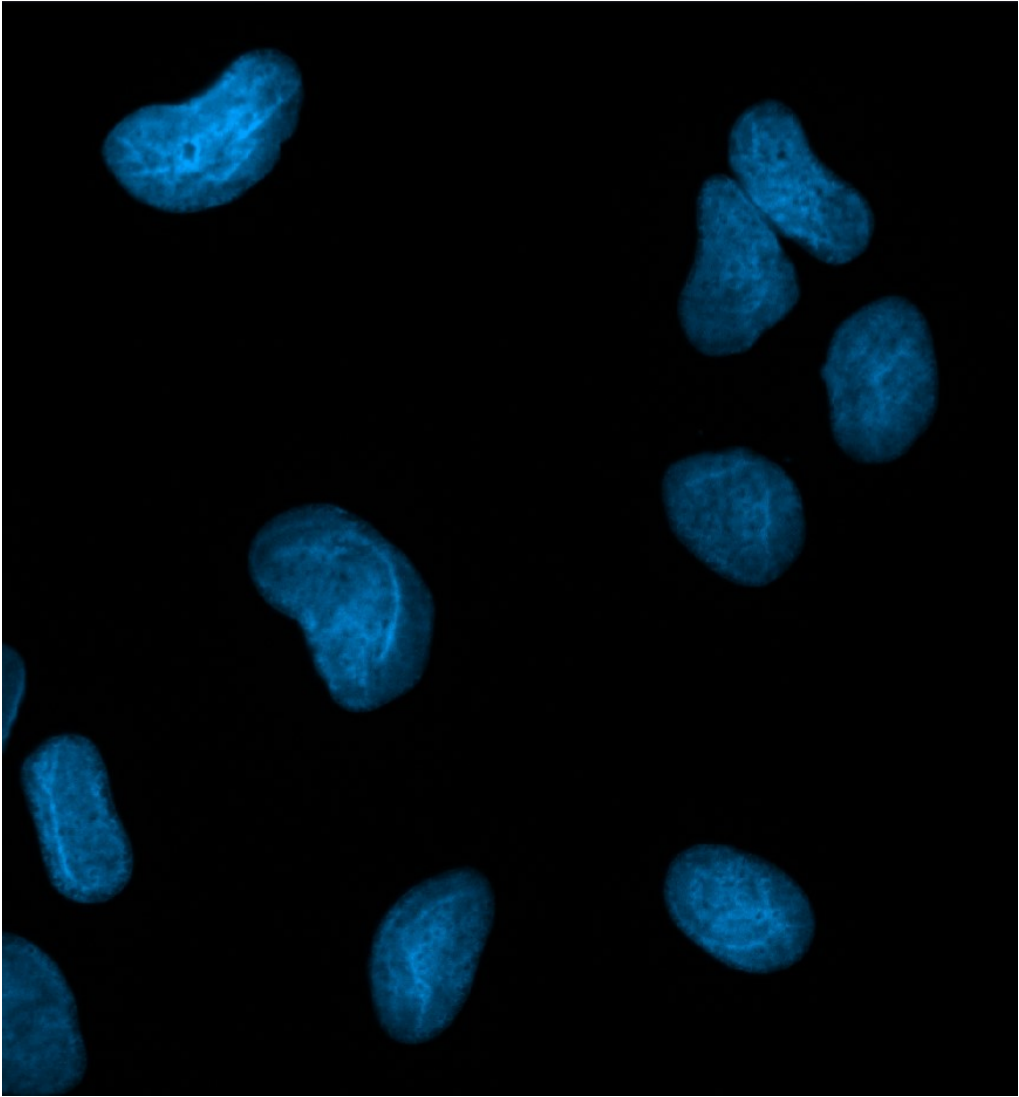

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 647
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	998 X 1069
Image Size (Scaled)	142.57µm X 152.71µm
Bit Depth	14 Bit
Image Center Position	X: 52.09mm, Y: 50.17mm
Stage Position	X: 52.09mm, Y: 50.17mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	50 Cy5 49 DAPI
Beam Splitter	660 395
Filter Excitation Wavelength	625 - 655 335 - 383
Filter Emission Wavelength	665 - 715 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	653 353
Emission Wavelength	668 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.30µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Goat Anti-Human IgG (H + L)

Cat ID# - SA 000058

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 998 X 1069

Image Size (Scaled) = 142.57µm X 152.71µm

Bit Depth = 14 Bit

Image Center Position = X: 52.09mm, Y: 50.17mm

Stage Position = X: 52.09mm, Y: 50.17mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

647 -

Reflector = 50 Cy5

Beam Splitter = 660

Filter Excitation Wavelength = 625 - 655

Filter Emission Wavelength = 665 - 715

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 653

Emission Wavelength = 668

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.30µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000058
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	3B1FC86835CE2925A29AA45AF74D0D196A75A2A7669DBAC0FF24666BC1BF9F2B
Method PDF SHA-256	F2D89ABEE121EF7556B0D7879369E2023600695C413626F89E469F32E6DE43E9