

IDENTIFYN RESOLUTION STEPDOWN

SECONDARY ANTIBODY MICROSCOPY EVIDENCE ASSESSMENT

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit, IgG (H + L)

WESTERN BLOT -> WIDEFIELD -> WIDEFIELD SIM / APOTOME -> CONFOCAL -> AIRYSCAN
SIM -> SIM²

7


ORDERED EVIDENCE TIERS

1

SUPPORTING CONTROL

PENDING

SCIENTIST REVIEW

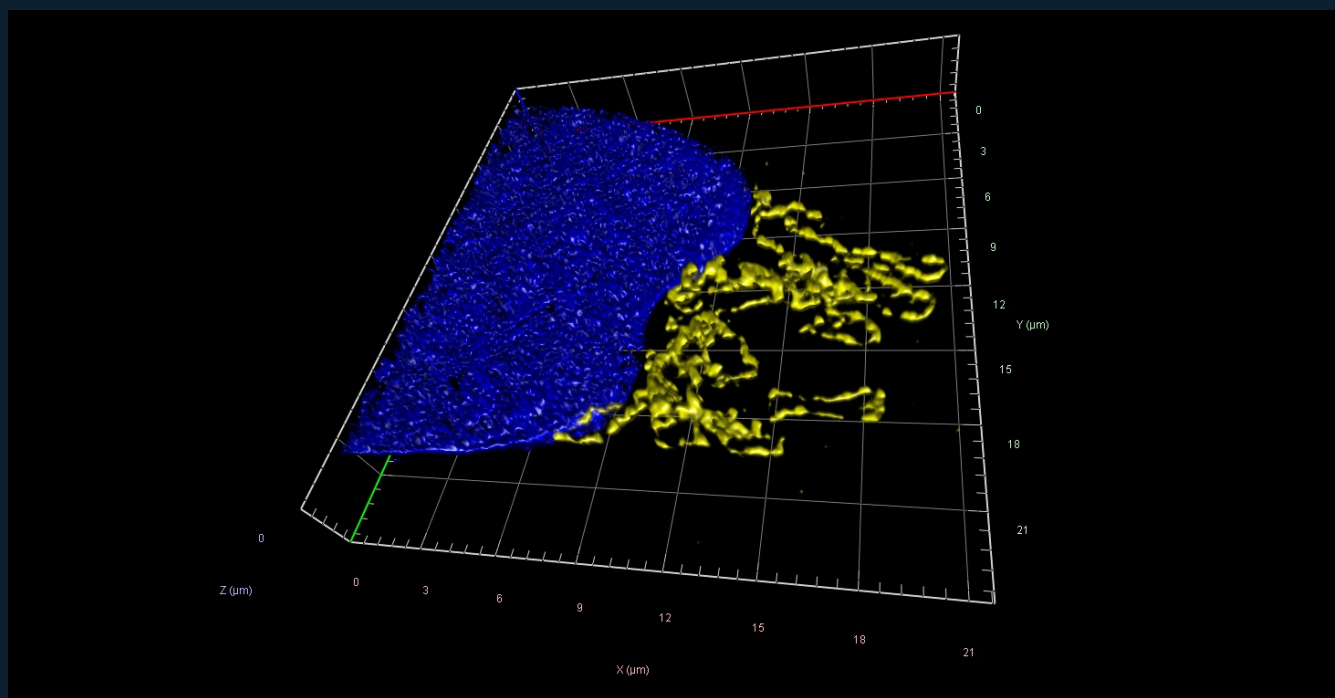

REPRESENTATIVE HIGHEST-AVAILABLE VIEW / SIM²

SA-000024 | PRODUCT-LEVEL EVIDENCE DOSSIER

SOURCE-GROUNDED MICROSCOPY + BIOLOGICAL REVIEW / IDENTIFYN EVIDENCE PRESERVED

CRO CAPABILITY DEMONSTRATED

Identifyn performed the experimental and imaging work documented below. BENNETT converts that operating history into a governed, source-bounded scientific record.

- Source-documented sample preparation and antibody/reagent workflow
- Western blot imaging
- Widefield fluorescence microscopy
- Widefield structured illumination / Apotome optical sectioning
- Confocal microscopy
- Airyscan acquisition and reconstruction workflow
- Structured illumination microscopy
- SIM² reconstruction workflow
- Preserved control experiment
- Z-stack / 3D acquisition
- Multi-channel imaging

SCIENTIFIC QUESTION

What source-bounded experimental and microscopy evidence is preserved for Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit, IgG (H + L), and what is the strongest interpretation allowed by the available modalities and controls?

EVIDENCE AVAILABLE

Western blot -> Widefield -> Widefield SIM — Apotome -> Confocal -> Airyscan -> SIM -> SIM² -> Control

BENNETT ASSESSMENT

The preserved SA-000024 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

BENNETT / MICROSCOPY AND BIOLOGICAL EVIDENCE REVIEW

SOURCE-GROUNDED RETROSPECTIVE REVIEW / PENDING SCIENTIST REVIEW

SOURCE STATE	SOURCE HASH VERIFIED / EVENT RECORDED
BENNETT EVIDENCE REVIEW	COMPLETED / SOURCE-BOUNDED
NATIVE IMAGE REPROCESSING	NOT PERFORMED
LITERATURE SOURCES	INDEPENDENTLY VERIFIED
SCIENTIST REVIEW	PENDING
PUBLICATION	NOT AUTHORIZED

OVERALL CALL / NIR-BOUNDED PARTIAL STEPDOWN

The preserved SA-000024 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

SETTINGS ASSESSMENT

The record preserves 7 ordered evidence tier(s). Structured source metadata values are reported by stage and remain tied to the original method objects.

BIOLOGICAL SUPPORT

The record supports retrospective, source-bounded microscopy review. Target identity, specificity, treatment effect, binding, interaction and mechanism are not inferred from presentation images alone.

CONTROL ASSESSMENT

A true preserved control record is present and appears last.

CLAIM CEILING

This draft supports historical capability documentation, source registration, modality ordering, settings review and bounded interpretation. It does not independently validate antibody specificity, quantify a population response, establish binding, interaction, causation or mechanism, or authorize publication.

PROCESSING DISCLOSURE

Retrospective BENNETT architectural evaluation of preserved Identifyn source evidence. BENNETT did not perform native-file reprocessing, new deterministic measurement or the historical image post-processing represented in this record.

REAGENTS AND ACQUIRED CHANNELS

Preparation reagents are reported separately from channels actually documented in the acquisition settings. Fluorophores mentioned in staining are not automatically treated as acquired channels.

EVIDENCE OBJECT	REAGENTS / FLUOROPHORES PRESENT	CHANNELS ACTUALLY ACQUIRED
Western blot	532	Not deterministically stated in the acquisition block
Widefield	532; 555; 750	2; AF532-49014; 515 - 545; 555 - 595; 528; 553
Widefield SIM — Apotome	405/DAPI; 532; 555	2; AF532-49014; AF405-49027; 515 - 545; 372 - 408; 555 - 595; 417 - 445; 528
Confocal	405/DAPI; 488; 532; 750	2; None; 488nm @0.7%; 405nm @0.5%; 528; 353; 553; 465
Airyscan	405/DAPI; 488; 532	2; None; 488nm @1.50%; 405nm @.08%; 528; 353; 553; 465
SIM	405/DAPI; 488; 532; 647	2; 488nm @5.0%; 642nm @2.0%; 488; 642; 655; 488.000000 (532); 642.000000
SIM ²	405/DAPI; 488; 532	2; 488nm @5.0%; 405nm @20.0%; 488; 405; 655; 515; 488.000000 (532)
Control	405/DAPI; 532; 555	2; AF532-49014; 49 DAPI; 515 - 545; 335 - 383; 555 - 595; 420 - 470; 528

MODALITY ASSESSMENT / PART 1 OF 8

Western blot / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: Azure Sapphire FL. Objective: not explicitly stated. Platform: Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s). Structured source metadata values: 11.

VISIBLE OBSERVATION

A preserved presentation image is available for Western blot. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Western blot, documents platform context as Azure Sapphire FL, and records detected dye/channel descriptors as 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Band position and displayed intensity do not independently establish analytical specificity, target identity, linearity, or treatment effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Western blot evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 2 OF 8

Widefield / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Imager.M1. Objective: Plan-Apochromat 63x/1.40 Oil M27. Platform: Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye: Acquisition: Channels: 2; Scaling (Per Pixel): 0.110µm X 0.110µm; Image size (Pixels): 1216 X 1028; Image Size (Pixels): 1216 X 1028; Bit Depth: 12 Bit. Detector/camera: Imaging Device: Axiocam MR R3. Timing: Exposure Time: 50ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @20.00%. Optical/scan: Effective NA: 1.4. Structured source metadata values: 21.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield, documents platform context as ZEISS Axio Imager.M1, and records detected dye/channel descriptors as 532; 555; 750.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

A presentation image supports gross localization and source-bounded co-occurrence, not direct binding, molecular colocalization, or native quantitative measurement.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 3 OF 8

Widefield SIM — Apotome / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm X 0.100µm; Image size (Pixels): 779 X 875; Image Size (Pixels): 779 X 875; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807 mono. Timing: Exposure Time: 80ms; 150ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @6.0%. Optical/scan: Effective NA: 1.25. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 41.

VISIBLE OBSERVATION

A preserved presentation image is available for Widefield SIM — Apotome. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Widefield SIM — Apotome, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Improved optical sectioning does not by itself establish reconstructed super-resolution or molecular interaction.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Widefield SIM — Apotome evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 4 OF 8

Confocal / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 57 Slices (2.8µm); Channels: 2; Scaling (Per Pixel): 0.071µm X 0.071µm X 0.050µm; Image size (Pixels): 1026 X 1026; Image Size (Pixels): 1026 X 1026; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; GaAsP-PMT; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.51µs; Frame Time: 5.07s. Illumination: Laser Wavelength: 488nm @0.7%; 405nm @0.5%. Optical/scan: Pinhole: 1.00AU / 45µm; 1.00 AU / 37µm; Scan Zoom: 1.4; LSM Scan Speed: 9; Scan Direction: Bidirectional; Averaging: 8; Effective NA: 1.4. Processing: Projection: View Angle 40°, Scale Z 41%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Confocal. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Confocal, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 750.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy and antibody-validation guidance that presentation evidence must remain bounded by modality, controls and source documentation.

INTERPRETIVE LIMIT

Optical sectioning and source-bounded 3D presentation do not establish binding, mechanism, or population-level effect.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Confocal evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 5 OF 8

Airyscan / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS LSM 900. Objective: with a Plan-Apochromat 63X/1.40 oil DIC M27. Platform: Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 81 Slices (4.0µm); Channels: 2; Scaling (Per Pixel): 35.29nm X 35.29nm X 50.00nm; Image size (Pixels): 1122 X 910; Image Size (Pixels): 1122 X 910; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Airyscan; Detector Type: GaAsP-PMT. Timing: Pixel Time: 0.55µs; Frame Time: 4.74s. Illumination: Laser Wavelength: 488nm @1.50%; 405nm @.08%. Optical/scan: Pinhole: 5.00AU / 226µm; 5.00AU / 181µm; Scan Zoom: 1.5; LSM Scan Speed: 7; Scan Direction: Bidirectional; Effective NA: 1.4. Processing: Projection: View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black; AiryScan Mode: Super Resolution 8.7 (3D, Auto); Super Resolution 7.9 (3D, Auto). Structured source metadata values: 46.

VISIBLE OBSERVATION

A preserved presentation image is available for Airyscan. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Airyscan, documents platform context as ZEISS LSM 900, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Displayed feature separation remains reconstruction dependent; no unsupported numerical resolution or molecular architecture claim is permitted.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Airyscan evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 6 OF 8

SIM / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye: Acquisition: Z-Stack: 50 Slices (2.94µm); Channels: 2; Scaling (Per Pixel): 31.30nm X 31.30nm X 60.00nm; Image size (Pixels): 1255 X 1315; Image Size (Pixels): 1255 X 1315; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 100ms; 50ms; Frame Time: 1.0s. Illumination: Laser Wavelength: 488nm @5.0%; 642nm @2.0%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 4; Effective NA: 1.46. Processing: Projection: View Angle 65°, Scale Z 107%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 54.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM, documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 532; 647.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Reconstruction-dependent structural separation requires native reconstruction QC before numerical resolution or molecular dimensions can be claimed.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 7 OF 8

SIM² / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Elyra 7. Objective: Plan-Apochromat 63X/1.46 oil Korr M27. Platform: Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:. Acquisition: Z-Stack: 47Slices (2.76µm); Channels: 2; Scaling (Per Pixel): 15.65nm X 15.65nm X 60.00nm; Image size (Pixels): 1384X 1516; Image Size (Pixels): 1384X 1516; Bit Depth: 16 Bit. Detector/camera: Imaging Device: Edge 4.2 SIM; Detector Type: Camera. Timing: Exposure Time: 100ms; 50ms; Frame Time: 1.0s. Illumination: Laser Wavelength: 488nm @5.0%; 405nm @20.0%. Optical/scan: Scan Zoom: 1.0; Scan Direction: Unidirectional; Averaging: 4; Effective NA: 1.46. Processing: Projection: View Angle 55°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected); Background: Black. Structured source metadata values: 62.

VISIBLE OBSERVATION

A preserved presentation image is available for SIM². BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as SIM², documents platform context as ZEISS Elyra 7, and records detected dye/channel descriptors as 405/DAPI; 488; 532.

LITERATURE CORRELATION

Consistent with peer-reviewed microscopy guidance that advanced-resolution outputs remain dependent on acquisition, reconstruction or localization quality control.

INTERPRETIVE LIMIT

Further reconstruction-dependent segmentation does not establish molecular dimensions, structure identity, or mechanism without independent support.

CLAIM CEILING

This tier supports a source-bounded statement that the historical SIM² evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

MODALITY ASSESSMENT / PART 8 OF 8

Control / SOURCE EVIDENCE PRESERVED

DOCUMENTED SETTINGS

Platform: ZEISS Axio Observer.Z1/7 Apotome. Objective: EC Plan-Neofluar 63X/1.25 oil M27. Platform: Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:. Acquisition: Channels: 2; Scaling (Per Pixel): 0.143µm X 0.143µm; Image size (Pixels): 1284 X 901; Image Size (Pixels): 1284 X 901; Bit Depth: 14 Bit. Detector/camera: Imaging Device: Axiocam 807 mono. Timing: Exposure Time: 200ms; 30ms. Illumination: Light source: X-Cite Xylis Lamp Intensity @25%; X-Cite Xylis Lamp Intensity @10%. Optical/scan: Effective NA: 1.25. Structured source metadata values: 32.

VISIBLE OBSERVATION

A preserved presentation image is available for Control. BENNETT records the visible source presentation without creating a new native measurement.

SOURCE-REPORTED RESULT

The source identifies this object as Control, documents platform context as ZEISS Axio Observer.Z1/7 Apotome, and records detected dye/channel descriptors as 405/DAPI; 532; 555.

LITERATURE CORRELATION

Consistent with peer-reviewed antibody-control guidance that a no-primary or related control bounds background but is not a complete specificity system.

INTERPRETIVE LIMIT

A control bounds the recorded comparison only; it does not establish universal specificity, zero background, zero cross-reactivity, or absence of free dye.

CLAIM CEILING

This tier supports a source-bounded statement that the historical Control evidence and documented settings are preserved for review; stronger biological or quantitative conclusions require compatible native data, controls and scientist review.

CROSS-MODALITY SYNTHESIS

The preserved SA-000024 record contains the applicable NIR sequence: Western blot; Widefield; Widefield SIM — Apotome; Confocal; Airyscan; SIM; SIM². The historical boundary reflects the validated application available for the specific conjugate, deployed platform and workflow, not a universal scientific impossibility.

Different acquisitions are not treated as a quantitative intensity ladder. Any continuity statement remains qualitative and source bounded.

BENNETT LITERATURE CORRELATION — SOURCES INDEPENDENTLY VERIFIED

BENNETT uses independently verified primary scientific sources to test whether the interpretation is scientifically consistent and properly bounded. Literature correlation does not retroactively validate the Identifyn dataset or replace scientist review.

- Antibody and microscopy evidence require orthogonal controls and modality-appropriate interpretation; presentation images alone do not establish analytical specificity.
- Advanced-resolution presentation depends on acquisition, reconstruction or localization quality control before numerical resolution, molecular dimensions or cluster identity can be claimed.

REFERENCES

1. Uhlen M, et al. A proposal for validation of antibodies. Nat Methods. 2016;13:823-827. doi:10.1038/nmeth.3995.
2. Gustafsson MGL. Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy. J Microsc. 2000;198:82-87. doi:10.1046/j.1365-2818.2000.00710.x.
3. Sage D, et al. Quantitative evaluation of software packages for single-molecule localization microscopy. Nat Methods. 2015;12:717-724. doi:10.1038/nmeth.3442.

RELEASE BOUNDARY

This assessment remains a draft pending scientist review. No publication is authorized. Any future native-file analysis must enter BENNETT as a new evidence snapshot with routed engine authority, native-data QA, methods, limitations, and an independently reviewed claim ceiling.

SOURCE RECONCILIATION / QA / PROVENANCE

SOURCE STATUS - PRESERVED / AUTO-RECONCILED

Historical Identifyn source material was subjected to BENNETT retrospective QA. Original source state and correction history remain preserved in provenance.

HIGH CONFIDENCE - AUTO-RECONCILED

Airyscan / Scaling (Per Pixel): 35.29nm X 35.29nm X 50.00nm -> 0.03529 μ m X 0.03529 μ m X 0.05000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1122 X 910; Image Size (Scaled)=39.59 μ m X 32.11 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.01%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM / Scaling (Per Pixel): 31.30nm X 31.30nm X 60.00nm -> 0.03130 μ m X 0.03130 μ m X 0.06000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1255 X 1315; Image Size (Scaled)=39.29 μ m X 41.46 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.72%.

HIGH CONFIDENCE - AUTO-RECONCILED

SIM² / Scaling (Per Pixel): 15.65nm X 15.65nm X 60.00nm -> 0.01565 μ m X 0.01565 μ m X 0.06000 μ m. Both numerical magnitude and unit were normalized, then checked against the reported field dimensions; no unit label was changed alone. Supporting evidence: Image Size (Pixels)=1384 X 1516; Image Size (Scaled)=21.66 μ m X 23.73 μ m; PASS - INTERNALLY CONSISTENT. Documented sampling reproduces the reported physical field dimensions. Maximum field-dimension error 0.02%.

BENNETT EVIDENCE REVIEW COMPLETE

INTERPRETATION TO EVIDENCE

BENNETT has completed the bounded retrospective assessment.
The following pages switch ownership to the preserved Identifyn record.



IDENTIFYN PRESERVED SOURCE EVIDENCE

ORIGINAL RECORD BEGINS

Images, source methods, acquisition settings and provenance follow in the locked Stepdown order. Scientific wording is preserved; missing modalities are not invented. A preserved control follows all available Stepdown tiers as supporting evidence.

PRODUCT-LEVEL STEPDOWN MAP


This dossier preserves the complete available evidence progression for SA-000024. Each tier remains tied to its original image, source method, instrument context and cryptographic identity. A preserved control is retained as supporting evidence and appears after every resolution tier.

CANONICAL STEPDOWN	EVIDENCE STAGE	INSTRUMENT	STATUS
1	Western blot	Azure Sapphire FL	PRESERVED
2	Widefield	ZEISS Axio Imager.M1	PRESERVED
3	Widefield SIM - Apotome	ZEISS Axio Observer.Z1/7 Apotome	PRESERVED
4	Confocal	ZEISS LSM 900	PRESERVED
5	Airyscan	ZEISS LSM 900	PRESERVED
6	SIM	ZEISS Elyra 7	PRESERVED
7	SIM ²	ZEISS Elyra 7	PRESERVED

INTERPRETIVE BOUNDARY

Ordering evidence by resolution tier does not make the modalities interchangeable. Each stage retains its acquisition environment and scientific limitations. This pilot organizes the record; it does not synthesize new biological claims.

ORIGINAL IMAGE EVIDENCE / WESTERN BLOT

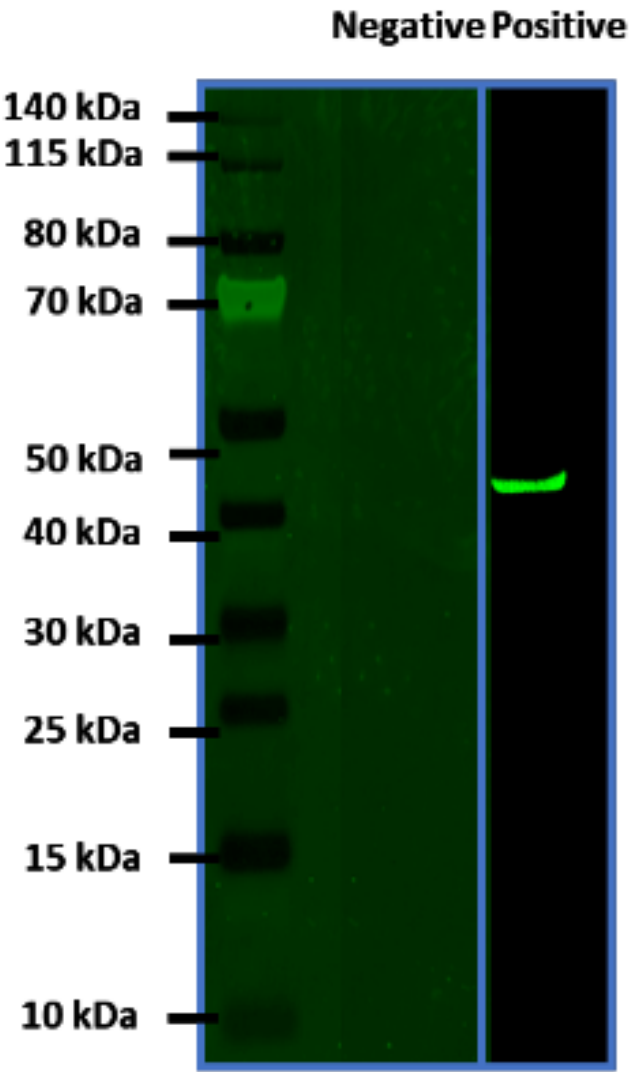

EVIDENCE TIER	Western blot
INSTRUMENT	Azure Sapphire FL
CELL MODEL	HeLa
DETECTED DYES	532
METADATA VALUES	11

STRUCTURED ACQUISITION METADATA / WESTERN BLOT

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Scanner	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Detection mode	Fluorescence
Laser	532nm
Emission Filter	572±28nm
Detector	Avalanche Photodiode
Intensity	4
Pixel size	100 µm
Scan Speed	High
Quality	High
Focus	Membrane (0 µm in Z-axis)
Image	Identifyn® LLC

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

HeLa cells were homogenized using an ultrasonication probe and centrifuged to pellet down the cell debris. The lysate supernatant was mixed with NuPAGE™ LDS Sample Buffer (4X), heated at 80oC for 10 min, and subjected to electrophoresis using NuPAGE™ 4 to 12%, Bis-Tris, 1.0-1.5 mm, Mini Protein Gels in a Mini Gel Tank. The Thermo Fisher's PageRuler™ Prestained Protein Ladder, 10 to 180 kDa, was used as a molecular weight marker.

Post electrophoresis, the proteins were transferred to a polyvinylidene difluoride (PVDF) membrane using iBlot™ 3 Starter Kit, PVDF. The transfer membrane was blocked with Blocker™ Casein for 1 hour and washed 5X with PBS. The blocked membrane was incubated with Actin Recombinant Monoclonal Antibody (JF47-01) (MA5-32530) for 2 hours, followed by 5X PBS washes. Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat anti-Rabbit IgG (H=L) was incubated for 1 hour and washed 5X with PBS. The membrane was dried for 10 minutes at room temperature before imaging.

Scanner = Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)

Scan Settings

Detection mode = Fluorescence

Laser = 532nm

Emission Filter = 572±28nm

Detector = Avalanche Photodiode

Intensity = 4

Pixel size = 100 µm

Scan Speed = High

Quality = High

Focus = Membrane (0 µm in Z-axis)

Post Processing

None

Image: Identifyn® LLC

Copyright 2026 GMD12 LLC.


SOURCE PROVENANCE

Catalog	SA-000024
Stage	Western blot
Instrument	Azure Sapphire FL

SOURCE PROVENANCE

Microscope configuration	Azure Sapphire FL with Customizable Laser and Filter Modules, Enabling UV- NIR and Z-Plane(s)
Structured source metadata values	11
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	42641AEF205ED7CADA0F42E548DEC27F9D39FA164CD99B45E5EAA057F5D56765
Method PDF SHA-256	2046020175A3746EA8E5B0BFD66C548EFB16F940DE10973B1480E124FE166107

ORIGINAL IMAGE EVIDENCE / WIDEFIELD

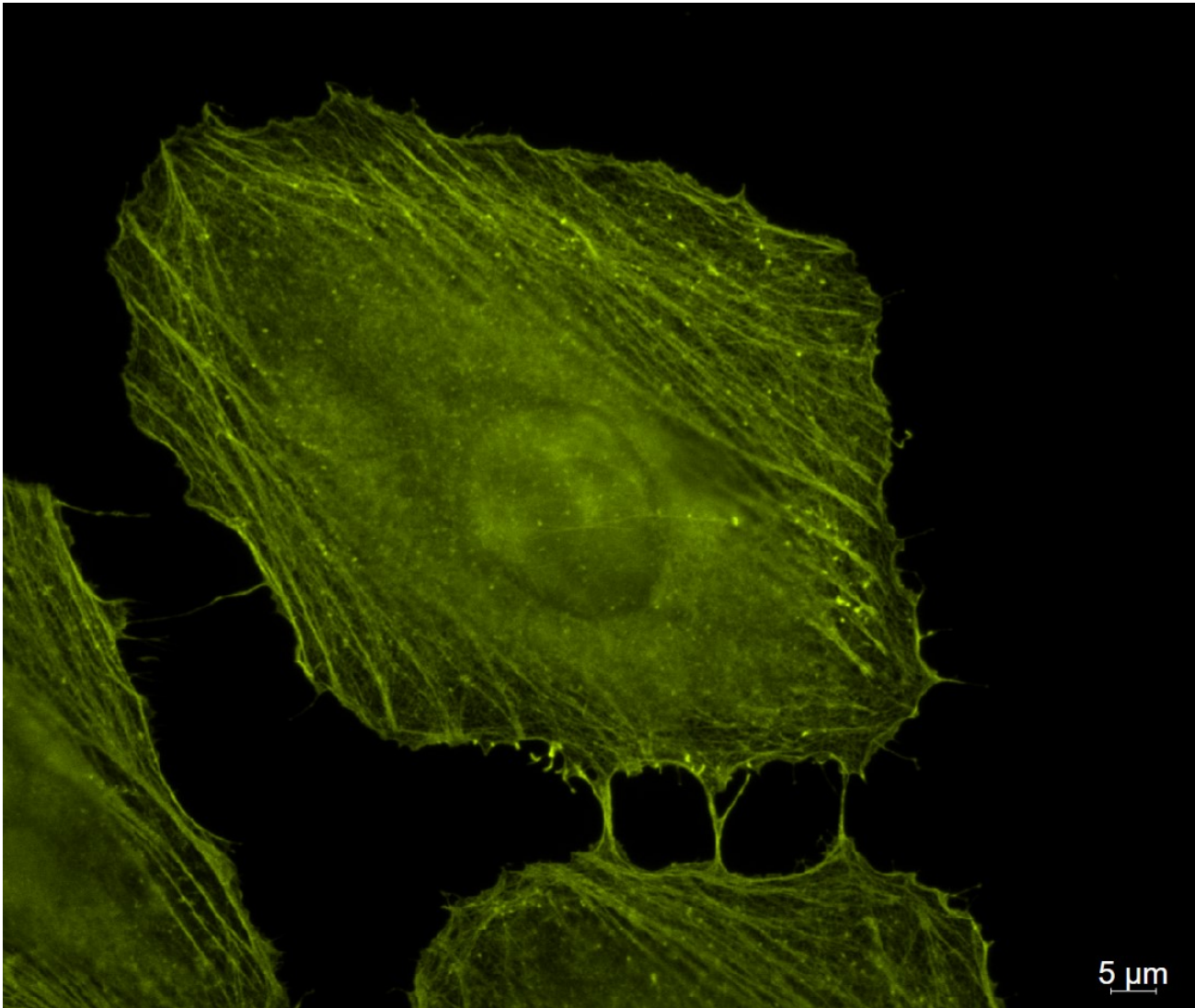

EVIDENCE TIER	Widefield
INSTRUMENT	ZEISS Axio Imager.M1
CELL MODEL	HeLa
DETECTED DYES	532; 555; 750
METADATA VALUES	21

STRUCTURED ACQUISITION METADATA / WIDEFIELD

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	Plan-Apochromat 63x/1.40 Oil M27
Channels	2
Scaling (Per Pixel)	0.110µm X 0.110µm
Image size (Pixels)	1216 X 1028
Image Size (Scaled)	133.18µm X 112.59µm
Bit Depth	12 Bit
Image Center Position	X: 0.00µm, Y: 0.00µm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	AF532-49014
Beam Splitter	550
Filter Excitation Wavelength	515 - 545
Filter Emission Wavelength	555 - 595
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @20.00%
Exposure Time	50ms
Excitation Wavelength	528
Emission Wavelength	553
Effective NA	1.4
Imaging Device	Axiocam MR R3
Camera Adapter	1X
Depth of Focus	0.86µm
Binning Mode	2,2

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The primary antibody, Actin Recombinant Rabbit Monoclonal Antibody (JF47-01), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 750 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using Prolong™ Diamond Antifade Mountant cat#36961.

Microscope

Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.110µm X 0.110µm

Image size (Pixels) = 1216 X 1028

Image Size (Scaled) = 133.18µm X 112.59µm

Bit Depth = 12 Bit

Image Center Position = X: 0.00µm, Y: 0.00µm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @20.00%

Exposure Time = 50ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.4

Imaging Device = Axiocam MR R3

Camera Adapter = 1X

Depth of Focus = 0.86µm

Binning Mode = 2,2

Post Processing

None

Image Correction

None


Image by E.F. Simon

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000024
Stage	Widefield
Instrument	ZEISS Axio Imager.M1
Microscope configuration	Zeiss Axio Imager.M1 Widefield, using a Plan-Apochromat 63x/1.40 Oil M27, Axiocam MR R3, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	21
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	4B816E02A7FD24F987F009AEBB5B8335835F6E93990B1F7681607DA2B595E319
Method PDF SHA-256	42E30E73462DF3922C0B99CD5027F451AC2DF968E1F4FA1F633C77131CB77A9F

ORIGINAL IMAGE EVIDENCE / WIDEFIELD SIM - APOTOME

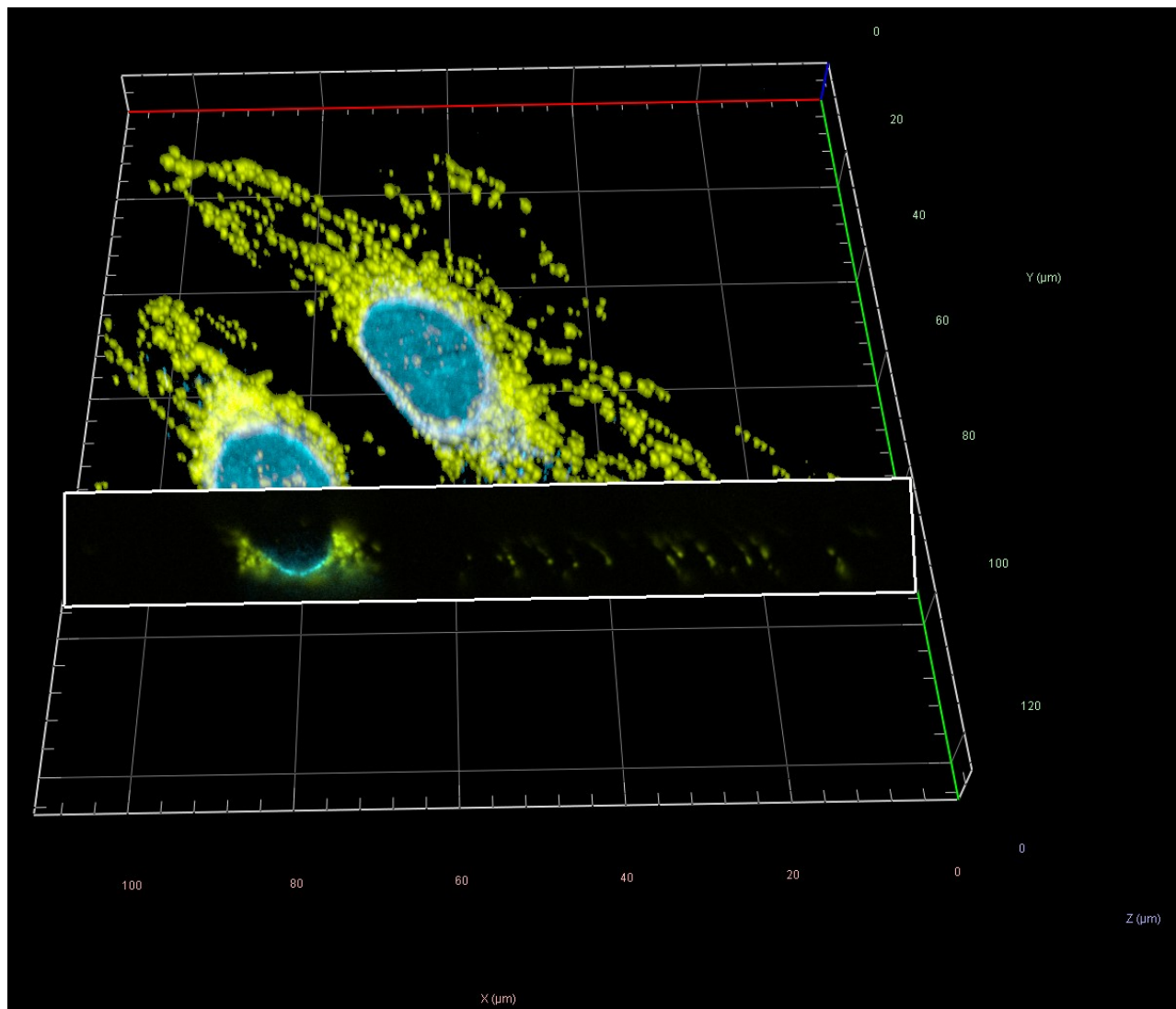

EVIDENCE TIER	Widefield SIM - Apotome
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	41

STRUCTURED ACQUISITION METADATA / WIDEFIELD SIM - APOTOME

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm X 0.100µm
Image Size (Pixels)	779 X 875
Image Size (Scaled)	111.29µm X 125.00µm
Bit Depth	14 Bit
Image Center Position	X: 94.52mm, Y: 46.81mm
Stage Position	: 94.52mm, Y: 46.81mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 AF405-49027
Beam Splitter	550 412
Filter Excitation Wavelength	515 - 545 372 - 408
Filter Emission Wavelength	555 - 595 417 - 445
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @6.0%
Exposure Time	80ms 150ms
Excitation Wavelength	528 401
Emission Wavelength	553 422
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.07µm 0.82µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.
3D Maximum Projection	Precise
Appearance	Threshold 13% (532), Threshold 13% (405), Ramp0% (532), Ramp 0% (DAPI), Maximum 100% all channels
Background	Black
Light	Brightness 100%, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)
Clipping Image Process	Clipping planes are introduced to pronounce regions of special interest.

FIELD	SOURCE-DOCUMENTED VALUE
X - Z	Active was selected, Textured Opaque was chosen to remove the obstruction from the image in front of the cut plane, and a white outline of the clipped perimeter was shown; Clip Back was Selected, and Clip Front was NOT selected.
Position of Clip	38%
Clip X	60°
Clip Y	0°

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Identifyn® SRM Alexa Fluor™ 405 Affinity Purified Donkey Anti-Mouse, F(ab')₂ (H + L)

Cat ID# - SA 000009

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Nuclear Pore Complex Monoclonal Antibody (39C7), was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 405 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Rabbit Recombinant Antibody (MTC02,2860R), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using.

Prolong™ Diamond Antifade Mountant cat#36961.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm X 0.100µm

Image Size (Pixels) = 779 X 875

Image Size (Scaled) = 111.29µm X 125.00µm

Bit Depth = 14 Bit

Image Center Position = X: 94.52mm, Y: 46.81mm

Stage Position = : 94.52mm, Y: 46.81mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

532 -

Reflector = AF532-49014
 Beam Splitter = 550
 Filter Excitation Wavelength = 515 - 545
 Filter Emission Wavelength = 555 - 595
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @6.0%
 Exposure Time = 80ms
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 1.07µm
 Binning Mode = 2,2

405 -

Reflector = AF405-49027
 Beam Splitter = 412
 Filter Excitation Wavelength = 372 - 408
 Filter Emission Wavelength = 417 - 445
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @6.0%
 Exposure Time = 150ms
 Excitation Wavelength = 401
 Emission Wavelength = 422
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.82µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.
 Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 13% (532), Threshold 13% (405), Ramp0% (532), Ramp 0% (DAPI), Maximum 100% all channels
 Background = Black
 Light = Brightness 100%, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)
 Clipping Image Process = Clipping planes are introduced to pronounce regions of special interest.

X - Z = Active was selected, Textured Opaque was chosen to remove the obstruction from the image in front of the cut plane, and a white outline of the clipped perimeter was shown; Clip Back was Selected, and Clip Front was NOT selected.

Position of Clip = 38%

Clip X = 60°

Clip Y = 0°

Image Corrections


Identifyn® SRM Alexa Fluor™ 405 was pseudo-colored in Zen Software from blue, the default color, to light blue to provide visual contrast between the 405/532 channels.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE	
Catalog	SA-000024
Stage	Widefield SIM - Apotome
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	41
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	AE1B31399503F177311B05A7FCBCEAC7A7EFBCB7F6C4489FFBB767D88B2748A8
Method PDF SHA-256	A6A1D61CB1AF3FC0026A12939DCAFF9076FC15129F48CA965D2BE30413287377

ORIGINAL IMAGE EVIDENCE / CONFOCAL

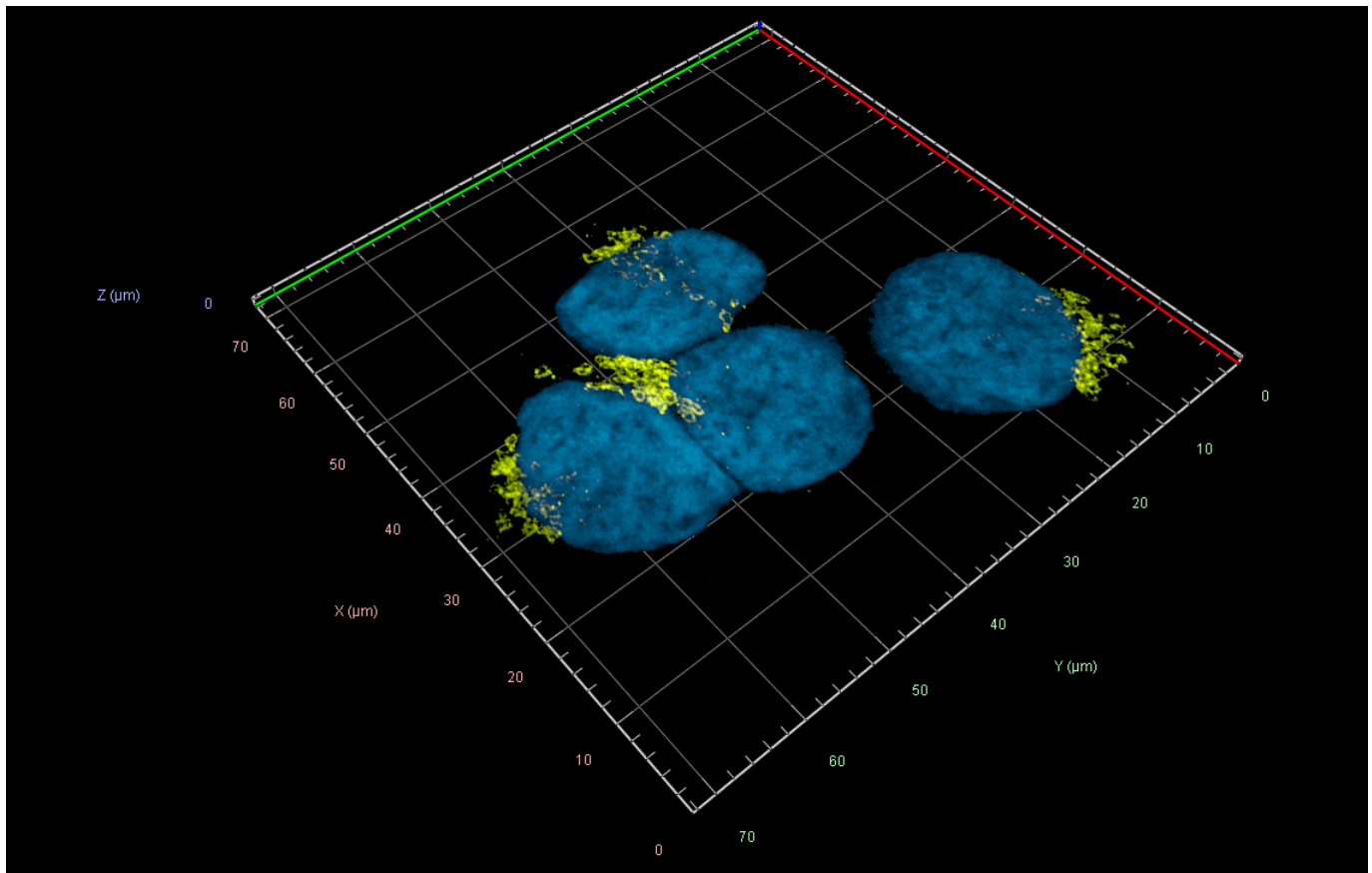

EVIDENCE TIER	Confocal
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 750
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / CONFOCAL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	57 Slices (2.8µm)
Channels	2
Scaling (Per Pixel)	0.071µm X 0.071µm X 0.050µm
Image Size (Pixels)	1026 X 1026
Image Size (Scaled)	72.44µm X 72.44µm
Bit Depth	16 Bit
Image Center Position	X: 76.85mm, Y: 50.47mm
Stage Position	X: 76.85mm, Y: 50.47mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	1.00AU / 45µm 1.00 AU / 37µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @0.7% 405nm @0.5%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.4
Rotation	0°
Sampling	1.00
Pixel Time	0.51µs
Frame Time	5.07s
LSM Scan Speed	9
Scan Direction	Bidirectional
Line Step	1
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Detection Wavelength	500 - 581 410-470
Imaging Device	Airyscan GaAsp-PMT
Detector Type	GaAsp-PMT
Detector Gain	750 V 650 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
Averaging	8
3D Maximum Projection	Precise
Appearance	Threshold 12% (532), Threshold 8% (DAPI), Ramp0% (532), Ramp 0% (DAPI), Maximum 100% all channels
Background	Black
Light	Brightness 87%, Enabled Tone Mapping
Projection	View Angle 40°, Scale Z 41%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack =57 Slices (2.8µm)

Channels = 2

Scaling (Per Pixel) = 0.071µm X 0.071µm X 0.050µm

Image Size (Pixels) = 1026 X 1026

Image Size (Scaled) = 72.44µm X 72.44µm

Bit Depth = 16 Bit

Image Center Position = X: 76.85mm, Y: 50.47mm

Stage Position = X: 76.85mm, Y: 50.47mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 1.00AU / 45µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @0.7%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.4
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.51µs
 Frame Time = 5.07s
 LSM Scan Speed = 9
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 8
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 500 - 581
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 750 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

DAPI -

Reflector - None
 Contrast Method - Fluorescence
 Pinhole = 1.00 AU / 37µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @0.5%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.4
 Rotation = 0°
 Sampling = 1.00
 Pixel Time = 0.51µs
 Frame Time = 5.07s
 LSM Scan Speed = 9
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging = 8
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 410-470
 Imaging Device = GaAsp-PMT
 Detector Type = GaAsp-PMT
 Detector Gain = 650 V
 Detector Offset = 0
 Detector Digital Gain = 1.0

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 12% (532), Threshold 8% (DAPI), Ramp 0% (532), Ramp 0% (DAPI), Maximum 100% all channels
 Background = Black
 Light = Brightness 87%, Enabled Tone Mapping
 Projection = View Angle 40°, Scale Z 41%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE


Catalog

SA-000024

SOURCE PROVENANCE

Stage	Confocal
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	402B6FBDBAF5C064547037C17C0EE1ED40D065101822866DC4280B34780ED7C8
Method PDF SHA-256	A9A2E8B64D294B50BC46FA631237077A0A650A004FDC66185AF45EE03E73EBE7

ORIGINAL IMAGE EVIDENCE / AIRYSCAN

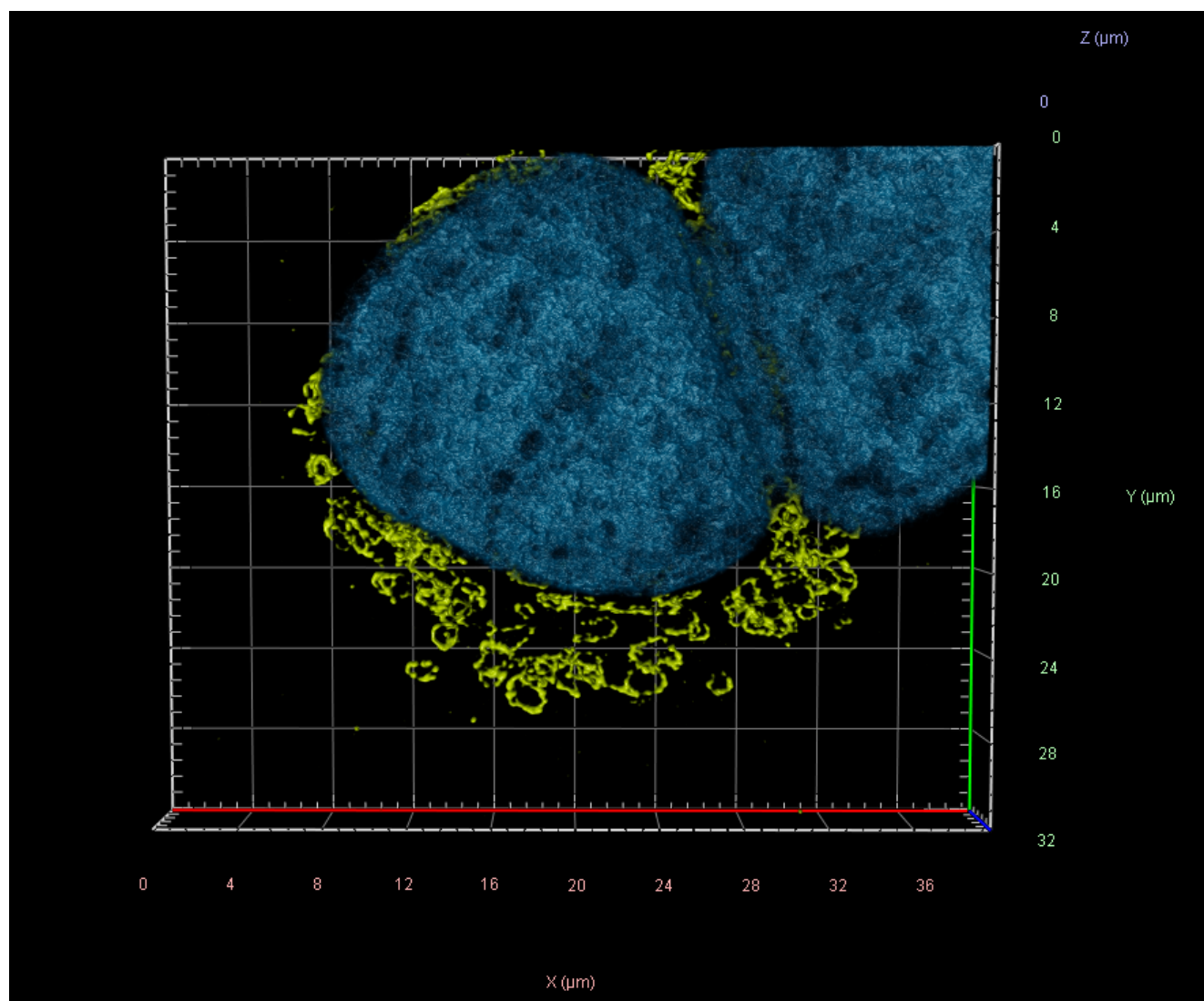

EVIDENCE TIER	Airyscan
INSTRUMENT	ZEISS LSM 900
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	46

STRUCTURED ACQUISITION METADATA / AIRYSCAN

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Objective	with a Plan-Apochromat 63X/1.40 oil DIC M27
Z-Stack	81 Slices (4.0µm)
Channels	2
Scaling (Per Pixel)	0.03529 µm X 0.03529 µm X 0.05000 µm
Image Size (Pixels)	1122 X 910
Image Size (Scaled)	39.59µm X 32.11µm
Bit Depth	16 Bit
Image Center Position	X: 78.10mm, Y: 46.58mm
Stage Position	X: 78.10mm, Y: 46.58mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Reflector	None
Contrast Method	Fluorescence
Pinhole	5.00AU / 226µm 5.00AU / 181µm
LMS MBS	MBS 405 + 488 + 561 + 640(T10/R90)
Laser Wavelength	488nm @1.50% 405nm @.08%
Laser Blanking	Enabled
Scan Mode	Frame
Scan Zoom	1.5
Rotation	0°
Sampling	2.00
Pixel Time	0.55µs
Frame Time	4.74s
LSM Scan Speed	7
Scan Direction	Bidirectional
Line Step	1
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.4
Detection Wavelength	505 - 578 400 - 462
Imaging Device	Airyscan
Detector Type	GaAsp-PMT
Detector Gain	850 V 800 V

FIELD	SOURCE-DOCUMENTED VALUE
Detector Offset	0
Detector Digital Gain	1.0
AiryScan Mode	Super Resolution 8.7 (3D, Auto) Super Resolution 7.9 (3D, Auto)
3D Volume Projection	Precise
Appearance	Threshold 12% (532), Threshold 0% (DAPI), Ramp 88% (532), Ramp 98% (DAPI), Maximum 100% all channels
Background	Black
Light	Brightness100%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods™

Identifyfyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyfyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:

Z Stack - (3D)

Z-Stack = 81 Slices (4.0µm)

Channels = 2

Scaling (Per Pixel) = 35.29nm X 35.29nm X 50.00nm

Image Size (Pixels) = 1122 X 910

Image Size (Scaled) = 39.59µm X 32.11µm

Bit Depth = 16 Bit

Image Center Position = X: 78.10mm, Y: 46.58mm

Stage Position = X: 78.10mm, Y: 46.58mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00AU / 226µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 488nm @1.50%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.5
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.55µs
 Frame Time = 4.74s
 LSM Scan Speed = 7
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 2
 Excitation Wavelength = 528
 Emission Wavelength = 553
 Effective NA = 1.4
 Detection Wavelength = 505 - 578
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 850 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 8.7 (3D, Auto)

DAPI -

Reflector = None
 Contrast Method = Fluorescence
 Pinhole = 5.00AU / 181µm
 LMS MBS = MBS 405 + 488 + 561 + 640(T10/R90)
 Laser Wavelength = 405nm @.08%
 Laser Blanking = Enabled
 Scan Mode = Frame
 Scan Zoom = 1.5
 Rotation = 0°
 Sampling = 2.00
 Pixel Time = 0.55µs
 Frame Time = 4.74s
 LSM Scan Speed = 7
 Scan Direction = Bidirectional
 Line Step = 1
 Averaging - 2
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.4
 Detection Wavelength = 400 - 462
 Imaging Device = Airyscan
 Detector Type = GaAsP-PMT
 Detector Gain = 800 V
 Detector Offset = 0
 Detector Digital Gain = 1.0
 AiryScan Mode = Super Resolution 7.9 (3D, Auto)

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 12% (532), Threshold 0% (DAPI), Ramp 88% (532), Ramp 98% (DAPI), Maximum 100% all channels
 Background = Black
 Light = Brightness 100%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
 Projection = View Angle 35°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by J.K. Bennett

Copyright 2026 GMD12 LLC.
 Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000024
Stage	Airyscan
Instrument	ZEISS LSM 900
Microscope configuration	Zeiss LSM 900 Confocal, Axio Observer.Z1/7, using with a Plan-Apochromat 63X/1.40 oil DIC M27 was used with the following parameters for each dye:
Structured source metadata values	46
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	09293C375E6040E8573199D3A8A246AFC02FA5B05012DA9F695424C5B4979594
Method PDF SHA-256	2F32AA9A325F0E72077434816B90617848C71BC929AC17CCFD75B55142ED5305

ORIGINAL IMAGE EVIDENCE / SIM

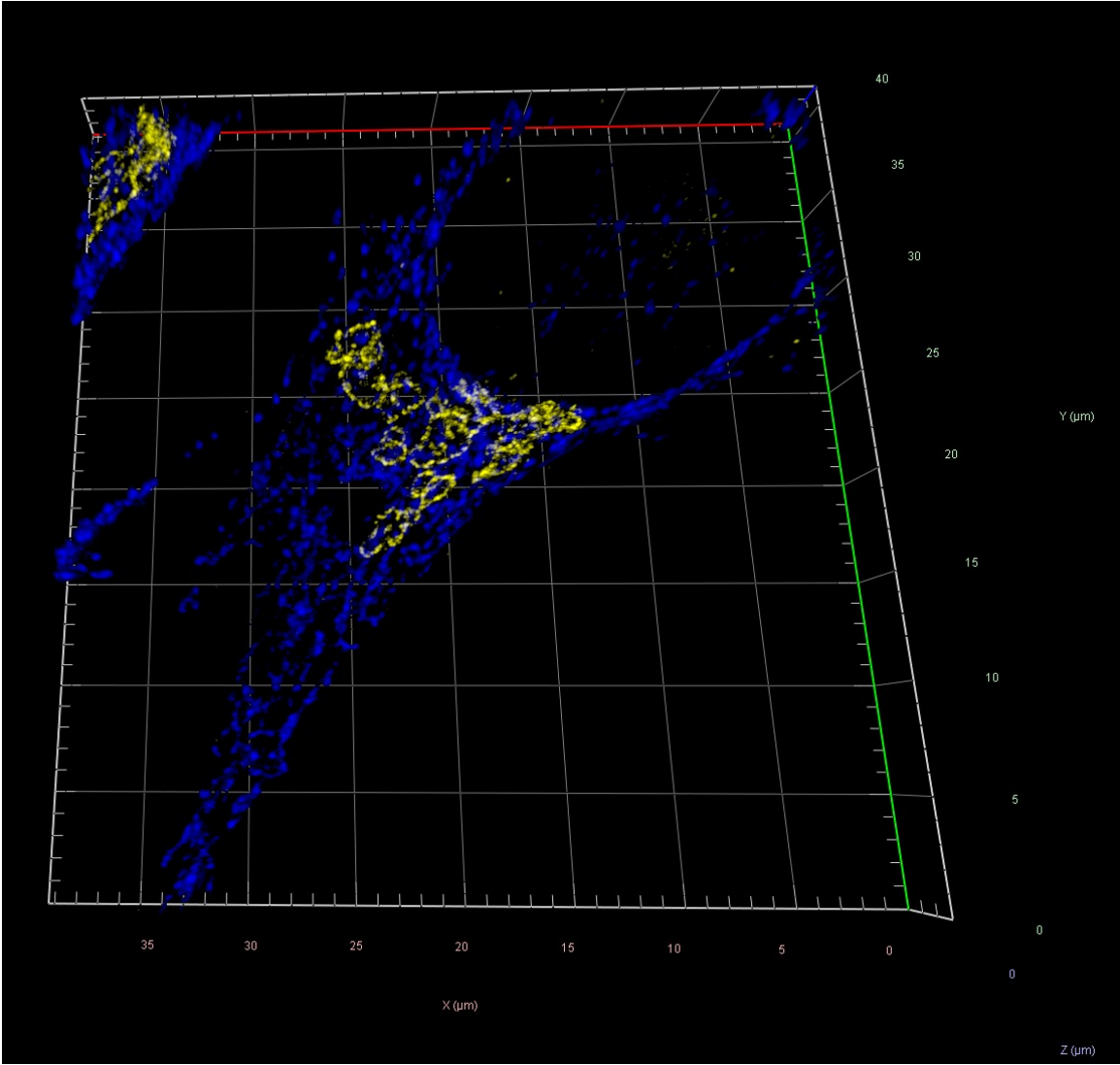

EVIDENCE TIER	SIM
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532; 647
METADATA VALUES	54

STRUCTURED ACQUISITION METADATA / SIM

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 647	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 640
Beam Splitter 532	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Z-Stack	50 Slices (2.94µm)
Channels	2
Scaling (Per Pixel)	0.03130 µm X 0.03130 µm X 0.06000 µm
Image Size (Pixels)	1255 X 1315
Image Size (Scaled)	39.29µm X 41.46µm
Bit Depth	16 Bit
Image Center Position	X: -6.62mm, Y: -4.37mm
Stage Position	X: -6.62mm, Y: -4.37mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	488nm @5.0% 642nm @2.0%
Scan Mode	Other
Scan Zoom	1.0
Rotation	0°
Frame Time	1.0s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	4
Excitation Wavelength	488 642
Emission Wavelength	655
Effective NA	1.46
Imaging Device	Edge 4.2 SIM
Exposure Time	100ms 50ms
Depth of Focus	0.93µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	488.000000 (532) 642.000000
Grating Period (µm)	27.500000 32.000000
Processing	3D
Auto Sharpness	Yes
Sharpness	11.487923 10.494319
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No
3D Maximum Projection	Precise
Appearance	Threshold 4% (647), Threshold 4% (532), Ramp 0% all channels, Maximum 100% all channels
Background	Black
Light	Brightness 114%, Enabled Tone Mapping
Projection	View Angle 65°, Scale Z 107%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Identifyn® SRM Alexa Fluor™ 647 Affinity Purified Donkey Anti-Mouse, IgG (H + L)

Cat ID# - SA 000059

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher, Fish Serum Blocker cat#37527X3. The first primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated overnight at 4°C at a dilution of 1:100. Cells were washed 3X in PBS. Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000.

The second primary antibody followed this, Mitochondria Antibody (MA5-12017), which was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS. Identifyn® 647 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 647 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 640

Beam Splitter 532 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 640

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Z Stack - (3D)

Z-Stack = 50 Slices (2.94µm)

Channels = 2

Scaling (Per Pixel) = 31.30nm X 31.30nm X 60.00nm

Image Size (Pixels) = 1255 X 1315

Image Size (Scaled) = 39.29µm X 41.46µm

Bit Depth = 16 Bit

Image Center Position = X: -6.62mm, Y: -4.37mm

Stage Position = X: -6.62mm, Y: -4.37mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Contrast Method = Fluorescence
 Laser Wavelength = 488nm @5.0%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 488
 Emission Wavelength = 655
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 100ms
 Depth of Focus = 0.93µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM Processing

Channel = 488.000000 (532)

Grating Period (µm): 27.500000

SIM

Processing: 3D
 Auto Sharpness: Yes
 Sharpness: 11.487923
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

647 -

Contrast Method = Fluorescence
 Laser Wavelength = 642nm @2.0%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 642
 Emission Wavelength = 655
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.93µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM Processing

Channel = 642.000000

Grating Period (µm): 32.000000

SIM

Processing: 3D
 Auto Sharpness: Yes
 Sharpness: 10.494319
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

3D Image Processing

3D Maximum Projection = Precise
 Appearance = Threshold 4% (647), Threshold 4% (532), Ramp 0% all channels, Maximum 100% all channels
 Background = Black
 Light = Brightness 114%, Enabled Tone Mapping
 Projection = View Angle 65°, Scale Z 107%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

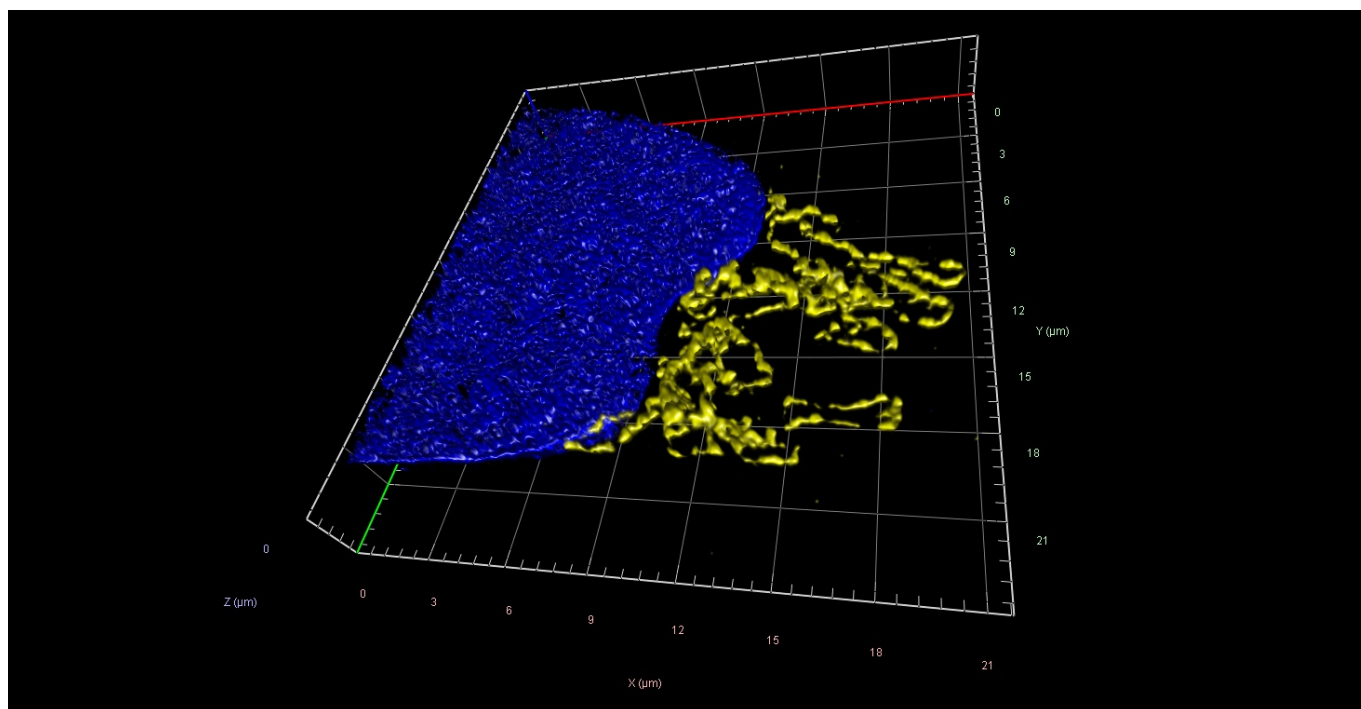
None

Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000024
Stage	SIM
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	54
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	705E51AD1C3A64EB9DB5A0F6B2DBE586D2056061F50D734824743A5B54DF9E8E
Method PDF SHA-256	3B7BA974D3DEB556791C5283AECFDDCD640BE09C662201D0C90049163620B7EC

ORIGINAL IMAGE EVIDENCE / SIM²

EVIDENCE TIER	SIM ²
INSTRUMENT	ZEISS Elyra 7
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 488; 532
METADATA VALUES	62

STRUCTURED ACQUISITION METADATA / SIM²

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Objective	Plan-Apochromat 63X/1.46 oil Korr M27
Filters	-2147483648
Beam Splitter 532	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Dual Camera Beam Splitter	SBS LP 640
Beam Splitter 405	Main Beam Splitter Non-Descanned: LBF 405/488/561/642
Calibration Marker Positions	Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm
Z-Stack	47Slices (2.76µm)
Channels	2
Scaling (Per Pixel)	0.01565 µm X 0.01565 µm X 0.06000 µm
Image Size (Pixels)	1384X 1516
Image Size (Scaled)	21.66µm X 23.73µm
Bit Depth	16 Bit
Image Center Position	X: -6.95mm, Y: -4.06mm
Stage Position	X: -6.95mm, Y: -4.06mm
ROI Center Offset	X: 0.00µm, Y: 0.00µm
Contrast Method	Fluorescence
Laser Wavelength	488nm @5.0% 405nm @20.0%
Scan Mode	Other
Scan Zoom	1.0
Rotation	0°
Frame Time	1.0s
Scan Direction	Unidirectional
Line Step	N/A
Averaging	4
Excitation Wavelength	488 405
Emission Wavelength	655 515
Effective NA	1.46
Imaging Device	Edge 4.2 SIM
Exposure Time	100ms 50ms
Depth of Focus	0.93µm 0.73µm
Binning Mode	1x1
Detector Type	Camera

FIELD	SOURCE-DOCUMENTED VALUE
Depth Offset	0
Detector Digital Gain	0.00
Channel	488.000000 (532) 405.000000
Grating Period (µm)	27.500000 32.000000
Processing	3D
Input SNR	Low Medium
Iterations	5 16
Regularization Weight	0.8000000 0.065000
Processing Sampling	4
Output Sampling	4
Filter	Median
Detrend	No
Sectioning	100.000000
Baseline	Yes
Scale to Raw Image	No
Experimental PSF	No
3D Volume Projection	Precise
Appearance	Threshold 5% (532), Threshold 6% (DAPI) Ramp 95% (532), Ramp 96% (DAPI), Maximum 100% both channels
Background	Black
Light	Brightness 215%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
Projection	View Angle 55°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

SOURCE METHOD

The text below is carried from the preserved Identifyfyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyfyn Standards for Optical Microscopy Methods™

Identifyfyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. The primary antibody, Golgi GOLPH2 Polyclonal Antibody, was incubated for 1 hour at room temperature at a dilution of 1:250. Cells were washed 3X in PBS, and Identifyfyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:

Filters = -2147483648

Beam Splitter 532 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 640

Beam Splitter 405 = Main Beam Splitter Non-Descanned: LBF 405/488/561/642

Dual Camera Beam Splitter = SBS LP 640

Calibration Marker Positions: Marker 1: 0.00µm, 0.00µm, Marker 2: 0.00µm, 0.00µm, Marker 3 0.00µm, 0.00µm

Z Stack - (3D)

Z-Stack = 47Slices (2.76µm)

Channels = 2

Scaling (Per Pixel) = 15.65nm X 15.65nm X 60.00nm

Image Size (Pixels) = 1384X 1516

Image Size (Scaled) = 21.66µm X 23.73µm

Bit Depth = 16 Bit

Image Center Position = X: -6.95mm, Y: -4.06mm

Stage Position =X: -6.95mm, Y: -4.06mm

ROI Center Offset = X: 0.00µm, Y: 0.00µm

532 -

Contrast Method = Fluorescence
 Laser Wavelength = 488nm @5.0%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 488
 Emission Wavelength = 655
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 100ms
 Depth of Focus = 0.93µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 488.000000 (532)

Grating Period (µm): 27.500000

SIM2

Processing: 3D
 Input SNR: Low
 Iterations: 5
 Regularization Weight: 0.8000000
 Processing Sampling: 4
 Output Sampling: 4
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

DAPI -

Contrast Method = Fluorescence
 Laser Wavelength = 405nm @20.0%
 Scan Mode = Other
 Scan Zoom = 1.0
 Rotation = 0°
 Frame Time = 1.0s
 Scan Direction = Unidirectional
 Line Step = N/A
 Averaging = 4
 Excitation Wavelength = 405
 Emission Wavelength = 515
 Effective NA = 1.46
 Imaging Device = Edge 4.2 SIM
 Exposure Time = 50ms
 Depth of Focus = 0.73µm
 Binning Mode = 1x1
 Detector Type = Camera
 Depth Offset = 0
 Detector Digital Gain = 0.00

Post Processing - SIM2 Processing

Channel = 405.000000
Grating Period (µm): 32.000000

SIM2

Processing: 3D
 Input SNR: Medium
 Iterations: 16
 Regularization Weight: 0.065000
 Processing Sampling: 4
 Output Sampling: 4
 Filter: Median
 Detrend: No
 Sectioning: 100.000000
 Baseline: Yes
 Scale to Raw Image: No
 Experimental PSF: No

3D Image Processing

3D Volume Projection = Precise
 Appearance = Threshold 5% (532), Threshold 6% (DAPI) Ramp 95% (532), Ramp 96% (DAPI), Maximum 100% both channels
 Background = Black
 Light = Brightness 215%, Azimuth 80°, Elongation -130°, Enabled Directional Lighting, Enabled Tone Mapping
 Projection = View Angle 55°, Scale Z 100%, Reverse Z-Axis NOT Selected, Stereo Anaglyph NOT Selected, Camera Separation 5% (deselected), Parallax Shift 0% (deselected)

Image Corrections

None


Image by P. Raut

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000024
Stage	SIM²
Instrument	ZEISS Elyra 7
Microscope configuration	Zeiss Elyra 7, Andor1, AxioObserver, using Plan-Apochromat 63X/1.46 oil Korr M27 was used with the following parameters for each dye:
Structured source metadata values	62
Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	DBA59D448438747D13188EB81D1D2B3D2FF8A3C23730511A123C4069E47771E3
Method PDF SHA-256	AABB0B5CC20ECE7F3459115671301C1FCBF7F61AABA61E6B595F097B2985DBED

ORIGINAL IMAGE EVIDENCE / CONTROL

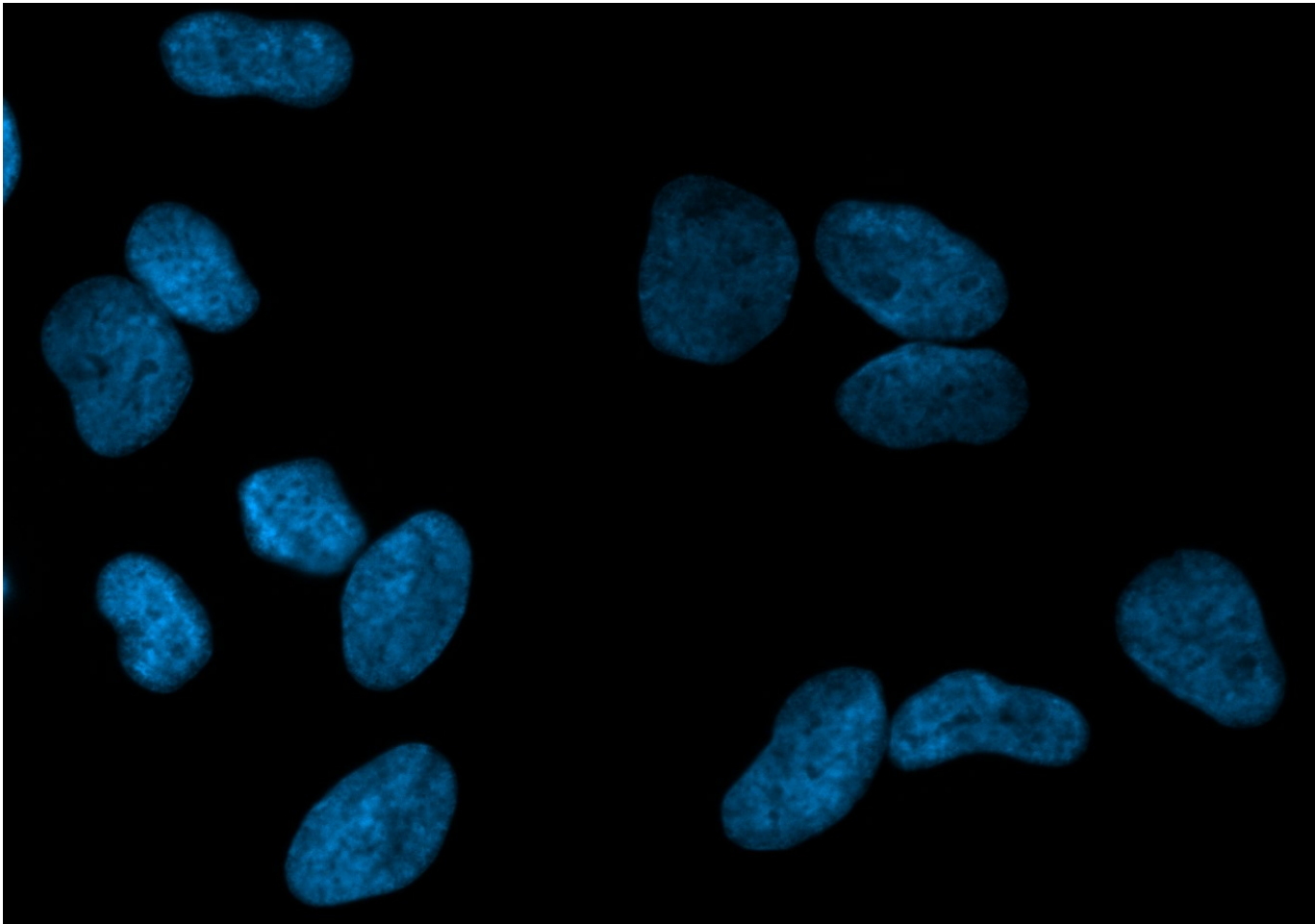

EVIDENCE TIER	Control
INSTRUMENT	ZEISS Axio Observer.Z1/7 Apotome
CELL MODEL	HeLa
DETECTED DYES	405/DAPI; 532; 555
METADATA VALUES	32

STRUCTURED ACQUISITION METADATA / CONTROL

The values below were deterministically extracted from the preserved Identifyn method PDF text. They remain source-documented metadata, not new BENNETT measurements.

FIELD	SOURCE-DOCUMENTED VALUE
Platform / configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Objective	EC Plan-Neofluar 63X/1.25 oil M27
Channels	2
Scaling (Per Pixel)	0.143µm X 0.143µm
Image Size (Pixels)	1284 X 901
Image Size (Scaled)	183.43µm X 128.71µm
Bit Depth	14 Bit
Image Center Position	X: 51.87mm, Y: 50.25mm
Stage Position	X: 51.87mm, Y: 50.25mm
ROI Center Offset	X: 1.14µm, Y: 0.00µm
Reflector	AF532-49014 49 DAPI
Beam Splitter	550 395
Filter Excitation Wavelength	515 - 545 335 - 383
Filter Emission Wavelength	555 - 595 420 - 470
Contrast Method	Fluorescence
Light source	X-Cite Xylis Lamp Intensity @25% X-Cite Xylis Lamp Intensity @10%
Exposure Time	200ms 30ms
Excitation Wavelength	528 353
Emission Wavelength	553 465
Effective NA	1.25
Imaging Device	Axiocam 807 mono
Camera Adapter	1X
Depth of Focus	1.07µm 0.90µm
Binning Mode	2,2
Apotome SIM Image Processing	Display Mode "Optical Sectioning" - enabled correction selected.

SOURCE METHOD

The text below is carried from the preserved Identifyn source record. Typography and pagination have been normalized for this preservation dossier; scientific wording has not been editorially rewritten.

Identifyn Standards for Optical Microscopy Methods™

Identifyn® SRM Alexa Fluor™ 532 Affinity Purified Goat Anti-Rabbit IgG (H + L)

Cat ID# - SA 000024

Hela Cells were grown on acid-etched #1.5H cover glass to ~50% confluency, followed by -20°C methanol fixation and blocked with Thermo Fisher Fish Serum Blocker cat#37527X3. In the absence of a primary antibody, the Identifyn® 532 secondary antibody was incubated at room temperature for 1 hour at a dilution of 1:1000. Cells were washed 5X in PBS and mounted using ProLong™ Diamond Antifade Mountant with DAPI cat#P36962.

Microscope

Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:

Single Plane (2D)

Channels = 2

Scaling (Per Pixel) = 0.143µm X 0.143µm

Image Size (Pixels) = 1284 X 901

Image Size (Scaled) = 183.43µm X 128.71µm

Bit Depth = 14 Bit

Image Center Position = X: 51.87mm, Y: 50.25mm

Stage Position = X: 51.87mm, Y: 50.25mm

ROI Center Offset = X: 1.14µm, Y: 0.00µm

532 -

Reflector = AF532-49014

Beam Splitter = 550

Filter Excitation Wavelength = 515 - 545

Filter Emission Wavelength = 555 - 595

Contrast Method = Fluorescence

Light source = X-Cite Xylis Lamp Intensity @25%

Exposure Time = 200ms

Excitation Wavelength = 528

Emission Wavelength = 553

Effective NA = 1.25

Imaging Device = Axiocam 807 mono

Camera Adapter = 1X

Depth of Focus = 1.07µm

Binning Mode = 2,2

DAPI -

Reflector = 49 DAPI
 Beam Splitter = 395
 Filter Excitation Wavelength = 335 - 383
 Filter Emission Wavelength = 420 - 470
 Contrast Method = Fluorescence
 Light source = X-Cite Xylis Lamp Intensity @10%
 Exposure Time = 30ms
 Excitation Wavelength = 353
 Emission Wavelength = 465
 Effective NA = 1.25
 Imaging Device = Axiocam 807 mono
 Camera Adapter = 1X
 Depth of Focus = 0.90µm
 Binning Mode = 2,2

Post Processing - SIM Processing

Apotome SIM Image Processing = Display Mode "Optical Sectioning" - enabled correction selected.

Normalization, Phase Corrector, Fourier Filter, and Deconvolution were NOT used.

Image Corrections

None

Control Summary -

In the absence of a primary antibody, Identifyn Standards for Optical Microscopy Methods™ were applied using standard staining methodology and image acquisition settings using the Zeiss Apotome SIM. With no image correction, the image demonstrates that the secondary antibody, minus a primary target, had no off-target binding and no significant presence of free dye within the sample.

The control data suggests the secondary antibody is highly specific and free of unconjugated (free) dye post-conjugation.

Image by J.K. Bennett

Copyright 2026 GMD12 LLC.

SOURCE PROVENANCE

Catalog	SA-000024
Stage	Control
Instrument	ZEISS Axio Observer.Z1/7 Apotome
Microscope configuration	Zeiss Axio Observer.Z1/7 Apotome SIM, using an EC Plan-Neofluar 63X/1.25 oil M27, Axiocam 807, and X-Cite Xylis Lamp was used with the following parameters for each dye:
Structured source metadata values	32

SOURCE PROVENANCE

Metadata source	PRESERVED_SOURCE_METHOD_PDF_TEXT
Cell model	HeLa
Image SHA-256	E0D052C44A42DF06F8782A1D85C690E471E0D13EB75DBBD023DA2BA865C7DA67
Method PDF SHA-256	F05D2F2166E5E1916912225DEF9145A7A63BBC7A47462CC34AE6692B44676836